

```
txtFileName.SetFocus
Exit Sub
ElseIf ULStat% <> 0 Then
Stop
```

# UNIVERSAL LIBRARY™

## *Data Acquisition and Control Programming Tools*

```
File Name, SetFocus
t Sub
ULStat% <> 0 Then
p
how many data points were collected
```

# Function Reference

```
Rate.C
LoChan
HiChan
Options
Gain
File
Total
PreTrig
Count
AcqStat
ShowRate
ShowLowChan
ShowHiChan
ShowOptions
ShowGain
ShowFile
Show
ULStat% <> 0 Then
ShowOptions.Caption =
ShowGain.Caption = Format$(Gain, "0")
ShowFile.Caption = FileName
ShowCount.Caption = Format$(CBCount, "0")
how many data points were col
the filename for the collected
data collection options
```



# Universal Library™

## Function Reference



**MEASUREMENT  
COMPUTING™**

Document Revision 7.9, September, 2008  
© Copyright 2008, Measurement Computing Corporation

**Your new Measurement Computing product comes with a fantastic extra —**

## **Management committed to your satisfaction!**

Refer to [www.mccdaq.com/execteam.html](http://www.mccdaq.com/execteam.html) for the names, titles, and contact information of each key executive at Measurement Computing.

Thank you for choosing a Measurement Computing product—and congratulations! You own the finest, and you can now enjoy the protection of the most comprehensive warranties and unmatched phone tech support. It's the embodiment of our mission:

- To provide PC-based data acquisition hardware and software that will save time and save money.

Simple installations minimize the time between setting up your system and actually making measurements. We offer quick and simple access to outstanding live FREE technical support to help integrate MCC products into a DAQ system.

**Limited Lifetime Warranty:** Most MCC products are covered by a limited lifetime warranty against defects in materials or workmanship for the life of the product, to the original purchaser, unless otherwise noted. Any products found to be defective in material or workmanship will be repaired, replaced with same or similar device, or refunded at MCC's discretion. For specific information, please refer to the terms and conditions of sale.

**Harsh Environment Warranty® Program:** Any Measurement Computing product that is damaged due to misuse, or any reason, may be eligible for replacement with the same or similar device for 50% of the current list price. I/O boards face some harsh environments, some harsher than the boards are designed to withstand. Contact MCC to determine your product's eligibility for this program

**30 Day Money-Back Guarantee:** Any Measurement Computing Corporation product may be returned within 30 days of purchase for a full refund of the price paid for the product being returned. If you are not satisfied, or chose the wrong product by mistake, you do not have to keep it.

*These warranties are in lieu of all other warranties, expressed or implied, including any implied warranty of merchantability or fitness for a particular application. The remedies provided herein are the buyer's sole and exclusive remedies. Neither Measurement Computing Corporation, nor its employees shall be liable for any direct or indirect, special, incidental or consequential damage arising from the use of its products, even if Measurement Computing Corporation has been notified in advance of the possibility of such damages.*

### **Licensing Information**

Each original copy of Universal Library is licensed for development use on one CPU at a time. It is theft to make copies of this program for simultaneous program development. If a customer creates an application using the Universal Library, they may distribute the necessary runtime files (Universal Library driver files) with their application royalty free. They may not distribute any files that give their customer the ability to develop applications using the Universal Library.

### **Trademark and Copyright Information**

TracerDAQ, Universal Library, *Harsh* Environment Warranty, Measurement Computing Corporation, and the Measurement Computing logo are either trademarks or registered trademarks of Measurement Computing Corporation.

Windows, Microsoft, and Visual Studio are either trademarks or registered trademarks of Microsoft Corporation

LabVIEW is a trademark of National Instruments.

CompactFlash is a registered trademark of SanDisk Corporation.

XBee and XBee-PRO are trademarks of MaxStream, Inc.

All other trademarks are the property of their respective owners.

Information furnished by Measurement Computing Corporation is believed to be accurate and reliable. However, no responsibility is assumed by Measurement Computing Corporation neither for its use; nor for any infringements of patents or other rights of third parties, which may result from its use. No license is granted by implication or otherwise under any patent or copyrights of Measurement Computing Corporation.

All rights reserved. No part of this publication may be reproduced, stored in a retrieval system, or transmitted, in any form by any means, electronic, mechanical, by photocopying, recording, or otherwise without the prior written permission of Measurement Computing Corporation.

### **Notice**

Measurement Computing Corporation does not authorize any Measurement Computing Corporation product for use in life support systems and/or devices without prior written consent from Measurement Computing Corporation. Life support devices/systems are devices or systems which, a) are intended for surgical implantation into the body, or b) support or sustain life and whose failure to perform can be reasonably expected to result in injury. Measurement Computing Corporation products are not designed with the components required, and are not subject to the testing required to ensure a level of reliability suitable for the treatment and diagnosis of people.

---

# Table of Contents

## Universal Library Functions (32-bit)

<b>1</b>	<b>Overview – Universal Library .....</b>	<b>13</b>
	Introduction .....	13
	Analog I/O functions.....	13
	Configuration functions .....	14
	Counter functions .....	14
	Data Logger functions.....	15
	Digital I/O functions.....	16
	Error handling functions .....	16
	Memory board functions .....	16
	Revision control functions .....	16
	Streamer file functions .....	17
	Synchronous I/O functions .....	17
	Temperature input functions .....	17
	Windows memory management functions .....	17
	Miscellaneous functions.....	18
	Universal Library example programs .....	18
<b>2</b>	<b>Analog I/O Functions .....</b>	<b>28</b>
	Introduction .....	28
	cbAConvertData() .....	29
	cbAConvertPretrigData().....	31
	cbACalibrateData() .....	33
	cbAIn() .....	34
	cbAInScan() .....	35
	cbALoadQueue().....	40
	cbAOut().....	41
	cbAOutScan() .....	42
	cbAPretrig() .....	44
	cbATrig() .....	47
	cbVIn() .....	48
	cbVOut().....	49
<b>3</b>	<b>Configuration Functions.....</b>	<b>50</b>
	Introduction .....	50
	cbGetConfig().....	51
	cbGetConfigString() .....	55
	cbGetSignal() .....	56
	cbSelectSignal() .....	57
	cbSetConfig() .....	60
	cbSetConfigString().....	63
	cbSetTrigger() .....	64
<b>4</b>	<b>Counter Functions .....</b>	<b>67</b>
	Introduction .....	67
	cbC7266Config().....	68
	cbC8254Config() .....	70
	cbC8536Config() .....	71
	cbC9513Config().....	72
	cbC8536Init() .....	74
	cbC9513Init() .....	75
	cbCClear().....	77
	cbCCConfigScan().....	78
	cbCFreqIn() .....	81
	cbCIn() .....	83
	cbCIn32() .....	84

cbCInScan()	85
cbCLoad()	87
cbCLoad32()	89
cbCStatus()	90
cbCStoreOnInt()	91
cbTimerOutStart()	92
cbTimerOutStop()	93
<b>5 Data Logger Functions</b>	<b>94</b>
Introduction	94
cbLogConvertFile()	95
cbLogGetAIChannelCount()	96
cbLogGetAllInfo()	97
cbLogGetCJCInfo()	98
cbLogGetDIOInfo()	99
cbLogGetFileInfo()	100
cbLogGetFileName()	101
cbLogGetPreferences()	102
cbLogGetSampleInfo()	103
cbLogReadAIChannels()	104
cbLogReadJCChannels()	106
cbLogReadDIOChannels()	108
cbLogReadTimeTags()	109
cbLogSetPreferences()	111
<b>6 Digital I/O Functions</b>	<b>112</b>
Introduction	112
cbDBitIn()	113
cbDBitOut()	114
cbDConfigBit()	115
cbDConfigPort()	116
cbDIn()	117
cbDInScan()	118
cbDOut()	120
cbDOutScan()	121
<b>7 Error Handling Functions</b>	<b>123</b>
Introduction	123
cbErrHandling()	124
cbGetErrMsg()	126
<b>8 Memory Board Functions</b>	<b>127</b>
Introduction	127
cbMemRead()	128
cbMemReadPretrig()	129
cbMemReset()	130
cbMemSetDTMode()	131
cbMemWrite()	132
<b>9 Revision Control Functions</b>	<b>133</b>
Introduction	133
cbDeclareRevision()	134
cbGetRevision()	135
<b>10 Streamer File Functions</b>	<b>136</b>
Introduction	136
cbFileAIInScan()	137
cbFileGetInfo()	139
cbFilePretrig()	140
cbFileRead()	142

<b>11 Synchronous I/O Functions.....</b>	<b>143</b>
Introduction .....	143
cbDaqInScan() .....	144
cbDaqOutScan() .....	147
cbDaqSetSetpoints() .....	149
cbDaqSetTrigger() .....	151
<b>12 Temperature Input Functions.....</b>	<b>155</b>
Introduction .....	155
cbTIn() .....	156
cbTInScan() .....	158
<b>13 Windows Memory Management Functions .....</b>	<b>161</b>
Introduction .....	161
cbWinBufAlloc() .....	162
cbWinBufAlloc32() .....	163
cbWinBufFree() .....	164
cbWinArrayToBuf() .....	165
cbWinBufToArray() .....	166
cbWinBufToArray32() .....	167
<b>14 Miscellaneous Functions.....</b>	<b>168</b>
Introduction .....	168
cbDeviceLogin() .....	169
cbDeviceLogout() .....	169
cbDisableEvent().....	170
cbEnableEvent() .....	171
User Callback function .....	173
cbFlashLED() .....	174
cbFromEngUnits() .....	175
cbGetBoardName() .....	176
cbGetStatus() .....	177
cbGetTCValues() .....	179
cbInByte().....	180
cbInWord() .....	181
cbOutByte() .....	182
cbOutWord() .....	183
cbRS485() .....	184
cbStopBackground() .....	185
cbToEngUnits() .....	186

## Universal Library for .NET Classes, Methods, and Properties

<b>15 UL for .NET Class Library Overview .....</b>	<b>188</b>
MccDaq namespace .....	188
MccDaq classes .....	188
MccBoard class .....	188
ErrorInfo class .....	189
MccService class .....	190
GlobalConfig class .....	190
DataLogger class .....	190
Analog I/O methods .....	190
Configuration methods and properties .....	192
Counter methods .....	194
Data Logger methods and properties.....	195
Digital I/O methods .....	196
Error Handling method and properties .....	196
Memory board methods .....	196
Revision control methods and properties .....	197
Streamer file methods .....	197
Synchronous I/O methods .....	197

Temperature input methods .....	197
Windows memory management methods .....	198
Miscellaneous methods, properties, and delegates .....	198
Universal Library for .NET example programs .....	199

## **16 Analog I/O Methods.....208**

Introduction .....	208
AConvertData() .....	209
AConvertPretrigData() .....	210
ACalibrateData() .....	212
Aln() .....	213
AlnScan() .....	214
ALoadQueue() .....	219
AOut() .....	220
AOutScan() .....	221
APretrig().....	223
ATrig().....	226
VIn() .....	227
VOut() .....	228

## **17 Configuration Methods and Properties .....229**

Introduction .....	229
BoardConfig property.....	230
BoardConfig.DACUpdate().....	230
BoardConfig.GetAdRetrigCount().....	230
BoardConfig.GetBaseAdr() .....	231
BoardConfig.GetBoardType().....	231
BoardConfig.GetCiNumDevs() .....	231
BoardConfig.GetClock() .....	232
BoardConfig.GetDACStartup() .....	232
BoardConfig.GetDACUpdateMode() .....	232
BoardConfig.GetDeviceID() .....	233
BoardConfig.GetDeviceNotes() .....	233
BoardConfig.GetDiNumDevs() .....	233
BoardConfig.GetDmaChan() .....	234
BoardConfig.GetDtBoard() .....	234
BoardConfig.GetIntLevel().....	234
BoardConfig.GetNumAdChans() .....	235
BoardConfig.GetNumDaChans().....	235
BoardConfig.GetNumExps().....	235
BoardConfig.GetNumIoPorts() .....	236
BoardConfig.GetPANID() .....	236
BoardConfig.GetRange().....	236
BoardConfig.GetRFChannel() .....	237
BoardConfig.GetRSS().....	237
BoardConfig.GetUsesExps() .....	237
BoardConfig.GetWaitState().....	238
BoardConfig.SetAdRetrigCount() .....	238
BoardConfig.SetBaseAdr() .....	238
BoardConfig.SetClock().....	239
BoardConfig.SetDmaChan() .....	239
BoardConfig.SetDACStartup() .....	239
BoardConfig.SetDACUpdateMode().....	240
BoardConfig.SetDeviceID() .....	241
BoardConfig.SetDeviceNotes() .....	241
BoardConfig.SetIntLevel() .....	241
BoardConfig.SetNumAdChans() .....	242
BoardConfig.SetPANID().....	242
BoardConfig.SetRange() .....	242
BoardConfig.SetRFChannel() .....	243
BoardConfig.SetWaitState() .....	243
BoardNum property.....	244
CtrConfig property.....	245



CtrConfig.GetCtrType() .....	245
DioConfig property .....	246
DioConfig.GetDInMask() .....	246
DioConfig.GetDOutMask() .....	247
DioConfig.GetConfig() .....	247
DioConfig.GetCurVal() .....	247
DioConfig.GetDevType() .....	248
DioConfig.GetNumBits() .....	248
ExpansionConfig property .....	249
ExpansionConfig.GetBoardType() .....	249
ExpansionConfig.GetCjcChan() .....	249
ExpansionConfig.GetMuxAdChan1() .....	250
ExpansionConfig.GetMuxAdChan2() .....	250
ExpansionConfig.GetNumExpChans() .....	250
ExpansionConfig.GetRange1() .....	251
ExpansionConfig.GetRange2() .....	251
ExpansionConfig.GetThermType() .....	251
ExpansionConfig.SetCjcChan() .....	253
ExpansionConfig.SetMuxAdChan1() .....	253
ExpansionConfig.SetMuxAdChan2() .....	254
ExpansionConfig.SetRange1() .....	254
ExpansionConfig.SetRange2() .....	254
ExpansionConfig.SetThermType() .....	255
GetSignal() .....	256
NumBoards property .....	258
NumExpBoards property .....	258
SelectSignal() .....	259
SetTrigger() .....	263
Version property .....	266

## **18 Counter Methods.....267**

Introduction .....	267
C7266Config() .....	268
C8254Config() .....	270
C8536Config() .....	271
C8536Init() .....	273
C9513Config() .....	274
C9513Init() .....	277
CClear() .....	279
CConfigScan() .....	280
CFreqIn() .....	283
CIn() .....	285
CIn32() .....	286
CInScan() .....	287
CLoad() .....	289
CLoad32() .....	291
CStatus() .....	292
CStoreOnInt() .....	293
TimerOutStart() .....	294
TimerOutStop() .....	295

## **19 Data Logger Methods and Property.....296**

Introduction .....	296
FileName property .....	297
ConvertFile() .....	298
GetAIChannelCount() .....	299
GetAllInfo() .....	300
GetCJCInfo() .....	301
GetDIOInfo() .....	302
GetFileInfo() .....	303
GetFileName() .....	304
GetPreferences() .....	305
GetSampleInfo() .....	306

ReadAIChannels()	307
ReadCJCChannels()	309
ReadDIOChannels()	311
ReadTimeTags()	312
SetPreferences()	314
<b>20 Digital I/O Methods</b>	<b>315</b>
Introduction	315
DBitIn()	316
DBitOut()	317
DConfigBit()	318
DConfigPort()	319
DIn()	320
DInScan()	321
DOut()	323
DOutScan()	324
<b>21 Error Handling Methods and Properties</b>	<b>326</b>
Introduction	326
ErrHandling()	327
Message property	329
<b>22 Memory Board Methods</b>	<b>330</b>
MemRead()	331
MemReadPretrig()	332
MemReset()	333
MemSetDTMode()	334
MemWrite()	335
<b>23 Revision Control Methods</b>	<b>336</b>
Introduction	336
DeclareRevision()	337
GetRevision()	338
<b>24 Streamer File Methods</b>	<b>339</b>
Introduction	339
FileAInScan()	340
FileGetInfo()	342
FilePretrig()	343
FileRead()	345
<b>25 Synchronous I/O Methods</b>	<b>347</b>
Introduction	347
DaqInScan()	348
DaqOutScan()	351
DaqSetSetpoints()	353
DaqSetTrigger()	355
<b>26 Temperature Input Methods</b>	<b>359</b>
Introduction	359
TIn()	360
TInScan()	363
<b>27 Windows Memory Management Methods</b>	<b>366</b>
Introduction	366
WinBufAlloc()	367
WinBufAlloc32()	368
WinBufFree()	369
WinArrayToBuf()	370
WinBufToArray()	371
WinBufToArray32()	372

## **28 Miscellaneous Methods, Properties, and Delegates .....374**

Introduction .....	374
BoardName property .....	375
DeviceLogin() .....	376
DeviceLogout() .....	376
HideLoginDialog() .....	377
DisableEvent() .....	378
EnableEvent() .....	379
EngArrayToWinBuf() .....	381
EventCallback delegate .....	383
FlashLED() .....	384
FromEngUnits() .....	385
GetBoardName() .....	386
GetStatus() .....	387
GetTCValues() .....	389
InByte() .....	391
InWord() .....	392
OutByte() .....	393
OutWord() .....	394
RS485() .....	395
StopBackground() .....	396
ToEngUnits() .....	397
WinBufToEngArray() .....	398

## **Appendix**

## **29 Error Codes.....401**

# **Universal Library Functions (32-bit)**

# Overview – Universal Library

## Introduction

This section contains a complete, detailed explanation of all 32-bit Universal Library functions. This chapter briefly explains each function, and provides you with a general idea of the capability of the Universal Library. We highly recommend that you refer to one of the many example programs provided. These programs present a "hands-on" explanation of the various functions, as well as providing you with a starting point from which to write your own programs.

## Analog I/O functions

These functions perform analog input or analog output.

Most PCI boards that support analog input and output scanning allow for simultaneous analog input and output scans. However, for most older boards, analog input scans (`cbAInScan()` and `cbAPretrig()`) cannot operate while an analog output scan (`cbAOutScan()`) is active.

- **cbAIn()** - Takes a single reading from an analog input channel (A/D).
- **cbAInScan()** - Repeatedly scans a range of analog input (A/D) channels. You can specify the channel range, the number of iterations, the sampling rate, and the A/D range. The data that is collected is stored in an array.
- **cbALoadQueue()** - Loads a series of channel/gain pairs into A/D board's queue. These channel/gains are used with all subsequent analog input functions.
- **cbAOut()** - Outputs a single value to an analog output (D/A).
- **cbAOutScan()** - Repeatedly scans a range of analog output (D/A) channels. You can specify the channel range, the number of iterations, and the rate. The data values from consecutive elements of an array are sent to each D/A channel in the scan.
- **cbAPretrig()** - Repeatedly scans a range of analog input (A/D) channels waiting for a trigger signal. When a trigger occurs, it returns the specified number of samples and points before the trigger occurred. You can specify the channel range, the sampling rate, and the A/D range. All of the data that is collected is stored in an array.
- **cbATrig()** - Reads the analog input and waits until it goes above or below a specified threshold. When the trigger condition is met, the current sample is returned.
- **cbAConvertData()** - Converts raw analog data into 12-bit A/D values. Each raw sample from analog input is a 16-bit value. For many 12-bit A/D boards, the raw data is a 16-bit value that contains a 12-bit A/D value and a 4-bit channel tag. This function is not intended for use with 16-bit A/D boards.

This conversion is handled automatically by the `cbAIn()` function. It can also be done automatically by the [cbAInScan\(\)](#) function with the `CONVERTDATA` option. In some cases though, it may be useful or necessary to collect the data and then do the conversion sometime later. The [cbAConvertData\(\)](#) function takes a buffer full of unconverted data and converts it.

- **cbACalibrateData()** - Calibrates analog data. Each raw sample from a board with software calibration factors that must be applied to the sample may be acquired and calibrated, then passed to an array. Alternatively, they can be acquired then passed to the array without calibration. When this second method is used, `cbACalibrateData()` may be used to apply the calibration factors to an array of data after the acquisition is complete. The only case where you would withhold calibration until after the acquisition run was complete is on slower CPUs, or when the processing time is at a premium. Applying calibration factors in real time on a per sample basis does eat up machine cycles.

To disable the automatic calibration so that you may apply the calibration later, specify the

`NOCALIBRATEDATA` option when collecting data with `cbAInScan()`.

- **cbAConvertPretrigData()** - Converts and re-orders pre-trigger data from data plus channel tags to separate the data and channel tags.

For devices with a hardware implementation of pretrigger, when data is collected with the `cbAPretrig()` function the same data conversion needs to be done as is performed by the `cbAConvertData()` function. There is a further complication because `cbAPretrig()` collects analog data into an array. It treats the array like a circular buffer. While it is waiting for the trigger to occur, it fills the array. When it gets to the end it resets to the start and begins again. When the trigger signal occurs it continues collecting data into the circular buffer until the requested number of samples have been collected. When the data acquisition is complete, all of the data is in the array but it is in the wrong order. The first element of the array does not contain the first data point. The data has to be rotated in the correct order.

This conversion can be done automatically by the `cbAPretrig()` function with the `CONVERTDATA` option. In some cases though, it may be useful or necessary to collect the data and then do the conversion sometime later. The `cbAConvertPretrigData()` function takes a buffer full of unconverted data, converts it, and arranges the data in the correct order.

- **cbAVIn()** - Reads an A/D input channel, and returns a voltage value.
- **cbVOut()** - Sets the value of a D/A output.

## Configuration functions

The configuration information for all boards is stored in the configuration file `CB.CFG`. This information is loaded from `CB.CFG` by all programs that use the library. The library includes the following functions to retrieve or change configuration options:

- **cbGetConfig()** - Returns the current value for a specified configuration option.
- **cbGetConfigString()** - Retrieves configuration or device information as a null-terminated string.
- **cbSetConfig()** - Sets the current value for a specified configuration option.
- **cbSetConfigString()** - Sets the configuration or device information as a null-terminated string.
- **cbGetSignal()** - Retrieves the configured auxiliary or DAQ Sync connection and polarity for the specified timing and control signal. This function is intended for advanced users.
- **cbSelectSignal()** - Configures timing and control signals to use specific auxiliary or DAQ Sync connections as a source or destination. This function is intended for advanced users.
- **cbSetTrigger()** - Sets up trigger parameters used with the `EXTTRIGGER` option for `cbAInScan()`.

## Counter functions

Counter functions load, read, and configure counters. There are five types of counter chips used in MCC counter boards: 8254's, 8536's, 7266's, 9513's, and generic event counters. Some of the counter commands only apply to one type of counter.

- **cbC7266Config()** - Selects the operating mode of an LS7266 counter.
- **cbC8254Config()** - Selects the operating mode of the 8254 counter.
- **cbC8536Config()** - Selects the operating mode of the 8536 counter.
- **cbC8536Init()** - Initializes and selects all of the chip-level features for a 8536 counter board. The options set by this command are associated with each counter chip, not the individual counters within it.
- **cbC9513Config()** - Sets the operating mode of the 9513 counter. This function sets all of the programmable options that are associated with a 9513 counter. It is similar in purpose to `cbC8254Config()` except that it is used with a 9513 counter.
- **cbC9513Init()** - Initializes and selects all of the chip level features for a 9513 counter board. The options set by this command are associated with each counter chip, not the individual counters within it.

- **cbCClear()** - Clears a scan counter value (sets it to zero).
- **cbCConfigScan()** - Configures a scan counter channel. `cbCConfigScan()` only works with counter boards that have counter scan capability.
- **cbCFreqIn()** - Measures the frequency of a signal by counting it for a specified period of time (`GateInterval`), and then converting the count to count/sec (Hz). This function only works with 9513 counters.
- **cbCIn()** - Reads a counter's current value as a 16-bit integer. (`cbCIn32()` is the preferred counter read function.)
- **cbCIn32()** - Reads a counter's current value as a 32-bit integer.
- **cbCInScan()** - Scans a range of scan counter channels, and stores the samples in an array.
- **cbCLoad()** - Loads a counter with an initial count value as a 16-bit integer. (`cbCLoad32()` is the preferred counter loading function.)
- **cbCLoad32()** - Loads a counter with a 32-bit integer initial value.
- **cbCStatus()** - Read the counter status of a counter. Returns various bits that indicate the current state of a counter; currently only applies to LS7266 counters.
- **cbCStoreOnInt()** - Installs an interrupt handler that stores the current count whenever an interrupt occurs. This function only works with 9513 counters.
- **cbTimerOutStart()** - Starts a timer square wave output.
- **cbTimerOutStop()** - Stops a timer square wave output.

## Data Logger functions

The data logger functions read and convert binary files logged by MCC hardware equipped with a data logger capability.

- **cbLogConvertFile()** - Converts a binary log file to a comma-separated values (.CSV) text file or another text file format that you specify.
- **cbLogGetAICChannelCount()** - Retrieves the total number of analog input channels logged in a binary file.
- **cbLogGetAllInfo()** - Retrieves the channel number and unit value of each analog input channel logged in a binary file.
- **cbLogGetCJCInfo()** - Retrieves the number of CJC temperature channels logged in a binary file.
- **cbLogGetDIOInfo()** - Retrieves the number of digital I/O channels logged in a binary file.
- **cbLogGetFileInfo()** - Retrieves the version level and byte size of a binary file.
- **cbLogGetFileName()** - Retrieves the name of the  $n^{\text{th}}$  file in the directory containing binary log files.
- **cbLogGetPreferences()** - Retrieves API preference settings for time stamped data, analog data, and CJC temperature data. Returns the default values unless changed using `cbLogSetPreferences()`.
- **cbLogGetSampleInfo()** - Retrieves the sample interval, sample count, and the date and time of the first data point contained in a binary file.
- **cbLogReadAChannels()** - Retrieves analog input data from a binary file, and stores the values in an array.
- **cbLogReadCJCChannels()** - Retrieves CJC temperature data from a binary file, and stores the values in an array.
- **cbLogReadDIOChannels()** - Retrieves digital I/O channel data from a binary file, and stores the values in an array.
- **cbLogReadTimeTags()** - Retrieves date and time values logged in a binary file. This function stores date values in the `DateTags` array, and time values in the `TimeTags` array.

- **cbLogSetPreferences()** - Sets preferences for returned time stamped data, analog temperature data, and CJC temperature data.

## Digital I/O functions

The digital I/O functions perform digital input and output operations on various types of digital I/O ports.

- **cbDBitIn()** - Reads a single bit from a digital input port.
- **cbDBitOut()** - Sets a single bit on a digital output port.
- **cbDConfigBit()** - Configures a specific digital bit as input or output.
- **cbDConfigPort()** - Selects whether a digital port is an input or an output.
- **cbDIn()** - Reads a specified digital input port.
- **cbDInScan()** - Reads a specified number of bytes or words from a digital input port at a specified rate.
- **cbDOut()** - Writes a byte to a digital output port.
- **cbDOutScan()** - Writes a series of bytes or words to a digital output port at a specified rate.

## Error handling functions

All library functions return error codes. The Universal Library includes two functions for handling errors. The different methods built into the functions for handling errors include stopping the program when an error occurs, and printing error messages versus error codes.

- **cbErrHandling()** - Sets the method of reporting and handling errors for all function calls.
- **cbGetErrMsg()** - Returns the error message associated with a specific error code.

## Memory board functions

The memory board functions read and write data to and from a memory board, and also set modes that control memory boards (MEGA-FIFO).

The most common use for the memory boards is to store large amounts of data from an A/D board via a DT-Connect cable between the two boards. To do this, use the `EXTMEMORY` option with [cbAInScan\(\)](#) or [cbAPretrig\(\)](#).

Once the data has been transferred to the memory board, you can use the memory functions to retrieve it.

- **cbMemSetDTMode()** - Sets DT-Connect mode on a memory board. Memory boards have a DT-Connect interface which can be used to transfer data through a cable between two boards rather than through the PC's system memory. The DT-Connect port on the memory board can be configured as either an input (from an A/D) or as an output (to a D/A). This function configures the port to one of these settings.
- **cbMemReset()** - Resets the memory board address. The memory board is organized as a sequential device. When data is transferred to the memory board, it is automatically put in the next address location. This function resets the current address to the location 0.
- **cbMemRead()** - Reads a specified number of points from a memory board starting at a specified address.
- **cbMemWrite()** - Writes a specified number of points to a memory board starting at a specified address.
- **cbMemReadPretrig()** - Reads data collected with `cbAPretrig()`. The `cbAPretrig()` function writes the pre-triggered data to the memory board in a scrambled order. This function unscrambles the data and returns it in the correct order.

## Revision control functions

As new revisions of the library are released, bugs from previous revisions are fixed and occasionally new functions are added. It is the manufacturers goal to preserve existing programs you have written and therefore



to never change the order or number of arguments in a function. However, sometimes it is not possible to achieve this goal.

The revision control function initializes the DLL so that the functions are interpreted according to the format of the revision you wrote and compiled your program in.

- **cbDeclareRevision()** - Declares the revision # of the Universal Library that your program was written with.
- **cbGetRevision()** - Returns the version number of the installed Universal Library.

## Streamer file functions

The streamer file functions explained below create, fill, and read streamer files.

- **cbFileAInScan()** - Transfer analog input data directly to file. Very similar to **cbAInScan()** except that the data is stored in a file instead of an array.
- **cbFilePretrig()** - Pre-triggered analog input to a file. Very similar to **cbAPretrig()** except that the data is stored in a file instead of an array.
- **cbFileGetInfo()** - Reads streamer file information on how much data is in the file, and the conditions under which it was collected (sampling rate, channels, etc.).
- **cbFileRead()** - Reads a selected number of data points from a streamer file into an array.

## Synchronous I/O functions

The synchronous I/O functions synchronously read, set, or write data from analog channels, counter channels, thermocouple channels, and digital ports.

- **cbDaqInScan()** – Scans analog, digital, temperature, and counter inputs synchronously, and stores the values in an array.
- **cbDaqOutScan()** – Outputs values synchronously to analog output channels and digital output ports.
- **cbDaqSetSetpoints()** – Configures up to 16 detection setpoints associated with the input channels within a scan group.
- **cbDaqSetTrigger()** – Selects a trigger source and sets up its parameters. This method starts or stops a synchronous data acquisition operation using [cbDaqInScan\(\)](#) with the `EXTTRIGGER` option.

## Temperature input functions

The temperature sensor functions convert a raw analog input from an EXP or other temperature sensor board to temperature.

- **cbTIn()** - Reads a channel from a digital input board, filters it (if specified), determines the cold junction compensation, linearizes and converts it to temperature.
- **cbTInScan()** - Scans a range of temperature inputs. Reads input temperatures from a range of channels, and returns the temperature values in an array.

## Windows memory management functions

The Windows memory management functions take care of allocating, freeing and copying to/from Windows global memory buffers.

- **cbWinBufAlloc()** - Allocates a Windows memory buffer.
- **cbWinBufAlloc32()** - Allocates a Windows global memory buffer for use with 32-bit scan functions, and returns a memory handle for the buffer.
- **cbWinBufFree()** - Frees a Windows buffer.

- **cbWinArrayToBuf()** - Copies data from an array to a Windows buffer.
- **cbWinBufToArray()** - Copies data from a Windows buffer to an array.
- **cbWinBufToArray32()** - Copies 32-bit data from a Windows global memory buffer into an array. This function is typically used to retrieve data from the buffer after executing an input scan function.

## Miscellaneous functions

These functions do not as a group fit into a single category. They get and set board information, convert units, manage events and background operations, and perform serial communication operations.

- **cbDeviceLogin()** – Opens a device session with a shared device.
- **cbDeviceLogout()** – Releases the device session with a shared device.
- **cbEnableEvent()** – Binds one or more event conditions to a user-defined callback function.
- **User Callback Function** – Defines the prototype for the user function for `cbEnableEvent()`. This defines the format for the user-defined handlers to be called when the events set up using `cbEnableEvent()` occurs.
- **cbDisableEvent()** – Disables one or more events set up with `cbEnableEvent()` and disconnects their user-defined handlers.
- **cbFlashLED()** – Causes the LED on a USB device to flash.
- **cbFromEngUnits()** – Converts a single precision voltage (or current) value in engineering units to an integer D/A count value for output to a D/A.
- **cbGetBoardName()** – Returns the name of a specified board.
- **cbGetStatus()** – Returns the status of a background operation. Once a background operation starts, your program needs to periodically check on its progress. This function returns the current status of the process.
- **cbGetTCValues()** – Converts raw thermocouple data gathered with `cbDagInScan()` to Celsius, Fahrenheit, or Kelvin.
- **cbInByte()** – Reads a byte from a hardware register on a board.
- **cbInWord()** – Reads a word from a hardware register on a board.
- **cbOutByte()** - Writes a byte to a hardware register on a board.
- **cbOutWord()** - Writes a word to a hardware register on a board.
- **cbRS485()** - Sets the transmit and receive buffers on an RS485 port.
- **cbStopBackground()** - Stop a background process. It is sometimes necessary to stop a background process even though the process has been set up to run continuously. This function stops a background process that is running. [cbStopBackground\(\)](#) should be executed after normal termination of all background functions in order to clear variables and flags.
- **cbToEngUnits()** - Converts an integer A/D count value to an equivalent single precision voltage (or current) value.

## Universal Library example programs

Universal Library contains many example programs to aid the user in learning and applying UL functions. We strongly recommend running appropriate example programs before attempting to use the functions.

Table 1 lists Universal Library example programs sorted by the program name. It includes their featured function calls, special aspects, and other function calls included in the program. All example programs include `cbDeclareRevision()` and `cbErrHandling()` functions. Table 2 lists the Universal Library example programs sorted by the function name.

**CWIN sample programs**

The CWIN sample program directory contains programs A101, A102 and A103 only.

Table 1. UL example programs – sorted by program name

Program Name	Featured UL Function Call	Notes	Other UL Function Calls
CInScan01	cbCInScan ()	Scans a range of counter input channels, and writes the data to an array. Board 0 must support counter scans.	cbWinBuffAlloc32 () cbWinBuffToArray32 () cbWinBufFree ()
CInScan02	cbCInScan () cbCConfigScan ()	Scans a counter input channel in decrement mode, and writes the data to an array. Board 0 must support counter scans.	cbWinBuffAlloc32 () cbWinBuffToArray32 () cbWinBufFree ()
DaqInScan01	cbDaqInScan ()	Synchronously scans analog input channels, digital input ports and counter input channels in the foreground. Board 0 must support synchronous input.	cbDConfigPort () cbCConfigScan ()
DaqInScan02	cbDaqInScan ()	Synchronously scans analog input channels, digital input ports, and counter input channels in the background. Board 0 must support synchronous input.	cbDConfigPort () cbCConfigScan () cbGetStatus () cbStopBackground ()
DaqInScan03	cbDaqInScan () cbGetTCValues ()	Synchronously scans analog input channels, digital input ports, and thermocouple input channels in the foreground. Board 0 must support synchronous input.	cbDConfigPort () cbCConfigScan ()
DaqOutScan01	cbDaqOutScan ()	Synchronously writes to an analog output channel and a digital output port in the background. Board 0 must support synchronous output.	cbDConfigPort ()
DaqSetSetpoints01	cbDaqSetSetpoints ()	Configures setpoints, adds the setpoint status to the scanlist, and performs asynchronous reads of the setpoint status. Board 0 must support cbDaqInScan ().	cbDaqInScan () cbDConfigPort () cbGetStatus () cbStopBackground ()
DaqSetTrigger01	cbDaqSetTrigger ()	Configures start and stop triggers. These triggers are used to initiate and terminate A/D conversion using cbDaqInScan () with the EXTTRIGGER option selected.	cbGetStatus () cbStopBackground ()

Program Name	Featured UL Function Call	Notes	Other UL Function Calls
TimerOut01	cbTimerOutputStart() cbTimerOutputStop()	Sends a frequency output to an output timer channel. Board 0 must have a timer output.	
ULAI01	cbAIn()		cbToEngUnits()
ULAI02	cbAInScan()	BACKGROUND mode	cbWinBufToArray() cbWinBufFree() cbWinBufAlloc()
ULAI03	cbAInScan()	BACKGROUND mode	cbGetStatus() cbStopBackground() cbWinBufToArray() cbWinBufFree() cbWinBufAlloc()
ULAI04	cbAConvertData()		cbAInScan() cbGetStatus() cbStopBackground() cbWinBufToArray() cbWinBufFree() cbWinBufAlloc()
ULAI05	cbAInScan()	with manual data conversion	cbGetStatus() cbStopBackground() cbWinBufToArray() cbWinBufFree() cbWinBufAlloc()
ULAI06	cbAInScan()	CONTINUOUS, BACKGROUND mode	cbAConvertData() cbGetStatus() cbStopBackground() cbWinBufToArray() cbWinBufFree() cbWinBufAlloc()
ULAI07	cbATrig()		cbFromEngUnits()
ULAI08	cbAPretrig()		cbWinBufToArray() cbWinBufFree() cbWinBufAlloc()
ULAI09	cbAConvertPretrigData()	BACKGROUND	cbAPretrig() cbGetStatus() cbStopBackground() cbWinBufToArray() cbWinBufFree() cbWinBufAlloc()
ULAI10	cbALoadQueue()		cbAInScan() cbWinBufToArray() cbWinBufFree() cbWinBufAlloc()
ULAI11	cbToEngUnits()		cbAIn()
ULAI12	cbAInScan()	EXTCLOCK mode	cbWinBufToArray() cbWinBufFree() cbWinBufAlloc()
ULAI13	cbAInScan()	Various sampling mode options	cbWinBufToArray() cbWinBufFree() cbWinBufAlloc()

Program Name	Featured UL Function Call	Notes	Other UL Function Calls
ULAI14	cbSetTrigger ()	with EXTTRIGGER selected	cbAInScan () cbFromEngUnits () cbWinBufToArray () cbWinBufFree () cbWinBufAlloc ()
ULAI001	cbAInScan () cbAOutScan ()	Concurrent analog input and analog output scans	cbGetStatus () cbStopBackground () cbWinArraytoBuf () cbWinBufToArray () cbWinBufFree () cbWinBufToAlloc ()
ULAO01	cbAOut ()		cbFromEngUnits ()
ULAO02	cbAOutScan ()		cbWinBufToArray () cbWinBufFree () cbWinBufAlloc ()
ULAO03	cbAOut () cbSetConfig ()	Demonstrates the difference between BIDACUPDATEMODE settings of UPDATEIMMEDIATE and UPDATEONCOMMAND. Board 0 must support BIDACUPDATEMODE settings, such as the PCI-DAC6700 Series boards.	cbFromEngUnits ()
ULCT01	cbC8254Config ()		cbCLoad () cbCIn ()
ULCT02	cbC9513Init () cbC9513Config ()		cbCLoad () cbCIn ()
ULCT03	cbCStoreOnInt ()		cbC9513Init cbC9513Config () cbCLoad () cbCIn ()
ULCT04	cbCFreqIn ()		cbC9513Init ()
ULCT05	cbC8536Init () cbC8536Config ()		cbCLoad () cbCIn ()
ULCT06	cbC7266Config ()		cbCLoad32 () cbCIn32 () cbCStatus ()
ULDI01	cbDIn ()		cbDConfigPort ()
ULDI02	cbDBitIn ()		cbDConfigPort ()
ULDI03	cbDInScan ()		cbDConfigPort () cbGetStatus () cbStopBackground () cbWinBufToArray () cbWinBufFree () cbWinBufAlloc ()
ULDI04	cbDIn ()	using the AUXPORT	
ULDI05	cbDBitIn ()	using the AUXPORT	
ULDI06	cbDConfigBit ()		cbDBitIn ()
ULDO01	cbDOut ()		cbDConfigPort ()
ULDO02	cbDBitOut ()		cbDOut () cbDConfigPort ()

Program Name	Featured UL Function Call	Notes	Other UL Function Calls
ULDO04	cbDOut ()	using the AUXPORT	
ULDO05	cbDBitOut ()	using the AUXPORT	cbDOut ()
ULEV01*	cbEnableEvent ()	using ONEXTERNALINTERRUPT	cbDisableEvent () cbDConfigPort () cbDIn ()
ULEV02*	cbEnableEvent ()	using ON_SCAN_ERROR, ON_DATA_AVAILABLE and ON_END_OF_AI_SCAN	cbAInScan () cbStopBackground () cbToEngUnits () cbWinBufAlloc () cbWinBufFree () cbWinBufToArray ()
ULEV03*	cbEnableEvent ()	using ON_SCAN_ERROR, ON_PRETRIGGER, and ON_END_OF_AI_SCAN	cbAPretrig () cbAConvertPretrigData cbDConfigPort () cbDOut () cbStopBackground () cbToEngUnits () cbWinBufAlloc () cbWinBufFree () cbWinBufToArray ()
ULEV04*	cbEnableEvent ()	using ON_END_OF_AO_SCAN	cbAOutScan () cbDConfigPort () cbDOut () cbFromEngUnits () cbStopBackground () cbWinBufAlloc () cbWinBufFree () cbWinBufToArray ()
ULFI01	cbFileAInScan ()		cbFileGetInfo ()
ULFI02	cbFileRead ()		cbFileAInScan () cbFileGetInfo ()
ULFI03	cbFilePretrig ()		cbFileGetInfo () cbFileRead ()
ULGT01	cbGetErrMsg ()		cbAIn ()
ULGT03	cbGetConfig ()		cbGetBoardName ()
ULGT04	cbGetBoardName ()		cbGetConfig ()
ULLOG01	cbLogGetFileName ()	Retrieves the name of a binary log file.	
ULLOG02	cbLogGetFileInfo () cbLogGetSampleInfo () cbLogGetAIChannelCount () cbLogGetCJCInfo () cbLogGetDIOInfo ()	Retrieves information about the analog data, CJC data, and digital I/O channel data contained in a binary log file.	cbLogGetFileName ()
ULLOG03	cbLogReadAIChannels () cbLogReadCJCChannels () cbLogReadDIOChannels () cbLogReadTimeTags ()	Retrieves the analog input data, CJC temperature data, digital I/O channel data, date values, and time values logged in a binary file, and writes the data to separate arrays.	cbLogGetFileName () cbLogGetSampleInfo () cbLogGetAIChannelCount () cbLogGetCJCInfo () cbLogGetDIOInfo () cbLogSetPreferences ()

Program Name	Featured UL Function Call	Notes	Other UL Function Calls
ULLOG04	cbLogConvertFile()	Converts a binary log file to a comma-separated values (.csv) text file or another text file format that you specify.	cbLogGetSampleInfo()
ULMBDI01	cbDIn()	Reads a digital input port on a MetraBus card	
ULMBDI02	cbDBitIn()	Reads the status of a single digital input bit from a MetraBus card	
ULMBDO01	cbDOut()	Writes a byte to a digital output port on a MetraBus card	
ULMBDO02	cbDBitOut()	Sets the state of a single digital output bit for a MetraBus card	
ULMM01	cbMemReadPretrig()		cbAPretrig()
ULMM02	cbMemRead() cbMemWrite()		
ULMM03	cbAInScan()	With the EXTMEMORY option	cbMemReset() cbMemRead()
ULTI01	cbTIn()		cbGetConfig()
ULTI02	cbTInScan()		cbGetConfig()
VIn01	cbVin()	Reads an A/D input channel.	
VOut01	cbVOut()	Writes to a D/A output channel.	
*Sample programs ULEV01, ULEV02, ULEV03 and ULEV04 are not available for the C Console.			

Table 2. UL Example Programs – sorted by function

UL Function Call	UL Example Program Name	Special Features / Notes
cbAConvertData()	ULAI04 ULA106	
cbAConvertPretrigData()	ULAI09 ULEV03*	
cbACalibrateData()	None	No example programs at this time
cbAIn()	ULAI01 ULGT01 ULAI11	
cbAInScan()	ULAI02 ULA110 ULAI03 ULA112 ULAI04 ULA113 ULAI05 ULA114 ULAI06 ULMM03 ULEV02*	BACKGROUND, BACKGROUND mode with manual data conversion CONTINUOUS BACKGROUND mode EXTCLOCK mode Various sampling mode options
cbALoadQueue()	ULAI10	
cbAOut()	ULAO01 ULAO03	ULAO03 demonstrates the difference between BIDACUPDATEMODE settings of UPDATEIMMEDIATE and UPDATEONCOMMAND. Board 0 must support BIDACUPDATEMODE settings, such as the PCI-DAC6700 Series.
cbAOutScan()	ULAO02 ULAIO01 ULEV04*	
cbAPretrig()	ULAI08 ULEV03* ULAI09 ULMM01 ULFI03	

UL Function Call	UL Example Program Name	Special Features / Notes
cbATrig()	ULAI07 ULMM01	
cbC7266Config()	ULCT06	
cbC8254Config()	ULCT01	
cbC8536Config()	ULCT05	
cbC8536Init()	ULCT05	
cbC9513Config()	ULCT02 ULCT03	
cbC9513Init()	ULCT02 ULCT04 ULCT03	
cbCClear()		
cbCConfigScan()	CInScan02	Demonstrates how to scan a counter input channel in decrement mode, and then write the data to an array. Board 0 must support counter scans.
cbCFreqIn()	ULCT04	
cbCIn()	ULCT01 ULCT05 ULCT02	
cbCIn32()	ULCT06 ULCT07	For ULCT07, board 0 must have an event counter, such as the miniLAB 1008 or USB-1208LS.
cbCInScan()	CInScan01 CInScan02	Demonstrates how to scan one or more counter input channels and then write the data to an array. Board 0 must support counter scans.
cbCLoad()	ULCT01 ULCT03 ULCT02 ULCT05	
cbCLoad32()	ULCT06	
cbCStoreOnInt()	ULCT03	
cbCStatus()	ULCT06	
cbDaqInScan()	DaqInScan01 DaqInScan02 DaqInScan03	Demonstrates how to synchronously scan analog, counter, and thermocouple input channels, and digital input ports. Board 0 must support synchronous input.
cbDaqOutScan()	DaqOutScan01	Demonstrates how to synchronously write to an analog output channel and a digital output port in the background. Board 0 must support synchronous output.
cbDaqSetSetpoints()	DaqSetSetpoints01	Demonstrates how to configure and use setpoints, including how to add the setpoint status to the scanlist and perform asynchronous reads of the setpoint status. Board 0 must support cbDaqInScan().
cbDaqSetTrigger()	DaqSetTrigger01	Demonstrates how to set up start and stop trigger events and display input channel data.
cbDBitIn()	ULDIO2 ULDI06 ULDIO5 ULMBDIO2	
cbDBitOut()	ULDO02 ULDO05 ULMBDO02	
cbDConfigBit()	ULDIO6	
cbDConfigPort()	ULDIO1 ULDO01 ULDIO2 ULDO02 ULDIO3 ULDO05 ULEV01* ULEV04* ULEV03*	



UL Function Call	UL Example Program Name	Special Features / Notes
cbDIn()	ULDI01 ULDI04 ULDI03 ULMBDI01 ULEV04*	
cbDInScan()	ULDI03	
cbDeclareRevision()	All samples	All example programs use this function
cbDOut()	ULDO01 ULDO05 ULDO02 ULMBDO01 ULDO04 ULMBDO02 ULEV03* ULEV04*	
cbDOutScan()	None	No example programs at this time
cbEnableEvent() cbDisableEvent()	ULEV01* ULEV03* ULEV02* ULEV04*	ON_EXTERNAL_INTERRUPT ON_DATA_AVAILABLE ON_PRETRIGGER ON_END_OF_AO_SCAN ON_SCAN_ERROR ON_END_OF_AI_SCAN
cbErrHandling()	All samples	All example programs use this function
cbFileAInScan()	ULFI01 ULFI02	
cbFilePretrig()	ULFI03	
cbFileRead()	ULFI02 ULFI03	
cbFromEngUnits()	ULAI01 ULAO03 ULAI07 ULEV04* ULAI14	
cbGetBoardName()	ULGT03 ULGT04	
cbGetConfig()	ULGT03 ULTI01 ULGT04 ULTI02	
cbGetErrMsg()	ULGT01	
cbGetRevision()	None	No example programs at this time
cbGetStatus()	ULAI03 – ULAI06 ULAI001 ULAI09 ULCT03 ULDI03	
cbGetTCValues()	DaqInScan03	Demonstrates how to retrieve analog, thermocouple, and digital data from a synchronous scan operation. Board 0 must support synchronous output.
cbInByte()	None	No example programs at this time
cbInWord()	None	No example programs at this time
cbLogConvertFile()	ULLOG04	Demonstrates how to convert a binary log file to a .CSV file.
cbLogGetAIChannelCount()	ULLOG02 ULLOG03	Demonstrates how to retrieve the number of analog channels contained in a binary log file.
cbLogGetAIInfo()	ULLOG03	Demonstrates how to retrieve information about the analog input data contained in a binary log file and then write the data to an array.
cbLogGetCJCInfo() cbLogGetDIOInfo()	ULLOG02 ULLOG03	Demonstrates how to retrieve information about the CJC temperature data and digital I/O channel data contained in a binary log file.
cbLogGetFileInfo()	ULLOG02	Demonstrates how to retrieve the version level and byte size of a binary log file.
cbLogGetFileName()	ULLOG01 – ULLOG03	Demonstrates how to retrieve a binary log file name.
cbLogGetPreferences()	None	No example programs at this time.

UL Function Call	UL Example Program Name	Special Features / Notes
cbLogGetSampleInfo()	ULLOG02 ULLOG03 ULLOG04	Demonstrates how to retrieve the sample interval, sample count, and the date and time of the first data point logged in a binary file.
cbLogReadAIChannels() cbLogReadCJCChannels() cbLogReadDIOChannels() cbLogReadTimeTags()	ULLOG03	Demonstrates how to retrieve analog input data, CJC data, DIO channel data, and date/time values contained in a binary log file, and store the data in separate arrays.
cbLogSetPreferences()	ULLOG03	Demonstrates how to store preference settings for time stamped data, analog data, and CJC temperature data.
cbMemRead()	ULMM01 – ULMM03	
cbMemReadPretrig()	ULMM01	
cbMemReset()	ULMM03	
cbMemSetDTMode()	None	No example programs at this time
cbMemWrite()	ULMM02	
cbOutByte()	None	No example programs at this time
cbOutWord()	None	No example programs at this time
cbRS485()	None	No example programs at this time
cbSetConfig()	ULAO03	Demonstrates the difference between BIDACUPDATEMODE settings of UPDATEIMMEDIATE and UPDATEONCOMMAND. Board 0 must support BIDACUPDATEMODE settings, such as the PCI-DAC6700 Series boards.
cbSetTrigger()	ULAI14	
cbStopBackground()	ULAI03 – ULAI06 ULAI09 ULAI001 ULCT03 ULDI03 ULEV02* – ULEV04*	Concurrent cbAInScan() and cbAOutScan()
cbTimerOutStart() cbTimerOutStop()	TimerOut01	Demonstrates how to send a frequency output to a timer output channel. Board 0 must have a timer output.
cbToEngUnits()	ULAI01 ULAI11 ULAI07 ULEV02* ULEV03*	
cbTIn()	ULTI01	
cbTInScan()	ULTI02	
cbVIn()	VIn01	Demonstrates how to read an A/D input channel.
cbVout()	VOut01	Demonstrates how to write to a D/A output channel.
cbWinArrayToBuf()	ULAI01 ULAO02 ULEV04*	
cbWinBufAlloc() cbWinBufFree() cbWinBufToArray()	ULAI01 – ULAI06 ULAI08 – ULAI10 ULAI12 – ULAI14 ULAO02 ULCT03 ULDI03 ULEV02* – ULEV04* (ULEV04:WinBufAlloc and WinBufFree only)	
cbWinArrayToBuff()	ULAI001 ULAO02 ULEV04*	

UL Function Call	UL Example Program Name	Special Features / Notes
cbWinBufAlloc32()	CInScan01 CInScan02	Demonstrates how to allocate a Windows global memory buffer for use with 32-bit scan functions.
cbWinBufToArray32()	CInScan01 CInScan02	Demonstrates how to copy 32-bit data from a Windows memory buffer into an array.
*Sample programs ULEV01, ULEV02, ULEV03 and ULEV04 are not available for the C Console.		

## Analog I/O Functions

### Introduction

The functions explained in this chapter handle analog input, analog output and analog data manipulation. To determine which of these functions are compatible with your hardware, refer to the *Universal Library User's Guide* (available in PDF format on our website at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf)).

Most of the functions in this section provide options that may not be compatible with your hardware. Again, you should refer to the *Universal Library User's Guide* to determine if the options you are considering using with a particular function are compatible with your hardware.

Table 3 below lists the constants you can use in the **Range** argument found in most of the functions explained in this chapter. These values are also used in the [cbALoadQueue\(\)](#) function's `GainArray` argument. Valid ranges for your hardware are listed in the *Universal Library User's Guide*.

Table 3. Range constants

UL settings	Value	UL settings	Value
BIP20VOLTS	±20 volts (V)	UNI10VOLTS	0 to 10 V
BIP10VOLTS	±10 V	UNI5VOLTS	0 to 5 V
BIP5VOLTS	±5 V	UNI2PT5VOLTS	0 to 2.5 V
BIP4VOLTS	±4 V	UNI2VOLTS	0 to 2 V
BIP2PT5VOLTS	±2.5 V	UNI1PT25VOLTS	0 to 1.25 V
BIP2VOLTS	±2 V	UNI1PT67VOLTS	0 to 1.67 V
BIP1PT25VOLTS	±1.25 V	UNI1VOLTS	0 to 1 V
BIP1VOLTS	±1 V	UNIPT5VOLTS	0 to 0.5 V
BIP1PT67VOLTS	±1.67 V	UNIPT25VOLTS	0 to 0.25 V
BIPPT625VOLTS	±0.625 V	UNIPT2VOLTS	0 to 0.2 V
BIPPT5VOLTS	±0.5 V	UNIPT1VOLTS	0 to 0.1 V
BIPPT25VOLTS	±0.25 V	UNIPT01VOLTS	0 to 0.01 V
BIPPT2VOLTS	±0.2 V	UNIPT02VOLTS	0 to 0.02 V
BIPPT1VOLTS	±0.1 V	MA4TO20	4 to 20 milliamperes (mA)
BIPPT05VOLTS	±0.05 V	MA2TO10	2 to 10 mA
BIPPT01VOLTS	±0.01 V	MA1TO5	1 to 5 mA
BIPPT005VOLTS	±0.005 V	MAPT5TO2PT5	0.5 to 2.5 mA
		MA0TO20	0 to 20 mA

## cbAConvertData()

### Changed R3.3 RW

Converts the raw data collected by [cbAInScan\(\)](#) into 12-bit A/D values. The `cbAInScan()` function can return either raw A/D data or converted data, depending on whether or not the `CONVERTDATA` option is used. For many 12-bit A/D boards, the raw data is a 16-bit value that contains a 12-bit A/D value and a 4-bit channel tag (refer to board-specific information in the *Universal Library User's Guide*). The data returned to `ADDData` consists of just the 12-bit A/D value. The data returned to `ChanTags` consists of just the channel numbers.

### Function prototype:

C/C++:	<code>int cbAConvertData(int BoardNum, long NumPoints, unsigned short ADDData[ ], unsigned short ChanTags[ ])</code>
Visual Basic:	<code>Function cbAConvertData(ByVal BoardNum&amp;, ByVal NumPoints&amp;, ADDData%, ChanTags%) As Long</code>
Delphi:	<code>function cbAConvertData(BoardNum:Integer; NumPoints:Longint; var ADDData:Word; var ChanTags:Word):Integer;</code>

### Arguments:

BoardNum	The board number used to collect the data. BoardNum may be 0 to 99. Refers to the number associated with the board used to collect the data when it was installed with the <i>InstaCal</i> configuration program.
NumPoints	Number of samples to convert
ADDData	Pointer or reference to the start of the data array
ChanTags	Pointer or reference to the start of the channel tag array.

### Returns:

[Error code](#) or 0 if no errors.

`ADDData` - converted data.

`ChanTags` - channel tags if available.

When collecting data using [cbAInScan\(\)](#) without the `CONVERTDATA` option, use this function to convert the data after it has been collected. There are cases where the `CONVERTDATA` option is not allowed. For example - if you are using both the `DMAIO` and `BACKGROUND` option with [cbAInScan\(\)](#) on some devices, the `CONVERTDATA` option is not allowed. In those cases this function should be used to convert the data after the data collection is complete.

For some boards, each raw data point consists of a 12-bit A/D value with a 4-bit channel number. This function pulls each data point apart and puts the A/D value into the `ADDData` array and the channel number into the `ChanTags` array.

### Notes:

#### 12-bit A/D boards

- The name of the array must match that used in [cbAInScan\(\)](#) or [cbWinBufToArray\(\)](#).
- Upon returning from `cbAConvertData()`, `ADDData` array contains only 12-bit A/D data.

**16-bit A/D boards**

This function is not for use with 16-bit A/D boards because 16-bit boards do not have channel tags. The argument `BoardNum` was added in revision 3.3 to prevent applying this function to 16-bit data. If you wrote your program for a 12-bit board then later upgrade to a 16-bit board, all you need change is the *InstaCal* configuration file. If this function is called for a 16-bit board, it is simply ignored, and no errors are generated.

## cbAConvertPretrigData()

### Changed R3.3 RW

For products with pretrigger implemented in hardware (most products), this function converts and aligns the raw data collected by [cbAPretrig\(\)](#). The [cbAPretrig\(\)](#) function can return either raw A/D data or converted data, depending on whether or not the CONVERTDATA option was used. The raw data as it is collected is not in the correct order. After the data collection is completed it must be rearranged into the correct order. This function correctly orders the data also, starting with the first pretrigger data point and ending with the last post-trigger point.

Change at revision 3.3 is to support multiple background tasks. It is now possible to run two boards with DMA or REP-INSW background convert-and-transfer features active, therefore, the convert function must know which board the data came from. The data value assigned to BoardNum should be assigned in the header file so it is easy to locate if a change is needed.

### Function prototype:

C/C++:	<code>int cbAConvertPretrigData(int BoardNum, long PretrigCount, long TotalCount, unsigned short ADData[], unsigned short ChanTags[])</code>
Visual Basic:	<code>Function cbAConvertPretrigData(ByVal BoardNum&amp;, ByVal PretrigCount&amp;, ByVal TotalCount&amp;, ADData%, ChanTags%) As Long</code>
Delphi:	<code>function cbAConvertPretrigData(BoardNum:Integer; PretrigCount:Longint; TotalCount:Longint; var ADData:Word; var ChanTags:Word):Integer;</code>

### Arguments:

BoardNum	The board number used to collect the data. BoardNum may be 0 to 99. Refers to the number associated with the board used to collect the data when it was installed with the <i>InstaCal</i> configuration program.
PretrigCount	Number of pre-trigger samples—must match the value returned by the PretrigCount argument in the <a href="#">cbAPretrig()</a> function
TotalCount	Total number of samples that were collected
ADData	Pointer to the data array—must match the array name used in the <a href="#">cbAPretrig()</a> function
ChanTags	Pointer to the channel tag array or a NULL pointer may be passed if using 16-bit boards or if channel tags are not desired (see the note regarding <a href="#">16-bit boards</a> on page 32).

### Returns:

[Error code](#) or 0 if no errors.

ADData - converted data.

When you collect data with [cbAPretrig\(\)](#) and you don't use the CONVERTDATA option, you must use this function to convert the data after it is collected. There are cases where the CONVERTDATA option is not allowed: for example, if you use the BACKGROUND option with [cbAPretrig\(\)](#) on some devices, the CONVERTDATA option is not allowed. In those cases this function should be used to convert the data after the data collection is complete.

**Notes:****12-bit A/D boards:**

- On some 12-bit boards, each raw data point consists of a 12-bit A/D value with a 4-bit channel number. This function pulls each data point apart and puts the A/D value into the ADData and the channel number into the ChanTags array.
- The name of the ADData array must match that used in [cbAInScan\(\)](#) or [cbWinBufToArray\(\)](#).
- Upon returning from `cbAConvertPretrigData()`, ADData array contains only 12-bit A/D data.

**16-bit A/D boards:**

This function is for use with 16-bit A/D boards only insofar as ordering the data. No channel tags are returned.

**Visual Basic programmers:**

After the data is collected with [cbAPretrig\(\)](#) it must be copied to an array with [cbWinBufToArray\(\)](#).

**IMPORTANT**

The entire array must be copied. This array includes the extra 512 samples needed by [cbAPretrig\(\)](#). Example code is given below.

```
Count& = 10000

Dim ADData%(Count& + 512)
Dim ChanTags%(Count& + 512)

cbAPretrig%(BoardNum, LowChan, HighChan, PretrigCount&, Count&...)
cbWinBufToArray%(MemHandle%, ADData%, Count& + 512)
cbAConvertPretrigData%(PretrigCount&, Count&, ADData%, ChanTags%)
```



## cbACalibrateData()

### New R3.3

Calibrates the raw data collected by [cbAInScan\(\)](#) from boards with real time software calibration when the real time calibration has been turned off. The [cbAInScan\(\)](#) function can return either raw A/D data or calibrated data, depending on whether or not the NOCALIBRATEDATA option was used.

#### Function prototype:

C/C++:	<code>int cbACalibrateData(int BoardNum, long NumPoints, int Range, unsigned ADData[ ])</code>
Visual Basic:	<code>Function cbACalibrateData(ByVal BoardNum&amp;, ByVal NumPoints&amp;, ByVal Range&amp;, ADData%) As Long</code>
Delphi:	<code>function cbACalibrateData(BoardNum:Integer; var NumPoints:Longint; Range:Integer; var ADData:Word):Integer;</code>

#### Arguments:

BoardNum	May be 0 to 99. Refers to the number associated with the board when it was installed using <i>InstaCal</i> .
NumPoints	Number of samples to convert
Range	The programmable gain/range used when the data was collected. See Table 3 on page 28 for valid values.
ADData	Pointer to data array.

#### Returns:

[Error code](#) or 0 if no errors.

ADData - converted data.

#### Notes:

When collecting data using [cbAInScan\(\)](#) with the NOCALIBRATEDATA option, use this function to calibrate the data once collected.

- The name of the array must match that used in [cbAInScan\(\)](#) or [cbWinBufToArray\(\)](#).
- Applying software calibration factors in real time on a per sample basis eats up machine cycles. If your CPU is slow, or if processing time is at a premium, do not calibrate until the acquisition run finishes. Turn off real time software calibration to save CPU time during high speed acquisitions by using the NOCALIBRATEDATA option to turn off real-time software calibration. After the acquisition is run, calibrate the data with [cbACalibrateData\(\)](#).

## cbAIn()

Reads an A/D input channel. This function reads the specified A/D channel from the specified board. If the specified A/D board has programmable gain then it sets the gain to the specified range. The raw A/D value is converted to an A/D value and returned to `DataValue`.

### Function prototype:

C/C++:	<code>int cbAIn(int BoardNum, int Channel, int Range, unsigned short *DataValue);</code>
Visual Basic:	<code>Function cbAIn(ByVal BoardNum&amp;, ByVal Channel&amp;, ByVal Range&amp;, DataValue%) As Long</code>
Delphi:	<code>function cbAIn(BoardNum:Integer; Channel:Integer; Range:Integer; var DataValue:Word):Integer;</code>

### Arguments:

BoardNum	The board number used to collect the data. BoardNum may be 0 to 99. Refers to the number associated with the board used to collect the data when it was installed with the <i>InstaCal</i> configuration program. The specified board must have an A/D.
Channel	A/D channel number. The maximum allowable channel depends on which type of A/D board is being used. For boards with both single ended and differential inputs, the maximum allowable channel number also depends on how the board is configured. For example, a CIO-DAS1600 has 8 channels for differential, 16 for single ended. Expansion boards are also supported by this function, so this argument can contain values up to 272. Refer to board-specific information in the <i>Universal Library User's Guide</i> for EXP boards if you are using an expansion board.
Range	A/D range code. If the selected A/D board does not have a programmable gain feature, this argument is ignored. If the A/D board does have programmable gain, set the Range argument to the desired A/D range. See Table 3 on page 28 for valid values.
DataValue	Pointer or reference to the data value.

### Returns:

[Error code](#) or 0 if no errors.

`DataValue` - Returns the value of the A/D sample.

## cbAInScan()

### Changed R3.3 ID

Scans a range of A/D channels and stores the samples in an array. `cbAInScan()` reads the specified number of A/D samples at the specified sampling rate from the specified range of A/D channels from the specified board. If the A/D board has programmable gain, then it sets the gain to the specified range. The collected data is returned to the data array.

**Changes: Revision 3.3 added a "no real time calibration" option.**

### Function prototype:

C/C++:	<code>int cbAInScan(int BoardNum, int LowChan, int HighChan, long Count, long *Rate, int Range, int MemHandle, int Options)</code>
Visual Basic:	<code>Function cbAInScan(ByVal BoardNum&amp;, ByVal LowChan&amp;, ByVal HighChan&amp;, ByVal Count&amp;, Rate&amp;, ByVal Range&amp;, ByVal MemHandle&amp;, ByVal Options&amp;) As Long</code>
Delphi:	<code>function cbAInScan(BoardNum:Integer; LowChan:Integer; HighChan:Integer; Count:Longint; var Rate:Longint; Range:Integer; MemHandle:Integer; Options:Integer) : Integer;</code>

### Arguments:

BoardNum	The board number used to collect the data. BoardNum may be 0 to 99. Refers to the number associated with the board used to collect the data when it was installed with the <i>InstaCal</i> configuration program. The specified board must have an A/D.
LowChan	The first A/D channel of scan. When <code>cbALoadQueue()</code> is used, the channel count is determined by the total number of entries in the channel gain queue, and LowChan is ignored.
HighChan	The last A/D channel of scan. When <code>cbALoadQueue()</code> is used, the channel count is determined by the total number of entries in the channel gain queue, and HighChan is ignored.  <b>Low / High Channel number:</b> The maximum allowable channel depends on which type of A/D board is being used. For boards that have both single ended and differential inputs the maximum allowable channel number also depends on how the board is configured. For example, a CIO-DAS1600 has 8 channels for differential, 16 for single ended.
Count	The number of A/D samples to collect. Specifies the total number of A/D samples that will be collected. If more than one channel is being sampled, the number of samples collected per channel is equal to $\text{Count} / (\text{HighChan} - \text{LowChan} + 1)$ .
Rate	The rate at which samples are acquired, in samples per second per channel. For example, if you sample four channels, 0-3, at a rate of 10,000 scans per second (10 kHz), the resulting A/D converter rate is 40 kHz: four channels at 10,000 samples per channel per second. This is different from some software where you specify the total A/D chip rate. In those systems, the per channel rate is equal to the A/D rate divided by the number of channels in a scan.  The channel count is determined by the LowChan and HighChan parameters. $\text{Channel Count} = (\text{HighChan} - \text{LowChan} + 1)$ .  When <code>cbALoadQueue</code> is used, the channel count is determined by the total number of entries in the channel gain queue. LowChan and HighChan are ignored.  Rate also returns the value of the actual rate set, which may be different from the requested rate because of pacer limitations.

Range	A/D range code. If the selected A/D board does not have a programmable range feature, this argument is ignored. Otherwise, set the <code>Range</code> argument to any range that is supported by the selected A/D board. See Table 3 on page 28 for valid values.
MemHandle	Handle for Windows buffer to store data in (Windows). This buffer must have been previously allocated with the <a href="#">cbWinBufAlloc()</a> function.
Options	Bit fields that control various options. This field may contain any combination of non-contradictory choices from the values listed in the "Options argument values" section below.

**Returns:**

[Error code](#) or 0 if no errors.

`Rate` - Actual sampling rate used.

`MemHandle` - Collected A/D data returned via the Windows buffer.

**Options argument values:**

**Transfer method options:** The following four options determine how data is transferred from the board to PC memory. If none of these options are specified (recommended), the optimum sampling mode is automatically chosen based on board type and sampling speed.

SINGLEIO	A/D transfers to memory are initiated by an interrupt. One interrupt per conversion. Rates attainable using <code>SINGLEIO</code> are PC-dependent and generally less than 10 kHz. Use the default method unless you have a reason to select a specific transfer method.
DMAIO	A/D transfers are initiated by a DMA request.
BLOCKIO	A/D transfers are handled in blocks (by <code>REP-INSW</code> for example).  <b>BLOCKIO is not recommended for slow acquisition rates:</b> If the rate of acquisition is very slow (for example less than 200 Hz) <code>BLOCKIO</code> is probably not the best choice for transfer mode. The reason for this is that status for the operation is not available until one packet of data has been collected (typically 512 samples). The implication is that if acquiring 100 samples at 100 Hz using <code>BLOCKIO</code> , the operation will not complete until 5.12 seconds has elapsed.
BURSTIO	Allows higher sampling rates for sample counts up to full FIFO. Data is collected into the local FIFO. Data transfers to the PC are held off until after the scan is complete. For <code>BACKGROUND</code> scans, the count and index returned by <a href="#">cbGetStatus()</a> remain 0 and the status equals <code>RUNNING</code> until the scan finishes. When the scan is complete and the data is retrieved, the count and index are updated and the status equals <code>IDLE</code> .  <code>BURSTIO</code> is the default mode for non-Continuous fast scans (aggregate sample rates above 1000 Hz) with sample counts up to full FIFO. To avoid the <code>BURSTIO</code> default, specify <code>BLOCKIO</code> .
BURSTMODE	Enables burst mode sampling. Scans from <code>LowChan</code> to <code>HighChan</code> are clocked at the maximum A/D rate in order to minimize channel to channel skew. Scans are initiated at the rate specified by <code>Rate</code> .  <code>BURSTMODE</code> is not recommended for use with the <code>SINGLEIO</code> option. If this combination is used, the <code>Count</code> value should be set as low as possible, preferably to the number of channels in the scan. Otherwise, overruns may occur.

CONVERTDATA	<p>If the CONVERTDATA option is used for 12-bit boards then the data that is returned to the buffer will automatically be converted to 12-bit A/D values. If CONVERTDATA is not used then the data from 12-bit A/D boards will be return unmodified (which, for some boards is 16-bit values that contain both a 12-bit A/D value and a 4 bit channel number). After the data collection is complete you can call <a href="#">cbAConvertData()</a> to convert the data after the fact. On some devices, CONVERTDATA may not be specified if you are using the BACKGROUND option and DMA transfers. This option is ignored for the 16-bit boards.</p>
BACKGROUND	<p>If the BACKGROUND option is not used then the cbAInScan() function will not return to your program until all of the requested data has been collected and returned to the buffer. When the BACKGROUND option is used, control will return immediately to the next line in your program and the data collection from the A/D into the buffer will continue in the background. Use <a href="#">cbGetStatus()</a> with AIFUNCTION to check on the status of the background operation. Alternatively, some boards support <a href="#">cbEnableEvent()</a> for event notification of changes in status of BACKGROUND scans. Use <a href="#">cbStopBackground()</a> with AIFUNCTION to terminate the background process before it has completed. cbStopBackground() should be executed after normal termination of all background functions in order to clear variables and flags.</p>
CONTINUOUS	<p>This option puts the function in an endless loop. Once it collects the required number of samples, it resets to the start of the buffer and begins again. The only way to stop this operation is with cbStopBackground(). Normally this option should be used in combination with BACKGROUND so that your program will regain control.</p> <p><b>Count argument settings in CONTINUOUS mode:</b> For some DAQ hardware, Count must be an integer multiple of the <i>packet size</i>. Packet size is the amount of data that a DAQ device transmits back to the PC's memory buffer during each data transfer. Packet size can differ among DAQ hardware, and can even differ on the same DAQ product depending on the transfer method.</p> <p>In some cases, the minimum value for the Count argument may change when the CONTINUOUS option is used. This can occur for several reasons; the most common is that in order to trigger an interrupt on boards with FIFOs, the circular buffer must occupy at least half the FIFO. Typical half-FIFO sizes are 256, 512 and 1024.</p> <p>Another reason for a minimum Count value is that the buffer in memory must be periodically transferred to the user buffer. If the buffer is too small, data will be overwritten during the transfer resulting in garbled data.</p> <p>Refer to the board-specific information in the <i>Universal Library User's Guide</i> (available on our web site at <a href="http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf">www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf</a>) for packet size information for your particular DAQ hardware.</p>
EXTCLOCK	<p>If this option is used, conversions will be controlled by the signal on the external clock input rather than by the internal pacer clock. Each conversion will be triggered on the appropriate edge of the clock input signal (refer to the board-specific information in the <i>Universal Library User's Guide</i>). In most cases, when this option is used the Rate argument is ignored. The sampling rate is dependent on the clock signal. Options for the board will default to a transfer mode that will allow the maximum conversion rate to be attained unless otherwise specified.</p> <p>In some cases, such as with the PCI-DAS4020/12, an approximation of the rate is used to determine the size of the packets to transfer from the board. Set the Rate argument to an approximate maximum value.</p>

	<p><b>SINGLEIO is recommended for slow external clock rates:</b> If the rate of the external clock is very slow (for example less than 200 Hz) and the board you are using supports <code>BLOCKIO</code>, you may want to include the <code>SINGLEIO</code> option. The reason for this is that the status for the operation is not available until one packet of data has been collected (typically 512 samples). The implication is that, if acquiring 100 samples at 100 Hz using <code>BLOCKIO</code> (the default for boards that support it if <code>EXTCLOCK</code> is used), the operation will not complete until 5.12 seconds has elapsed.</p>
<code>EXTMEMORY</code>	<p>Causes the command to send the data to a connected memory board via the DT-Connect interface rather than returning the data to the buffer. Data for each call to this function will be appended unless <a href="#">cbMemReset()</a> is called. The data should be unloaded with the <a href="#">cbMemRead()</a> function before collecting new data. When <code>EXTMEMORY</code> option is used, the <code>MemHandle</code> argument can be set to null or 0. <code>CONTINUOUS</code> option cannot be used with <code>EXTMEMORY</code>. Do not use <code>EXTMEMORY</code> and <code>DTCONNECT</code> together. The transfer modes <code>DMAIO</code>, <code>SINGLEIO</code>, <code>BLOCKIO</code> and <code>BURSTIO</code> have no meaning when used with this option.</p>
<code>EXTTRIGGER</code>	<p>If this option is specified, the sampling will not begin until the trigger condition is met. On many boards, this trigger condition is programmable (see <a href="#">cbSetTrigger()</a> on page 60 and board-specific information for details) and can be programmed for rising or falling edge or an analog level.</p> <p>On other boards, only 'polled gate' triggering is supported. In this case, assuming active high operation, data acquisition will commence immediately if the trigger input is high. If the trigger input is low, acquisition will be held off until it goes high. Acquisition will then continue until <code>NumPoints</code> samples have been taken regardless of the state of the trigger input. For "polled gate" triggering, this option is most useful if the signal is a pulse with a very low duty cycle (trigger signal in TTL low state most of the time) so that triggering will be held off until the occurrence of the pulse.</p>
<code>NOTODINTS</code>	<p>If this option is specified, the system's time-of-day interrupts are disabled for the duration of the scan. These interrupts are used to update the system's real time clock and are also used by various other programs. These interrupts can limit the maximum sampling speed of some boards - particularly the PCM-DAS08. If the interrupts are turned off using this option, the real-time clock will fall behind by the length of time that the scan takes.</p>
<code>NOCALIBRATEDATA</code>	<p>Turns off real-time software calibration for boards which are software calibrated. This is done by applying calibration factors to the data on a sample by sample basis as it is acquired. Examples are the PCM-DAS16/330 and PCM-DAS16x/12. Turning off software calibration saves CPU time during a high speed acquisition run. This may be required if your processor is less than a 150 MHz Pentium and you desire an acquisition speed in excess of 200 kHz. These numbers may not apply to your system. Only trial will tell for sure. DO NOT use this option if you do not have to. If this option is used, the data must be calibrated after the acquisition run with the <a href="#">cbACalibrateData()</a> function.</p>
<code>DTCONNECT</code>	<p>All A/D values will be sent to the A/D board's DT-Connect port. This option is incorporated into the <code>EXTMEMORY</code> option. Use <code>DTCONNECT</code> only if the external board is not supported by Universal Library.</p>
<code>RETRIGMODE</code>	<p>Re-arms the trigger after a trigger event is performed. With this mode, the scan begins when a trigger event occurs. When the scan completes, the trigger is re-armed to acquire the next batch of data. You can specify the number of samples in the scan for each trigger event (described below). The <code>RETRIGMODE</code> option can be used with the <code>CONTINUOUS</code> option to continue arming the trigger until <code>cbStopBackground()</code> is called.</p>

You can specify the number of samples to acquire with each trigger event. This is the trigger count. Use the `ConfigItem` option `BIADTRIGCOUNT` with `cbSetConfig()` to set the trigger count. If you specify a trigger count that is either zero or greater than the value of the `cbAInScan()` `Count` argument, the trigger count will be set to the value of the `Count` argument.

Specify the `CONTINUOUS` option with the trigger count set to zero to fill the buffer with `Count` samples, re-arm the trigger, and refill the buffer upon the next trigger.

**Caution!** You will generate an error if you specify a total A/D rate beyond the capability of the board. For example, if you specify `LowChan = 0`, `HighChan = 7` (8 channels total), and `Rate = 20,000`, and you are using a CIO-DAS16/JR, you will get an error — you have specified a total rate of  $8 \times 20,000 = 160,000$ , but the CIO-DAS16/JR is capable of converting only 120,000 samples per second. The maximum sampling rate depends on the A/D board that is being used. It is also dependent on the sampling mode options.

**Important**

In order to understand the functions, you must read the board-specific information found in the *Universal Library User's Guide* (available on our web site at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf)). The example programs should be examined and run prior to attempting any programming of your own. Following this advice will save you hours of frustration, and possibly time wasted holding for technical support.

This note, which appears elsewhere, is especially applicable to this function. Now is the time to read the board-specific information for your board that is contained in the *Universal Library User's Guide*. We suggest that you make a copy of this information for reference as you read this manual and examine the example programs.

## cbALoadQueue()

Loads the A/D board's channel/gain queue. This function only works with A/D boards that have channel/gain queue hardware.

Some products do not support channel / gain queue, and some that do support it are limited on the order of elements, number of elements, and gain values that can be included, etc. Please refer to the device-specific information in the *Universal Library User's Guide* (available on our web site at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf)) to find details for your particular product.

### Function prototype:

C/C++:	<code>int cbALoadQueue(int BoardNum, short ChanArray[], short GainArray[], int Count)</code>
Visual Basic:	<code>Function cbALoadQueue(ByVal BoardNum&amp;, ChanArray%, GainArray%, ByVal Count&amp;) As Long</code>
Delphi:	<code>function cbALoadQueue(BoardNum:Integer; var ChanArray:SmallInt; var GainArray:SmallInt; Count:LongInt):Integer;</code>

### Arguments:

BoardNum	The board number used to collect the data. BoardNum may be 0 to 99. Refers to the number associated with the board used to collect the data when it was installed with the <i>InstaCal</i> configuration program. The specified board must have an A/D and a channel/gain queue.
ChanArray	Array containing channel values. This array should contain all of the channels that will be loaded into the channel gain queue.
GainArray	Array containing A/D range values. This array should contain each of the A/D ranges that will be loaded into the channel gain queue.
Count	Number of elements in ChanArray and GainArray or 0 to disable channel/gain queue. Specifies the total number of channel/gain pairs that will be loaded into the queue. ChanArray and GainArray should contain at least Count elements. Set Count = 0 to disable the board's channel/gain queue. The maximum value is specific to the queue size of the A/D boards channel gain queue.

### Returns:

[Error code](#) or 0 if no errors.

### Notes:

Normally the [cbAInScan\(\)](#) function scans a fixed range of channels (from LowChan to HighChan) at a fixed A/D range. If you load the channel gain queue with this function then all subsequent calls to `cbAInScan()` will cycle through the channel/range pairs that you have loaded into the queue.



## cbAOut()

Sets the value of a D/A output.

### Function prototype:

C/C++:	<code>int cbAOut(int BoardNum, int Channel, int Range, unsigned short DataValue)</code>
Visual Basic:	<code>Function cbAOut(ByVal BoardNum&amp;, ByVal Channel&amp;, ByVal Range&amp;, ByVal DataValue%) As Long</code>
Delphi:	<code>function cbAOut(BoardNum:Integer; Channel:Integer; Range:Integer; DataValue:Word):Integer;</code>

### Arguments:

BoardNum	The board number used to collect the data. BoardNum may be 0 to 99. Refers to the number associated with the board used to collect the data when it was installed with the <i>InstaCal</i> configuration program. The specified board must have a D/A.
Channel	D/A channel number. The maximum allowable channel depends on which type of D/A board is being used.
Range	D/A range code. The output range of the D/A channel can be set to any of those supported by the board. If the D/A board does not have programmable ranges then this argument will be ignored. See Table 3 on page 28 for valid values.
DataValue	Value to set D/A to. Must be in the range 0 - N where N is the value $2^{\text{Resolution}} - 1$ of the converter  Exception: Using 16-bit boards with Basic range is -32768 to 32767. Refer to the discussion of Basic signed integers in the "16-bit values using a signed integer data type" section in the "Universal Library Description & Use" chapter of the <i>Universal Library User's Guide</i> (available on our web site at <a href="http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf">www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf</a> .)

### Returns:

[Error code](#) or 0 if no errors

### Notes:

**"Simultaneous Update" or "Zero Power-Up" boards:** If you set the simultaneous update jumper for simultaneous operation, use [cbAOutScan\(\)](#) for simultaneous update of multiple channels. `cbAOut()` always writes the D/A data then reads the D/A, which causes the D/A output to be updated.

## cbAOutScan()

Outputs values to a range of D/A channels. This function can be used for paced analog output on hardware that supports paced output. It can also be used to update all analog outputs at the same time when the `SIMULTANEOUS` option is used.

### Function prototype:

C/C++:	<pre>int cbAOutScan(int BoardNum, int LowChan, int HighChan, long NumPoints, long *Rate, int Range, int MemHandle, int Options)</pre>
Visual Basic:	<pre>Function cbAOutScan(ByVal BoardNum&amp;, ByVal LowChan&amp;, ByVal HighChan&amp;, ByVal NumPoints&amp;, Rate&amp;, ByVal Range&amp;, ByVal MemHandle&amp;, ByVal Options&amp;) As Long</pre>
Delphi:	<pre>function cbAOutScan(BoardNum:Integer; LowChan:Integer; HighChan:Integer; NumPoints:Longint; var Rate:Longint; Range:Integer; MemHandle:Integer; Options:Integer):Integer;</pre>

### Arguments:

BoardNum	Refers to the board number associated with the board when it was installed with the <i>InstaCal</i> configuration program. The specified board must have a D/A. BoardNum may be 0 to 99.
LowChan	First D/A channel of scan.
HighChan	Last D/A channel of scan. LowChan/HighChan - The maximum allowable channel depends on which type of D/A board is being used.
NumPoints	Number of D/A values to output. Specifies the total number of D/A values that will be output. Most D/A boards do not support timed outputs. For these boards, set the count to the number of channels in the scan.
Rate	Sample rate in scans per second. For many D/A boards the Rate is ignored and can be set to <code>NOTUSED</code> . For D/A boards with trigger and transfer methods which allow fast output rates, such as the CIO-DAC04/12-HS, Rate should be set to the D/A output rate (in scans/sec). This argument returns the value of the actual rate set. This value may be different from the user specified rate due to pacer limitations.  If supported, this is the rate at which scans are triggered. If you are updating 4 channels, 0-3, then specifying a rate of 10,000 scans per second (10 kHz) will result in the D/A converter rates of 10 kHz (one D/A per channel). The data transfer rate is 40,000 words per second; 4 channels * 10,000 updates per scan.  The maximum update rate depends on the D/A board that is being used. It is also dependent on the sampling mode options.
Range	D/A range code. The output range of the D/A channel can be set to any of those supported by the board. If the D/A board does not have a programmable gain, this argument is ignored. See Table 3 on page 28 for valid values.
MemHandle	Handle for Windows buffer from which data will be output. This buffer must have been previously allocated with the <a href="#">cbWinBufAlloc()</a> function and data values loaded (perhaps using <a href="#">cbWinArrayToBuf()</a> ).
Options	Bit fields that control various options. This field may contain any combination of non-contradictory choices from the values listed in the "Options argument values" section on page 43.

**Returns:**

[Error code](#) or 0 if no errors.

Rate - Actual sampling rate used.

**Options argument values:**

CONTINUOUS	This option may only be used with boards which support interrupt, DMA or REP-INSW transfer methods. This option puts the function in an endless loop. Once it outputs the specified number (NumPoints) of D/A values, it resets to the start of the buffer and begins again. The only way to stop this operation is by calling <a href="#">cbStopBackground()</a> with AOFUNCTION. This option should only be used in combination with BACKGROUND so that your program can regain control.
BACKGROUND	This option may only be used with boards which support interrupt, DMA or REP-INSW transfer methods. When this option is used, the D/A operations will begin running in the background and control will immediately return to the next line of your program. Use <a href="#">cbGetStatus()</a> with AOFUNCTION to check the status of background operation. Alternatively, some boards support <a href="#">EnableEvent()</a> for event notification of changes in status of BACKGROUND scans. Use <a href="#">cbStopBackground()</a> with AOFUNCTION to terminate background operations before they are completed. <a href="#">cbStopBackground()</a> should be executed after normal termination of all background functions in order to clear variables and flags.
SIMULTANEOUS	When this option is used (if the board supports it and the appropriate switches are set on the board) all of the D/A voltages will be updated simultaneously when the last D/A in the scan is updated. This generally means that all the D/A values will be written to the board, then a read of a D/A address causes all D/As to be updated with new values simultaneously.
EXTCLOCK	<p>If this option is specified, conversions will be paced by the signal on the external clock input rather than by the internal pacer clock. Each conversion will be triggered on the appropriate edge of the clock input signal (refer to board-specific information contained in the <i>Universal Library Users Guide</i>).</p> <p>When this option is used the Rate argument is ignored. The sampling rate is dependent on the clock signal. Options for the board default to transfer types that allow the maximum conversion rate to be attained unless otherwise specified.</p>
EXTTRIGGER	If this option is specified, the sampling will not begin until the trigger condition is met. On many boards, this trigger condition is programmable (see <a href="#">cbSetTrigger()</a> on page 60 and board-specific information for details).
ADCCLOCKTRIG	Triggers a data output operation when the ADC clock starts.
ADCCLOCK	Paces the data output operation using the ADC clock.

**Caution!** You will generate an error if you specify a total D/A rate beyond the capability of the board. For example: If you specify LowChan = 0 and HighChan = 3 (4 channels total) and Rate = 100,000, and you are using a cSBX-DDA04, you will get an error. You have specified a total rate of 4\*100,000 = 400,000. The cSBX-DDA04 is rated to 330,000 updates per second. The maximum update rate depends on the D/A board that is being used. It is also dependent on the sampling mode options.

## cbAPretrig()

Waits for a trigger to occur and then returns a specified number of analog samples before and after the trigger occurred. If only 'polled gate' triggering is supported, the trigger input line (refer to the user's manual for the board) must be at TTL low before this function is called, or a TRIGSTATE error will occur. The trigger occurs when the trigger condition is met. Refer to [cbSetTrigger\(\)](#) on page 60 for details.

### Function prototype:

C/C++:	<pre>int cbAPretrig(int BoardNum, int LowChan, int HighChan, long *PretrigCount, long *TotalCount, long *Rate, int Range, int MemHandle, int Options)</pre>
Visual Basic:	<pre>Function cbAPretrig(ByVal BoardNum&amp;, ByVal LowChan&amp;, ByVal HighChan&amp;, PretrigCount&amp;, TotalCount&amp;, Rate&amp;, ByVal Range&amp;, ByVal MemHandle&amp;, ByVal Options&amp; ) As Long</pre>
Delphi:	<pre>function cbAPretrig(BoardNum:Integer; LowChan:Integer; HighChan:Integer; var PretrigCount:Longint; var TotalCount:Longint; var Rate:Longint; Range:Integer; MemHandle:Integer; Options:Integer):Integer;</pre>

### Arguments:

BoardNum	Refers to the board number associated with the board when it was installed with the <i>InstaCal</i> configuration program. The specified board must have an A/D. BoardNum may be 0 to 99.
LowChan	First A/D channel of scan.
HighChan	Last A/D channel of scan.  <b>LowChan/HighChan:</b> The maximum allowable channel depends on which type of A/D board is being used. For boards with both single ended and differential inputs, the maximum allowable channel number also depends on how the board is configured (e.g., 8 channels for differential inputs, 16 for single ended inputs).
PretrigCount	Number of pre-trigger A/D samples to collect. Specifies the number of samples to collect before the trigger occurs.  For products using a hardware implementation of pretrigger (most products), PretrigCount must be less than (TotalCount - 512). For these devices, if the trigger occurs too early, fewer than the requested number of pre-trigger samples will be collected, and a TOOFEW error will occur. The PretrigCount will be set to indicate how many samples were actually collected. The post trigger samples will still be collected.  For software implementations of pretrigger, PretrigCount must be less than TotalCount. For these devices, triggers that occur before the requested number of pre-trigger samples are collected are ignored. See board-specific information.
TotalCount	Total number of A/D samples to collect. Specifies the total number of samples that will be collected and stored in the buffer.  For products using a hardware implementation of pretrigger (most products), TotalCount must be greater than or equal to the PretrigCount + 512. If the trigger occurs too early, fewer than the requested number of samples will be collected, and a TOOFEW error will occur. The TotalCount will be set to indicate how many samples were actually collected.  For software implementations of pretrigger, TotalCount must be greater than PretrigCount. For these devices, triggers that occur before the requested number of pre-trigger samples are collected are ignored. See board-specific information.

	TotalCount must be evenly divisible by the number of channels being scanned. If it is not, this function will adjust the number (down) to the next valid value and return that value to the TotalCount argument.
	PretrigCount must also be evenly divisible by the number of channels being scanned. If it is not, this function will adjust the number (up) to the next valid value and return that value to the PretrigCount argument.
Rate	Sample rate in scans per second.
Range	A/D Range code. If the selected A/D board does not have a programmable gain feature, this argument is ignored. Otherwise, set to any range that is supported by the selected A/D board. See Table 3 on page 28 for valid values.
MemHandle	Handle for Windows buffer to store data. This buffer must have been previously allocated with the <a href="#">cbWinBufAlloc()</a> function.  For hardware trigger types, the buffer referenced by MemHandle must be big enough to hold at least TotalCount + 512 integers.
Options	Bit fields that control various options. This field may contain any combination of non-contradictory choices from the values listed in the "Options argument values" section below.

**Returns:**

[Error code](#) or 0 if no errors

PretrigCount - Number of pre-trigger samples

TotalCount - Total number of samples collected

Rate - Actual sampling rate

MemHandle - Collected A/D data returned via the Windows buffer

**Options argument values:**

CONVERTDATA	For hardware trigger types, the data is collected into a "circular" buffer. When the data collection is complete, the data is in the wrong order. If you use the CONVERTDATA option, the data is automatically rotated into the correct order (and converted to 12-bit values if required) when the data acquisition is complete. Otherwise, call <a href="#">cbAConvertPretrigData()</a> to rotate the data. You cannot use the CONVERTDATA option in combination with the BACKGROUND option for this function. The CONVERTDATA option is not required for software triggered types.
BACKGROUND	If the BACKGROUND option is not used, the cbAPretrig() function will not return to your program until all of the requested data has been collected and returned to the buffer. When the BACKGROUND option is used, control returns immediately to the next line in your program, and the data collection from the A/D into the buffer will continue in the background. Use <a href="#">cbGetStatus()</a> with AIFUNCTION to check on the status of the background operation. Alternatively, some boards support <a href="#">cbEnableEvent()</a> for event notification of changes in status of BACKGROUND scans. Use <a href="#">cbStopBackground()</a> with AIFUNCTION to terminate the background process before it has completed.  Call cbStopBackground() after normal termination of all background functions to clear variables and flags. For hardware trigger types, you cannot use the CONVERTDATA option in combination with the BACKGROUND option for this function. To correctly order and parse the data, use <a href="#">cbAConvertPretrigData()</a> after the function completes.
EXTCLOCK	This option is available only for boards that have separate inputs for external pacer and external trigger. See your hardware manual or board-specific information.

## EXTMEMORY

Causes this function to send the data to a connected memory board via the DT-Connect interface rather than returning the data to the buffer. If you use this option to send the data to a MEGA-FIFO memory board, then you must use [cbMemReadPretrig\(\)](#) to later read the pre-trigger data from the memory board. If you use [cbMemRead\(\)](#), the data will NOT be in the correct order.

Every time this option is used, it overwrites any data already stored in the memory board. All data should be read from the board (with [cbMemReadPretrig\(\)](#)) before collecting any new data. When this option is used, the `MemHandle` argument is ignored. The MEGA-FIFO memory must be fully populated in order to use the `cbAPretrig()` function with the `EXTMEMORY` option.

## DTCONNECT

When `DTCONNECT` option is used with this function the data from ALL A/D conversions is sent out the DT-Connect interface. While this function is waiting for a trigger to occur, it will send data out the DT-Connect interface continuously. If you have a Measurement Computing memory board plugged into the DT-Connect interface then you should use `EXTMEMORY` option rather than this option.

**Important:**

For hardware trigger types, the buffer referenced by `MemHandle` must be big enough to hold at least `TotalCount + 512` integers.

## cbATrig()

Waits for a specified analog input channel to go above or below a specified value. `cbATrig` continuously reads the specified channel and compares its value to `TrigValue`. Depending on whether `TrigType` is set to `TRIGABOVE` or `TRIGBELOW`, it waits for the first A/D sample that is above or below `TrigValue`. The first sample that meets the trigger criteria is returned to `DataValue`.

### Function prototype:

**C/C++:** `int cbATrig(int BoardNum, int Channel, int TrigType, int TrigValue, int Range, unsigned short *DataValue)`

**Visual Basic:** `Function cbATrig(ByVal BoardNum&, ByVal Channel&, ByVal TrigType&, ByVal TrigValue%, ByVal Range&, DataValue%) As Long`

**Delphi:** `function cbATrig (BoardNum:Integer; Channel:Integer; TrigType:Integer; TrigValue:Word; Range:Integer; var DataValue:Word):Integer;`

### Arguments:

BoardNum	Refers to the board number associated with the board when it was installed with the <i>InstaCal</i> configuration program. The specified board must have an A/D. BoardNum may be 0 to 99.
Channel	A/D channel number. The maximum allowable channel depends on which type of A/D board is being used. For boards with both single ended and differential inputs, the maximum allowable channel number also depends on how the board is configured. For example a CIO-DAS1600 has 8 channels for differential inputs and 16 channels for single ended inputs.
TrigType	TRIGABOVE or TRIGBELOW. Specifies whether to wait for the analog input to be ABOVE or BELOW the specified trigger value.
TrigValue	The threshold value that all A/D values are compared to. Must be in the range 0 - 4095 for 12-bit A/D boards, or 0-65,535 for 16-bit A/D boards. Refer to your BASIC manual for information on signed BASIC integer data types.
Range	Gain code. If the selected A/D board does not have a programmable gain feature, this argument is ignored. Otherwise, set to any range that is supported by the selected A/D board. See Table 3 on page 28 for valid values.
DataValue	Returns the value of the first A/D sample to meet the trigger criteria.

### Returns:

[Error code](#) or 0 if no errors

`DataValue` - Value of the first A/D sample to match the trigger criteria.

### Notes:

Pressing **Ctrl-C** will not terminate the wait for an analog trigger that meets the specified condition. There are only two ways to terminate this call: satisfy the trigger condition or reset the computer.

**Caution!** Use caution when using this function in Windows programs. All active windows will lock on the screen until the trigger condition is satisfied. The keyboard and mouse activity will also lock until the trigger condition is satisfied.

## cbVIn()

Reads an A/D input channel, and returns a voltage value. If the specified A/D board has programmable gain, then this function sets the gain to the specified range. The voltage value is returned to `DataValue`.

### Function prototype:

C/C++:	<code>int cbVIn(int BoardNum, int Channel, int Range, float *DataValue, int Options);</code>
Visual Basic:	<code>Function cbVIn(ByVal BoardNum&amp;, ByVal Channel&amp;, ByVal Range&amp;, DataValue!, ByVal Options&amp;) As Long</code>
Delphi:	<code>function cbVIn(BoardNum:Integer; Channel:Integer; Range:Integer; Var DataValue:Single; Options:Integer):Integer;</code>

### Arguments:

<code>BoardNum</code>	The board number associated with the board used to collect the data when it was installed with the <i>InstaCal</i> configuration program. <code>BoardNum</code> may be 0 to 99. The specified board must have an A/D.
<code>Channel</code>	A/D channel number. The maximum allowable channel depends on which type of A/D board is being used. For boards with both single-ended and differential inputs, the maximum allowable channel number also depends on how the board is configured.
<code>Range</code>	A/D range code. If the board has a programmable gain, it will be set according to this argument value. Keep in mind that some A/D boards have a programmable gain feature, and others set the gain via switches on the board. In either case, the range that the board is configured for must be passed to this function. Refer to Table 3 on page 28 for a list of valid range settings.
<code>DataValue</code>	A pointer or reference to the data value.
<code>Options</code>	Reserved for future use.

### Returns:

[Error code](#) or 0 if no errors.

`DataValue` - Returns the value in volts of the A/D sample.

### Options argument values:

<code>Default</code>	Reserved for future use.
----------------------	--------------------------



## cbVOut()

Sets the value of a D/A output.

### Function prototype:

C/C++:	<code>int cbVOut(int BoardNum, int Channel, int Range, float DataValue, int Options);</code>
Visual Basic:	<code>Function cbVOut(ByVal BoardNum&amp;, ByVal Channel&amp;, ByVal Range&amp;, ByVal DataValue!, ByVal Options&amp;) As Long</code>
Delphi:	<code>function cbVOut(BoardNum:Integer; Channel:Integer; Range:Integer; DataValue:Single; Options:Integer):Integer;</code>

### Arguments:

BoardNum	The board number associated with the board used to collect the data when it was installed with the <i>InstaCal</i> configuration program. BoardNum may be 0 to 99. The specified board must have an D/A.
Channel	The D/A channel number. The maximum allowable channel depends on which type of D/A board is being used.
Range	The D/A range code. If the board has a programmable gain, it will be set according to this argument value. The output range of the D/A channel can be set to any of those supported by the board. Keep in mind that some D/A boards have programmable voltage ranges, and others set the voltage range via switches on the board. In either case, the range that the board is configured for must be passed to this function.
DataValue	The voltage value to be written.
Options	Reserved for future use.

### Returns:

[Error code](#) or 0 if no errors.

### Options argument values:

Default	Reserved for future use.
---------	--------------------------

---

# Configuration Functions

## Introduction

This section covers Universal Library functions that retrieve or change configuration options on a board. The configuration information for all boards is stored in the configuration file CB.CFG. This information is loaded from CB.CFG by all programs that use the library.

To determine which of these functions are compatible with your hardware, refer to the *Universal Library User's Guide* (available in PDF format on our website at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf)).

## cbGetConfig()

Returns a configuration option for a board. The configuration information for all boards is stored in the CB.CFG file. This information is loaded from CB.CFG by all programs that use the library. You can change the current configuration within a running program with the [cbSetConfig\(\)](#) function. The `cbGetConfig()` function returns the current configuration information.

### Function prototype:

C/C++:	<code>int cbGetConfig(int InfoType, int BoardNum, int DevNum, int ConfigItem, int *ConfigVal)</code>
Visual Basic:	<code>Function cbGetConfig(ByVal InfoType&amp;, ByVal BoardNum&amp;, ByVal DevNum&amp;, ByVal ConfigItem&amp;, ConfigVal&amp;) As Long</code>
Delphi:	<code>function cbGetConfig(InfoType:Integer; BoardNum:Integer; DevNum:Integer; ConfigItem:Integer; var ConfigVal:Integer):Integer;</code>

### Arguments:

InfoType	The configuration information for each board is grouped into different categories. This argument specifies which category you want. Set it to one of the constants listed in the "InfoType argument values" section below.
BoardNum	Refers to the board number associated with a board when it was installed. BoardNum may be 0 to 99.
DevNum	Selects a particular device. If InfoType = DIGITALINFO, then DevNum specifies which of the board's digital devices you want information on. If InfoType = COUNTERINFO, then DevNum specifies which of the board's counter devices you want information from.
ConfigItem	Specifies which configuration item you wish to retrieve. Set it in conjunction with the InfoType argument using the table in the "ConfigItem argument values" section on page 52.
ConfigVal	The specified configuration item is returned to this variable.

### Returns:

[Error code](#) or 0 if no errors.

ConfigVal - returns the value of the specified configuration item here.

### InfoType argument values:

GLOBALINFO	Information about the configuration file.
BOARDINFO	General information about a board.
DIGITALINFO	Information about a digital device.
COUNTERINFO	Information about a counter device.
EXPANSIONINFO	Information about an expansion device.
MISCINFO	One of the miscellaneous options for the board.

**ConfigItem argument values:**

Valid ConfigItem constant settings for each InfoType constant are as follows:

InfoType	ConfigItem	Description
GLOBALINFO	GIVERSION	CB.CFG file format - used by the library to determine compatibility.
	GINUMBOARDS	Maximum number of installable boards
	GINUMEXPBOARDS	Maximum number of expansion boards allowed to be installed.
BOARDINFO	BIADCSETTLETIME	ADC settling time
	BIBASEADR	Base address of the board
	BIBOARDTYPE	Returns a unique number in the range of 0 to 8000 Hex describing the board type installed.
	BICIDEVNUM	Index into counter information for the first device.
	BICINUMDEVS	Number of counter devices
	BICLOCK	Clock frequency in MHz (40, 10, 8, 6, 5, 4, 3, 2, 1) or 0 for not supported.
	BIDACSTARTUP	Returns the setting of a DAC board's configuration register STARTUP bit. Refer to the "Notes" section for the cbSetConfig() method on page 61 for more information.
	BIDACUPDATEMODE	Setting of the update mode for a digital-to-analog converter (DAC). Refer to the "Notes" section on page 53 for more information.
	BIDIDEVNUM	Index into digital information for the first device
	BIDINUMDEVS	Number of digital devices
	BIDMACHAN	DMA channel. 0, 1 or 3.
	BIDTBOARD	Board number of the connected DT board
	BIFACTORYID	The factory serial number of a USB device, or the MAC address of a WEB device.
	BIHIDELOGINDLG	Enables or disables the Device Login dialog. Set to a nonzero value to disable the dialog. When disabled, the cbDeviceLogin() function must be used to log in to a device session.
	BIINTLEVEL	Interrupt level. 0 for none, or 1 - 15
	BINETIOTIMEOUT	The amount of time (in milliseconds) to wait for a WEB device to acknowledge a command or query sent to the device over a network connection. If no acknowledgement is received in this time a timeout occurs.
	BINUMADCHANS	Number of A/D channels
	BINUMDACHANS	Number of D/A channels
	BINUMIOPORTS	Number of IO ports used by board
	BIPANID	Personal Area Network (PAN) identifier for a USB device that supports wireless communication.
	BIRANGE	Selected voltage range. For switch-selectable gains only. If the selected A/D board does not have a programmable gain feature, this argument returns the range as defined by the installed InstaCal settings. If InstaCal and the board are installed correctly, the returned range will correspond to the input range as set via the switches on the board. Refer to board-specific information for a list of the A/D ranges supported by each board.
	BIRFCHANNEL	Returns the RF channel number used to transmit/receive data by a USB device that supports wireless communication.
	BIRSS	Returns the received signal strength in dBm of a remote device.
	BISERIALNUM	Returns the serial number assigned by a user to a USB device in InstaCal. This ConfigItem does not return the factory serial number.
	BIWAITSTATE	Setting of Wait State jumper. 1 = enabled, 0 = disabled
	BIUSEEXPS	Supports expansion boards TRUE/FALSE

InfoType	ConfigItem	Description
DIGITALINFO	DIDEVTYPE	Device Type - AUXPORT, FIRSTPORTA etc.
	DICONFIG	Current configuration INPUT or OUTPUT
	DINUMBITS	Number of bits in the port
	DICURVAL	Current value of outputs
	DIINMASK	Returns the bit configuration of the specified port. Any bits that return a value of 1 are configured for input. Refer to the " <a href="#">Notes</a> " section below for more information.
	DIOUTMASK	Returns the bit configuration of the specified port. Any bits that return a value of 1 are configured for output. Refer to the " <a href="#">Notes</a> " section below for more information.
COUNTERINFO	CICTRTYPE	Counter chip type, where 1 = 8254, 2 = 9513, 3 = 8536, 4 = 7266, 5 = event counter, 6 = scan counter, and 7 = timer counter.
EXPANSIONINFO	XIBOARDTYPE	Board type (refer to the "BoardType Codes" topic in the <i>Universal Library User's Guide</i> )
	XIMUXADCHAN1	A/D channel EXP board is connected to
	XIMUXADCHAN2	2nd A/D channel EXP board is connected to
	XIRANGE1	Range (gain) of low 16 channels
	XIRANGE2	Range (gain) of high 16 channels
	XICJCCHAN	A/D channel that CJC is connected to
	XITHERMTYPE	Sensor type. Use one of the sensor types listed below: J = 1 K = 2 T = 3 E = 4 R = 5 S = 6 B = 7 Platinum .00392 = 257 Platinum .00391 = 258 Platinum .00385 = 259 Copper .00427 = 260 Nickel/Iron .00581 = 261 Nickel/Iron .00527 = 262
	XINUMEXPCHANS	Number of channels on expansion board
	XIPARENTBOARD	Board number of parent A/D board

**Notes:**

- Use the DIINMASK and DIOUTMASK options to determine if an AUXPORT is configurable. Execute cbGetConfig() twice to the same port—once using DIINMASK and once using DIOUTMASK. If both of the ConfigVal arguments returned have input and output bits that overlap, the port is not configurable.

You can determine overlapping bits by *Anding* both arguments: For example, the PCI-DAS08 has seven bits of digital I/O (four outputs and three inputs). For this board, the ConfigVal returned by DIINMASK is always 7 (0000 0111), while the ConfigVal argument returned by DIOUTMASK is always 15 (0000 1111). When you *And* both ConfigVal arguments together, you get a non-zero number (7). Any non-zero number indicates that input and output bits overlap for the specified port, and the port is a non-configurable AUXPORT.

- Use the BIDACUPDATEMODE option to check the update mode for a DAC board.

With ConfigItem set to BIDACUPDATEMODE, if ConfigVal returns 0, the DAC update mode is immediate. Values written with [cbAOut\(\)](#) are automatically output by the DAC channels.

With ConfigItem set to BIDACUPDATEMODE, if ConfigVal returns 1, the DAC update mode is set to *on command*. Values written with [cbAOut\(\)](#) are not output by the DAC channels until a [cbSetConfig\(\)](#)

call is made with its ConfigItem argument set to BIDACUPDATECMD.

- Use the BIDACSTARTUP option (ConfigItem argument) Returns 0 if startup bit is disabled, or 1 to if startup bit is enabled to determine if the DAC values before the board was last powered down are stored.

Refer to the "[Notes](#)" section for cbSetConfig() on page 61 for more information.

To store the current DAC values as start-up values, call cbSetConfig() with a value of 1 for the BIDACSTARTUP value. Then, call [cbAOut\(\)](#) or [cbAOutScan\(\)](#) for each channel, and call cbSetConfig() again with a value of 0 for the BIDACSTARTUP value.

**Example:**

```
cbSetConfig(BOARDINFO, boardNumber, 0, BIDACSTARTUP, 1);
for (int i =1; i <8; i++)
{
    cbAOut(boardNumber, i, BIP5VOLTS, DACValue[i]);
}
cbSetConfig(BOARDINFO, boardNumber, 0, BIDACSTARTUP, 0);
```

To store the DAC's last settings, call cbSetConfig() with a BIDACSTARTUP value of 1. Leave this bit turned on until the application exits. The next time the board is powered up, it restores the values last written to the DACs.

## cbGetConfigString()

Retrieves configuration or device information as a null-terminated string.

### Function prototype:

**C/C++:**                    `int cbGetConfigString(int InfoType, int BoardNum, int ItemIndex, int ConfigItem, char *ConfigVal, int* maxConfigLen)`

**Visual Basic:**            `Function cbGetConfigString(ByVal InfoType, ByVal BoardNum&, ByVal ItemIndex&, ByVal ConfigItem&, ByVal ConfigVal$, ByRef maxConfigLen&) As Long`

**Delphi:**                    `function cbGetConfigString(InfoType:Integer; BoardNum:Integer; ItemIndex:Integer; ConfigItem:Integer; ConfigVal:PChar; var maxConfigLen:Integer):Integer;`

### Arguments:

InfoType	The configuration information for each board is grouped into different categories. This argument specifies which category you want. Always set this argument to BOARDINFO.
BoardNum	Refers to the board number associated with a board when it was installed. BoardNum may be 0 to 99.
ItemIndex	The location in the device memory (specified by ConfigItem) at which to start reading.
ConfigItem	Specifies the type of information (or memory area) to read from the device. Set it to one of the constants listed in the "ConfigItem argument values" section below.
ConfigVal	Pointer to a user-allocated buffer where the configuration string is copied.
maxConfigLen	Pointer to the value holding the maximum number of bytes to be read from the device into ConfigVal.

### Returns:

Error code or 0 if no errors.

maxConfigLen	The number of bytes read from the device into ConfigVal.
ConfigVal	The string read from the device.

### ConfigItem argument values:

BIDEVNOTES	Reads up to maxConfigLen characters/bytes from the device notes memory, starting at the location defined by ItemIndex. Currently supported only for WLS Series devices.
BIFACTORYID	Reads the MAC address of a WEB device.
BINODEID	Reads up to maxConfigLen character/bytes from the string identifier memory. Note that ItemIndex is not used for this ConfigItem.

## cbGetSignal()

Retrieves the configured Auxiliary or DAQ Sync connection and polarity for the specified timing and control signal.

This function is intended for advanced users. Except for the SYNC\_CLK input, you can easily view the settings for the timing and control signals using *InstaCal*.

**Note:** This function is not supported by all board types.

### Function prototype:

C/C++: `int cbGetSignal(int BoardNum, int Direction, int Signal, int Index, int* Connection, int* Polarity)`

Visual Basic: `Function cbGetSignal(ByVal BoardNum&, ByVal Direction&, ByVal Signal&, ByVal Index&, ByRef Connection, ByRef Polarity) As Long`

Delphi: `function cbGetSignal(BoardNum:Integer; Direction:Integer; Signal:Integer; Index:Integer; var Connection:Integer; var Polarity:Integer):Integer;`

### Arguments:

BoardNum	Refers to the board number associated with the A/D board when it was installed. The specified board must have configurable signal inputs and outputs.
Direction	Specifies whether retrieving the source (SIGNAL_IN) or destination (SIGNAL_OUT) of the specified signal.
Signal	Signal type whose connection is to be retrieved. See <a href="#">cbSelectSignal()</a> on page 57 for valid signal types.
Index	Used to indicate which connection to reference when there is more than one connection associated with the output Signal type. When querying output signals, increment this value until BADINDEX is returned or 0 is returned via the Connection parameter to determine all the output Connections for the specified output Signal. The first Connection is indexed by 0. For input signals (Direction=SIGNAL_IN), this should always be set to 0.
Connection	The specified connection is returned through this variable. This is set to 0 if no connection is associated with the Signal, or if the Index is set to an invalid value.
Polarity	Holds the polarity for the associated Signal and Connection.  For output Signals assigned an AUXOUT Connection, the return value is either INVERTED or NONINVERTED. For Signal settings of ADC_CONVERT, DAC_UPDATE, ADC_TB_SRC and DAC_TB_SRC input signals, either POSITIVEEDGE or NEGATIVEEDGE are returned. All other signals return 0.

### Returns:

[Error code](#) or 0 if no errors.

### Notes:

Timing and control configuration information can be viewed and edited inside InstaCal. Do the following:

1. Run InstaCal.
2. Click on the board and press the **Configure...** button or menu item. If the board supports DAQ Sync and Auxiliary Input/Output signal connections, a button labeled **Advanced Timing & Control Configuration** displays.
3. Press this button to open a display for viewing and modifying the above timing and control signals.



## cbSelectSignal()

Configures timing and control signals to use specific Auxiliary or DAQ Sync connections as a source or destination.

This function is intended for advanced users. Except for the SYNC\_CLK input, you can easily configure all the timing and control signals using *InstaCal*.

Note: This function is not supported by all board types. Please refer to board specific information for details.

### Function prototype:

**C/C++:** `int cbSelectSignal(int BoardNum, int Direction, int Signal, int Connection, int Polarity);`

**Visual Basic:** `Function cbSelectSignal(ByVal BoardNum&, ByVal Direction&, ByVal Signal&, ByVal Connection&, ByVal Polarity&) as Long`

**Delphi:** `Function cbSelectSignal(BoardNum:Integer; Direction:Integer; Signal:Integer; Connection:Integer; Polarity:Integer):Integer; StdCall;`

### Arguments:

BoardNum	Refers to the board number associated with the A/D board when it was installed. The specified board must have configurable signal inputs and outputs.
Direction	Direction of the specified signal type to be assigned a connector pin. For most signal types, this should be either <code>SIGNAL_IN</code> or <code>SIGNAL_OUT</code> . For the <code>SYNC_CLK</code> , <code>ADC_TB_SRC</code> and <code>DAC_TB_SRC</code> signals, the external source can also be disabled by specifying <code>DISABLED(=0)</code> such that it is neither input nor output. Set it in conjunction with the <code>Signal</code> , <code>Connection</code> , and <code>Polarity</code> arguments using the tables in the "Direction argument values" on page 58.
Signal	Signal type to be associated with a connector pin. Set it to one of the constants listed in the "Signal argument values" section below.
Connection	Designates the connector pin to associate the signal type and direction. Since individual pin selection is not allowed for the DAQ-Sync connectors, all DAQ-Sync pin connections are referred to as <code>DS_CONNECTOR</code> . The <code>AUXIN</code> and <code>AUXOUT</code> settings match their corresponding hardware pin names.
Polarity	<code>ADC_TB_SRC</code> and <code>DAC_TB_SRC</code> input signals ( <code>SIGNAL_IN</code> ) can be set for either rising edge ( <code>POSITIVEEDGE</code> ) or falling edge ( <code>NEGATIVEEDGE</code> ) signals. The <code>AUXOUT</code> connections can be set to <code>INVERTED</code> or <code>NONINVERTED</code> from their internal polarity.

### Returns:

[Error code](#) or 0 if no errors.

### Signal argument values:

<code>ADC_CONVERT</code>	A/D conversion pulse or clock.
<code>ADC_GATE</code>	External gate for A/D conversions.
<code>ADC_SCANCLK</code>	A/D channel scan signal.
<code>ADC_SCAN_STOP</code>	A/D scan completion signal.
<code>ADC_SSH</code>	A/D simultaneous sample and hold signal.
<code>ADC_STARTSCAN</code>	Start of A/D channel-scan sequence signal.
<code>ADC_START_TRIG</code>	A/D scan start trigger.

ADC_STOP_TRIG	A/D stop- or pre- trigger.
ADC_TB_SRC	A/D pacer timebase source.
CTR1_CLK	CTR1 clock source.
CTR2_CLK	CTR2 clock source.
DAC_START_TRIG	D/A start trigger.
DAC_TB_SRC	D/A pacer timebase source.
DAC_UPDATE	D/A update signal.
DGND	Digital ground.
SYNC_CLK	STC timebase signal.

**Direction argument values:**

Valid input (Direction= SIGNAL\_IN) settings include:

Signal	Connection	Polarity
ADC_CONVERT	AUXIN0..AUXIN5 DS_CONNECTOR	POSITIVEEDGE or NEGATIVEEDGE
ADC_GATE	AUXIN0..AUXIN5	See <a href="#">cbSetTrigger()</a> .
ADC_START_TRIG	AUXIN0..AUXIN5 DS_CONNECTOR	See <a href="#">cbSetTrigger()</a> .
ADC_STOP_TRIG	AUXIN0..AUXIN5 DS_CONNECTOR	See <a href="#">cbSetTrigger()</a>
ADC_TB_SRC	AUXIN0..AUXIN5	POSITIVEEDGE or NEGATIVEEDGE
DAC_START_TRIG	AUXIN0..AUXIN5 DS_CONNECTOR	Not assigned here.
DAC_TB_SRC	AUXIN0..AUXIN5	POSITIVEEDGE or NEGATIVEEDGE
DAC_UPDATE	AUXIN0..AUXIN5 DS_CONNECTOR	POSITIVEEDGE or NEGATIVEEDGE
SYNC_CLK	DS_CONNECTOR	Not assigned here.

Valid output (Direction= SIGNAL\_OUT) settings include:

Signal	Connection	Polarity
ADC_CONVERT	AUXOUT0..AUXOUT2 DS_CONNECTOR	INVERTED* or NONINVERTED
ADC_SCANCLK	AUXOUT0..AUXOUT2	
ADC_SCAN_STOP	AUXOUT0..AUXOUT2	
ADC_SSH	AUXOUT0..AUXOUT2	
ADC_STARTSCAN	AUXOUT0..AUXOUT2	
ADC_START_TRIG	AUXOUT0..AUXOUT2 DS_CONNECTOR	
ADC_STOP_TRIG	AUXOUT0..AUXOUT2 DS_CONNECTOR	
CTR1_CLK	AUXOUT0..AUXOUT2	
CTR2_CLK	AUXOUT0..AUXOUT2	
DAC_START_TRIG	AUXOUT0..AUXOUT2 DS_CONNECTOR	
DAC_UPDATE	AUXOUT0..AUXOUT2 DS_CONNECTOR	
DGND	AUXOUT0..AUXOUT2	Not assigned here.
SYNC_CLK	DS_CONNECTOR	Not assigned here.

\* INVERTED is only valid for Auxiliary Output (AUXOUT) connections.

Valid disabled settings (Direction=DISABLED):

Signal	Connection	Polarity
ADC_TB_SRC	Not assigned here.	Not assigned here.
DAC_TB_SRC		
SYNC_CLK		

**Notes:**

- You can view and edit the above timing and control configuration information from *InstaCal*. Open *InstaCal*, click on the board, and press the "Configure..." button or menu item. If the board supports DAQ Sync and Auxiliary Input/Output signal connections, a button labeled "Advanced Timing & Control Configuration" displays. Press that button to open a display for viewing and modifying the above timing and control signals.
- Except for the ADC\_TB\_SRC, DAC\_TB\_SRC and SYNC\_CLK signals, selecting an input signal connection does not necessarily activate it. However, assigning an output signal to a connection does activate the signal upon performing the respective operation. For instance, when running an EXTCLOCK cbAInScan(), ADC\_CONVERT SIGNAL\_IN selects the connection to use as an external clock to pace the A/D conversions; if cbAInScan() is run without setting the EXTCLOCK option, however, the selected connection is not activated and the signal at that connection is ignored. In both cases, the ADC\_CONVERT signal is output via the connection(s) selected for the ADC\_CONVERT SIGNAL\_OUT. Since there are no scan options for enabling the Timebase Source and the SYNC\_CLK, selecting an input for the A/D or D/A Timebase Source, or SYNC\_CLK does activate the input source for the next respective operations.
- Multiple input signals can be mapped to the same AUXINn connection by successive calls to cbSelectSignal; however, only one connection can be mapped to each input signal. If another connection had already been assigned to an input signal, the former selection is de-assigned and the new connection is assigned.
- Only one output signal can be mapped to the same AUXOUTn connection; however, multiple connections can be mapped to the same output signal by successive calls to cbSelectSignal. If an output signal had already been assigned to a connection, then the former output signal is de-assigned and the new output signal is assigned to the connection. Note that there are at most MAX\_CONNECTIONS (=4) connections that can be assigned to each output signal.
- When selecting DS\_CONNECTOR for a signal, only one direction per signal type can be defined at a given time. Attempting to assign both directions of a signal to the DS\_CONNECTOR results in only the latest selection being applied. If the signal type had formerly been assigned an input direction from the DS\_CONNECTOR, assigning the output direction for that signal type results in the input signal being reassigned to its default connection.

Default Input Signal Connections	Input signal	Default connection
	ADC_CONVERT	AUXIN0
	ADC_GATE	AUXIN5
	ADC_START_TRIG	<b>AUXIN1</b>
	ADC_STOP_TRIG	AUXIN2
	DAC_UPDATE	AUXIN3
	DAC_START_TRIG	AUXIN3

- ADC\_TB\_SRC and DAC\_TB\_SRC are intended to synchronize the timebase of the analog input and output pacers across two or more boards. Internal calculations of sampling and update rates assume that the external timebase has the same frequency as its internal clock. Adjust sample rates to compensate for differences in clock frequencies.

For instance, if the external timebase has a frequency of 10 MHz on a board that has a internal clock frequency of 40 MHz, the scan function samples or updates at a rate of about 1/4 the rate entered. However, while compensating for differences in external timebase and internal clock frequency, if the rate entered results in an invalid pacer count, the function returns a BADRATE error.

## cbSetConfig()

Sets a configuration option for a board. The configuration information for all boards is stored in the CB.CFG file. All programs that use the library read this file. You can use this function to override the configuration information stored in the CB.CFG file.

### Function prototype:

C/C++:	<code>int cbSetConfig(int InfoType, int BoardNum, int DevNum, int ConfigItem, int ConfigVal)</code>
Visual Basic:	<code>Function cbSetConfig(ByVal InfoType&amp;, ByVal BoardNum&amp;, ByVal DevNum&amp;, ByVal ConfigItem&amp;, ByVal ConfigVal&amp;) As Long</code>
Delphi:	<code>function cbSetConfig(InfoType:Integer; BoardNum:Integer; DevNum:Integer; ConfigItem:Integer; ConfigVal:Integer):Integer;</code>

### Arguments:

InfoType	The configuration information for each board is grouped into different categories. InfoType specifies which category you want. Set it to one of the constants listed in the "InfoType argument values" section below.
BoardNum	Refers to the board number associated with a board when it was installed. BoardNum may be 0 to 99.
DevNum	Selects a particular device. If InfoType = DIGITALINFO, then DevNum specifies which of the board's digital devices you want to set information on. If InfoType = COUNTERINFO then DevNum specifies which of the board's counter devices you want to set information on.
ConfigItem	Specifies which configuration item you wish to set. Set it in conjunction with the InfoType argument using the table under "ConfigItem argument values" on page 61.
ConfigVal	The value to set the specified configuration item to.

### Returns:

[Error code](#) or 0 if no errors.

### InfoType argument values:

BOARDINFO	General information about a board.
DIGITALINFO	Information about a digital device.
COUNTERINFO	Information about a counter device.
EXPANSIONINFO	Information about an expansion device.
MISCINFO	One of the miscellaneous options for the board.

**ConfigItem argument values:**

InfoType	ConfigItem	Description
BOARDINFO	BIADCSETTLETIME	ADC settling time
	BIADTRIGCOUNT	Trigger count
	BIBASEADR	Base address of the board
	BICALOUTPUT	Sets the voltage for the CAL pin on supported USB devices.
	BICLOCK	Clock frequency in MHz (1, 4, 6 or 10)
	BIDACSTARTUP	Sets the board's configuration register STARTUP bit to 0 or 1 to enable/disable the storing of digital-to-analog converter (DAC) startup values. Each time the board is powered up, the stored values are written to the DACs. Refer to the " <a href="#">Notes</a> " section below for more information.
	BIDACUPDATECMD	Updates all analog output channels. When ConfigItem is set to BIDACUPDATECMD, the DevNum and ConfigVal arguments are not used and can be set to 0. Refer to the " <a href="#">Notes</a> " section below for more information.
	BIDACUPDATEMODE	Sets the update mode for a digital-to-analog converter (DAC). Use this setting in conjunction with one of these ConfigVal settings: UPDATEIMMEDIATE UPDATEONCOMMAND Refer to the " <a href="#">Notes</a> " section below for more information.
	BIDIDEBOUNCESTATE	The state of the digital inputs when debounce timing is set.
	BIDIDEBOUNCETIME	Sets the debounce time of digital inputs
	BIDMACHAN	DMA channel
	BIINTLEVEL	Interrupt level
	BIHIDELOGINDLG	Enables or disables the Device Login dialog. Set to a nonzero value to disable the dialog. When disabled, the cbDeviceLogin() function must be used to log in to a device session.
	BINETIOTIMEOUT	Sets the amount of time (in milliseconds) to wait for a WEB device to acknowledge a command or query sent to the device over a network connection. If no acknowledgement is received in this time a timeout occurs.
	BINUMADCHANS	Number of A/D channels
	BIPANID	Sets the Personal Area Network (PAN) identifier of a USB device that supports wireless communication.
	BIRANGE	Selected voltage range
	BIRFCHANNEL	Sets the RF channel number used to transmit/receive data by a USB device that supports wireless communication.
	BIRSS	The received signal strength in dBm of a remote device.
	BISRCADPACER	Outputs the A/D pacer signal to the SYNC pin on supported USB devices.
	BIWAITSTATE	Sets the Wait State jumper
EXPANSIONINFO	XIMUXADCHAN1	A/D channel board is connect to
	XIMUXADCHAN2	2nd A/D channel board is connected to
	XIRANGE1	Range (gain) of low 16 channels
	XIRANGE2	Range (gain) of high 16 channels
	XICJCCHAN	A/D channel that CJC is connected to
	XITHERMTYPE	Thermocouple type

**Notes:**

Use the BIDACSTARTUP option (ConfigItem argument) to store either the current DAC values, or the DAC values before the board was last powered down.

- To store the current DAC values as start-up values, call `cbSetConfig()` with a value of 1 for the `BIDACSTARTUP` value. Then, call [cbAOut\(\)](#) or [cbAOutScan\(\)](#) for each channel (), and call `cbSetConfig()` again with a value of 0 for the `BIDACSTARTUP` value.

**Example:**

```
cbSetConfig(BOARDINFO, boardNumber, 0, BIDACSTARTUP, 1);
    for (int i =1; i <8; i++)
    {
        cbAOut(boardNumber, i, BIP5VOLTS, DACValue[i]);
    }
cbSetConfig(BOARDINFO, boardNumber, 0, BIDACSTARTUP, 0);
```

- To store the DAC's last settings, call `cbSetConfig()` with a `BIDACSTARTUP` value of 1. Leave this bit turned on until the application exits. The next time the board is powered up, it restores the values last written to the DACs.

Use the `BIDACUPDATEMODE` option (`ConfigItem` argument) to set the update mode for a DAC board.

- With `ConfigItem` set to `BIDACUPDATEMODE`, and `ConfigVal` set to 0, the DAC update mode is *immediate*. Values written with [cbAOut\(\)](#) or [cbAOutScan\(\)](#) are automatically output by the DAC channels
- With `ConfigItem` set to `BIDACUPDATEMODE` and `ConfigVal` set to 1, the DAC update mode is *on command*. Values written with `cbAOut()` or `cbAOutScan()` are not output by the DAC channels until another `cbSetConfig()` call is made with `ConfigItem` set to `BIDACUPDATECMD`.

## cbSetConfigString()

Sets the configuration or device information as a null-terminated string.

### Function prototype:

**C/C++:**                    `int cbSetConfigString(int InfoType, int BoardNum, int ItemIndex, int ConfigItem, char *ConfigVal, int* maxConfigLen)`

**Visual Basic:**           `Function cbSetConfigString(ByVal InfoType, ByVal BoardNum&, ByVal ItemIndex&, ByVal ConfigItem&, ByVal ConfigVal$, ByRef maxConfigLen&) As Long`

**Delphi:**                   `function cbSetConfigString(InfoType:Integer; BoardNum:Integer; ItemIndex:Integer; ConfigItem:Integer; ConfigVal:PChar; var maxConfigLen:Integer):Integer;`

### Arguments:

InfoType	The configuration information for each board is grouped into different categories. This argument specifies which category you want. Always set this argument to BOARDINFO.
BoardNum	Refers to the board number associated with a board when it was installed. BoardNum may be 0 to 99.
ItemIndex	The location in the device memory (specified by ConfigItem) at which to start writing.
ConfigItem	The type of information (or memory area) to write to the device. Set it to one of the constants listed in the "ConfigItem argument values" section below.
ConfigVal	Pointer to the user-allocated buffer containing the string to copy to the device's memory.
maxConfigLen	Pointer to the value specifying the number of bytes to be written to the device from ConfigVal.

### Returns:

[Error code](#) or 0 if no errors.

maxConfigLen              The number of bytes written to the device.

### ConfigItem argument values:

BIDEVNOTES	Writes up to maxConfigLen characters/bytes from the ConfigVal buffer to the device notes memory, beginning at the location defined by ItemIndex. Currently supported only for WLS Series devices.
BINODEID	Writes up to maxConfigLen characters/bytes from the ConfigVal buffer to the string identifier memory on the device. Note that ItemIndex is not used for this ConfigItem.

## cbSetTrigger()

Selects the trigger source and sets up its parameters. This trigger is used to initiate analog to digital conversions using the following Universal Library functions:

- [cbAInScan\(\)](#), if the EXTTRIGGER option is selected.
- [cbDInScan\(\)](#), if the EXTTRIGGER option is selected.
- [cbCInScan\(\)](#), if the EXTTRIGGER option is selected.
- [cbAPretrig\(\)](#)
- [cbFilePretrig\(\)](#)

### Function prototype:

C/C++:	<code>int cbSetTrigger(int BoardNum, int TrigType, unsigned short LowThreshold, unsigned short HighThreshold);</code>
Visual Basic:	<code>Function cbSetTrigger(ByVal BoardNum&amp;, ByVal TrigType&amp;, ByVal LowThreshold%, ByVal HighThreshold%) As Long</code>
Delphi:	<code>Function cbSetTrigger(BoardNum:Integer; TrigType:Integer; LowThreshold:Word; HighThreshold:Word):Integer;</code>

### Arguments:

BoardNum	Specifies the board number associated with the board when it was installed with the configuration program. The board must have the software selectable triggering source and/or options. BoardNum may be 0 to 99.
TrigType	Specifies the type of triggering based on the external trigger source. Set it to one of the constants in the "TrigType argument values" section on page 65.
LowThreshold	Selects the low threshold used when the trigger input is analog. The range depends upon the resolution of the trigger circuitry. Must be 0 to 255 for 8-bit trigger circuits, 0 to 4095 for 12-bit trigger circuits, and 0 to 65535 for 16-bit trigger circuits. Refer to the " <a href="#">Notes</a> " section on page 65.
HighThreshold	Selects the high threshold used when the trigger input is analog. The range depends upon the resolution of the trigger circuitry. Must be 0 to 255 for 8-bit trigger circuits, 0 to 4095 for 12-bit trigger circuits, and 0 to 65535 for 16-bit trigger circuits. Refer to the " <a href="#">Notes</a> " section on page 65.

### Returns:

[Error code](#) or 0 if no errors.



**TrigType argument values:**

Trigger Source	TrigType	Explanation
Analog	GATE_NEG_HYS	AD conversions are enabled when the external analog trigger input is more positive than HighThreshold. AD conversions are disabled when the external analog trigger input more negative than LowThreshold. Hysteresis is the level between LowThreshold and HighThreshold.
	GATE_POS_HYS	AD conversions are enabled when the external analog trigger input is more negative than LowThreshold. AD conversions are disabled when the external analog trigger input is more positive than HighThreshold. Hysteresis is the level between LowThreshold and HighThreshold.
	GATE_ABOVE	AD conversions are enabled as long as the external analog trigger input is more positive than HighThreshold.
	GATE_BELOW	AD conversions are enabled as long as the external analog trigger input is more negative than LowThreshold.
	TRIG_ABOVE	AD conversions are enabled when the external analog trigger makes a transition from below HighThreshold to above. Once conversions are enabled, the external trigger is ignored.
	TRIG_BELOW	AD conversions are enabled when the external analog trigger input makes a transition from above LowThreshold to below. Once conversions are enabled, the external trigger is ignored.
	GATE_IN_WINDOW	AD conversions are enabled as long as the external analog trigger is inside the region defined by LowThreshold and HighThreshold.
	GATE_OUT_WINDOW	AD conversions are enabled as long as the external analog trigger is outside the region defined by LowThreshold and HighThreshold.
Digital	GATE_HIGH	AD conversions are enabled as long as the external digital trigger input is 5 V (logic HIGH or 1).
	GATE_LOW	AD conversions are enabled as long as the external digital trigger input is 0 V (logic LOW or 0).
	TRIG_HIGH	AD conversions are enabled when the external digital trigger is 5 V (logic HIGH or 1). Once conversions are enabled, the external trigger is ignored.
	TRIG_LOW	AD conversions are enabled when the external digital trigger is 0 V (logic LOW or 0). Once conversions are enabled, the external trigger is ignored.
	TRIG_POS_EDGE	AD conversions are enabled when the external digital trigger makes a transition from 0 V to 5 V (logic LOW to HIGH). Once conversions are enabled, the external trigger is ignored.
	TRIG_NEG_EDGE	AD conversions are enabled when the external digital trigger makes a transition from 5 V to 0 V (logic HIGH to LOW). Once conversions are enabled, the external trigger is ignored.

**Notes:**

The threshold value must be within the range of the analog trigger circuit associated with the board. Refer to the board-specific information in the *Universal Library User's Guide*. For example, on the PCI-DAS1602/16, the analog trigger circuit handles  $\pm 10$  V. A value of 0 corresponds to -10 V, whereas a value of 65535 corresponds to +10 V.

Since Visual Basic does not support unsigned integer types, the thresholds range from -32768 to 32767 for 16-bit boards, instead of from 0 to 65535. In this case, the unsigned value of 65535 corresponds to a value of -1, 65534 corresponds to -2, ..., 32768 corresponds to -32768.

For most boards that support analog triggering, you can pass the required trigger voltage level and the appropriate Range to cbFromEngUnits/FromEngUnits to calculate the HighThreshold and LowThreshold values.

For some boards (refer to the "Analog Input Boards" chapter in the Universal Library User's Guide), you must manually calculate the threshold by first calculating the least significant bit (LSB) for a particular range for the trigger resolution of your hardware. You then use the LSB to find the threshold in counts based on an analog voltage trigger threshold.

To calculate the threshold, do the following:

1. Calculate the LSB by dividing the full scale range (FSR) by  $2^{\text{resolution}}$ . FSR is the entire span from  $-FS$  to  $+FS$  of your hardware for a particular range. For example, the full scale range of  $\pm 10$  V is 20 V.
2. Calculate how many times you need to add the LSB calculate in step 1 to the negative full scale ( $-FS$ ) to reach the trigger threshold value.

The maximum threshold value is  $2^{\text{resolution}} - 1$ . The formula is shown here:

$$\text{Abs}(-FS - \text{threshold in volts}) \div (\text{LSB}) = \text{threshold in counts}$$

Here are two examples that use this formula—one for 8-bit trigger resolution and one for 12-bit trigger resolution.

- 8-bit example using the  $\pm 10$  V range with a  $-5$  V threshold:

**Calculate LSB:**  $\text{LSB} = 20 \div 2^8 = 20 \div 256 = .078125$

**Calculate threshold:**  $\text{Abs}(-10 - (-5)) \div .078125 = 5 \div .078125 = 64$  (round this result if it is not an integer). A count of 64 translates to a voltage threshold of  $-5.0$  V.

- 12-bit example using the  $\pm 10$  V range with a  $+1$  V threshold:

**Calculate LSB:**  $\text{LSB} = 20 \div 2^{12} = 20 \div 4096 = .00488$

**Calculate threshold:**  $\text{Abs}(-10 - 1) \div .00488 = 11 \div .00488 = 2254$  (rounded from 2254.1). A count of 2254 translates to a voltage threshold of  $0.99952$  V.

---

# Counter Functions

## Introduction

This section covers Universal Library functions that load, read, and configure counters. There are five types of counter chips used in MCC counter boards: 8254s, 8536s, 7266s, 9513s and generic event counters. Some of the counter commands only apply to one type of counter.

## cbC7266Config()

Configures 7266 counter for desired operation. This function can only be used with boards that contain a 7266 counter chip (Quadrature Encoder boards). For more information, refer to the LS7266R1 data sheet in the accompanying ls7266r1.pdf file located in the *Documents* subdirectory where you installed UL (C:\Program files\Measurement Computing\DAQ by default).

This data sheet is also available on our web site at [www.mccdaq.com/PDFmanuals/LS7266R1.pdf](http://www.mccdaq.com/PDFmanuals/LS7266R1.pdf)

### Function prototype:

C/C++:	<pre>int cbC7266Config(int BoardNum, int CounterNum, int Quadrature, int CountingMode, int DataEncoding, int IndexMode, int InvertIndex, int FlagPins, int Gating)</pre>
Visual Basic:	<pre>Function cbC7266Config(ByVal BoardNum&amp;, ByVal CounterNum&amp;, ByVal Quadrature&amp;, ByVal CountingMode&amp;, ByVal DataEncoding&amp;, ByVal IndexMode&amp;, ByVal InvertIndex&amp;, ByVal FlagPins&amp;, ByVal Gating&amp;) As Long</pre>
Delphi:	<pre>function cbC7266Config(BoardNum:Integer; CounterNum:Integer; Quadrature:Integer; CountingMode:Integer; DataEncoding:Integer; IndexMode:Integer; InvertIndex:Integer; FlagPins:Integer; Gating:Integer):Integer;</pre>

### Arguments:

BoardNum	Refers to the board number associated with the board when it was installed with the <i>InstaCal</i> configuration program. The specified board must have an LS7266 counter. BoardNum may be 0 to 99.
CounterNum	Counter Number (1 - n), where n is the number of counters on the board.
Quadrature	Selects the resolution multiplier for quadrature input, or disables quadrature input (NO_QUAD) so that the counters can be used as standard TTL counters. NO_QUAD, X1_QUAD, X2_QUAD, or X4_QUAD.
CountingMode	Selects operating mode for the counter. NORMAL_MODE, RANGE_LIMIT, NO_RECYCLE, MODULO_N. Set it to one of the constants in the "CountingMode argument values" section on page 69.
DataEncoding	Selects the format of the data that is returned by the counter - either Binary or BCD format. BCD_ENCODING, BINARY_ENCODING.
IndexMode	Selects which action will be taken when the Index signal is received. The IndexMode must be set to INDEX_DISABLED whenever a Quadrature is set to NON_QUAD or when Gate is set to ENABLED. Set it to one of the constants in the "IndexMode argument values" section on page 69.
InvertIndex	Selects the polarity of the Index signal. If set to DISABLED the Index signal is assumed to be positive polarity. If set to ENABLED the Index signal is assumed to be negative polarity.
FlagPins	Selects which signals will be routed to the FLG1 and FLG2 pins. Set it to one of the constants in the "FlagPins argument values" section on page 69.
Gating	If gating is set to ENABLED, then the channel INDEX input is routed to the RCNTR pin of the LS7266 chip, and is used as a gating signal for the counter. Whenever Gating = ENABLED the IndexMode must be set to INDEX_DISABLED.

### Returns:

[Error code](#) or 0 if no error occurs

**CountingMode argument values:**

NORMAL_MODE	Each counter operates as a 24-bit counter that rolls over to 0 when the maximum count is reached.
RANGE_LIMIT	In range limit count mode, an upper and lower limit is set, mimicking limit switches in the mechanical counterpart. The upper limit is set by loading the PRESET register with the <a href="#">cbCLoad()</a> function after the counter has been configured. The lower limit is always 0. When counting up, the counter freezes whenever the count reaches the value that was loaded into the PRESET register. When counting down, the counter freezes at 0. In either case the counting is resumed only when the count direction is reversed.
NO_RECYCLE	In non-recycle mode the counter is disabled whenever a count overflow or underflow takes place. The counter is re-enabled when a reset or load operation is performed on the counter.
MODULO_N	In modulo-n mode, an upper limit is set by loading the PRESET register with a maximum count. Whenever counting up, when the maximum count is reached, the counter will roll-over to 0 and continue counting up. Likewise when counting down, whenever the count reaches 0, it will roll over to the maximum count (in the PRESET register) and continue counting down.

**IndexMode argument values:**

INDEX_DISABLED	The Index signal is ignored.
LOAD_CTR	The channel INDEX input is routed to the LCNTR pin of the LS7266 counter chip. The counter is loaded whenever the signal occurs.
LOAD_OUT_LATCH	The channel INDEX input is routed to the LCNTR pin of the LS7266 counter chip. The current count is latched whenever the signal occurs. When this mode is selected, the <a href="#">cbCIn()</a> function will return the same count value each time it is called until the Index signal occurs.
RESET_CTR	The channel INDEX input is routed to the RCNTR pin of the LS7266 counter chip. The counter is reset whenever the signal occurs.

**FlagPins argument values:**

CARRY_BORROW	FLG1 pin is CARRY output, FLG2 is BORROW output.
COMPARE_BORROW	FLG1 pin is COMPARE output, FLG2 is BORROW output.
CARRYBORROW_UPDOWN	FLG1 pin is CARRY/BORROW output, FLG2 is UP/DOWN signal.
INDEX_ERROR	FLG1 pin is INDEX output, FLG2 is error output.

## cbC8254Config()

Configures 8254 counter for desired operation. This function can only be used with 8254 counters. For more information, refer to the 82C54 data sheet in the accompanying 82C54.pdf file located in the *Documents* subdirectory where you installed UL (C:\Program Files\Measurement Computing\DAQ by default). This data sheet is also available on our web site at [www.mccdaq.com/PDFmanuals/82C54.pdf](http://www.mccdaq.com/PDFmanuals/82C54.pdf)

### Function prototype:

C/C++:	int cbC8254Config(int BoardNum, int CounterNum, int Config)
Visual Basic:	Function cbC8254Config(ByVal BoardNum&, ByVal CounterNum&, ByVal Config&) As Long
Delphi:	function cbC8254Config(BoardNum:Integer; CounterNum:Integer; Config:Integer):Integer;

### Arguments:

BoardNum	Refers to the number associated with the board when it was installed with the <i>InstaCal</i> configuration program. Board must have an 82C54 installed. BoardNum may be 0 to 99.
CounterNum	Selects one of the counter channels. An 8254 has 3 counters. The value may be 1 - n, where n is the number of 8254 counters on the board (see board-specific information in the <i>Universal Library User's Guide</i> ).
Config	Refer to the 8254 data sheet for a detailed description of each of the configurations. Set it to one of the constants in the "Config argument values" section below.

### Returns:

[Error code](#) or 0 if no errors

### Config argument values:

HIGHONLASTCOUNT	Output of counter (OUT N) transitions from low to high on terminal count and remains high until reset. See Mode 0 on 8254 data sheet in accompanying 82C54.pdf file located in the <i>Documents</i> subdirectory where you installed UL (C:\Program Files\Measurement Computing\DAQ by default).
ONESHOT	Output of counter (OUT N) transitions from high to low on rising edge of GATE N, then back to high on terminal count. See mode 1 on 8254 data sheet in the 82C54.pdf file located in the <i>Documents</i> subdirectory where you installed UL (C:\Program Files\Measurement Computing\DAQ by default).
RATEGENERATOR	Output of counter (OUT N) pulses low for one clock cycle on terminal count, reloads counter and recycles. See mode 2 on 8254 data sheet in the 82C54.pdf file in the <i>Documents</i> subdirectory where you installed UL (C:\Program Files\Measurement Computing\DAQ by default).
SQUAREWAVE	Output of counter (OUT N) is high for count < 1/2 terminal count then low until terminal count, whereupon it recycles. This mode generates a square wave. See mode 3 on the 8254 data sheet in accompanying 82C54.pdf file located in the <i>Documents</i> subdirectory where you installed UL.
SOFTWARESTROBE	Output of counter (OUT N) pulses low for one clock cycle on terminal count. Count starts after counter is loaded. See mode 4 on 8254 data sheet in the 82C54.pdf file located in the <i>Documents</i> subdirectory where you installed UL (C:\Program Files\Measurement Computing\DAQ by default).
HARDWARESTROBE	Output of counter (OUT N) pulses low for one clock cycle on terminal count. Count starts on rising edge at GATE N input. See mode 5 on 8254 data sheet in accompanying 82C54.pdf file located in the <i>Documents</i> subdirectory where you installed UL (C:\Program Files\Measurement Computing\DAQ by default).

## cbC8536Config()

Configures 8536 counter for desired operation. This function can only be used with 8536 counters. For more information, refer to the *Zilog 8536* product specification. The document is available on our web site at [www.mccdaq.com/PDFmanuals/Z8536.pdf](http://www.mccdaq.com/PDFmanuals/Z8536.pdf).

### Function prototype:

**C/C++:**                    `int cbC8536Config(int BoardNum, int CounterNum, int OutputControl, int RecycleMode, int TrigType)`

**Visual Basic:**           `Function cbC8536Config(ByVal BoardNum&, ByVal CounterNum&, ByVal OutputControl&, ByVal RecycleMode&, ByVal TrigType&) As Long`

**Delphi:**                    `function cbC8536Config(BoardNum:Integer; CounterNum:Integer; OutputControl:Integer; RecycleMode:Integer; TrigType:Integer):Integer;`

### Arguments:

BoardNum	Refers to the board number associated with the board when it was installed with the <i>InstaCal</i> configuration program. The board must have an 8536. BoardNum may be 0 to 99.
CounterNum	Selects one of the counter channels. An 8536 has 3 counters. The value may be 1, 2 or 3. INT32 Series boards have two chips installed, so the CounterNum value may be 1 to 6.
OutputControl	Specifies the action of the output signal. Set it to one of the constants in the "OutputControl argument values" section below.
RecycleMode	If set to RECYCLE (as opposed to ONETIME), the counter automatically reloads to the starting count every time it reaches 0, then counting continues.
TrigType	Specifies the trigger type. Set it to one of the constants in the "TrigType argument values" section below.

### Returns:

[Error code](#) or 0 if no errors

OutputControl argument values:

HIGHPULSEONTC	Output transitions from low to high for one clock pulse on the terminal count.
TOGGLEONTC	Output changes state on the terminal count.
HIGHUNTILTC	Output transition to high at the start of counting, then goes low on the terminal count.

### TrigType argument values:

HW_START_TRIG	The first trigger on the counter's trigger input initiates loading of the initial count. Counting proceeds from the initial count.
HW_RETRIG	Every trigger on the counter's trigger input initiates loading of the initial count. Counting proceeds from the initial count.
SW_START_TRIG	The <a href="#">cbCLoad()</a> function initiates loading of the initial count. Counting proceeds from the initial count.

## cbC9513Config()

Sets all of the configurable options of a 9513 counter. For more information, refer to the AM9513A data sheet in the 9513A.pdf file located in the *Documents* subdirectory where you installed UL (C:\Program files\Measurement Computing\DAQ by default).

The data sheet is also available on our web site at [www.mccdaq.com/PDFmanuals/9513A.pdf](http://www.mccdaq.com/PDFmanuals/9513A.pdf)

### Function prototype:

**C/C++:**                   int cbC9513Config(int BoardNum, int CounterNum, int GateControl, int CounterEdge, int CountSource, int SpecialGate, int Reload, int RecycleMode, int BCDMode, int CountDirection, int OutputControl);

**Visual Basic:**           Function cbC9513Config(ByVal BoardNum&, ByVal CounterNum&, ByVal GateControl&, ByVal CounterEdge&, ByVal CountSource&, ByVal SpecialGate&, ByVal Reload&, ByVal RecycleMode&, ByVal BCDMode&, ByVal CountDirection&, ByVal OutputControl&) As Long

**Delphi:**                   function cbC9513Config(BoardNum:Integer; CounterNum:Integer; GateControl:Integer; CounterEdge:Integer; CountSource:Integer; SpecialGate:Integer; Reload:Integer; RecycleMode:Integer; BCDMode:Integer; CountDirection:Integer; OutputControl:Integer):Integer;

### Arguments:

BoardNum	Refers to the board number associated with the board when it was installed with the <i>InstaCal</i> configuration program. The specified board must have a 9513 counter. BoardNum may be 0 to 99.
CounterNum	Counter number (1 - n) where n is the number of counters on the board. For example, a CIO-CTR5 has 5, a CIO-CTR10 has 10, etc. See board specific info.
GateControl	Sets the gating response for level, edge, etc. Set it to one of the constants in the "GateControl argument values" section on page 73.
CounterEdge	Which edge to count. Referred to as "Source Edge" in 9513 data book. Can be set to POSITIVEEDGE (count on rising edge) or NEGATIVEEDGE (count on falling edge).
CountSource	Each counter may be set to count from one of 16 internal or external sources. Set it to one of the constants in the "CountSource argument values" section on page 73.
SpecialGate	Special gate may be enabled or disabled (CBENABLED or CBDISABLED: in Visual Basic or Delphi).
Reload	Reload the counter from the load register (Reload = LOADREG) or alternately load from the load register, then the hold register (Reload = LOADANDHOLDREG).
RecycleMode	Execute once (RecycleMode = ONETIME) or reload and recycle (RecycleMode = RECYCLE).
BCDMode	Counter may operate in <i>binary coded decimal</i> count (ENABLED) or <i>binary</i> count (DISABLED) (CBENABLED or CBDISABLED in Visual Basic or Delphi).
CountDirection	AM9513 may count up (COUNTUP) or down (COUNTDOWN).
OutputControl	The type of output desired. Set it to one of the constants in the "OutputControl argument values" section on page 73.

### Returns:

[Error code](#) or 0 if no errors



**GateControl argument values:**

NOGATE	No gating
AHLTCPREVCTR	Active high TCN -1
AHLNEXTGATE	Active High Level GATE N + 1
AHLPREVGATE	Active High Level GATE N – 1
AHLGATE	Active High Level GATE N
ALLGATE	Active Low Level GATE N
AHEGATE	Active High Edge GATE N
ALEGATE	Active Low Edge GATE N

**CountSource argument values:**

TCPREVCTR	TCN - 1 (Terminal count of previous counter)
CTRINPUT1	SRC 1 (Counter Input 1)
CTRINPUT2	SRC 2 (Counter Input 2)
CTRINPUT3	SRC 3 (Counter Input 3)
CTRINPUT4	SRC 4 (Counter Input 4)
CTRINPUT5	SRC 5 (Counter Input 5)
GATE1	GATE1
GATE2	GATE2
GATE3	GATE3
GATE4	GATE4
GATE5	GATE 5
FREQ1	F1
FREQ2	F2
FREQ3	F3
FREQ4	F4
FREQ5	F5
ALWAYSLOW	Inactive, Output Low

**OutputControl argument values:**

HIGHPULSEONTC	High pulse on Terminal Count
TOGGLEONTC	TC Toggled
DISCONNECTED	Inactive, Output High Impedance
LOWPULSEONTC	Active Low Terminal Count Pulse
3, 6, 7	(numeric values) Illegal

**Notes:**

The information provided here and in the [cbC9513Init\(\)](#) data sheet will only help you understand how Universal Library syntax corresponds to information in the 9513 data sheet. It is not a substitute for the data sheet. You cannot program and use a 9513 without this data sheet.

Refer to the accompanying 9513A.pdf file located in the *Documents* subdirectory where you installed UL (C:\Program Files\Measurement Computing\DAQ by default). The data sheet is also available on our web site at [www.mccdaq.com/PDFmanuals/9513A.pdf](http://www.mccdaq.com/PDFmanuals/9513A.pdf).

## cbC8536Init()

Initializes the counter linking features of an 8536 counter chip. Counters 1 and 2 must be linked before enabling the counters.

Refer to the *Zilog 8536* product specification for a description of the hardware affected by this mode. The document is available on our web site at [www.mccdaq.com/PDFmanuals/Z8536.pdf](http://www.mccdaq.com/PDFmanuals/Z8536.pdf).

### Function prototype:

C/C++:	<code>int cbC8536Init(int BoardNum, int ChipNum, int CtrlOutput)</code>
Visual Basic:	<code>Function cbC8536Init(ByVal BoardNum&amp;, ByVal ChipNum&amp;, ByVal CtrlOutput&amp;) As Long</code>
Delphi:	<code>function cbC8536Init(BoardNum:Integer; ChipNum:Integer; CtrlOutput:Integer):Integer;</code>

### Arguments:

BoardNum	Refers to the board number associated with the board when it was installed with the <i>InstaCal</i> configuration program. The specified board must have an 8536. BoardNum may be 0 to 99.
ChipNum	Selects one of the 8536 chips on the board, 1 to <i>n</i> .
CtrlOutput	Specifies how counter 1 is to be linked to counter 2, if at all. Set it to one of the constants in the "CtrlOutput argument values" section below.

### Returns:

[Error code](#) or 0 if no errors.

### CtrlOutput argument values:

NOTLINKED	Counter 1 is not connected to any other counters inputs.
GATECTR2	Output of counter 1 is connected to the GATE of counter #2.
TRIGCTR2	Output of counter 1 is connected to the trigger of counter #2.
INCTR2	Output of counter 1 is connected to counter #2 clock input.

## cbC9513Init()

Initializes all of the chip level features of a 9513 counter chip. This function can only be used with 9513 counters. For more information, refer to the AM9513A data sheet in the 9513A.pdf file located in the *Documents* subdirectory where you installed UL (C:\Program Files\Measurement Computing\DAQ by default).

This data sheet is also available on our web site at [www.mccdaq.com/PDFmanuals/9513A.pdf](http://www.mccdaq.com/PDFmanuals/9513A.pdf).

### Function prototype:

C/C++:	<code>int cbC9513Init(int BoardNum, int ChipNum, int FOutDivider, int FOutSource, int Compare1, int Compare2, int TimeOfDay)</code>
Visual Basic:	<code>Function cbC9513Init(ByVal BoardNum&amp;, ByVal ChipNum&amp;, ByVal FOutDivider&amp;, ByVal FOutSource&amp;, ByVal Compare1&amp;, ByVal Compare2&amp;, ByVal TimeOfDay&amp;) As Long</code>
Delphi:	<code>function cbC9513Init(BoardNum:Integer; ChipNum:Integer; FOutDivider:Integer; FOutSource:Integer; Compare1:Integer; Compare2:Integer; TimeOfDay:Integer):Integer;</code>

### Arguments:

BoardNum	Refers to the board number associated with the board when it was installed with the <i>InstaCal</i> configuration program. The specified board must have a 9513 counter. BoardNum may be 0 to 99.
ChipNum	Specifies which 9513 chip is to be initialized. For a CTR05 board this should be set to 1. For a CTR10 board it should be either 1 or 2, and for a CTR20 it should be 1-4.
FOutDivider	F-Out divider (0-15). If set to 0, FoutDivider is the rate of FoutSource divided by 16. If set to a number between 1 and 15, FoutDivider is the rate of FoutSource divided by FoutDivider.
FOutSource	Specifies source of the signal for F-Out signal. Set it to one of the constants in the "FOutSource argument values" section on page 76.
Compare1	Compare1 ENABLED or Compare1 DISABLED (CBENABLED or CBDISABLED in Visual Basic or Delphi).
Compare2	Compare2 ENABLED or Compare2 DISABLED. (CBENABLED or CBDISABLED in Visual Basic or Delphi).
TimeOfDay	TimeOfDay ENABLED or TimeOfDay DISABLED. (CBENABLED or CBDISABLED in Visual Basic or Delphi). The options for this argument are listed in the "TimeOfDay argument values" section on page 76.

### Returns:

[Error code](#) or 0 if no errors

**FOutSource argument values:**

FOutSource	9513 Data Sheet Equivalent
CTRINPUT1	SRC 1 (Counter Input 1)
CTRINPUT2	SRC 2 (Counter Input 2)
CTRINPUT3	SRC 3 (Counter Input 3)
CTRINPUT4	SRC 4 (Counter Input 4)
CTRINPUT5	SRC 5 (Counter Input 5)
GATE1	GATE1
GATE2	GATE2
GATE3	GATE3
GATE4	GATE4
GATE5	GATE5
FREQ1	F1
FREQ2	F2
FREQ3	F3
FREQ4	F4
FREQ5	F5

**TimeOfDay argument values:**

TimeOfDay	9513 Data Sheet Equivalent
CBDISABLED	TOD Disabled
1	TOD Enabled / 5 Input
2	TOD Enabled / 6 Input
3	TOD Enabled / 10 Input
No arguments for:	9513 data sheet equivalent
0 (FOUT on)	FOUT Gate
0 (Data bus matches board)	Data Bus Width
1 (Disable Increment)	Data Pointer Control
1 (BCD Scaling)	Scalar Control

**Notes:**

The information provided here and in [cbC9513Config\(\)](#) will help you understand how the Universal Library syntax corresponds to the 9513 data sheet, but is not a substitute for the data sheet. You cannot program and use a 9513 without this data sheet.

Refer to the accompanying 9513A.pdf file located in the *Documents* subdirectory where you installed UL (C:\Program files\Measurement Computing\DAQ by default). The data sheet is also available on our web site at [www.mccdaq.com/PDFmanuals/9513A.pdf](http://www.mccdaq.com/PDFmanuals/9513A.pdf).

## cbCClear()

Clears a scan counter value (sets it to zero). This function only works with counter boards that have counter scan capability.

### Function prototype:

C/C++:	<code>int cbCClear(int BoardNum, int CounterNum)</code>
Visual Basic:	<code>Function cbCClear(ByVal BoardNum&amp;, ByVal CounterNum&amp;) As Long</code>
Delphi:	<code>function cbCClear(BoardNum:Integer; CounterNum:Integer):Integer;</code>

### Arguments:

BoardNum	The board number associated with the board when it was installed with the <i>InstaCal</i> configuration program. The specified board must have a counter. BoardNum may be 0 to 99.
CounterNum	The counter to clear. <b>Note:</b> This argument is zero-based (the first counter number to clear is "0").

### Returns:

[Error code](#) or 0 if no errors

## cbCConfigScan()

Configures a counter channel. This function only works with counter boards that have counter scan capability.

### Function prototype:

**C/C++:**                    `int cbCConfigScan(int BoardNum, short CounterNum, int Mode, int DebounceTime, int DebounceTrigger, int EdgeDetection, int TickSize, int MapCounter)`

**Visual Basic:**           `Function cbCConfigScan(ByVal BoardNum&, ByVal CounterNum&, ByVal Mode&, ByVal DebounceTime&, ByVal DebounceTrigger&, ByVal EdgeDetection&, ByVal TickSize&, ByVal MapCounter&) As Long`

**Delphi:**                   `function cbCConfigScan(BoardNum:Integer; CounterNum:SmallInt; Mode:Integer; DebounceTime:Integer; DebounceTrigger:Integer; EdgeDetection:Integer; TickSize:Integer; MapCounter:Integer):Integer;`

### Arguments:

**BoardNum**                    The board number used to collect the data. BoardNum may be 0 to 99. Refers to the number associated with the board used to collect the data when it was installed with the *InstaCal* configuration program. The specified board must have a scan counter.

**CounterNum**                The counter to set up.

**Note:** This argument is zero-based (the first counter number to set up is "0").

**Mode**                        Bit fields that control various options. This field may contain any combination of non-contradictory choices in the "Mode argument values" section on page 79.

**DebounceTime**              Used to bypass the debounce mode, or to set a channel's comparator output to one of 16 debounce times. Debounce is used to eliminate switch-induced transients typically associated with electro-mechanical devices including relays, proximity switches, and encoders. The choices are:

```
CTR_DEBOUNCE500ns
CTR_DEBOUNCE1500ns
CTR_DEBOUNCE3500ns
CTR_DEBOUNCE7500ns
CTR_DEBOUNCE15500ns
CTR_DEBOUNCE31500ns
CTR_DEBOUNCE63500ns
CTR_DEBOUNCE127500ns
CTR_DEBOUNCE100us
CTR_DEBOUNCE300us
CTR_DEBOUNCE700us
CTR_DEBOUNCE1500us
CTR_DEBOUNCE3100us
CTR_DEBOUNCE6300us
CTR_DEBOUNCE12700us
CTR_DEBOUNCE25500us
CTR_DEBOUNCE_NONE
```

DebounceMode	<p>Sets the mode of the debounce module to CTR_TRIGGER_AFTER_STABLE or to CTR_TRIGGER_BEFORE_STABLE.</p> <p>CTR_TRIGGER_AFTER_STABLE: This mode rejects glitches, and only passes state transitions after a specified period of stability (the debounce time). This mode is used with electro-mechanical devices like encoders and mechanical switches to reject switch bounce and disturbances due to a vibrating encoder that is not otherwise moving. The debounce time should be set short enough to accept the desired input pulse but longer than the period of the undesired disturbance.</p> <p>CTR_TRIGGER_BEFORE_STABLE: Use this mode when the input signal has groups of glitches and each group is to be counted as one. The trigger before stable mode will recognize and count the first glitch within a group but reject the subsequent glitches within the group if the debounce time is set accordingly. In this case the debounce time should be set to encompass one entire group of glitches.</p>
EdgeDetection	<p>Selects whether to detect rising edge or falling edge. Choices are: CTR_RISING_EDGE and CTR_FALLING_EDGE.</p> <p>If a counter is configured for CTR_FALLING_EDGE, calling cbCIn() or cbCIn32() for that counter will result in a BADCOUNTERMODE error.</p>
TickSize	Reserved.
MapCounter	Used to select the mapped counter. A mapped counter is one of the counter input channels other than CounterNum that can participate with the input signal of the counter defined by CounterNum by gating the counter or decrementing the counter.

**Returns:**

[Error code](#) or 0 if no errors

**Mode argument values:**

CLEAR_ON_READ	The counter counts up and is cleared at the beginning of every sample. By default, the counter counts up and only clears the counter at the start of a new scan command.
STOP_AT_MAX	The counter will stop at the top of its count. For the cbCIn32() function, the top of the count depends on whether the BIT_32 option is used. If it is, the top of the count is FFFFFFFF hex. If not, the top of the count is FFFF hex. By default, the counter counts upward and rolls over on the 32-bit boundary.
DECREMENT_ON	<p>Allows the mapped channel to decrement the counter. With this option, the main counter will increment the counter, and the mapped counter can be used to decrement the counter. By default, the counter decrement option is set to "off."</p> <p>This mode is not compatible with cbCIn() or cbCIn32(). If a counter is configured for DECREMENT_ON, calling cbCIn() or cbCIn32() for that counter will result in a BADCOUNTERMODE error.</p>
GATING_ON	<p>Selects gating "on." When "on", the counter is enabled when the mapped channel to gate the counter is high. When the mapped counter is low, the counter is disabled but holds the count value. By default, the counter gating option is set to "off."</p> <p>This mode is not compatible with cbCIn() or cbCIn32(). If a counter is configured for GATING_ON, calling cbCIn() or cbCIn32() for that counter will result in a BADCOUNTERMODE error.</p>

LATCH_ON_MAP	<p>Causes the count to be latched by the signal on the mapped counter. By default, the count is latched by the internal "start of scan" signal, so the count is updated each time it's read.</p> <p>This mode is not compatible with <code>cbCIn()</code> or <code>cbCIn32()</code>. If a counter is configured for <code>LATCH_ON_MAP</code>, calling <code>cbCIn()</code> or <code>cbCIn32()</code> for that counter will result in a <code>BADCOUNTERMODE</code> error.</p>
BIT_32	<p>Selects a 32-bit counter. This mode affects only <code>cbCIn32()</code> and <code>cbCIn()</code>, and only when the counter is configured for <code>STOP_AT_MAX</code>. Recommended for use only with <code>cbCIn32()</code>. (Using the <code>BIT_32</code> option with <code>cbCIn()</code> is not very useful, since the value returned by <code>cbCIn()</code> is only 16 bits. The effect is that the value returned by <code>cbCIn()</code> rolls over at 64k 65,535 times before stopping.)</p>
ENCODER	<p>Sets the specified counter to encoder mode.</p>
ENCODER_MODE_X1	<p>Sets the encoder measurement mode to X1.</p>
ENCODER_MODE_X2	<p>Sets the encoder measurement mode to X2.</p>
ENCODER_MODE_X4	<p>Sets the encoder measurement mode to X4.</p>
LATCH_ON_Z	<p>Selects the Encoder Z mapped signal to latch the counter outputs. This allows the user to know the exact counter value when an edge is present on another counter.</p>
CLEAR_ON_Z_ON	<p>Selects "clear on Z" on. The counter is cleared on the rising edge of the mapped (Z) counter. By default, the "clear on Z" option is off, and the counter is not cleared.</p>



## cbCFreqIn()

Measures the frequency of a signal. This function is only used with 9513 counters. This function uses internal counters #4 and #5.

### Function prototype:

**C/C++:**                    `int cbCFreqIn(int BoardNum, int SigSource, int GateInterval, unsigned short *Count, long *Freq)`

**Visual Basic:**           `Function cbCFreqIn(ByVal BoardNum&, ByVal SigSource&, ByVal GateInterval&, Count%, Freq&) As Long`

**Delphi:**                    `function cbCFreqIn(BoardNum:Integer; SigSource:Integer; GateInterval:Integer; var Count:Word; var Freq:Longint):Integer;`

### Arguments:

BoardNum	The board number associated with the board when it was installed with the <i>InstaCal</i> configuration program. The specified board must have a 9513 counter. BoardNum may be 0 to 99.
SigSource	Specifies the source of the signal from which the frequency is calculated. The signal to be measured is routed internally from the source specified by SigSource to the clock input of counter 5. On boards with more than one 9513 chip, there is more than one counter 5. Which counter 5 is used is also determined by SigSource. Set it to one of the constants in the "SigSource argument values" section on page 82.  The value of SigSource determines which chip will be used. CTRINPUT6 through CTRINPUT10, FREQ6 through FREQ10 and GATE6 through GATE9 indicate chip two will be used. The signal to be measured must be present at the chip two input specified by SigSource. Also, the gating connection from counter 4 output to counter 5 gate must be made between counters 4 and 5 of this chip (see below). Refer to board-specific information to determine valid values for your board.
GateInterval	Gating interval in milliseconds (must be > 0). Specifies the time (in milliseconds) that the counter will be counting. The optimum GateInterval depends on the frequency of the measured signal. The counter can count up to 65535. If the gating interval is too low, the count will be too low and the resolution of the frequency measurement will be poor. For example, if the count changes from 1 to 2, the measured frequency doubles. If the gating interval is too long, then the counter overflows and a FREQOVERFLOW error occurs.  The cbCFreqIn function does not return until the GateInterval has expired. There is no background option. Under Windows, this means that window activity will stop for the duration of the call. Adjust the GateInterval so this does not pose a problem to your user interface.
Count	The raw count is returned here.
Freq	The measured frequency in Hz is returned here.

### Returns:

[Error code](#) or 0 if no errors.

Count - Count that frequency calculation based on returned here.

Freq - Measured frequency in Hz returned here.

**SigSource argument values:**

One 9513 chip (Chip 1 used):

- CTRINPUT1 through CTRINPUT5
- GATE1 through GATE4
- FREQ1 through FREQ5

Two 9513 chips (Chip 1 or Chip 2 used):

- CTRINPUT1 through CTRINPUT10
- GATE1 through GATE9 (excluding gate 5)
- FREQ1 through FREQ10

Four 9513 chips (Chips 1- 4 may be used):

- CTRINPUT1 through CTRINPUT20
- GATE1 through GATE19 (excluding gates 5, 10 & 15)
- FREQ1 through FREQ20

**Notes:**

- This function requires an electrical connection between counter 4 output and counter 5 gate. This connection must be made between counters 4 and 5 *on the chip determined by SigSource*.
- [cbC9513Init\(\)](#) must be called for each ChipNum that will be used by this function. The values of FOutDivider, FOutSource, Compare1, Compare2, and TimeOfDay are irrelevant to this function and may be any value shown in the cbC9513Init() function description.
- If you select an external clock source for the counters, the GateInterval, Count, and Freq settings are only valid if the external source is 1 MHz. Otherwise, you need to scale the values according to the frequency of the external clock source. For example, for an external clock source of 2 MHz, increase your GateInterval setting by a factor of 2, and also double the Count and Freq values returned when analyzing your results.

## cbCIn()

Reads the current count from a counter channel.

### Function prototype:

C/C++:	<code>int cbCIn(int BoardNum, int CounterNum, unsigned short *Count)</code>
Visual Basic:	<code>Function cbCIn (ByVal BoardNum&amp;, ByVal CounterNum&amp;, Count%) As Long</code>
Delphi:	<code>function cbCIn(BoardNum:Integer; CounterNum:Integer; var Count:Word):Integer;</code>

### Arguments:

BoardNum	The board number associated with the board when it was installed with the <i>InstaCal</i> configuration program. The specified board must have a counter. BoardNum may be 0 to 99.
CounterNum	The counter to read the current count from. Valid values are 1 to 20, up to the number of counters on the board.
Count	Counter value returned here. See the " <a href="#">Notes</a> " section below.

### Returns:

[Error code](#) or 0 if no errors.

### Notes:

**Count** - The range of counter values returned are: 0 to 65,535 for C or PASCAL languages. Refer to your BASIC manual for information on BASIC integer data types. -32,768 to 32,767 for BASIC languages. BASIC reads counters as:

- 65535 reads as -1
- 32768 reads as -32768
- 32767 reads as 32767
- 2 reads as 2
- 0 reads as 0

**cbCIn() vs. cbCIn32():** Although the `cbCIn()` and [cbCIn32\(\)](#) functions perform the same operation, `cbCIn32()` is the preferred function to use.

The only difference between the two is that `cbCIn()` returns a 16-bit count value and `cbCIn32()` returns a 32-bit value. Both `cbCIn()` and `cbCIn32()` can be used, but `cbCIn32()` is required whenever you need to read count values greater than 16-bits (counts > 65535).

## cbCIn32()

Reads the current count from a counter and returns it as a 32-bit integer.

### Function prototype:

C/C++:	<code>int cbCIn32(int BoardNum, int CounterNum, unsigned long *Count)</code>
Visual Basic:	<code>Function cbCIn32(ByVal BoardNum&amp;, ByVal CounterNum&amp;, Count&amp;) As Long</code>
Delphi:	<code>function cbCIn32(BoardNum:Integer; CounterNum:Integer; var Count:Longint):Integer;</code>

### Arguments:

BoardNum	The board number associated with the board when it was installed with the <i>InstaCal</i> configuration program. The specified board must have a counter. BoardNum may be 0 to 99.
CounterNum	The counter to read current count from. Valid values are 1 to N, where N is the number of counters on the board.
Count	Current count value from selected counter is returned here.

### Returns:

[Error code](#) or 0 if no error occurs.

### Notes:

**cbCIn() vs. cbCIn32():** Although the [cbCIn\(\)](#) and `cbCIn32()` functions perform the same operation, `cbCIn32()` is the preferred function to use.

The only difference between the two is that `cbCIn()` returns a 16-bit count value and `cbCIn32()` returns a 32-bit value. Both `cbCIn()` and `cbCIn32()` can be used, but `cbCIn32()` is required whenever you need to read count values greater than 16-bits (counts > 65535).

## cbCInScan()

Scans a range of scan counter channels, and stores the samples in an array.

### Function prototype:

C/C++:	<code>int cbCInScan(int BoardNum, int FirstCtr, int LastCtr, long Count, long *Rate, int MemHandle, int Options)</code>
Visual Basic:	<code>Function cbCInScan(ByVal BoardNum&amp;, ByVal FirstCtr&amp;, ByVal LastCtr&amp;, ByVal Count&amp;, Rate&amp;, ByVal MemHandle&amp;, ByVal Options&amp;) As Long</code>
Delphi:	<code>function cbCInScan(BoardNum:Integer; FirstCtr:Integer; LastCtr:Integer; Count:Longint; var Rate:Longint; MemHandle:Integer; Options:Integer):Integer;</code>

### Arguments:

BoardNum	The board number used to collect the data. BoardNum may be 0 to 99. Refers to the number associated with the board used to collect the data when it was installed with the <i>InstaCal</i> configuration program. The specified board must have a counter with scan capability.
FirstCtr	First counter channel of the scan.  This argument is zero-based, so the first counter number is "0".
LastCtr	Last counter channel of the scan.  This argument is zero-based, so the first counter number is "0".  The maximum allowable channel for both FirstCtr and LastCtr depends on how many scan counters are available on the Measurement Computing device in use.
Count	The total number of counter samples to collect. If more than one channel is being sampled then the number of samples collected per channel is equal to $\text{Count} / (\text{LastCtr} - \text{FirstCtr} + 1)$ .
Rate	The rate at which samples are taken in samples per second.  Rate also returns the value of the actual rate set, which may be different from the requested rate because of pacer limitations.
MemHandle	The handle for the Windows buffer to store data. This buffer must have been previously allocated with the <a href="#">cbWinBufAlloc32()</a> function.
Options	Bit fields that control various options. This field may contain any combination of non-contradictory choices in the "Options argument values" section below.

### Returns:

[Error code](#) or 0 if no errors

Rate – the actual sampling rate used.

MemHandle – the collected counter data returned via the Windows buffer.

### Options argument values:

BACKGROUND	When the BACKGROUND option is used, control returns immediately to the next line in your program and the data collection from the counters into the buffer continues in the background. If the BACKGROUND option is not used, the <code>cbCInScan()</code> function does not return to your program until all of the requested data has been collected and returned to the buffer.
------------	--

	<p>Use <code>cbGetStatus()</code> with <code>CTRFUNCTION</code> to check on the status of the background operation. Use <code>cbStopBackground()</code> with <code>CTRFUNCTION</code> to terminate the background process before it has completed. Execute <code>cbStopBackground()</code> after normal termination of all background functions in order to clear variables and flags.</p>
CONTINUOUS	<p>This option puts the function in an endless loop. Once it collects the required number of samples, it resets to the start of the buffer and begins again. The only way to stop this operation is by using <code>cbStopBackground()</code> with <code>CTRFUNCTION</code>. Normally, you should use this option with <code>BACKGROUND</code> so that your program regains control.</p>
EXTTRIGGER	<p>If this option is specified, sampling does not begin until the trigger condition is met. You can set the trigger condition to rising edge, falling edge, or the level of the digital trigger input with the <code>cbSetTrigger()</code> function. Refer to board-specific information in the <i>UL User's Guide</i>.</p>
EXTCLOCK	<p>If this option is specified, conversions will be controlled by the signal on the external clock input rather than by the internal pacer clock. Each conversion will be triggered on the appropriate edge of the clock input signal (refer to board-specific information in the <i>UL User's Guide</i>). When this option is used the <code>Rate</code> argument is ignored. The sampling rate is dependent on the clock signal. Options for the board will default to a transfer mode that will allow the maximum conversion rate to be attained unless otherwise specified.</p>

## cbCLoad()

Loads the specified counter's LOAD, HOLD, ALARM, COUNT, PRESET or PRESCALER register with a count. When loading a counter with a starting value, it is never loaded directly into the counter's count register. Rather, it is loaded into the load or hold register. From there, the counter, after being enabled, loads the count from the appropriate register, generally on the first valid pulse.

### Function prototype:

C/C++:	<code>int cbCLoad(int BoardNum, int RegNum, unsigned LoadValue)</code>
Visual Basic:	<code>Function cbCLoad(ByVal BoardNum&amp;, ByVal RegNum&amp;, ByVal LoadValue&amp;) As Long</code>
Delphi:	<code>function cbCLoad(BoardNum:Integer; RegNum:Integer; LoadValue:Word):Integer;</code>

### Arguments:

BoardNum	The board number associated with the board when it was installed with the <i>InstaCal</i> configuration program. The specified board must have a counter. BoardNum may be 0 to 99.
----------	--

**RegNum** The register to load the count to. Set it to one of the constants in the "RegNum argument values" section below.

LoadValue	The value to be loaded. Must be between 0 and $2^{\text{resolution}} - 1$ of the counter. For example, a 16-bit counter is $2^{16} - 1$ , or 65,535. Refer to the discussion of Basic signed integers in the "16-bit values using a signed integer data type" section in the "Universal Library Description & Use" chapter of the <i>Universal Library User's Guide</i> (available on our web site at <a href="http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf">www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf</a> ).
-----------	---

### Returns:

[Error code](#) or 0 if no errors.

### RegNum argument values:

LOADREG1 .. 20	Load registers 1 through 20. This may span several chips.
HOLDREG1 .. 20	Hold registers 1 through 20. This may span several chips. (9513 only)
ALARM1CHIP1	Alarm register 1 of the first counter chip. (9513 only)
ALARM2CHIP1	Alarm register 2 of the first counter chip. (9513 only)
ALARM1CHIP2	Alarm register 1 of the second counter chip. (9513 only)
ALARM2CHIP2	Alarm register 2 of the second counter chip. (9513 only)
ALARM1CHIP3	Alarm register 1 of the third counter chip. (9513 only)
ALARM2CHIP3	Alarm register 2 of the third counter chip. (9513 only)
ALARM1CHIP4	Alarm register 1 of the four counter chip. (9513 only)
ALARM2CHIP4	Alarm register 2 of the four counter chip. (9513 only)
COUNT1 .. 4	Current Count (LS7266 only)
PRESET1 .. 4	Preset register (LS7266 only)
PRESCALER1 .. 4	Prescaler register (LS7266 only)

**Notes:**

You cannot load a count-down-only counter with less than 2.

**Counter types:** There are several counter types supported. Please refer to the counter chip's data sheet for the registers that are available.

**cbCLoad() vs. cbCLoad32():** Although the `cbCLoad()` and [cbCLoad32\(\)](#) functions perform the same operation, `cbCLoad32()` is the preferred function to use.

The only difference between the two is that `cbCLoad()` loads a 16-bit count value, and `cbCLoad32()` loads a 32-bit value. The only time you need to use `cbCLoad32()` is to load counts that are larger than 32-bits (counts > 65535).



## cbCLoad32()

Loads the specified counter's COUNT, PRESET, or PRESCALER register with a count.

### Function prototype:

C/C++:	<code>int cbCLoad32(int BoardNum, int RegNum, unsigned long LoadValue)</code>
Visual Basic:	<code>Function cbCLoad32(ByVal BoardNum&amp;, ByVal RegNum&amp;, ByVal LoadValue&amp;) As Long</code>
Delphi:	<code>function cbCLoad32(BoardNum:Integer; RegNum:Integer; LoadValue:Longint):Integer;</code>

### Arguments:

BoardNum	Refers to the board number associated with the board when it was installed with the <i>InstaCal</i> configuration program. BoardNum may be 0 to 99.
RegNum	The register to load the value into. Set it to one of the constants in the "RegNum argument values" section below.

### Returns:

[Error code](#) or 0 if no error occurs.

### RegNum argument values:

LOADREG1 .. 20	Load registers 1 through 20. This may span several chips.
HOLDREG1 .. 20	Hold registers 1 through 20. This may span several chips. (9513 only)
ALARM1CHIP1	Alarm register 1 of the first counter chip. (9513 only)
ALARM2CHIP1	Alarm register 2 of the first counter chip. (9513 only)
ALARM1CHIP2	Alarm register 1 of the second counter chip. (9513 only)
ALARM2CHIP2	Alarm register 2 of the second counter chip. (9513 only)
ALARM1CHIP3	Alarm register 1 of the third counter chip. (9513 only)
ALARM2CHIP3	Alarm register 2 of the third counter chip. (9513 only)
ALARM1CHIP4	Alarm register 1 of the four counter chip. (9513 only)
ALARM2CHIP4	Alarm register 2 of the four counter chip. (9513 only)
COUNT1 .. 4	Current Count (LS7266 only)
PRESET1 .. 4	Preset register (LS7266 only)
PRESCALER1 .. 4	Prescaler register (LS7266 only)

### Notes:

**cbCLoad() vs. cbCLoad32():** Although the `cbCLoad()` and `cbCLoad32()` functions perform the same operation, `cbCLoad32()` is the preferred function to use.

The only difference between the two is that `cbCLoad()` loads a 16-bit count value, and `cbCLoad32()` loads a 32-bit value. The only time you need to use `cbCLoad32()` is to load counts that are larger than 32-bits (counts > 65535).

## cbCStatus()

Returns status information about the specified counter (7266 counters only). For more information, see the LS7261 data sheet in the LS7266R1pdf file located in the *Documents* subdirectory where you installed UL (C:\Program files\Measurement Computing\DAQ by default). This data sheet is also available on our web site at [www.mccdaq.com/PDFmanuals/LS7266R1.pdf](http://www.mccdaq.com/PDFmanuals/LS7266R1.pdf).

### Function prototype:

C/C++:	<code>int cbCStatus(int BoardNum, int CounterNum, unsigned long *StatusBits)</code>
Visual Basic:	<code>Function cbCStatus(ByVal BoardNum&amp;, ByVal CounterNum&amp;, StatusBits&amp;) As Long</code>
Delphi:	<code>function cbCStatus(BoardNum:Integer; CounterNum:Integer; var StatusBits:Longint):Integer;</code>

### Arguments:

BoardNum	The board number associated with the board when it was installed with the <i>InstaCal</i> configuration program. The specified board must have an LS7266 counter. BoardNum may be 0 to 99.
CounterNum	The number of the counter whose status bits you want to read. Valid values are 1 to N, where N is the number of counters on the board.
StatusBits	Current status from selected counter is returned here. The status consists of individual bits that indicate various conditions within the counter. Set it to one of the constants in the "StatusBits argument values" section below.

### Returns:

[Error code](#) or 0 if no error occurs.

### StatusBits argument values:

C_UNDERFLOW	Set to 1 whenever the count decrements past 0. Is cleared to 0 whenever <code>cbCStatus()</code> is called.
C_OVERFLOW	Set to 1 whenever the count increments past it's upper limit. Is cleared to 0 whenever <code>cbCStatus()</code> is called.
C_COMPARE	Set to 1 whenever the count matches the preset register. Is cleared to 0 whenever <code>cbCStatus()</code> is called.
C_SIGN	Set to 1 when the MSB of the count is 1. Is cleared to 0 whenever the MSB of the count is set to 0.
C_ERROR	Set to 1 whenever an error occurs due to excessive noise on the input. Is cleared to 0 by calling <a href="#">cbC7266Config()</a> set to 1 when index is valid. Is cleared to 0 when index is not valid.
C_UP_DOWN	Set to 1 when counting up. Is cleared to 0 when counting down
C_INDEX	Set to 1 when index is valid. Is cleared to 0 when index is not valid.

## cbCStoreOnInt()

### Changed R4.0 RW

Installs an interrupt handler that will store the current count whenever an interrupt occurs. This function can only be used with 9513 counters. This function will continue to operate in the background until either `IntCount` has been satisfied or [cbStopBackground\(\)](#) with `CTRFUNCTION` is called.

### Function prototype:

**C/C++:** `int cbCStoreOnInt(int BoardNum, int IntCount, short CntrControl[], int MemHandle)`

**Visual Basic:** `Function cbCStoreOnInt(ByVal BoardNum&, ByVal IntCount&, CntrControl%, ByVal MemHandle&) As Long`

**Delphi:** `function cbCStoreOnInt(BoardNum:Integer; IntCount:Integer; var CntrControl:SmallInt; MemHandle:Integer):Integer;`

### Arguments:

<code>BoardNum</code>	The board number associated with the board when it was installed with the <i>InstaCal</i> configuration program. The specified board must have a 9513 counter. <code>BoardNum</code> may be 0 to 99.
<code>IntCount</code>	The counters will be read every time an interrupt occurs until <code>IntCount</code> number of interrupts have occurred. If <code>IntCount</code> is = 0 then the function will run until <code>cbStopBackground()</code> is called. (refer to <code>MemHandle</code> ).
<code>CntrControl</code>	The array should have an element for each counter on the board. (5 elements for CTR-05 board, 10 elements for a CTR-10, etc.). Each element corresponds to a counter channel. Each element should be set to either <code>CBDISABLED</code> or <code>CBENABLED</code> . All channels that are set to <code>CBENABLED</code> will be read when an interrupt occurs.
<code>MemHandle</code>	Handle for Windows buffer. If <code>IntCount</code> is non-zero, the buffer referenced by <code>MemHandle</code> must be of sufficient size to hold ( <code>IntCount * Number of Counters</code> ) points.

### Returns:

[Error code](#) or 0 if no errors.

### Notes:

**New functionality:** If the Library Revision is set to 4.0 or greater, the following code changes are required:

- If `IntCount` is non-zero, the buffer referenced by `MemHandle` must be able to hold (`IntCount * Number of Counters`) points.  
For example, if you set `IntCount` to 100 for a CTR-05 board, you must allocate the size of the buffer to be  $(100 * 5) = 500$ . This new functionality keeps the user application from having to move the data out of the buffer for every interrupt, before it is overwritten. Now, for each interrupt, the counter values will be stored in adjacent memory locations in the buffer.

#### Allocate the proper buffer size for non-zero `IntCount` settings

Specifying `IntCount` as a non-zero value and failing to allocate the proper sized buffer results in a runtime error. There is no way for the Universal Library to determine if the buffer has been allocated with the proper size.

- If `IntCount` = 0, the functionality is unchanged.

## cbTimerOutStart()

Starts a timer square wave output. Use `cbTimerOutStop()` to stop the output.

### Function prototype:

C/C++:	<code>int cbTimerOutStart(int BoardNum, int TimerNum, double *Frequency)</code>
Visual Basic:	<code>Function cbTimerOutStart(ByVal BoardNum&amp;, ByVal TimerNum&amp;, Frequency#) As Long</code>
Delphi:	<code>function cbTimerOutStart(BoardNum:Integer; TimerNum:Integer; var Frequency:Double):Integer;</code>

### Arguments:

BoardNum	The board number associated with the board when it was installed with the <i>InstaCal</i> configuration program. The specified board must have a timer-type counter. BoardNum may be 0 to 99.
TimerNum	The timer to output the square wave from. Valid values are zero (0) up to the number of timers – 1 on the board.
Frequency	The desired square wave frequency. The timers clock will be divided down by integer values to produce the frequency. The actual frequency output will be returned. Valid values are dependant on the timer's clock and the timer resolution.

### Returns:

[Error code](#) or 0 if no errors

Frequency – the actual frequency set.

## cbTimerOutStop()

Stops a timer square wave output. Use `cbTimerOutStart()` to start the output.

### Function prototype:

C/C++:	<code>int cbTimerOutStop(int BoardNum, int TimerNum)</code>
Visual Basic:	<code>Function cbTimerOutStop(ByVal BoardNum&amp;, ByVal TimerNum&amp;) As Long</code>
Delphi:	<code>function cbTimerOutStop(BoardNum:Integer; TimerNum:Integer):Integer;</code>

### Arguments:

BoardNum	The number associated with the board when it was installed with the <i>InstaCal</i> configuration program. The specified board must have a timer-type counter. BoardNum may be 0 to 99.
TimerNum	The timer to stop. Valid values are zero up to the number of timers on the board – 1.

### Returns:

[Error code](#) or 0 if no errors

---

# Data Logger Functions

## Introduction

This section covers Universal Library functions that read and convert data logged to a binary file from a data acquisition product equipped with data logging functionality. The data is typically logged to a CompactFlash® memory card which may then be inserted into a media reader for reading and conversion using these functions.

Data is stored in a binary file. The data may consist of analog data, digital I/O data, time stamped data and information about the device configuration. You can use the data logger functions to read this information, apply conversions to the data, and convert the files to a comma separated values (.CSV) text file or other specified text file formats.

## cbLogConvertFile()

Converts a binary log file to a comma-separated values (.CSV) text file or another text file format that you specify.

### Function prototype:

```
C/C++:      int cbLogConvertFile(char* srcFile, char* destFile, int startSample,
              int count, int delimiter)
```

### Arguments:

srcFile	The name and path of the binary file to read.
destFile	The name and destination path of the converted file. Use the file extension of the file type that you want to create.
startSample	The index number of the first sample to read.
count	The number of samples to read.
delimiter	Specifies the character used between fields in the converted file. Set to one of the Delimiter constants. Choices are:  0 = DELIMITER_COMMA 1 = DELIMITER_SEMICOLON 2 = DELIMITER_SPACE 3 = DELIMITER_TAB

### Returns:

[Error code](#) or 0 if no errors.

### Notes:

- Timestamp data is stored according to the `TimeZone` and `TimeFormat` arguments. Refer to [cbLogSetPreferences\(\)](#) on page 111.
- Time stamps in the converted file may be in either 12-hour or 24-hour format based on the setting of the `TimeFormat` argument. Time stamps can optionally be converted to local time based on the setting of the `TimeZone` argument.
- AI temperature data is returned according to the `Units` preference. Refer to [cbLogSetPreferences\(\)](#) on page 111.
- The `Units` preference is only applied to the AI data if the data was logged as temperature data. Refer to [cbLogGetAIInfo\(\)](#) on page 97. This value is ignored if the AI data was logged as raw data.
- The `Units` preference is always applied to CJC data, since it is always logged as temperature data.
- If the `destFile` argument ends with a .CSV extension, the `delimiter` argument must be set to `DELIMITER_COMMA`. Otherwise, an `INVALIDDELIMITER` error is returned.
- You can open a comma-separated values text file (.CSV) directly in Microsoft Excel. Text files with extensions other than .CSV can only be imported into Excel.

## cbLogGetAIChannelCount()

Retrieves the total number of analog input channels logged in a binary file.

### Function prototype:

C/C++: `int cbLogGetAIChannelCount(char* Filename, int* AICount)`

### Arguments:

Filename	The name of the file to retrieve the information from.
AICount	The number of analog input channels logged in the binary file.

### Returns:

[Error code](#) or 0 if no errors.

AICount – Returns the number of analog input channels logged in the binary file.



## cbLogGetAllInfo()

Retrieves the channel number and unit value of each analog input channel logged in a binary file.

### Function prototype:

C/C++:                      int cbLogGetAllInfo(char\* Filename, int\* ChannelNumbers, int\* Units)

### Arguments:

Filename	The name of the file to retrieve the information.
ChannelNumbers	An array that contains the analog input channel numbers logged in the file.
Units	An array that contains the unit values set for the device in <i>InstaCal</i> for each analog input channel logged in the file.

### Returns:

[Error code](#) or 0 if no errors.

ChannelNumbers – Returns the analog input channel numbers logged in the binary file.

Units – Returns the unit values set for the device in *InstaCal* for each analog input channel logged in the binary file. Returned values include:

0 = UNITS\_TEMPERATURE

1 = UNITS\_RAW

## cbLogGetCJCInfo()

Retrieves the number of CJC temperature channels logged in a binary file.

### Function prototype:

C/C++: `int cbLogGetCJCInfo(char* Filename, int* CJCCcount)`

### Arguments:

Filename                      The name of the file to retrieve the information from.

CJCCcount                    The number of CJC temperature channels logged in the file.

### Returns:

[Error code](#) or 0 if no errors.

CJCCcount – Returns the number of CJC channels logged in the binary file.

## cbLogGetDIOInfo()

Retrieves the number of digital I/O channels logged in a binary file.

### Function prototype:

C/C++: `int cbLogGetDIOInfo(char* Filename, int* DIOCount)`

### Arguments:

Filename	The name of the file to retrieve the information from.
DIOCount	The number of digital I/O channels logged in the binary file.

### Returns:

[Error code](#) or 0 if no errors.

DIOCount – Returns the number of digital I/O channels logged in the binary file.

## cbLogGetFileInfo()

Retrieves the version level and byte size of a binary file.

### Function prototype:

```
C/C++:          int cbLogGetFileInfo(char* Filename, int* Version, int* Size)
```

### Arguments:

Filename	The name of the file to retrieve the information from.
Version	The version level of the binary file.
Size	The size in bytes of the binary file.

### Returns:

[Error code](#) or 0 if no errors.

Version – Returns the version level of the binary file.

Size – Returns the size in bytes of the binary file.

## cbLogGetFileName()

Retrieves the name of the  $n^{\text{th}}$  file in the directory containing binary log files.

### Function prototype:

C/C++: `int cbLogGetFileName(int FileNumber, char* Path, char*Filename)`

### Arguments:

FileNumber	Index of the file whose name you want to return. Specify one of the following:  The number ( $n$ ) that represents the location of the file in the directory (where $n = 0$ , 1, 2, and so on), or  GETFIRST – get the first file in the directory, or  GETNEXT – get the next file in the directory, based on the current index.  This parameter is the index of the file in the directory, and is not part of the filename.
Path	The full path to the directory containing the binary file. The path must be NULL terminated and cannot be longer than 256 characters.
Filename	A NULL terminated string containing the full path to the file.

### Returns:

[Error code](#) or 0 if no errors.

Filename – Returns a NULL terminated string containing the full path to the file.

### Notes:

Set FileNumber to GETFIRST to access the first binary file in a directory. Subsequent calls with FileNumber = GETNEXT returns each successive file in the directory. When you call the function after accessing the last file in the directory, the function returns the error code [NOMOREFILES](#).

## cbLogGetPreferences()

Retrieves API preference settings for time stamped data, analog temperature data, and CJC temperature data. Returns the default values unless changed using [cbLogSetPreferences\(\)](#).

### Function prototype:

```
C/C++:      int cbLogGetPreferences(int* TimeFormat, int* TimeZone,
                    int* Units)
```

### Arguments:

TimeFormat	<p>The time format to apply to time stamped data. Set to one of the <code>TimeFormat</code> constants. Choices are:</p> <p>0 = <code>TIMEFORMAT_12HOUR</code> - for example 2:32:51PM.</p> <p>1 = <code>TIMEFORMAT_24HOUR</code> - for example 14:32:51.</p>
TimeZone	<p>The time zone to store time stamped data. Set to one of the <code>TimeZone</code> constants. Choices are:</p> <p>0 = <code>TIMEZONE_LOCAL</code>. Converts time stamped data to the local time zone on your computer.</p> <p>1 = <code>TIMEZONE_GMT</code>. Leaves time stamped data in Greenwich Mean Time.</p>
Units	<p>The unit to use for temperature data. Set to one of the <code>Units</code> constants. Choices are:</p> <p>0 = <code>FAHRENHEIT</code></p> <p>1 = <code>CELSIUS</code></p> <p>2 = <code>KELVIN</code></p>

### Returns:

[Error code](#) or 0 if no errors.

`TimeFormat` – Returns the format to apply to time stamped data from API functions that return time data.

`TimeZone` – Returns the time zone to apply to time stamped data from API functions that return time data.

`Units` – Returns the unit to use when converting temperature data from API functions that return temperature data.

## cbLogGetSampleInfo()

Retrieves the sample interval, sample count, and the date and time of the first data point contained in a binary file.

### Function prototype:

```
C/C++:      int cbLogGetSampleInfo(char* Filename, int* SampleInterval,
              int* SampleCount, int* StartDate, int* StartTime)
```

### Arguments:

Filename	The name of the file to retrieve sample information from.
SampleInterval	The time interval, in seconds, between samples.
SampleCount	The number of samples contained in the file.
StartDate	The date when the first data point was logged in the file. Date values are packed in the following format:  Byte 0:        day Byte 1:        month Byte 2 - 3:    year
StartTime	The time when the first data point was logged in the file. Time values are packed in the following format:  Byte 0:        seconds Byte 1:        minutes Byte 2:        hours Byte 3:        0xff = 24hour format 0x0 = AM 0x1 = PM

### Returns:

[Error code](#) or 0 if no errors.

SampleInterval – Returns the time interval, in seconds, between samples.

SampleCount – Returns the number of samples in the file.

StartDate – Returns the date when the first data point was logged in the file.

StartTime – Returns the time when the first data point was logged in the file.

### Notes:

Time stamped data is returned according to the `TimeZone` and `TimeFormat` preferences. Refer to [cbLogSetPreferences\(\)](#) on page 111.

## cbLogReadAIChannels()

Retrieves analog input data from a binary file, and stores the values in an array.

### Function prototype:

```
C/C++:      int cbLogReadAIChannels(char* Filename, int StartSample,
              int Count, float* AIChannels)
```

### Arguments:

Filename	The name of the file to retrieve the information from.
StartSample	The first sample to read from the binary file.
Count	The number of samples to read from the binary file.
AIChannels	Receives the analog input values.

### Returns:

[Error code](#) or 0 if no errors.

AIChannels – Returns the analog input values logged in the file.

### Notes:

The units of the analog input data that is returned is set by the value of the `Units` preference. Refer to [cbLogSetPreferences\(\)](#) on page 111.


The `units` preference is only applied if the logged data is temperature data. This value is ignored if the data logged is raw.

### Analog array:

The user is responsible for allocating the size of the analog array, and ensuring that it is large enough to hold the data that will be returned. You can calculate the array allocation using the `SampleCount` value from [cbLogGetSampleInfo\(\)](#), and the `AICount` value from [cbLogGetAIChannelCount\(\)](#):

```
float* aiChannels = new float[SampleCount * AICount];
```

The figure below shows the layout of the analog array, and how the elements should be indexed.

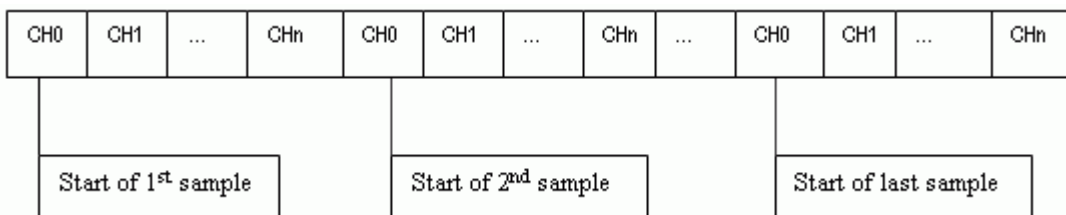

Where `n` is (`AICount - 1`). `CH0 - CHn` refer to the channels in the array, not the input channels of the device.

For example, assume that all of the even number input channels are logged. The analog array channels are mapped as shown here:

Array Channel	Device Input Channel
0	0
1	2
2	4
3	6

Use the following code fragment to access the elements of the analog array:



```
for (i=0; i<numberOfSamples; i++)  
{  
    for (j=0; j<numberOfAIChannels; j++)  
    {  
        a = analogArray[(i *numberOfAIChannels) + j];  
    }  
}
```

where

the `numberOfSamples` is set by the `SampleCount` value from [cbLogGetSampleInfo\(\)](#).

the `numberOfAIChannels` is set by the `AICount` value from [cbLogGetAIChannelCount\(\)](#).

## cbLogReadCJCChannels()

Retrieves CJC temperature data from a binary file, and stores the values in an array.

### Function prototype:

```
C/C++:      int cbLogReadCJCChannels(char* Filename, int StartSample, int Count,
            float* CJCChannels)
```

### Arguments:

Filename	The name of the file to retrieve the information from.
StartSample	The first sample to read from the binary file.
Count	The number of samples to read from the binary file.
CJCChannels	Receives the CJC temperature values.

### Returns:

[Error code](#) or 0 if no errors.

CJCChannels – Returns the CJC temperature values logged in the file.

### Notes:


The temperature scale of the CJC data that is returned is set by the value of the `Units` preference. Refer to [cbLogSetPreferences\(\)](#) on page 111.

### CJC array:

The user is responsible for allocating the size of the CJC array, and ensuring that it is large enough to hold the data that will be returned. You can calculate the array allocation using the `SampleCount` value from [cbLogGetSampleInfo\(\)](#), and the `CJCCount` value from [cbLogGetCJCInfo\(\)](#):

```
float* CJCChannels = new float[SampleCount * CJCCount];
```

The figure below shows the layout of the CJC array, and how the elements should be indexed.

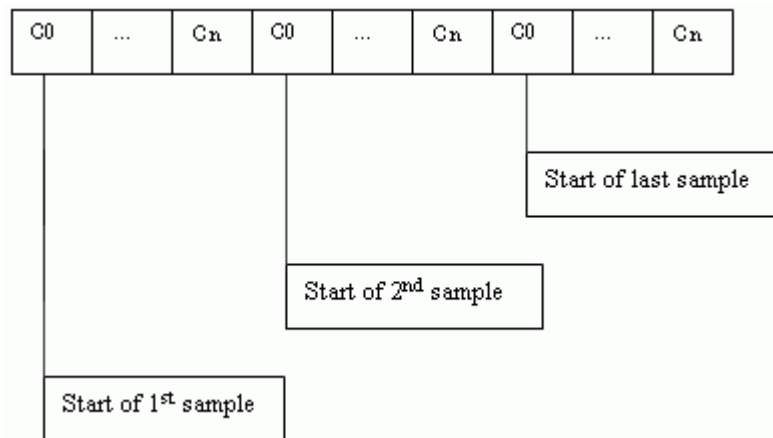

Where  $n$  is  $(CJCCount - 1)$

Use the following code fragment to access the elements of the CJC array.

```
for (i=0; i<numberOfSamples; i++)
{
    for (j=0; j<numberOfCJCChannels; j++)
    {
        c = cjcArray[(i * numberOfCJCChannels) + j];
    }
}
```

where

the `numberOfSamples` is set by the `SampleCount` value from [cbLogGetSampleInfo\(\)](#).

the `numberOfCJCChannels` is set by the `CJCCount` value from [cbLogGetCJCInfo\(\)](#).

## cbLogReadDIOChannels()

Retrieves digital I/O channel data from a binary file, and stores the values in an array.

### Function prototype:

```
C/C++:      int cbLogReadDIOChannels(char* Filename, int StartSample, int Count,
              int* DIOChannels)
```

### Arguments:

Filename	The name of the file to retrieve the information from.
StartSample	The first sample to read from the binary file.
Count	The number of samples to read from the binary file.
DIOChannels	Receives the DIO input values.

### Returns:

[Error code](#) or 0 if no errors.


**DIOChannels** – Returns the DIO channel values logged in the file. Each element of the array contains the value of one bit from a digital channel.

### DIO array:

The user is responsible for allocating the size of the DIO array, and ensuring that it is large enough to hold the data that will be returned. You can calculate the array allocation using the `SampleCount` value from `cbLogGetSampleInfo()` and the `DIOCount` value from `cbLogGetDIOInfo()`:

```
int* DIOChannels = new int[SampleCount * DIOCount];
```

The figure below shows the layout of the DIO array, and how the elements should be indexed.

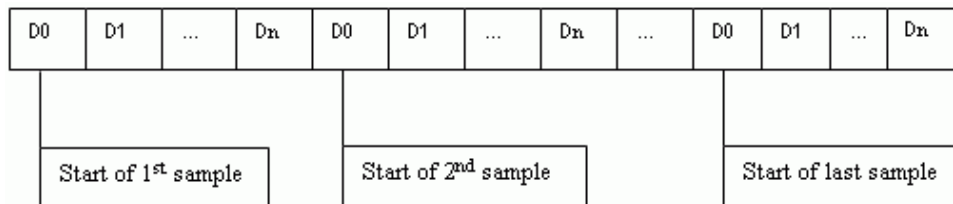

where  $n$  is  $(\text{DIOCount} - 1)$

Use the following code fragment to access the elements of the DIO array:

```
for (i=0; i<numberOfSamples; i++)
{
    for (j=0; j<numberOfDIOChannels; j++)
    {
        d = dioArray[(i * numberOfDIOChannels) + j];
    }
}
```

where

`numberOfSamples` is set by the `SampleCount` value from [cbLogGetSampleInfo\(\)](#)

`numberOfDIOChannels` is set by the `DIOCount` value from [cbLogGetDIOInfo\(\)](#)

## cbLogReadTimeTags()

Retrieves date and time values logged in a binary file. This function stores date values in the DateTags array, and time values in the TimeTags array.

### Function prototype:

```
C/C++:      int cbLogReadTimeTags(char* Filename, int StartSample,
              int Count, int* DateTags, int* TimeTags)
```

### Arguments:

Filename	The name of the file to retrieve the information from.
StartSample	The first sample to read from the binary file.
Count	The number of samples to read from the binary file.
DateTags	Receives the date value for each sample logged in the file. The dates are packed in the following format:  Byte 0:        day Byte 1:        month Byte 2 -3:     year
TimeTags	Receives the time value for each sample logged in the file. The times are packed in the following format:  Byte 0:        seconds Byte 1:        minutes Byte 2:        hours Byte 3:        0xff = 24hour format 0x0 = AM 0x1 = PM

### Returns:

[Error code](#) or 0 if no errors.

DateTags – Returns the date values for each sample logged in the file.

TimeTags – Returns the time values for each sample logged in the file.

### Notes:

Time stamped data is returned according to the `TimeZone` preference value and the `TimeFormat` preference value. Refer to [cbLogSetPreferences\(\)](#) on page 111.

Time stamped data are logged in the file if *InstaCal* is configured to do so. If time stamps are not logged, the `TimeTags` and `DateTags` arrays are filled with values calculated from the file header information.

### Array size:

The user is responsible for allocating the size of the `DateTags` and `TimeTags` arrays, and ensuring that they are large enough to hold the data that is returned. You can calculate the array allocation using the `SampleCount` value from [cbLogGetSampleInfo\(\)](#) on page 102.

```
int*   dates  = new int[SampleCount];
int*   times  = new int[SampleCount];
```

**DateTags array**

The figure below shows the layout of the DateTags array, and how the elements should be indexed.

D0	D1	D2	...	Dn
----	----	----	-----	----

where: n is (SampleCount - 1)

Each sample has only one date. Use the following code fragment to access the elements of the DateTags array:

```
for (i=0; i<numberOfSamples; i++)
{
    d = DateTagsArray[i];
}
```

**TimeTags array**

The figure below shows the layout of the TimeTags array, and how the elements should be indexed.

T0	T1	T2	...	Tn
----	----	----	-----	----

where: n is (SampleCount - 1)

Each sample has only one time stamp. Use the following code fragment to access the elements of the TimeTags array:

```
for (i=0; i<numberOfSamples; i++)
{
    t = TimeTagsArray[i];
}
```

## cbLogSetPreferences()

Sets preferences for returned time stamped data, analog temperature data, and CJC temperature data.

### Function prototype:

C/C++: `int cbLogSetPreferences(int TimeFormat, int TimeZone, int Units)`

### Arguments:

TimeFormat	Specifies the time format to apply when returning time stamped data (when using <a href="#">cbLogReadTimeTags()</a> for example). Set to one of the TimeFormat constants. Choices are:  0 = TIMEFORMAT_12HOUR - for example 2:32:51PM (default). 1 = TIMEFORMAT_24HOUR - for example 14:32:51.
TimeZone	Specifies whether to convert time stamped data that is returned (when using <code>cbLogReadTimeTags()</code> for example) to the local time zone or to return the time stamps as they are stored in the file (in the GMT time zone). Set to one of the TimeZone constants. Choices are:  0 = TIMEZONE_LOCAL. Converts timestamp data to the local time zone on your computer (default). 1 = TIMEZONE_GMT. Leaves time stamped data in Greenwich Mean Time.
Units	Specifies whether to convert temperature data returned (when using <a href="#">cbLogReadAIChannels()</a> for example) to Fahrenheit or Kelvin, or return temperature data as they are stored in the file (in Celsius units).  Set to one of the Units constants. Choices are:  0 = FAHRENHEIT (Default) 1 = CELSIUS 2 = KELVIN  This value is ignored if raw data is logged.

### Returns:

[Error code](#) or 0 if no errors.

### Notes:

- The TimeFormat and TimeZone preferences are applied to all time data returned using API functions that return time data.
- The Units preference specifies the temperature scale that the API applies when reading and converting analog temperature and CJC data.

---

## Digital I/O Functions

### Introduction

Use the functions explained in this chapter to read and set digital values. Most digital ports are configurable, while some others are non-configurable. Some types of hardware allow readback of the values that output ports are set to on configurable port types. Devices using 8255 chips for digital I/O are one example. For these devices, input functions such as `cbDIn()` are valid for ports configured as output.



## cbDBitIn()

Reads the state of a single digital input bit.

This function treats all of the DIO ports of a particular type on a board as a single port. It lets you read the state of any individual bit within this port.

Note that for some port types—such as 8255 ports—if the port is configured for `DIGITALOUT`, this function provides readback of the last output value.

### Function prototype:

C/C++:	<code>int cbDBitIn(int BoardNum, int PortType, int BitNum, unsigned short *BitValue)</code>
Visual Basic:	<code>Function cbDBitIn Lib (ByVal BoardNum&amp;, ByVal PortType&amp;, ByVal BitNum&amp;, BitValue%) As Long</code>
Delphi:	<code>function cbDBitIn(BoardNum:Integer; PortType:Integer; BitNum:Integer; var BitValue:Word):Integer;</code>

### Arguments:

BoardNum	The number associated with the board when it was installed with the <i>InstaCal</i> configuration program. BoardNum may be 0 to 99.
PortType	There are three general types of digital ports—ports that are programmable as input or output, ports that are fixed input or output, and ports for which each bit may be programmed as input or output. For the first of these types, set PortType to <code>FIRSTPORTA</code> . For the latter two types, set PortType to <code>AUXPORT</code> . Some boards have both types of digital ports (DAS1600). Set PortType to either <code>FIRSTPORTA</code> or <code>AUXPORT</code> , depending on which digital inputs you wish to read.
BitNum	Specifies the bit number within the single large port.
BitValue	Place holder for return value of bit. Value will be 0 or 1. A 0 indicates a logic low reading, a 1 indicates a logic high reading. Logic high does not necessarily mean 5V. See the board manual for chip input specifications.

### Returns:

[Error code](#) or 0 if no errors.

BitValue - value (0 or 1) of specified bit returned here.

## cbDBitOut()

Sets the state of a single digital output bit. This function treats all of the DIO ports of a particular type on a board as a single very large port. It lets you set the state of any individual bit within this large port. If the port type is not AUXPORT, you **must** use [cbDConfigPort\(\)](#) to configure the port for output first. If the port type is AUXPORT, you **may** need to use [cbDConfigBit\(\)](#) or [cbDConfigPort\(\)](#) to configure the bit for output first. Refer to the board-specific information in the *Universal Library User's Guide* (available on our web site at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf)) to determine if AUXPORT should be configured for your hardware.

### Function prototype:

C/C++:	<pre>int cbDBitOut(int BoardNum, int PortType, int BitNum, unsigned short BitValue)</pre>
Visual Basic:	<pre>Function cbDBitOut(ByVal BoardNum&amp;, ByVal PortType&amp;, ByVal BitNum&amp;, ByVal BitValue%) As Long</pre>
Delphi:	<pre>function cbDBitOut(BoardNum:Integer; PortType:Integer; BitNum:Integer; BitValue:Word):Integer;</pre>

### Arguments:

BoardNum	The number associated with the board when it was installed with the <i>InstaCal</i> configuration program. BoardNum may be 0 to 99.
PortType	There are three general types of digital ports - ports that are programmable as input or output, ports that are fixed input or output and ports for which each bit may be programmed as input or output. For the first of these types, set PortType to FIRSTPORTA. For the latter two types, set PortType to AUXPORT. Some boards have both types of digital ports (DAS1600). Set PortType to either FIRSTPORTA or AUXPORT depending on which digital port you wish to write to.
BitNum	Specifies the bit number within the single large port. The specified bit must be in a port that is currently configured as an output.
BitValue	The value to set the bit to. Value will be 0 or 1. A 0 indicates a logic low output, a 1 indicates a logic high output. Logic high does not necessarily mean 5V. See the board manual for chip specifications.

### Returns:

[Error code](#) or 0 if no errors.

## cbDConfigBit()

Configures a specific digital bit as Input or Output. This function treats all DIO ports of the AUXPORT type on a board as a single port. This function is NOT supported by 8255 type DIO ports. Refer to the board-specific information for details.

### Function prototype:

C/C++:	<code>int cbDConfigBit(int BoardNum, int PortType, int BitNum, int Direction)</code>
Visual Basic:	<code>Function cbDConfigBit(ByVal BoardNum&amp;, ByVal PortType&amp;, ByVal BitNum&amp;, ByVal Direction&amp;) As Long</code>
Delphi:	<code>function cbDConfigBit(Boardnum:Integer; PortType:Integer; BitNum:Integer; Direction:Integer) :Integer;</code>

### Arguments:

BoardNum	The number associated with the board when it was installed with the <i>InstaCal</i> configuration program. BoardNum may be 0 to 99.
PortType	The port (AUXPORT) whose bits are to be configured. The port specified must be bitwise configurable. Check the board-specific information in the <i>Universal Library User's Guide</i> for details.
BitNum	The bit number to configure as input or output. See board-specific information for details.
Direction	DIGITALOUT or DIGITALIN configures the specified bit for output or input, respectively.

### Returns:

[Error code](#) or 0 if no errors.

## cbDConfigPort()

Configures a digital port as input or output. This function is for use with ports that may be programmed as input or output, such as those on the 82C55 chips and 8536 chips. Refer to the Zilog 8536 manual for details of chip operation. Also refer to the 82C55 data sheet, which is available on our web site at [www.mccdaq.com/PDFmanuals/82C55A.pdf](http://www.mccdaq.com/PDFmanuals/82C55A.pdf).

### Function prototype:

C/C++:	<code>int cbDConfigPort(int BoardNum, int PortNum, int Direction)</code>
Visual Basic:	<code>Function cbDConfigPort (ByVal BoardNum&amp;, ByVal PortNum&amp;, ByVal Direction&amp;) As Long</code>
Delphi:	<code>function cbDConfigPort (Boardnum:Integer; PortNum:Integer; Direction:Integer) :Integer;</code>

### Arguments:

BoardNum	The number associated with the board when it was installed with the <i>InstaCal</i> configuration program. BoardNum may be 0 to 99.
PortNum	The specified port must be configurable. For most boards, AUXPORT is not configurable. Check the board-specific information in the <i>Universal Library User's Guide</i> for details.
Direction	DIGITALOUT or DIGITALIN configures entire eight or four bit port for output or input.

### Returns:

[Error code](#) or 0 if no errors.

### Notes:

When used on ports within an 8255 chip, this function will reset all ports on that chip configured for output to a zero state. This means that if you set an output value on FIRSTPORTA and then change the configuration on FIRSTPORTB from OUTPUT to INPUT, the output value at FIRSTPORTA will be all zeros. You can, however, set the configuration on SECONDPORTX without affecting the value at FIRSTPORTA. For this reason, this function is usually called at the beginning of the program for each port requiring configuration.

## cbDIn()

Reads a digital input port. Note that for some port types, such as 8255 ports, if the port is configured for DIGITALOUT, this function will provide readback of the last output value.

### Function prototype:

C/C++:	<code>int cbDIn(int BoardNum, int PortNum, unsigned short *DataValue)</code>
Visual Basic:	<code>Function cbDIn(ByVal BoardNum&amp;, ByVal PortNum&amp;, DataValue%) As Long</code>
Delphi:	<code>function cbDIn(BoardNum:Integer; PortNum:Integer; var DataValue:Word):Integer; StdCall;</code>

### Arguments:

BoardNum	The number associated with the board when it was installed with the <i>InstaCal</i> configuration program. BoardNum may be 0 to 99.
PortNum	Specifies which digital I/O port to read. Some hardware does allow readback of the state of the output using this function. Check the board-specific information in the <i>Universal Library User's Guide</i> .
DataValue	Digital input value returned here.

### Returns:

[Error code](#) or 0 if no errors.

DataValue - Digital input value returned here.

### Notes:

The size of the ports vary. If it is an eight bit port then the returned value will be in the range 0 - 255. If it is a four bit port the value will be in the range 0 - 15.

Refer to the example programs and the board-specific information contained in the *Universal Library User's Guide* (available on our web site at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf)) for clarification of valid PortNum values.

## cbDInScan()

Multiple reads of digital input port of a high speed digital port on a board with a pacer clock such as the CIO-PDMA16.

### Function prototype:

**C/C++:** `int cbDInScan(int BoardNum, int PortNum, long Count, long *Rate, int MemHandle, int Options)`

**Visual Basic:** `Function cbDInScan(ByVal BoardNum&, ByVal PortNum&, ByVal Count&, Rate&, ByVal MemHandle&, ByVal Options&) As Long`

**Delphi:** `function cbDInScan(BoardNum:Integer; PortNum:Integer; Count:Longint; var Rate:Longint; MemHandle:Integer; Options:Integer):Integer;`

### Arguments:

BoardNum	The number associated with the board when it was installed with the <i>InstaCal</i> configuration program. BoardNum may be 0 to 99.
PortNum	Specifies which digital I/O port to read (usually FIRSTPORTA or FIRSTPORTB).
Count	The number of times to read digital input.
Rate	Number of times per second (Hz) to read the port. The actual sampling rate in some cases will vary a small amount from the requested rate. The actual rate will be returned to the Rate argument.
MemHandle	Handle for Windows buffer to store data. This buffer must have been previously allocated with the <a href="#">cbWinBufAlloc()</a> function.
Options	Bit fields that control various options. Refer to the constants in the "Options argument values" section below.

### Returns:

[Error code](#) or 0 if no errors.

Rate - actual sampling rate returned.

MemHandle - digital input value returned via the allocated Windows buffer.

### Options argument values:

BACKGROUND	<p>If the BACKGROUND option is not used then the <code>cbDInScan()</code> function will not return to your program until all of the requested data has been collected and returned to MemHandle.</p> <p>When the BACKGROUND option is used, control will return immediately to the next line in your program and the transfer from the digital input port to MemHandle will continue in the background. Use <a href="#">cbGetStatus()</a> with DIFUNCTION to check on the status of the background operation. Use <a href="#">cbStopBackground()</a> with DIFUNCTION to terminate the background process before it has completed.</p>
CONTINUOUS	<p>This option puts the function in an endless loop. Once it transfers the required number of bytes it resets to the start of DataBuffer and begins again. The only way to stop this operation is with <a href="#">cbStopBackground()</a> with DIFUNCTION. Normally this option should be used in combination with BACKGROUND so that your program will regain control.</p>

EXTCLOCK	If this option is used then transfers will be controlled by the signal on the trigger input line rather than by the internal pacer clock. Each transfer will be triggered on the appropriate edge of the trigger input signal (refer to board-specific information contained in the <i>Universal Library Users Guide</i> ). When this option is used the <code>Rate</code> argument is ignored. The transfer rate is dependent on the trigger signal.
EXTTRIGGER	If this option is used then the scan will not begin until the signal on the trigger input line meets the trigger criteria.
WORDXFER	Normally this function reads a single (byte) port. If <code>WORDXFER</code> is specified then it will read two adjacent ports on each read and store the value of both ports together as the low and high byte of a single array element in the buffer. When <code>WORDXFER</code> is used, it is generally required to set <code>PortNum</code> to <code>FIRSTPORTA</code> .

**Notes:**

**Transfer method** - May not be specified. DMA is used.

## cbDOut()

Writes a byte to a digital output port. If the port type is not AUXPORT, you **must** use [cbDConfigPort\(\)](#) to configure the port for output first. If the port type is AUXPORT, you **may** need to use [cbDConfigPort\(\)](#) to configure the port for output first. Check the board specific information in the *Universal Library User's Guide* (available on our web site at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf)) to determine if AUXPORT should be configured for your hardware.

### Function prototype:

C/C++:	<code>int cbDOut(int BoardNum, int PortNum, unsigned short DataValue)</code>
Visual Basic:	<code>Function cbDOut(ByVal BoardNum&amp;, ByVal PortNum&amp;, ByVal DataValue%) As Long</code>
Delphi:	<code>function cbDOut(BoardNum:Integer; PortNum:Integer; DataValue:Word):Integer;</code>

### Arguments:

BoardNum	The number associated with the board when it was installed with the <i>InstaCal</i> configuration program. BoardNum may be 0 to 99.
PortNum	There are three general types of digital ports - ports that are programmable as input or output, ports that are fixed input or output and ports for which each bit may be programmed as input or output. For the first of these types, set PortNum to FIRSTPORTA. For the latter two types, set PortNum to AUXPORT. Some boards have both types of digital ports (for example the DAS1600 Series). Set PortNum to either FIRSTPORTA or AUXPORT depending on the digital port you want to set.
DataValue	Digital input value to be written.

### Returns:

[Error code](#) or 0 if no errors.

### Notes:

The size of the ports vary. If it is an eight bit port then the output value should be in the range 0 – 255. If it is a four bit port the value should be in the range 0 – 15. Refer to the example programs and the board-specific information in the *Universal Library User's Guide* for clarification of valid PortNum values.



## cbDOutScan()

Performs multiple writes to a digital output port of a high speed digital port on a board with a pacer clock, such as the CIO-PDMA16 or CIO-PMA32.

### Function prototype:

**C/C++:**                    `int cbDOutScan(int BoardNum, int PortNum, long Count, long *Rate, int MemHandle, int Options)`

**Visual Basic:**           `Function cbDOutScan(ByVal BoardNum&, ByVal PortNum&, ByVal Count&, Rate&, ByVal MemHandle&, ByVal Options&) As Long`

**Delphi:**                    `function cbDOutScan(BoardNum:Integer; PortNum:Integer; Count:Longint; var Rate:Longint; MemHandle:Integer; Options:Integer):Integer;`

### Arguments:

BoardNum	The number associated with the board when it was installed with the <i>InstaCal</i> configuration program. BoardNum may be 0 to 99.
PortNum	Specifies which digital I/O port to write (usually FIRSTPORTA or FIRSTPORTB). The specified port must be configured as an output.
Count	The number of times to write digital output.
Rate	Number of times per second (Hz) to write to the port. The actual update rate in some cases will vary a small amount from the requested rate. The actual rate will be returned to the Rate argument.
MemHandle	Handle for Windows buffer to store data. This buffer must have been previously allocated with the <a href="#">cbWinBufAlloc()</a> function.
Options	Bit fields that control various options. Refer to the constants in the "Options argument values" section below.

### Returns:

[Error code](#) or 0 if no errors.

Rate - actual sampling rate returned.

### Options argument values:

BACKGROUND	<p>If the BACKGROUND option is not used then the cbDOutScan() function will not return to your program until all of the requested data has been output.</p> <p>When the BACKGROUND option is used, control returns immediately to the next line in your program and the transfer to the digital output port from MemHandle will continue in the background. Use <a href="#">cbGetStatus()</a> with DOFUNCTION to check on the status of the background operation. Use <a href="#">cbStopBackground()</a> with DOFUNCTION to terminate the background process before it has completed.</p>
CONTINUOUS	<p>This option puts the function in an endless loop. Once it transfers the required number of bytes it resets to the start of the buffer and begins again. The only way to stop this operation is with <a href="#">cbStopBackground()</a> with DOFUNCTION. Normally this option should be used in combination with BACKGROUND so that your program will regain control.</p>
EXTCLOCK	<p>If this option is used then transfers will be controlled by the signal on the trigger input line rather than by the internal pacer clock. Each transfer will be triggered on the appropriate edge of the trigger input signal (refer to board-specific information contained in the <i>Universal Library Users Guide</i>). When this option is used the Rate argument is ignored. The transfer rate is dependent on the trigger signal.</p>

WORDXFER	Normally this function writes a single (byte) port. If WORDXFER is specified then it will write two adjacent ports as the low and high byte of a single array element in the buffer. When WORDXFER is used, it is generally required to set PortNum to FIRSTPORTA.
NONSTREAMEDIO	<p>When this option is used, you can output non-streamed data to a specific DAC output channel.</p> <p>To load the data output buffer into the device's internal output FIFO, the aggregate size of the data output buffer must be <math>\leq</math> the size of the internal data output FIFO in the device. Once the sample data are transferred or downloaded to the device, the device is responsible for outputting the data. You can't make any changes to the output buffer once the output begins.</p> <p>With NONSTREAMEDIO mode, you do not have to periodically feed output data through the program to the device for the data output to continue. However, the size of the buffer is limited.</p>
ADCCLOCKTRIG	Triggers a data output operation when the ADC clock starts.
ADCCLOCK	Paces the data output operation using the ADC clock.

**Notes:**

- BYTEXFER is the default option. Make sure you are using an array when your data is arranged in bytes. Use the WORDXFER option for word array transfers.
- NONSTREAMEDIO can only be used with the number of samples (Count) set equal to the size of the FIFO or less.
- Transfer method may not be specified. DMA is used.

---

# Error Handling Functions

## Introduction

Use the functions explained in this chapter to get information from error codes returned by other UL functions. Most library functions return error codes. The different methods built in to the functions for handling errors include stopping the program when an error occurs, and printing error messages versus error codes.

## cbErrHandling()

Sets the error handling for all subsequent function calls. Most functions return error codes after each call. In addition, other error handling features are built into the library. This function controls those features. If the Universal Library cannot find the configuration file CB.CFG, it always terminates the program, regardless of the `cbErrHandling()` setting.

### Function prototype:

C/C++:	<code>int cbErrHandling(int ErrReporting, int ErrHandling)</code>
Visual Basic:	<code>Function cbErrHandling(ByVal ErrReporting&amp;, ByVal ErrHandling&amp;) As Long</code>
Delphi:	<code>function cbErrHandling(ErrReporting:Integer; ErrHandling:Integer):Integer;</code>

### Arguments:

<code>ErrReporting</code>	This argument controls when the library will print error messages on the screen. The default is <code>DONTPRINT</code> . Set it to one of the constants in the "ErrReporting argument values" section below.
<code>ErrHandling</code>	This argument specifies what class of error will cause the program to halt. The default is <code>DONTSTOP</code> . Set it to one of the constants in the "ErrHandling argument values" section below.

### Returns:

Always returns 0.

### ErrReporting argument values:

<code>DONTPRINT</code>	Errors will not generate a message to the screen. In that case your program must always check the returned error code after each library call to determine if an error occurred.
<code>PRINTWARNINGS</code>	Only warning errors will generate a message to the screen. Your program will have to check for fatal errors.
<code>PRINTFATAL</code>	Only fatal errors will generate a message to the screen. Your program must check for warning errors.
<code>PRINTALL</code>	All errors will generate a message to the screen.

### ErrHandling argument values:

<code>DONTSTOP</code>	The program will always continue executing when an error occurs.
<code>STOPFATAL</code>	The program will halt if a "fatal" error occurs.
<code>STOPALL</code>	Will stop whenever any error occurs. If you are running in an Integrated Development Environment (IDE) then when errors occur, the environment may be shut down along with the program. If your IDE behaves this way, (QuickBasic and VisualBasic do), then set <code>ErrHandling</code> to <code>DONTSTOP</code> . Refer to " <a href="#">Error Codes</a> " on page 401 for a complete list of error codes and their associated messages.

### Notes:

Warnings vs. Fatal Errors: All errors that can occur are classified as either "warnings" or "fatal":

- Errors that can occur in normal operation in a bug free program (disk is full, too few samples before trigger occurred) are classified as "warnings".
- All other errors indicate a more serious problem and are classified as "fatal".

STOPALL is not intended for 32-bit C console programs: Do not use the the STOPALL option in 32-bit C console applications. Instead, use other methods to end the program, such as checking the return value of the function.

## cbGetErrMsg()

Returns the error message associated with an error code. Each function returns an error code. An error code that is not equal to 0 indicates that an error occurred. Call this function to convert the returned error code to a descriptive error message.

### Function prototype:

C/C++:	<code>int cbGetErrMsg(int ErrCode, char ErrMsg[ERRSTRLEN])</code>
Visual Basic:	<code>Function cbGetErrMsg(ByVal ErrCode&amp;, ByVal ErrMsg\$) As Long</code>
Delphi:	<code>function cbGetErrMsg(ErrCode:Integer; ErrMsg:PChar):Integer;</code>

### Arguments:

<code>ErrCode</code>	Error code that is returned by any function in library.
<code>ErrMsg</code>	Error message returned here. The <code>ErrMsg</code> variable must be pre-allocated to be at least as large as <code>ERRSTRLEN</code> . This size is guaranteed to be large enough to hold the longest error message.

### Returns:

[Error code](#) or 0 if no errors.

\*ErrMsg - error message string is returned here.

### Notes:

See also [cbErrHandling\(\)](#) on page 124 for an alternate method of handling errors.

---

# Memory Board Functions

## Introduction

Use the functions explained in this chapter to read and write data to and from a memory board, and also set modes that control memory boards (MEGA-FIFO).

The most common use for the memory boards is to store large amounts of data from an A/D board via a DT-Connect cable to a memory board. To do this, use the `EXTMEMORY` option with [cbAInScan\(\)](#) or [cbAPretrig\(\)](#). Once the data has been transferred to the memory board, use the memory functions to retrieve it.

## cbMemRead()

Reads data from a memory board into an array.

### Function prototype:

C/C++:	<code>int cbMemRead(int BoardNum, unsigned short DataBuffer[], long FirstPoint, long Count)</code>
Visual Basic:	<code>Function cbMemRead(ByVal BoardNum&amp;, DataBuffer%, ByVal FirstPoint&amp;, ByVal Count&amp;) As Long</code>
Delphi:	<code>function cbMemRead(BoardNum:Integer; var DataBuffer:Word; FirstPoint:Longint; Count:Longint):Integer;</code>

### Arguments:

BoardNum	The number associated with the board when it was installed with the <i>InstaCal</i> configuration program. BoardNum may be 0 to 99.
DataBuffer	Pointer to the data array
FirstPoint	Index of first point to read, or FROMHERE. Use FirstPoint to specify the first point to read. For example, to read data sample numbers 200 through 250, set FirstPoint = 200 and Count = 50.
Count	Number of data points (words) to read

### Returns:

[Error code](#) or 0 if no errors.

DataBuffer - data read from memory board.

### Notes:

When reading a large amount of data from the board in small chunks, set FirstPoint to FROMHERE to read each successive chunk. Using FROMHERE speeds up a cbMemRead() operation when working with large amounts of data.

For example, to read 300,000 points in 100,000 point chunks, the calls would look like this:

```
cbMemRead (0, DataBuffer, 0, 100000)
cbMemRead (0, DataBuffer, FROMHERE, 1000000)
cbMemRead (0, DataBuffer, FROMHERE, 1000000)
```

**DT-Connect Conflicts** - The cbMemRead() function can not be called while a DT-Connect transfer is in progress. For example, if you start collecting A/D data to the memory board in the background (by calling [cbAInScan\(\)](#) with the DTCONNECT + BACKGROUND options) you can not call cbMemRead() until the cbAInScan() has completed. If you do you will get a DTACTIVE error.



## cbMemReadPretrig()

Reads pre-trigger data collected with the [cbAPretrig\(\)](#) function from a memory board, and re-arranges the data in the correct order (pre-trigger data first, then post-trigger data). This function can only be used to retrieve data that was collected with the [cbAPretrig\(\)](#) function with EXTMEMORY set in the options argument. After each [cbAPretrig\(\)](#) call, all data must be unloaded from the memory board with this function. If any more data is sent to the memory board then the pre-trigger data will be lost.

### Function prototype:

**C/C++:**                    `int cbMemReadPretrig(int BoardNum, unsigned short DataBuffer[], long FirstPoint, long Count)`

**Visual Basic:**           `Function cbMemReadPretrig(ByVal BoardNum&, DataBuffer%, ByVal FirstPoint&, ByVal Count&) As Long`

**Delphi:**                    `function cbMemReadPretrig(BoardNum:Integer; var DataBuffer:Word; FirstPoint:Longint; Count:Longint):Integer;`

### Arguments:

BoardNum	The number associated with the board when it was installed with the <i>InstaCal</i> configuration program. BoardNum may be 0 to 99.
DataBuffer	The pointer to the data array.
FirstPoint	Index of first point to read, or FROMHERE. Use FirstPoint to specify the first point to read. For example, to read data sample numbers 200 through 250, set FirstPoint = 200 and Count = 50.
Count	Number of data samples (words) to read

### Returns:

[Error code](#) or 0 if no errors.

DataBuffer - data read from memory board.

### Notes:

When reading a large amount of data from the board in small chunks, set FirstPoint to FROMHERE to read each successive chunk. Using FROMHERE speeds up a [cbMemRead\(\)](#) operation when working with large amounts of data. For example, to read 300,000 points in 100,000 chunks the calls would look like this:

```
cbMemReadPretrig(0, DataBuffer, 0, 100000)
cbMemReadPretrig(0, DataBuffer, FROMHERE, 1000000)
cbMemReadPretrig(0, DataBuffer, FROMHERE, 1000000)
```

**DT-Connect Conflicts** - The [cbMemReadPretrig\(\)](#) function can not be called while a DT-Connect transfer is in progress. For example, if you start collecting A/D data to the memory board in the background (by calling [cbAInScan\(\)](#) with the DTCONNECT + BACKGROUND options), you can not call [cbMemReadPretrig\(\)](#) until the [cbAInScan\(\)](#) has completed. If you do you will get a DTACTIVE error.

## cbMemReset()

Resets the memory board pointer to the start of the data. The memory boards are sequential devices. They contain a counter which points to the 'current' word in memory. Every time a word is read or written this counter increments to the next word.

### Function prototype:

C/C++:	<code>int cbMemReset(int BoardNum)</code>
Visual Basic:	<code>Function cbMemReset(ByVal BoardNum&amp;) As Long</code>
Delphi:	<code>function cbMemReset(BoardNum:Integer):Integer;</code>

### Arguments:

BoardNum	The number associated with the board when it was installed with the <i>InstaCal</i> configuration program. BoardNum may be 0 to 99.
----------	---

### Returns:

[Error code](#) or 0 if no errors.

### Notes:

This function is used to reset the counter back to the start of the memory. Between successive calls to [cbAInScan\(\)](#), you should call this function so that the second `cbAInScan()` overwrites the data from the first call. Otherwise, the data from the first `cbAInScan()` will be followed by the data from the second `cbAInScan()` in the memory on the card.

Likewise, anytime you call [cbMemRead\(\)](#) or [cbMemWrite\(\)](#) it will leave the counter pointing to the next memory location after the data that you read or wrote. Call `cbMemReset()` to reset back to the start of the memory buffer before the next call to `cbAInScan()`.

## cbMemSetDTMode()

Sets the DT-Connect Mode of a memory board.

### Function prototype:

C/C++:	<code>int cbMemSetDTMode(int BoardNum, int Mode)</code>
Visual Basic:	<code>Function cbMemSetDTMode(ByVal BoardNum&amp;, ByVal Mode&amp;) As Long</code>
Delphi:	<code>function cbMemSetDTMode(BoardNum:Integer; Mode:Integer):Integer;</code>

### Arguments:

BoardNum	The number associated with the board when it was installed with the <i>InstaCal</i> configuration program. BoardNum may be 0 to 99.
Mode	Must be set to either DTIN or DTOUT. Set the Mode on the memory board to DTIN to transfer data from an A/D board to the memory board. Set Mode = DTOUT to transfer data from a memory board to a D/A board.

### Returns:

[Error code](#) or 0 if no errors.

### Notes:

- This command only controls the direction of data transfer between the memory board and its parent board that is connected to it via a DT-Connect cable.
- If you are using the EXTMEMORY option, do not use `cbMemSetDTMode()`, as the memory board mode is already set with EXTMEMORY. Only use `cbMemSetDTMode()` when the parent board is not supported by the Universal Library.

## cbMemWrite()

Writes data from an array to the memory card.

### Function prototype:

C/C++:	<code>int cbMemWrite(int BoardNum, unsigned short DataBuffer[], long FirstPoint, long Count);</code>
Visual Basic:	<code>Function cbMemWrite(ByVal BoardNum&amp;, DataBuffer%, ByVal FirstPoint&amp;, ByVal Count&amp;) As Long</code>
Delphi:	<code>function cbMemWrite(BoardNum:Integer; var DataBuffer:Word; FirstPoint:Longint; Count:Longint):Integer;</code>

### Arguments:

BoardNum	The number associated with the board when it was installed with the <i>InstaCal</i> configuration program. BoardNum may be 0 to 99.
DataBuffer	Pointer to the data array.
FirstPoint	Index of first point to write, or FROMHERE. Use FirstPoint to specify the first point to write data to. For example, to write to location numbers 200 through 250, set FirstPoint = 200 and Count = 50.
Count	Number of data points (words) to write.

### Returns:

[Error code](#) or 0 if no errors.

### Notes:

To write a large amount of data to the board in small chunks, set FirstPoint to FROMHERE to write each successive chunk. For example, to write 300,000 points in 100,000 point chunks:

```
cbMemWrite(0, DataBuffer, 0, 100000)
cbMemWrite(0, DataBuffer, FROMHERE, 100000)
cbMemWrite(0, DataBuffer, FROMHERE, 100000)
```

**DT-Connect Conflicts** - The `cbMemWrite()` function cannot be called while a DT-Connect transfer is in progress. For example, if you start collecting A/D data to the memory board in the background (by calling [cbAInScan\(\)](#) with the DTCONNECT + BACKGROUND options). You cannot call `cbMemWrite()` until the `cbAInScan()` is complete. Doing so will generate a DTACTIVE error.

---

## Revision Control Functions

### Introduction

Use the functions explained in this chapter to initialize the Universal Library DLL so that the functions are interpreted according to the format of the revision that you wrote and compiled your program in. As new revisions of the library are released, bugs from previous revisions are fixed and occasionally new functions are added. It is Measurement Computing's goal to preserve the existing programs that you have written, and therefore to never change the order or number of arguments in a function. However, it is not always possible to achieve this goal.

## cbDeclareRevision()

### New R3.3 ID

Initializes the Universal Library with the revision number of the library used to write your program. Must be the first Universal Library function to be called by your program.

### Function prototype:

C/C++:	<code>int cbDeclareRevision(float* RevNum);</code>
Visual Basic:	<code>Function cbDeclareRevision(RevNum!) As Long</code>
Delphi:	<code>Function cbDeclareRevision(var RevNum:single):Integer;</code>

### Arguments:

RevNum	Revision number of the Universal Library to interpret function arguments.
--------	---

**Default setting:** Any program using the 32-bit library and not containing this line of code will be defaulted to revision 5.4 argument assignments.

### Returns:

[Error Code](#) or 0 if no errors.

### Notes:

As new revisions of the library are released, bugs from previous revisions are fixed and occasionally new functions are added. It is Measurement Computing's goal to preserve existing programs you have written, and therefore to never change the order or number of arguments in a function. Sometimes this is not possible, as in the changes from revision 3.2 to 3.3. In revision 3.3, we added support for multiple background tasks, a feature that users have requested.

Allowing multiple background tasks required adding the argument `BoardNum` to several functions. Doing so would have meant that programs written for version 3.2 would not run with 3.3 if they called those functions. If not for the new `cbDeclareRevision()` function, the programs would have had to be rewritten in each line where the affected functions are used, and the program recompiled.

The revision control function initializes the DLL so that the functions are interpreted according to the format of the revision you wrote and used to compile your program. This function is new in revision 3.3. To take advantage of it, the function must be added to your program and the program recompiled.

The function works by interpreting the UL function call from your program and filling in any arguments needed to run with the new revision. For example, the function `cbAConvertData()` had the argument `BoardNum` added in Revision 3.3.

The two revisions of the function look like this:

#### Rev 3.2

```
int cbAConvertData(long NumPoints, unsigned ADData[], int ChanTags[])
```

#### Rev 3.3

```
int cbAConvertData(int BoardNum, long NumPoints, unsigned ADData[], int
ChanTags[])
```

If your program has declared you are running code written for revision 3.2, and you call this function, the argument `BoardNum` is ignored. If you want the benefits afforded by `BoardNum`, you must rewrite your program with the new argument and declare revision 3.3 (or higher) in `cbDeclareRevision()`.

If a revision less than 3.2 is declared, revision 3.2 is assumed.

## cbGetRevision()

Gets the revision level of Universal Library DLL and the VXD.

### Function prototype:

C/C++:	<code>int cbGetRevision(float* DLLRevNum, float* VXDRevNum);</code>
Visual Basic:	<code>Function cbGetRevision(DLLRevNum!, VXDRevNum!) As Long</code>
Delphi:	<code>function cbGetRevision(var DLLRevNum:Single; var VXDRevNum:Single):Integer;</code>

### Arguments:

DLLRevNum	Place holder for the revision number of Library DLL.
VXDRevNum	Place holder for the revision number of Library VXD.

### Returns:

DLLRevNum - Revision number of the Library DLL

VXDRevNum - Revision number of the Library VXD

[Error Code](#) if revision levels of VXD and DLL are incompatible.

---

## Streamer File Functions

### Introduction

Use the streamer file functions explained in this chapter to create, fill, and read streamer files.



## cbFileAInScan()

Scans a range of A/D channels and stores the samples in a disk file. `cbFileAInScan` reads the specified number of A/D samples at the specified sampling rate from the specified range of A/D channels from the specified board. If the A/D board has programmable gain, it sets the gain to the specified range. The collected data is returned to a file in binary format. Use [cbFileRead\(\)](#) to load data from that file into an array. See board-specific information to determine if this function is supported on your board.

### Function prototype:

C/C++:	<code>int cbFileAInScan(int BoardNum, int LowChan, int HighChan, long Count, long *Rate, int Range, char *FileName, unsigned Options)</code>
Visual Basic:	<code>Function cbFileAInScan(ByVal BoardNum&amp;, ByVal LowChan&amp;, ByVal HighChan&amp;, ByVal Count&amp;, Rate&amp;, ByVal Range&amp;, ByVal FileName\$, ByVal Options&amp;) As Long</code>
Delphi:	<code>function cbFileAInScan(BoardNum:Integer; LowChan:Integer; HighChan:Integer; Count:Longint; var Rate:Longint; Range:Integer; FileName:PChar; Options:Integer):Integer;</code>

### Arguments:

BoardNum	The number associated with the board when it was installed with the <i>InstaCal</i> configuration program. The specified board must have an A/D. BoardNum may be 0 to 99.
LowChan	First A/D channel of scan
HighChan	Last A/D channel of scan  The maximum allowable channel depends on which type of A/D board is being used. For boards with both single ended and differential inputs, the maximum allowable channel number also depends on how the board is configured (for example, eight channels for differential, 16 for single ended).
Count	Specifies the total number of A/D samples that will be collected. If more than one channel is being sampled, the number of samples collected per channel is equal to $\text{Count} / (\text{HighChan} - \text{LowChan} + 1)$ .
Rate	Sample rate in samples per second (Hz) per channel. The maximum sampling rate depends on the A/D board that is being used (see Rate explanation <a href="#">cbAInScan()</a> ).
Range	If the selected A/D board does not have a programmable range feature, this argument is ignored. Otherwise set the Range argument to any range that is supported by the selected A/D board. Refer to board specific information for a list of the supported A/D ranges of each board.
FileName	The name of the file in which to store the data. If the file doesn't exist, it will be created.
Options	Bit fields that control various options. Refer to the constants in the "Options argument values" section on page 138.

### Returns:

[Error code](#) or 0 if no errors.

Rate = actual sampling rate.

**Options argument values:**

EXTCLOCK	If this option is used, conversions are controlled by the signal on the trigger input line rather than by the internal pacer clock. Each conversion is triggered on the appropriate edge of the trigger input signal (refer to board-specific information contained in the <i>Universal Library Users Guide</i> ). Additionally, the <code>Rate</code> argument is ignored. The sampling rate is dependent on the trigger signal.
EXTTRIGGER	<p>If this option is specified, the sampling does not begin until the trigger condition is met.</p> <p>On many boards, this trigger condition is programmable (refer to the <a href="#">cbSetTrigger()</a> function and board-specific information for details) and can be programmed for rising or falling edge or an analog level.</p> <p>On other boards, only 'polled gate' triggering is supported. Assuming active high operation, data acquisition commences immediately if the trigger input is high. If the trigger input is low, acquisition is held off until it goes high. Acquisition continues until <code>NumPoints</code>&amp; samples are taken, regardless of the state of the trigger input. For 'polled gate' triggering, this option is most useful if the signal is a pulse with a very low duty cycle (trigger signal in TTL low state most of the time) to hold off triggering until the pulse occurs.</p>
DTCONNECT	Samples are sent to the DT-Connect port if the board is equipped with one.

**Notes:**

[OVERRUN Error](#) - (Error code 29): This error indicates that the data was not written to the file as fast as the data was sampled. Consequently some data was lost. The value returned from [cbFileGetInfo\(\)](#) in `TotalCount` is the number of points that were successfully collected.

**Important**

In order to understand the functions, read the board-specific information contained in the *Universal Library User's Guide* (available on our web site at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf)).

We also urge you to examine and run one or more of the example programs supplied prior to attempting any programming of your own. Following this advice may save you hours of frustration, and wasted time.

This note, which appears elsewhere, is especially applicable to this function. Now is the time to read the board specific information for your board. We suggest that you make a copy of that page to refer to as you read this manual and examine the example programs.

## cbFileGetInfo()

Returns information about a streamer file. When [cbFileAInScan\(\)](#) or [cbFilePretrig\(\)](#) fills the streamer file, information is stored about how the data was collected (sample rate, channels sampled etc.). This function returns that information. Refer to board-specific information in the *Universal Library User's Guide* to determine if your board supports `cbFileAInScan()` and/or `cbFilePretrig()`.

### Function prototype:

C/C++:	<code>int cbFileGetInfo(char *FileName, short *LowChan, short *HighChan, long *PretrigCount, long *TotalCount, long *Rate, int *Range)</code>
Visual Basic:	<code>Function cbFileGetInfo(ByVal FileName\$, LowChan%, HighChan%, PretrigCount&amp;, TotalCount&amp;, Rate&amp;, Range&amp;) As Long</code>
Delphi:	<code>function cbFileGetInfo(FileName:PChar; var LowChan:SmallInt; var HighChan:SmallInt; var PretrigCount:Longint; var TotalCount:Longint; var Rate:Longint; var Range:LongInt):Integer;</code>

### Arguments:

FileName	Name of streamer file.
LowChan	Variable to return LowChan to.
HighChan	Variable to return HighChan to.
PretrigCount	Variable to return PretrigCount to.
TotalCount	Variable to return TotalCount to.
Rate	Variable to return sampling rate to.
Range	Variable to return A/D range code to.

### Returns:

[Error code](#) or 0 if no errors.

LowChan - low A/D channel of scan.

HighChan - high A/D channel of scan.

TotalCount - total number of points collected.

PretrigCount - number of pre-trigger points collected.

Rate - sampling rate when data was collected.

Range - Range of A/D when data was collected .

## cbFilePretrig()

Scan a range of channels continuously while waiting for a trigger. Once the trigger occurs, return the specified number of samples including the specified number of pre-trigger samples to a disk file. This function waits for a trigger signal to occur on the Trigger Input. Once the trigger occurs, it returns the specified number (TotalCount) of A/D samples including the specified number of pre-trigger points. It collects the data at the specified sampling rate (Rate) from the specified range (LowChan-HighChan) of A/D channels from the specified board. If the A/D board has programmable gain then it sets the gain to the specified range. The collected data is returned to a file. See board-specific info to determine if this function is supported by your board.

### Function prototype:

C/C++:	<pre>int cbFilePretrig(int BoardNum, int LowChan, int HighChan, long *PretrigCount, long *TotalCount, long *Rate, int Range, char *FileName, unsigned Options)</pre>
Visual Basic:	<pre>Function cbFilePretrig(ByVal BoardNum&amp;, ByVal LowChan&amp;, ByVal HighChan&amp;, PretrigCount&amp;, TotalCount&amp;, Rate&amp;, ByVal Range&amp;, ByVal FileName\$, ByVal Options&amp;) As Long</pre>
Delphi:	<pre>function cbFilePretrig(BoardNum:Integer; LowChan:Integer; HighChan:Integer; var PretrigCount:Longint; var TotalCount:Longint; var Rate:Longint; Range:Integer; FileName:PChar; Options:Integer):Integer;</pre>

### Arguments:

BoardNum	The number associated with the board when it was installed with the <i>InstaCal</i> configuration program. The specified board must have an A/D and pretrigger capability. BoardNum may be 0 to 99.
LowChan	First A/D channel of scan
HighChan	Last A/D channel of scan  The maximum allowable channel depends on which type of A/D board is being used. For boards that have both single ended and differential inputs the maximum allowable channel number also depends on how the board is configured. Refer to board-specific information for the maximum number of channels allowed in differential and single ended modes.
PretrigCount	Specifies the number of samples before the trigger that will be returned. PretrigCount must be less than 16000 and PretrigCount must also be less than TotalCount - 512.  If the trigger occurs too early, then fewer than the requested number of pre-trigger samples will be collected. In that case a TOOFEW error will occur. The PretrigCount will be set to indicate how many samples were collected and the post trigger samples will still be collected.
TotalCount	Specifies the total number of samples that will be collected and stored in the file. TotalCount must be greater than or equal to PretrigCount + 512. If the trigger occurs too early then fewer than the requested number of samples will be collected. In that case a TOOFEW error will occur. The TotalCount will be set to indicate how many samples were actually collected.

Rate	Sample rate in samples per second (Hz) per channel. The maximum sampling rate depends on the A/D board that is being used. This is the rate at which scans are triggered. If you are sampling 4 channels, 0 - 3, then specifying a rate of 10,000 scans per second (10 kHz) will result in the A/D converter rate of 40 kHz: 4 channels at 10,000 samples per channel per second. This is different from some software where you specify the total A/D chip rate. In those systems, the per channel rate is equal to the A/D rate divided by the number of channels in a scan. This argument also returns the value of the actual set. This may be different from the requested rate because of pacer limitations.
Range	If the selected A/D board does not have a programmable range feature, this argument is ignored. Otherwise, set the Range argument to any range that is supported by the selected A/D board. Refer to board specific information for a list of the supported A/D ranges of each board.
FileName	The name of the file in which to store the data. If the file doesn't exist, it will be created.
Options	Bit fields that control various options. Refer to the constants in the "Options argument values" section below.

#### Returns:

[Error code](#) or 0 if no errors.

PretrigCount - actual number of pre-trigger samples collected.

TotalCount - actual number of samples collected.

Rate = actual sampling rate.

#### Options argument values:

EXTCLOCK	If this option is used then conversions will be controlled by the signal on the trigger input line rather than by the internal pacer clock. Each conversion will be triggered on the appropriate edge of the trigger input signal (refer to board-specific information contained in the <i>Universal Library Users Guide</i> ). When this option is used the Rate argument is ignored. The sampling rate is dependent on the trigger signal.
DTCONNECT	Samples are sent to the DT-Connect port if the board is equipped with one.

#### Notes:

[OVERRUN Error](#) - (Error code 29): This error indicates that the data was not written to the file as fast as the data was sampled. Consequently some data was lost. The value in TotalCount will be the number of points that were successfully collected.

## cbFileRead()

Reads data from a streamer file. When [cbFileAInScan\(\)](#) or [cbFilePretrig\(\)](#) fills the streamer file, this function returns the content of that file. Refer to information on your board in the *Universal Library User's Guide* to determine if your board supports `cbFileAInScan()` and/or `cbFilePreTrig()`.

### Function prototype:

C/C++:	<code>int cbFileRead(char *FileName, long FirstPoint, long *NumPoints, int *DataBuffer)</code>
Visual Basic:	<code>Function cbFileRead(ByVal FileName\$, ByVal FirstPoint&amp;, NumPoints &amp;, DataBuffer%) As Long</code>
Delphi:	<code>function cbFileRead(FileName:PChar; FirstPoint:Longint; var NumPoints:Longint; var DataBuffer:Word):Integer;</code>

### Arguments:

<code>FileName</code>	Name of streamer file
<code>FirstPoint</code>	Index of first point to read
<code>NumPoints</code>	Number of points to read from file
<code>DataBuffer</code>	Pointer to data buffer that data will be read into.

### Returns:

[Error code](#) or 0 if no errors.

`DataBuffer` - data read from file.

`NumPoints` - number of points actually read.

`NumPoints` may be less than the requested number of points if an error occurs.

### Notes:

**Data format:** The data is returned as 16-bits. The 16-bits may represent 12-bits of analog, 12-bits of analog plus 4 bits of channel, or 16-bits of analog. Use [cbAConvertData\(\)](#) to correctly load the data into an array.

**Loading portions of files:** The file may contain much more data than can fit in `DataBuffer`. In those cases use `NumPoints` and `FirstPoint` to read a selected piece of the file into `DataBuffer`. Call [cbFileGetInfo\(\)](#) first to find out how many points are in the file.

---

## Synchronous I/O Functions

### Introduction

Use the functions discussed in this chapter to synchronously read and write data from analog channels, counter channels, thermocouple channels, and digital ports. These functions can be used with hardware equipped with synchronous input and output capability.

## cbDaqInScan()

Scans analog, digital, counter, and temperature input channels synchronously, and stores the samples in an array. This function only works with boards that support synchronous input.

### Function prototype:

**C/C++:** `int cbDaqInScan(int BoardNum, short ChanArray[], short ChanTypeArray[], short GainArray[], int ChanCount, long* Rate, long *PretrigCount, long *TotalCount, int MemHandle, int Options);`

**Visual Basic:** `Function cbDaqInScan(ByVal BoardNum&, ChanArray%, ChanTypeArray%, GainArray%, ByVal ChanCount&, CBRate&, PretrigCount&, CBCount&, ByVal MemHandle&, ByVal Options&) As Long`

**Delphi:** `function cbDaqInScan(BoardNum:Integer; var ChanArray:SmallInt; var ChanTypeArray:SmallInt; var GainArray:SmallInt; ChanCount:Integer; var Rate:LongInt; var PretrigCount:LongInt; var TotalCount:LongInt; MemHandle:Integer; Options:Integer):Integer;`

### Arguments:

BoardNum	The board number used to collect the data. BoardNum may be 0 to 99. Refers to the number associated with the board used to collect the data when it was installed with the <i>InstaCal</i> configuration program. The board must support synchronous input.
ChanArray	Array containing channel values. Valid channel values are analog input channels, digital ports, counter input channels, and temperature input channels of the device.
ChanTypeArray	Array containing channel types. Each element of this array defines the type of the corresponding element in the ChanArray. Set to one of the constants in the "ChanTypeArray argument values" section on page 145.
GainArray	Array containing A/D range codes. If the corresponding element in the ChanArray is not an analog input channel, the range code for this channel is ignored.
ChanCount	Number of elements in each of the three arrays - ChanArray, ChanTypeArray and GainArray.
Rate	The sample rate at which samples are acquired, in samples per second per channel.  Rate also returns the value of the actual rate set, which may be different from the requested rate because of pacer limitations.
PretrigCount	Sets the number of pre-trigger samples to collect. Specifies the number of samples to collect before the trigger occurs. This function won't run in pre-trigger mode if PreTrigCount is set to zero. PreTrigCount is ignored if the EXTTRIGGER option is not specified.  PreTrigCount also returns the value of the actual pre-trigger count set, which may be different from the requested pre-trigger count because pre-trigger count must be a multiple of ChanCount.  PreTrigCount must be evenly divisible by the number of channels being scanned (ChanCount). If it is not, this function adjusts the number (up) to the next valid value, and returns that value to the PreTrigCount argument.
TotalCount	Total number of samples to collect. Specifies the total number of samples to collect and store in the buffer. TotalCount must be greater than PretrigCount.  TotalCount also returns the value of the actual total count set, which may be different from the requested total count, because total count must be a multiple of ChanCount.



	TotalCount must be evenly divisible by the number of channels being scanned (ChanCount). If it is not, this function adjusts the number (down) to the next valid value, and returns that value to the TotalCount argument.
MemHandle	Handle for the Windows buffer to store data in. This buffer must have been previously allocated with the <a href="#">cbWinBufAlloc()</a> function.
Options	Bit fields that control various options. This field may contain any combination of non-contradictory choices in the "Options argument values" section below.

**ChanTypeArray argument values:**

ANALOG	Analog input channel.
DIGITAL8	8-bit digital input port.
DIGITAL16	16-bit digital input port. (FIRSTPORTA only)
CTR16	16-bit counter.
CTR32LOW	Lower 16-bits of a 32-bit counter.
CTR32HIGH	Upper 16-bits of a 32-bit counter.
CJC	CJC channel.
TC	Thermocouple channel.

The [cbGetTCValues\(\)](#) function can be used to convert raw thermocouple data to data on a temperature scale (Celsius, Fahrenheit or Kelvin). **Note:** If at least one TC channel is listed in the channel array, and averaging is enabled for that channel, the averaging will be applied to all of the channels listed in the channel array.

SETPOINTSTATUS	The setpoint status register. This is a bitfield indicating the state of each of the setpoints. A "1" indicates that the setpoint criteria has been met.
----------------	--

**ChanTypeArray flag values:**

SETPOINT_ENABLE	Enables a setpoint. When this option is specified, it must be OR'ed with the ChanTypeArray argument values.
-----------------	---

You set the setpoint criteria with the [cbDagSetSetpoints\(\)](#) function. The number of channels set with the SETPOINT\_ENABLE flag must match the number of setpoints set by the [cbDagSetSetpoints\(\)](#) function's SetpointCount argument.

**Options argument values:**

BACKGROUND	When the BACKGROUND option is used, control returns immediately to the next line in your program and the data collection into the buffer continues in the background. If the BACKGROUND option is not used, the <a href="#">cbDagInScan()</a> function does not return control to your program until all of the requested data has been collected and returned to the buffer.
------------	---

Use [cbGetStatus\(\)](#) with DAQIFUNCTION to check on the status of the background operation. Use [cbStopBackground\(\)](#) with DAQIFUNCTION to terminate the background process before it has completed. Execute [cbStopBackground\(\)](#) after normal termination of all background functions, in order to clear variables and flags.

CONTINUOUS	This option puts the function in an endless loop. Once it collects the required number of samples, it resets to the start of the buffer and begins again. The only way to stop this operation is by using <a href="#">cbStopBackground()</a> with DAQIFUNCTION. Normally this option should be used in combination with BACKGROUND so that your program will regain control.
------------	--

EXTCLOCK	If this option is used, conversions will be controlled by the signal on the external clock input rather than by the internal pacer clock. Each conversion will be triggered on the appropriate edge of the clock input signal. When this option is used the <code>Rate</code> argument is ignored. The sampling rate is dependent on the clock signal. Options for the board will default to a transfer mode that will allow the maximum conversion rate to be attained unless otherwise specified.
EXTTRIGGER	If this option is specified, the sampling will not begin until the trigger condition is met (refer to the <a href="#">cbDaqSetTrigger()</a> function).

**Returns:**

[Error code](#) or 0 if no errors

`Rate` – Actual sampling rate used.

`PreTrigCount` – Actual pre-trigger count used.

`TotalCount` – Actual total count used.

`MemHandle` - Collected data returned via the Windows buffer.

## cbDagOutScan()

Outputs values synchronously to analog output channels and digital output ports. This function only works with boards that support synchronous output.

### Function prototype:

**C/C++:**                    `int cbDagOutScan(int BoardNum, short ChanArray[], short ChanTypeArray[], short GainArray[], int ChanCount, long* Rate, long Count, int MemHandle, int Options);`

**Visual Basic:**           `Function cbDagOutScan(ByVal BoardNum%, ChanArray%, ChanTypeArray%, GainArray%, ByVal ChanCount&, CBRate&, ByVal CBCount&, ByVal MemHandle&, ByVal Options&) As Long`

**Delphi:**                   `function cbDagOutScan(BoardNum:Integer; var ChanArray:SmallInt; var ChanTypeArray:SmallInt; var GainArray:SmallInt; ChanCount:Integer; var Rate:LongInt; Count:LongInt; MemHandle:Integer; Options:Integer):Integer;`

### Arguments:

BoardNum	The board number used to collect the data. BoardNum may be 0 to 99. Refers to the number associated with the board used to collect the data when it was installed with the <i>InstaCal</i> configuration program. The board must support synchronous output.
ChanArray	Array containing channel values. Valid channel values are analog output channels and digital ports.
ChanTypeArray	Array containing channel types. Each element of this array defines the type of the corresponding element in the ChanArray. Choices are:  ANALOG                    Analog output channel.  DIGITAL16                16-bit digital output port. (FIRSTPORTA only)
GainArray	Array containing D/A range codes. If the corresponding element in the ChanArray is not an analog output channel, the range code for this channel is ignored. If the board does not have programmable gain, this parameter is ignored, and therefore can be set to null.
ChanCount	Number of elements in each of the three arrays - ChanArray, ChanTypeArray and GainArray.
Rate	Sample rate in scans per second. Rate also returns the value of the actual rate set, which may be different from the requested rate because of pacer limitations.
Count	Sets the total number of values to output.Count also returns the value of the actual count set, which may be different from the requested total count because count must be a multiple of the channel count.
MemHandle	Handle for the Windows buffer from which data will be output. This buffer must have been previously allocated with the <a href="#">cbWinBufAlloc()</a> function and data values loaded (for example using <a href="#">cbWinArrayToBuf()</a> ).
Options	Bit fields that control various options. This field may contain any combination of non-contradictory choices in the "Options argument values" section on page 148.

**Options argument values :**

BACKGROUND	When this option is used the output operations will begin running in the background and control will immediately return to the next line of your program. Use <a href="#">cbGetStatus()</a> with the DAQOFUNCTION option to check the status of background operation. Use <a href="#">cbStopBackground()</a> with DAQOFUNCTION to terminate background operations before they are completed. Execute <code>cbStopBackground()</code> with DAQOFUNCTION after normal termination of all background functions in order to clear variables and flags.
CONTINUOUS	This option puts the function in an endless loop. Once it outputs the specified number (Count) of output values, it resets to the start of the buffer and begins again. The only way to stop this operation is by calling <a href="#">cbStopBackground()</a> with DAQOFUNCTION. This option should only be used in combination with BACKGROUND so that your program can regain control.
EXTCLOCK	<p>If this option is used, conversions will be paced by the signal on the external clock input rather than by the internal pacer clock. Each conversion will be triggered on the appropriate edge of the clock input signal.</p> <p>When this option is used, the Rate argument is ignored. The sampling rate is dependent on the clock signal. Options for the board will default to transfer types that allow the maximum conversion rate to be attained unless otherwise specified.</p>
ADCCLOCKTRIG	If this option is used, the data output operation will be triggered upon the start of the ADC clock.
ADCCLOCK	When this option is specified, the data output operation will be paced by the ADC clock.
NONSTREAMEDIO	<p>This option allows non-streamed data output to be generated to a specified output channel.</p> <p>In this mode, the aggregate size of data output buffer must be less than or equal to the size of the internal data output FIFO on the Measurement Computing device. This allows the data output buffer to be loaded into the device's internal output FIFO.</p> <p>Once the sample updates are transferred (or downloaded) to the device, the device is responsible for outputting the data. While the size is limited, and the output buffer cannot be changed once the output is started, this mode has the advantage being able to continue data output without having to periodically feed output data through the program to the device.</p>

**Returns:**

[Error code](#) or 0 if no errors

Rate – Actual sampling rate used.

## cbDagSetSetpoints()

Configures up to 16 detection setpoints associated with the input channels within a scan group. This function only works with boards that support synchronous input.

### Function Prototype:

**C/C++:**

```
int cbDagSetSetpoints(int BoardNum, float *LimitAArray, float
*LimitBArray, float *reserved, int *SetpointFlagsArray, int
*SetpointOutputArray, float *Output1Array, float *Output2Array,
float *OutputMask1Array, float *OutputMask2Array, int
SetpointCount);
```

**Visual Basic:**

```
Function cbDagSetSetpoints(ByVal BoardNum&, LimitAArray!,
LimitBArray!, Reserved!, SetpointFlagsArray&, SetpointOutputArray&,
Output1Array!, Output2Array!, OutputMask1Array!, OutputMask2Array!,
ByVal SetpointCount&) As Long
```

**Delphi:**

```
function cbDagSetSetpoints(BoardNum:Integer; var LimitAArray:Single;
var LimitBArray:Single; var Reserved:Single; var
SetpointFlagsArray:Integer; var SetpointOutputArray:Integer; var
Output1Array:Single; var Output2Array:Single; var
OutputMask1Array:Single; var OutputMask2Array:Single;
SetpointCount:Integer):Integer;
```

### Arguments:

BoardNum	The board number used to collect the data. BoardNum may be 0 to 99. Refers to the number associated with the board used to collect the data when it was installed with the <i>InstaCal</i> configuration program. The specified board must support synchronous input.
LimitAArray	Array containing the limit A values for the input channels used for the setpoint. Limit A specifies a value used to determine if the setpoint criteria are met.
LimitBArray	Array containing the limit B values for the input channels used for the setpoint. Limit B specifies a value used to determine if the setpoint criteria are met.
Reserved	Reserved for future use.
SetpointFlagsArray	Array containing the setpoint flags. Set to one of the constants in the "SetpointFlagsArray argument values" section below.
SetpointOutputArray	Array containing output sources. Set to one of the constants in the "SetpointOutputArray argument values" section on page 150.
Output1Array	Array containing the values for the output channels used for the setpoint.
Output2Array	Array containing the values for the output channels used for the setpoint.
OutputMask1Array	Array containing the output masks for output value 1 – for FIRSTPORTC only.
OutputMask2Array	Array containing the output masks for output value 2 – for FIRSTPORTC only.
SetpointCount	Number of setpoints to configure (0 -16). Set to 0 to disable the setpoints.

### SetpointFlagsArray argument values:

Flag	Description
SF_EQUAL_LIMITA	Setpoint criteria: The input channel = limit A.
SF_LESSTHAN_LIMITA	Setpoint criteria: The input channel < limit A.
SF_GREATERTHAN_LIMITB	Setpoint criteria: The input channel > limit B.

SF_INSIDE_LIMITS	Setpoint criteria: The input channel > limit A and < limit B.
SF_OUTSIDE_LIMITS	Setpoint criteria: The input channel < limit A and > limit B.
SF_HYSTERESIS	Setpoint criteria: If the input channel > limit A then output value 1. If the input channel < limit B then output value 2.
SF_UPDATEON_TRUEONLY	If the criteria is met then output value 1.
SF_UPDATEON_TRUEANDFALSE	If the criteria is met then output value 1, else output value 2.

**SetpointOutputArray argument values:**

Output Source	Description
SO_NONE	Perform no outputs.
SO_FIRSTPORTC	Output to FIRSTPORTC when the criteria is met.
SO_DAC0	Output to DAC0 when the criteria is met. You must have a device with DAC0.
SO_DAC1	Output to DAC1 when the criteria is met. You must have a device with DAC1.
SO_DAC2	Output to DAC2 when the criteria is met. You must have a device with DAC2.
SO_DAC3	Output to DAC3 when the criteria is met. You must have a device with DAC3.
SO_TMR0	Output to timer 0 when the criteria is met.
SO_TMR1	Output to timer 1 when the criteria is met.

**Returns:**

Error code or 0 if no errors

## cbDaqSetTrigger()

Selects the trigger source and sets up its parameters. This trigger is used to initiate or terminate an acquisition using the [cbDaqInScan\(\)](#) function if the EXTTRIGGER option is selected. This function only works with boards that support synchronous output.

### Function prototype:

C/C++:	<pre>int cbDaqSetTrigger(int BoardNum, int TrigSource, int TrigSense, int TrigChan, int ChanType, int Gain, float Level, float Variance, int TrigEvent);</pre>
Visual Basic:	<pre>Function cbDaqSetTrigger(ByVal BoardNum&amp;, ByVal TrigSource&amp;, ByVal TrigSense&amp;, ByVal TrigChan&amp;, ByVal ChanType&amp;, ByVal Gain&amp;, ByVal Level!, ByVal Variance!, ByVal TrigEvent&amp;) As Long</pre>
Delphi:	<pre>function cbDaqSetTrigger(BoardNum:Integer; TrigSource:Integer; TrigSense:Integer; TrigChan:Integer; ChanType:Integer; Gain:Integer; Level:Single; Variance:Single; TrigEvent:Integer):Integer;</pre>

### Arguments:

BoardNum	The board number used to collect the data. BoardNum may be 0 to 99. Refers to the number associated with the board used to collect the data when it was installed with the <i>InstaCal</i> configuration program. The board must support synchronous output.
TrigSource	Specifies the type of triggering based on the external trigger source. Set to one of the constants specified in "TrigSource argument values" section on page 152.
TrigSense	Specifies the trigger sensitivity. The trigger sensitivity normally defines the way in which a trigger event is detected based upon the characteristics of the trigger input signal. However, it often defines the way in which the trigger input signal(s) should be compared to the trigger level parameter value. Set to one of the constants specified in "TrigSense argument values" section on page 152.
TrigChan	Specifies the trigger channel. The trigger channel must be a configured channel in the channel array (refer to <a href="#">cbDaqInScan()</a> ).
ChanType	Specifies the channel type and should match the channel type setting for the trigger channel configured using the <a href="#">cbDaqInScan()</a> function.
Gain	Specifies the trigger channel gain code. If the device has programmable gain, this argument should match the gain code setting when the channel is configured using the <a href="#">cbDaqInScan()</a> function. The Gain parameter is ignored if TrigChan is not an analog channel.
Level	<p>A single precision floating point value which represents, in engineering units, the level at or around which the trigger event should be detected.</p> <p>This option is used for trigger types that depend on an input channel comparison to detect the start trigger or stop trigger event.</p> <p>The actual level at which the trigger event is detected depends upon trigger sensing and variability. Refer to <a href="#">Trigger levels</a> on page 153 for more information.</p>
Variance	<p>A single-precision floating point value which represents, in engineering units, the amount that the trigger event can vary from the Level parameter.</p> <p>While the TrigSense parameter indicates the direction of the input signal relative to the Level parameter, the Variance parameter specifies the degree to which the input signal can vary relative to the Level parameter.</p>

TrigEvent	<p>Specifies the trigger event type. Valid values indicate either a start trigger event (START_EVENT) or a stop trigger event (STOP_EVENT).</p> <p>START_EVENT: The start trigger event defines the conditions under which post-trigger acquisition data collection should be initiated or triggered. The start trigger event can vary in complexity from starting immediately, to starting on complex channel value definitions.</p> <p>STOP_EVENT: The stop trigger event signals the current data acquisition process to terminate. The stop event can be as simple as that of a scan count, or as complex as involving a channel value level condition.</p>
-----------	---

**TrigSource argument values:**

TRIG_IMMEDIATE	Start trigger event only. Acquisition begins immediately upon invocation the <a href="#">cbDaqInScan</a> function. No pre-trigger data acquisition is possible with this trigger type.
TRIG_EXTTTL	Start trigger event only. Acquisition begins on the selectable edge of an external TTL signal. No pre-trigger data acquisition is possible with this trigger type.
TRIG_ANALOGHW	Start trigger event only. Data acquisition begins upon a selectable criteria of the input signal (above level, below level, rising edge, etc.) TrigChan must be defined as the first channel in the channel scan group. No pre-trigger data acquisition is possible with this trigger type.
TRIG_ANALOGSW	Post-trigger data acquisition begins upon a selectable criteria of the input signal (above level, below level, rising edge, etc.)
TRIG_DIGPATTERN	Post-trigger data acquisition begins upon receiving a specified digital pattern on the specified digital port.
TRIG_COUNTER	Post-trigger data acquisition begins upon detection of specified counter criteria.
TRIG_SCANCOUNT	Stop trigger event only. Stops collecting post-trigger data when the specified number of post-trigger scans are completed.

**TrigSense argument values:**

RISING_EDGE:	Triggers when the signal goes from low to high (TTL trigger) or rises through a specified level (hardware analog, software analog, and counter).
FALLING_EDGE:	Triggers when the signal goes from high to low (TTL trigger) or falls through a specified level (hardware analog, software analog, and counter).
ABOVE_LEVEL:	Triggers when the signal is above a specified level (hardware analog, software analog, counter, and digital pattern).
BELOW_LEVEL:	Triggers when the signal is below a specified level (hardware analog, software analog, counter, and digital pattern).
EQ_LEVEL:	Triggers when the signal equals a specified level (hardware analog, software analog, counter, and digital pattern).
NE_LEVEL:	Triggers when the signal does not equal a specified level (hardware analog, software analog, counter, and digital pattern).

**Returns:**

[Error code](#) or 0 if no errors



**Notes:**

**Trigger levels:** The actual level at which the trigger event is detected depends upon trigger sensing and variability. The various ranges of possible values for the `Level` parameter based on the trigger source are listed here:

TRIG_ANALOG_HW:	The voltage used to define the trigger level. Trigger detection is performed in hardware.
TRIG_ANALOG_SW:	The voltage used to define the trigger level. Trigger detection is performed in software.
TRIG_DIGPATTERN:	Sets the bit pattern for the digital channel trigger. Choices are: 0.0 (no bits set): 255.0 (all bits set) for 8-bit digital ports. 0.0 (no bits set): 65,535.0 (all bits set) for 16-bit digital ports.
TRIG_COUNTER:	Selects either Pulse or Totalize counter values (0.0 – 65,535).
TRIG_IMMEDIATE:	Ignored
TRIG_SCANCOUNT:	Ignored

**Trigger start and stop criteria:** The table below lists the trigger start and stop criteria based on the selected trigger type and sensitivity.

Table 4. Interaction of trigger variance with trigger level and trigger sensitivities

Trigger Start/Stop Source (TrigSource)	Trigger Sensitivity (TrigSense)	Trigger Start/Stop Criteria
TRIG_ANALOGHW (Start trigger event only)	RISING_EDGE	Triggers when the signal value $< (\text{Level} - \text{Variance})$ Then, the signal value $> \text{Level}$
	FALLING_EDGE	Triggers when the signal value $> (\text{Level} + \text{Variance})$ Then, the signal value $< \text{Level}$
	ABOVE_LEVEL	Triggers when the signal value $> (\text{Level})$
	BELOW_LEVEL	Triggers when the signal value $< (\text{Level})$
TRIG_ANALOGSW	RISING_EDGE	Triggers/Stops when the signal value $< (\text{Level} - \text{Variance})$ Then, the signal value $> \text{Level}$
	FALLING_EDGE	Triggers/Stops when the signal value $> (\text{Level} + \text{Variance})$ Then, the signal value $< \text{Level}$
	ABOVE_LEVEL	Triggers/Stops when the signal value $> (\text{Level})$
	BELOW_LEVEL	Triggers/Stops when the signal value $< (\text{Level})$
	EQ_LEVEL	Triggers/Stops when the $(\text{Level} - \text{Variance}) < \text{signal value} < (\text{Level} + \text{Variance})$
	NE_LEVEL	Triggers/Stops when the signal value $< (\text{Level} - \text{Variance})$ OR when the signal value $> (\text{Level} + \text{Variance})$
TRIG_DIGPATTERN	ABOVE_LEVEL	Triggers/Stops when the $(\text{digital port value AND (bitwise) Variance}) > (\text{Level AND (bitwise) Variance})$
	BELOW_LEVEL	Triggers/Stops when the $(\text{digital port value AND (bitwise) Variance}) < (\text{Level AND (bitwise) Variance})$
	EQ_LEVEL	Triggers/Stops when the $(\text{digital port value AND (bitwise) Variance}) = (\text{Level AND (bitwise) Variance})$
	NE_LEVEL	Triggers/Stops when the $(\text{digital port value AND (bitwise) Variance}) \neq (\text{Level AND (bitwise) Variance})$

Trigger Start/Stop Source (TrigSource)	Trigger Sensitivity (TrigSense)	Trigger Start/Stop Criteria
TRIG_COUNTER	RISING_EDGE	Triggers/Stops when the counter channel $< (\text{Level} - \text{Variance})$ Then, the counter channel $> \text{Level}$
	FALLING_EDGE	Triggers/Stops when counter channel $> (\text{Level} + \text{Variance})$ Then, the counter channel $< \text{Level}$
	ABOVE_LEVEL	Triggers/Stops when counter channel $> (\text{Level} - \text{Variance})$
	BELOW_LEVEL	Triggers/Stops when counter channel $< (\text{Level} + \text{Variance})$
	EQ_LEVEL	Triggers/Stops when the $(\text{Level} - \text{Variance}) < \text{counter channel} < (\text{Level} + \text{Variance})$
	NE_LEVEL	Triggers/Stops when the counter channel $< (\text{Level} - \text{Variance})$ OR when the counter channel $> (\text{Level} + \text{Variance})$

---

# Temperature Input Functions

## Introduction

Use the functions discussed in this chapter to convert a raw analog input from an EXP board, or other temperature sensor board, to temperature.

## cbTIn()

### Changed R3.3 ID

Reads an analog input channel, linearizes it according to the selected temperature sensor type, and returns the temperature in degrees. The CJC channel, the gain, and sensor type, are read from the *InstaCal* configuration file. They should be set by running the *InstaCal* configuration program.

### Function prototype:

C/C++:	<code>int cbTIn(int BoardNum, int Chan, int Scale, float *TempVal, int Options)</code>
Visual Basic:	<code>Function cbTIn(ByVal BoardNum&amp;, ByVal Chan&amp;, ByVal Scale&amp;, TempVal!, ByVal Options&amp;) As Long</code>
Delphi:	<code>function cbTIn(BoardNum:Integer; Chan:Integer; Scale:Integer; var TempValue:Single; Options:Integer):Integer;</code>

### Arguments:

BoardNum	The number associated with the board when it was installed with the <i>InstaCal</i> configuration program. BoardNum may be 0 to 99.
Chan	Input channel to read.
Scale	Specifies the temperature scale that the input will be converted to. Choices are CELSIUS, FAHRENHEIT and KELVIN.
TempVal	The temperature in degrees is returned here. Thermocouple resolution is approximately 0.25 °C, depending on scale, range and thermocouple type. RTD resolution is 0.1 °C.
Options	Bit fields that control various options. Refer to the constants in the "Options argument values" section below.

### Returns:

[Error code](#) or 0 if no errors.

\*TempVal - Temperature returned here.

### Options argument values:

FILTER	When selected, a smoothing function is applied to temperature readings, very much like the electrical smoothing inherent in all hand held temperature sensor instruments. This is the default. When selected, 10 samples are read from the specified channel and averaged. The average is the reading returned. Averaging removes normally distributed signal line noise.
NOFILTER	If you use the NOFILTER option, then the readings will not be smoothed and you will see a scattering of readings around a mean.

### Notes:

**Using CIO-EXP boards:** For CIO-EXP boards, the channel number is calculated using the following formula, where:

- ADChan is the A/D channel that is connected to the multiplexer.
  - MuxChan is a number ranging from 0 to 15 that specifies the channel number on a particular bank of the multiplexer board
- $$\text{Chan} = (\text{ADChan} * 16) + (16 + \text{MuxChan})$$

For example, you have an EXP16 connected to a CIO-DAS08 via the CIO-DAS08 channel 0. (Remember that DAS08 channels are numbered 0, 1, 2, 3, 4, 5, 6 & 7). If you connect a thermocouple to channel 5 of the EXP16, the value for Chan would be  $(0 * 16) + (16 + 5) = 0 + 21 = 21$ .

**Using 6K-EXP boards:** For 6K-EXP boards, the channel number (Chan) is calculated using one of the following formulas, where:

- ADChan is the A/D channel that is connected to the multiplexer.
- MuxChan is a number ranging from 0 to 15 that specifies the channel number on a particular bank of the multiplexer board.
- If the A/D board has 16 or less single-ended channels:  

$$\text{Chan} = (\text{ADChan} * 16) + (16 + \text{MuxChan})$$

For example, you have a 6K-EXP16 connected to a PCI-DAS6052 via the a PCI-DAS6052 channel 0. If you connect a thermocouple to channel 5 of the 6K-EXP16, the value for Chan would be  $(0 * 16) + (16 + 5) = 0 + 21 = 21$ .

- If the A/D board has 64 single-ended channels and the A/D multiplexer channel is less than or equal to 7:  

$$\text{Chan} = (\text{ADChan} * 16) + (64 + \text{MuxChan})$$

For example, you have a 6K-EXP16 connected to a PCI-DAS6031 via the a PCI-DAS6031 channel 7. If you connect a thermocouple to channel 5 of the 6K-EXP16, the value for Chan would be  $(7 * 16) + (64 + 5) = 112 + 69 = 181$ .

- If the A/D board has 64 single-ended channels and the A/D multiplexer channel is  $\geq 31$ :  

$$\text{Chan} = (\text{ADChan} * 16 - 320) + \text{MuxChan}$$

For example, you have a 6K-EXP16 connected to a PCI-DAS6031 via the PCI-DAS6031 channel 32. If you connect a thermocouple to channel 5 of the 6K-EXP16, the value for Chan would be  $(32 * 16 - 320) + 5 = 192 + 5 = 197$ .

**CJC Channel:** The CJC channel is set in the *InstaCal* install program. If you have multiple EXP boards, Universal Library will apply the CJC reading to the linearization formula in the following manner:

- If you have chosen a CJC channel for the EXP board that the channel you are reading is on, it will use the CJC temp reading from that channel.
- If you left the CJC channel for the EXP board that the channel you are reading is on to NOT SET, the library will use the CJC reading from the next lower EXP board with a CJC channel selected.

For example: You have four CIO-EXP16 boards connected to a CIO-DAS08 on channel 0, 1, 2 and 3. You choose CIO-EXP16 #1 (connected to CIO-DAS08 channel 0) to have its CJC read on CIO-DAS08 channel 7, AND, you leave the CIO-EXP16's 2, 3 and 4 CJC channels to NOT SET. Result: The CIO-EXP boards will all use the CJC reading from CIO-EXP16 #1, connected to channel 7 for linearization. As you can see, it is important to keep the CIO-EXP boards in the same case and out of any breezes to ensure a clean CJC reading.

#### Important

For an EXP board connected to an A/D board that does not have programmable gain (DAS08, DAS16, DAS16F), the A/D board range is read from the configuration file (cb.cfg). In most cases, set hardware-selectable ranges to  $\pm 5$  V for thermocouples, and to 0 to 10 V for RTDs. Refer to the board-specific information in the *Universal Library User's Guide* (available on our web site at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf)) or in the user manual for your board. If the board has programmable RTDs gains, the cbTIn() function sets the appropriate A/D range.

**Specific Errors:** If an OutOfRange or OpenConnection error occurs, the value returned is -9999.0. If a NotReady error occurs, the value returned is -9000.

## cbTInScan()

### Changed R3.3 ID

Reads a range of channels from an analog input board, linearizes them according to temperature sensor type, and returns the temperatures to an array in degrees. The CJC channel, the gain, and temperature sensor type are read from the configuration file. Use the *InstaCal* configuration program to change any of these options.

### Function prototype:

**C/C++:**                    `int cbTInScan(int BoardNum, int LowChan, int HighChan, int Scale, float DataBuffer[], int Options)`

**Visual Basic:**           `Function cbTInScan(ByVal BoardNum&, ByVal LowChan&, ByVal HighChan&, ByVal Scale&, DataBuffer!, ByVal Options&) As Long`

**Delphi:**                   `function cbTInScan(BoardNum:Integer; LowChan:Integer; HighChan:Integer; Scale:Integer; var DataBuffer:Single; Options:Integer):Integer;`

### Arguments:

BoardNum	The number associated with the board when it was installed with the <i>InstaCal</i> configuration program. BoardNum may be 0 to 99.
LowChan	Low mux channel of scan.
HighChan	High mux channel of scan.
Scale	Specifies the temperature scale that the input will be converted to. Choices are CELSIUS, FAHRENHEIT and KELVIN.
DataBuffer	The temperature is returned in degrees. Each element in the array corresponds to a channel in the scan. DataBuffer must be at least large enough to hold HighChan - LowChan + 1 temperature values. Thermocouple resolution is approximately 0.25° C, depending on scale, range and thermocouple type. RTD resolution is 0.1 °C.
Options	Bit fields that control various options. Refer to the constants in the "Options argument values" section below.

### Returns:

[Error code](#) or 0 if no errors.

DataBuffer[] - Temperature values in degrees are returned here for each channel in scan.

### Options argument values:

FILTER	When selected, a smoothing function is applied to temperature readings, very much like the electrical smoothing inherent in all hand held temperature sensor instruments. This is the default. When selected, 10 samples are read and averaged on each channel. The average is the reading returned. Averaging removes normally distributed signal line noise.
NOFILTER	If you use the NOFILTER option then the readings will not be smoothed, and you will see a scattering of readings around a mean.

### Notes:

**Using EXP boards:** For EXP boards, these channel numbers are calculated using the following formula:

- ADChan = A/D channel that is connected to the multiplexer

- MuxChan is a number ranging from 0 to 15 that specifies the channel number on a particular bank of the multiplexer board

$$\text{Chan} = (\text{ADChan} * 16) + (16 + \text{MuxChan})$$

For example, you have an EXP16 connected to a CIO-DAS08 via the CIO-DAS08 channel 0. (Remember, DAS08 channels are numbered 0, 1, 2, 3, 4, 5, 6 & 7). If you connect thermocouples to channels 5, 6, and 7 of the EXP16, the value for LowChan would be  $(0 * 16) + (16 + 5) = 0 + 21 = 21$ , and the value for HighChan would be  $(0 * 16) + (16 + 7) = 0 + 23 = 23$ .

#### Important

For an EXP board connected to an A/D board that does not have programmable gain (DAS08, DAS16, DAS16F), the A/D board range is read from the configuration file (cb.cfg). In most cases, set hardware-selectable ranges to  $\pm 5$  V for thermocouples, and to 0 to 10 V for RTDs. Refer to the board-specific information in the *Universal Library User's Guide* (available on our web site at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf)) or in the user manual for your board. If the board has programmable RTDs gains, the cbTIn() function sets the appropriate A/D range.

**Using 6K-EXP boards:** For 6K-EXP boards, the channel number is calculated using one of the following formulas, where:

- ADChan is the A/D channel that is connected to the multiplexer.
- MuxChan is a number ranging from 0 to 15 that specifies the channel number (Chan) on a particular bank of the multiplexer board.
- If the A/D board has 16 or less single-ended channels:

$$\text{Chan} = (\text{ADChan} * 16) + (16 + \text{MuxChan})$$

For example, you have a 6K-EXP16 connected to a PCI-DAS6052 via the a PCI-DAS6052 channel 0. If you connect a thermocouple to channels 5, 6, and 7 of the 6K-EXP16, the value for LowChan would be  $(0 * 16) + (16 + 5) = 0 + 21 = 21$ , and the value for HighChan would be  $(0 * 16) + (16 + 7) = 0 + 23 = 23$ .

- If the A/D board has 64 single-ended channels and the A/D multiplexer channel is less than or equal to 7:

$$\text{Chan} = (\text{ADChan} * 16) + (64 + \text{MuxChan})$$

For example, you have a 6K-EXP16 connected to a PCI-DAS6031 via the a PCI-DAS6031 channel 7. If you connect a thermocouple to channels 5, 6, and 7 of the 6K-EXP16, the value for LowChan would be  $(7 * 16) + (64 + 5) = 112 + 69 = 181$ , and the value for HighChan would be  $(7 * 16) + (64 + 7) = 112 + 71 = 183$ .

- If the A/D board has 64 single-ended channels and the A/D multiplexer channel is greater than or equal to 32:

$$\text{Chan} = (\text{ADChan} * 16 - 320) + \text{MuxChan}$$

For example, you have a 6K-EXP16 connected to a PCI-DAS6031 via the PCI-DAS6031 channel 32. If you connect a thermocouple to channels 5, 6, and 7 of the 6K-EXP16, the value for LowChan would be  $(32 * 16 - 320) + 5 = 192 + 5 = 197$ , and the value for HighChan would be  $(32 * 16 - 320) + 7 = 192 + 7 = 199$ .

**CJC Channel:** The CJC channel is set in the InstaCal install program. If you have multiple EXP boards, Universal Library will apply the CJC reading to the linearization formula in the following manner:

- First, if you have chosen a CJC channel for the EXP board that the channel you are reading is on, it will use the CJC temp reading from that channel.
- Second, if you have left the CJC channel for the EXP board that the channel you are reading is on to NOT SET, the library will use the CJC reading from the next lower EXP board with a CJC channel selected.

For example: You have four CIO-EXP16 boards connected to a CIO-DAS08 on channel 0, 1, 2 and 3. You choose CIO-EXP16 #1 (connected to CIO-DAS08 channel 0) to have its CJC read on CIO-DAS08 channel 7, AND, you leave the CIO-EXP16's 2, 3 and 4 CJC channels to NOT SET. Result: The CIO-EXP boards will

all use the CJC reading from CIO-EXP16 #1, connected to channel 7 for linearization. As you can see, it is important to keep the CIO-EXP boards in the same case and out of any breezes to ensure a clean CJC reading.

**Important**

In order to understand the functions, refer to the board-specific information in the *Universal Library User's Guide* (available on our web site at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf)) and also in the Readme files installed with the Universal Library. We also urge you to examine and run one or more of the example programs supplied prior to attempting any programming of your own. Following this advice may save you hours of frustration, and wasted time.

This note, which appears elsewhere, is especially applicable to this function. Now is the time to read the board specific information for your board. We suggest that you make a copy of that page to refer to as you read this manual and examine the example programs.

**Specific Errors:** For most boards, if an OUTOFRANGE or OPENCONNECTION error occurs, the value in the array element associated with the channel causing the error returned will be -9999.0 (Refer to the board-specific information in the *Universal Library User's Guide*).



---

# Windows Memory Management Functions

## Introduction

Use the functions explained in this chapter when you run the Windows version of the library. These functions allocate, free and copy to/from Windows global memory buffers.

## cbWinBufAlloc()

Allocates a Windows global memory buffer which can be used with the scan functions and returns a memory handle for it.

### Function prototype:

C/C++:	<code>int cbWinBufAlloc(long NumPoints)</code>
Visual Basic:	<code>Function cbWinBufAlloc(ByVal NumPoints&amp;) As Long</code>
Delphi:	<code>function cbWinBufAlloc(NumPoints:Longint):Integer;</code>

### Arguments:

<code>NumPoints</code>	The size of the buffer to allocate. Specifies how many data points (16-bit integers, NOT bytes) can be stored in the buffer.
------------------------	--

### Returns:

0 if the buffer could not be allocated, or a non-zero integer handle to the buffer.

### Notes:

Unlike most other functions in the library, this function does not return an error code. It returns a Windows global memory handle which can then be passed to the scan functions in the library. If an error occurs the handle will come back as 0 to indicate that the buffer was not allocated.

## cbWinBufAlloc32()

Allocates a Windows global memory buffer for use with 32-bit scan functions, and returns a memory handle for the buffer.

### Function prototype:

C/C++:	<code>int cbWinBufAlloc32(long NumPoints)</code>
Visual Basic:	<code>Function cbWinBufAlloc32(ByVal NumPoints&amp;) As Long</code>
Delphi:	<code>function cbWinBufAlloc32(NumPoints:Longint):Integer;</code>

### Arguments:

<code>NumPoints</code>	The size of the buffer to allocate. Specifies how many data points (32-bit integers, NOT bytes) can be stored in the buffer.
------------------------	--

### Returns:

0 if a buffer could not be allocated, or a non-zero integer handle to the buffer.

### Notes:

Unlike most other functions in the library, this function does not return an error code. It returns a Windows global memory handle which can then be passed to the scan functions in the library. If an error occurs, the handle will come back as 0 to indicate that the buffer was not allocated.

## cbWinBufFree()

Frees a Windows global memory buffer which was previously allocated with the [cbWinBufAlloc\(\)](#) or [cbWinBufAlloc32\(\)](#) function.

### Function prototype:

C/C++:	<code>int cbWinBufFree(int MemHandle)</code>
Visual Basic:	<code>Function cbWinBufFree(ByVal MemHandle&amp;) As Long</code>
Delphi:	<code>function cbWinBufFree (MemHandle:Integer) :Integer;</code>

### Arguments:

MemHandle	A Windows memory handle. This must be a memory handle that was returned by <code>cbWinBufAlloc()</code> or <code>cbWinBufAlloc32()</code> when the buffer was allocated.
-----------	--

### Returns:

[Error code](#) or zero if no errors.

## cbWinArrayToBuf()

Copies data from an array into a Windows memory buffer.

### Function prototype:

C/C++:	<code>int cbWinArrayToBuf(unsigned short *dataArray, int MemHandle, long FirstPoint, long Count)</code>
Visual Basic:	<code>Function cbWinArrayToBuf(DataArray%, ByVal MemHandle&amp;, ByVal FirstPoint&amp;, ByVal Count&amp;) As Long</code>
Delphi:	<code>function cbWinArrayToBuf(var DataArray:Word; MemHandle:Integer; FirstPoint:Longint; Count:Longint):Integer;</code>

### Arguments:

dataArray	The array containing the data to be copied.
MemHandle	This must be a memory handle that was returned by <a href="#">cbWinBufAlloc()</a> when the buffer was allocated. The data will be copied into this buffer.
FirstPoint	Index of first point in memory buffer where data will be copied to.
Count	Number of data points to copy.

### Returns:

[Error code](#) or zero if no errors.

### Notes:

This function copies data from an array to a Windows global memory buffer. This would typically be used to initialize the buffer with data before doing an output scan. Using the `FirstPoint` and `Count` arguments it is possible to fill a portion of the buffer. This can be useful if you want to send new data to the buffer after a BACKGROUND+CONTINUOUS output scan has been started – for example, during circular buffering.

Although this function is available to both Windows C and Delphi programs, it is not necessary, since you can manipulate the memory buffer directly by casting the `MemHandle` returned from `cbWinBufAlloc()` to the appropriate type. This method avoids having to copy the data from an array to a memory buffer. The following example illustrates this method:

```
long Count= 1000;
unsigned short *dataArray=NULL;
int MemHandle = 0;

/*allocate the buffer and cast it to an unsigned short*/
MemHandle = cbWinBufAlloc(Count);
dataArray = (unsigned short*)MemHandle;

/*calculate and store the waveform*/
for(int i=0; i<Count; ++i)
    dataArray[i] = 2047*(1.0 + sin(6.2832*i/Count));

/*output the waveform*/
cbAOutScan(.....,MemHandle,...);

/*free the buffer and NULL the pointer*/
cbWinBufFree(MemHandle);
dataArray = NULL;
```

## cbWinBufToArray()

Copies data from a Windows memory buffer into an array.

### Function prototype:

C/C++:	<code>int cbWinBufToArray(int MemHandle, unsigned short*DataArray, long FirstPoint, long Count)</code>
Visual Basic:	<code>Function cbWinBufToArray(ByVal MemHandle&amp;, DataArray%, ByVal FirstPoint&amp;, ByVal Count&amp;) As Long</code>
Delphi:	<code>function cbWinBufToArray(MemHandle:Integer; var DataArray:Word; FirstPoint:Longint; Count:Longint):Integer;</code>

### Arguments:

MemHandle	This must be a memory handle that was returned by <a href="#">cbWinBufAlloc()</a> when the buffer was allocated. The buffer should contain the data that you want to copy.
DataArray	The array that the data will be copied to.
FirstPoint	Index of first point in memory buffer that data will be copied from.
Count	Number of data points to copy.

### Returns:

[Error code](#) or zero if no errors.

### Notes:

This function copies data from a Windows global memory buffer to an array. This would typically be used to retrieve data from the buffer after executing an input scan function.

Using the `FirstPoint` and `Count` argument it is possible to copy only a portion of the buffer to the array. This can be useful if you want foreground code to manipulate previously collected data while a `BACKGROUND` scan continues to collect new data.

Although this function is available to both Windows C and Delphi programs, it is not necessary, since it is possible to manipulate the memory buffer directly by casting the `MemHandle` returned from `cbWinBufAlloc()` to the appropriate type. This method avoids having to copy the data from the memory buffer to an array. Refer to the following example:

```
/*declare and initialize the variables*/
long Count=1000;
unsigned short *DataArray=NULL;
int MemHandle=0;

/*allocate the buffer and cast it to a pointer to an unsigned short*/
MemHandle = cbWinBufAlloc(Count);
DataArray = (unsigned short*)MemHandle;

/*Scan the waveform data*/
cbAInScan(.....,MemHandle,...);

/*print the results*/
for(int i=0; i<Count; ++i)
    printf("Data[%d]=%d\n", i, DataArray[i]);

/*free the buffer and NULL the pointer*/
cbWinBufFree(MemHandle);
DataArray = NULL;
```

## cbWinBufToArray32()

Copies 32-bit data from a Windows global memory buffer into an array. This function is typically used to retrieve data from the buffer after executing an input scan function.

### Function prototype:

C/C++:	<code>int cbWinBufToArray32(int MemHandle, unsigned long* DataArray, long FirstPoint, long Count)</code>
Visual Basic:	<code>Function cbWinBufToArray32(ByVal MemHandle&amp;, DataArray%, ByVal FirstPoint&amp;, ByVal Count&amp;) As Long</code>
Delphi:	<code>function cbWinBufToArray32(MemHandle:Integer; var DataArray:Longint; FirstPoint:Longint; Count:Longint):Integer;</code>

### Arguments:

MemHandle	The memory handle that was returned by <a href="#">cbWinBufAlloc32()</a> when the buffer was allocated. The buffer should contain the data that you want to copy.
DataArray	The array that the data is copied to.
FirstPoint	The index of the first point in the memory buffer that data is copied from.
Count	The number of data points to copy.

### Returns:

[Error code](#) or 0 if no errors

### Notes:

You can copy only a portion of the buffer to the array using the `FirstPoint` and `Count` argument. This is useful if you want foreground code to manipulate previously collected data while a BACKGROUND scan continues to collect new data.

Although this function is available to both Windows C and Delphi programs, it is not necessary, since you can manipulate the memory buffer directly by casting the `MemHandle` returned from `cbWinBufAlloc32()` to the appropriate type. This method avoids having to copy the data from the memory buffer to an array. Refer to the following example:

```
/*declare and initialize the variables*/
long Count = 1000;
unsigned short *DataArray = NULL;
int MemHandle = 0;

/*allocate the buffer and cast it to a pointer to an unsigned long*/
MemHandle = cbWinBufAlloc32(Count);
DataArray = (unsigned long*)MemHandle;

/*scan in the data*/
cbCInScan(.....,MemHandle,...);

/*print the results*/
for(int i=0; i<Count; ++i)
    printf("Data[%d]=%d\n", i, DataArray[i]);

/*free the buffer and NULL the pointer*/
cbWinBufFree(MemHandle);
DataArray = NULL;
```

---

## Miscellaneous Functions

### Introduction

The functions explained in this chapter do not as a group fit into a single category. They get and set board information, convert units, manage events and background operations, and perform serial communication operations.



## cbDeviceLogin()

Opens a device session with a shared device.

### Function prototype:

C/C++:	<code>int cbDeviceLogin(int BoardNum, char* UserName, char* Password);</code>
Visual Basic:	<code>Function cbDeviceLogin (ByVal BoardNum&amp;, UserName\$, Password\$) As Long</code>

### Arguments:

BoardNum	The number associated with the board when it was installed with the <i>InstaCal</i> configuration program. BoardNum may be 0 to 99.
UserName	A null-terminated string that identifies the user name used to log in to a device session.
Password	A null-terminated string that identifies the password used to log in to a device session.

### Returns:

[Error code](#) or 0 if no errors.

### Notes:

If the user name or password is invalid, the function returns `INVALIDLOGIN`.

If the session is already opened by another user, the function returns `SESSIONINUSE`.

## cbDeviceLogout()

Releases the device session with a shared device.

### Function prototype:

C/C++:	<code>int cbDeviceLogout(int BoardNum);</code>
Visual Basic:	<code>Function cbDeviceLogout (ByVal BoardNum) As Long</code>

### Arguments:

BoardNum	The number associated with the board when it was installed with the <i>InstaCal</i> configuration program. BoardNum may be 0 to 99.
----------	---

### Returns:

[Error code](#) or 0 if no errors.

## cbDisableEvent()

Disables one or more event conditions and disconnects their user-defined handlers.

### Function prototype:

C/C++:	<code>int cbDisableEvent(int BoardNum, unsigned EventType)</code>
Visual Basic:	<code>Function cbDisableEvent(ByVal BoardNum&amp;, ByVal EventType&amp;) as Long</code>
Delphi:	<code>Function cbDisableEvent(BoardNum:Integer; EventType:Integer):Integer;StdCall</code>

### Arguments:

BoardNum	The board number used to indicate which device's event handling will be disabled. BoardNum may be 0 to 99. Refers to the number associated with the board when it was installed with the <i>InstaCal</i> configuration program.
EventType	Specifies one or more event conditions to disable. More than one event type can be specified by bitwise OR'ing the event types. Note that specifying an event that has not been enabled is benign and will not cause any errors. Refer to " <a href="#">EventType argument values</a> " on page 171 for valid EventType settings.  To disable all events in a single call, use ALL_EVENT_TYPES.

### Returns:

[Error code](#) or 0 if no errors.

### Notes:

For most event types, this function cannot be called while any background operations ([cbAInScan\(\)](#), [cbAPretrig\(\)](#), or [cbAOutScan\(\)](#)) are active. Perform a [cbStopBackground\(\)](#) before calling `cbDisableEvent()`. However, for ON\_EXTERNAL\_INTERRUPT events, you can call `cbDisableEvent()` while the board is actively generating events.

### Important

In order to understand the functions, refer to the board-specific information in the *Universal Library User's Guide* (available on our web site at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf)) and also in the Readme files installed with the Universal Library.

We also urge you to examine and run one or more of the example programs supplied prior to attempting any programming of your own. Following this advice may save you hours of frustration, and wasted time. This note, which appears elsewhere, is especially applicable to this function. Now is the time to read the board specific information for your board. We suggest that you make a copy of that page to refer to as you read this manual and examine the example programs.

## cbEnableEvent()

Binds one or more event conditions to a user-defined callback function. Upon detection of an event condition, the user-defined function is invoked with board- and event-specific data. Detection of event conditions occurs in response to interrupts. Typically, this function is used in conjunction with interrupt driven processes such as [cbAInScan\(\)](#), [cbAPretrig\(\)](#), or [cbAOutScan\(\)](#).

### Function prototype:

C/C++:	<code>int cbEnableEvent(int BoardNum, unsigned EventType, unsigned EventParam, void* CallbackFunc, void* UserData)</code>
Visual Basic:	<code>Function cbEnableEvent (ByVal BoardNum&amp;, ByVal EventType&amp;, ByVal EventParam&amp;, ByVal CallbackFunc&amp;, ByRef UserData as Any) as Long</code>
Delphi:	<code>Function cbEnableEvent (BoardNum:Integer; EventType:Integer; EventParam:Integer; CallbackFunc:Pointer; UserData:Pointer):Integer;StdCall</code>

### Arguments:

BoardNum	The board number used to indicate which device will generate the event conditions. BoardNum may be 0 to 99. Refers to the number associated with the board when it was installed with the <i>InstaCal</i> configuration program.
EventType	Specifies one or more event conditions that will be bound to the user-defined callback function. More than one event type can be specified by bitwise OR'ing the event types. Refer to the constants in the "EventType argument values" section below.
EventParam	Additional data required to specify some event conditions such as an ON_DATA_AVAILABLE event or ON_EXTERNAL_INTERRUPT event.  For ON_DATA_AVAILABLE events, EventParam is used to determine the minimum number of samples to acquire during an analog input scan before generating the event. For ON_EXTERNAL_INTERRUPT events, EventParam is used to latch digital bits on supported hardware by setting it to one of the constants in the "EventParam argument values" section on page 172. Most event conditions ignore this value.
CallbackFunc	The address of or pointer to the user-defined callback function to handle the above event type(s). This function must be defined using the standard call (__stdcall) calling convention. Consequently, Visual Basic programs must define their callback functions in standard modules(.bas) and cannot be object methods. C++ programs can define this callback function as either a global function or as a static member function of a class (note that static members do NOT have access to instance specific data).  Refer to the <a href="#">"User Callback function"</a> on page 172 for proper function syntax.
UserData	The address of or pointer to user-defined data that will be passed to the user-defined callback function. This parameter is NOT dereferenced by the library or its drivers; as a consequence, a NULL pointer can be supplied.

### Returns:

[Error code](#) or 0 if no errors.

### EventType argument values:

ON_SCAN_ERROR	Generates an event upon detection of a driver error during BACKGROUND input and output scans. This includes OVERRUN, UNDERRUN, and TOOFEW errors.
---------------	---

ON_EXTERNAL_INTERRUPT	For some digital and counter boards, generates an event upon detection of a pulse at the External Interrupt pin.
ON_PRETRIGGER	For <code>cbAPretrig()</code> , generates an event upon detection of the first trigger.
ON_DATA_AVAILABLE	Generates an event whenever the number of samples acquired during an analog input scan increases by <code>EventParam</code> samples or more. Note that for <code>BLOCKIO</code> scans, events will be generated on packet transfers; for example, even if <code>EventParam</code> is set to 1, events will only be generated every packet-size worth of data (256 samples for the PCI-DAS1602) for aggregate rates greater than 1 kHz for the default <code>cbAInScan()</code> mode.  For <a href="#">cbAPretrig()</a> , the first event is not generated until a minimum of <code>EventParam</code> samples after the pretrigger.
ON_END_OF_AI_SCAN	Generates an event upon completion or fatal error of a <a href="#">cbAInScan()</a> or <code>cbAPretrig()</code> . This event is NOT generated when scans are aborted using <a href="#">cbStopBackground()</a> .
ON_END_OF_AO_SCAN	Generates an event upon completion or fatal error of a <a href="#">cbAOutScan()</a> . This event is not generated when scans are aborted using <code>cbStopBackground()</code> .

**EventParam argument values:**

LATCH_DI	Returns the data that was latched in at the most recent interrupt edge.
LATCH_DO	Latches out the data most recently written to the hardware.

**Notes:**

- This function cannot be called while any background operations ([cbAInScan\(\)](#), [cbAPretrig\(\)](#), or [cbAOutScan\(\)](#)) are active. If a background operation is in progress when `cbEnableEvent()` is called, the function returns an `ALREADYACTIVE` error. Perform a [cbStopBackground\(\)](#) before calling `cbEnableEvent()`.
- Events can be generated no faster than the user callback function can handle them. If an event type becomes multi-signaled before the event handler returns, events are merged. The event handler is called once per event type and is supplied with the event data corresponding to the latest event. In addition, if more than one event type becomes signaled, the event handler for each event type is called in the same order in which they are listed above.
- Events are generated while handling board-generated interrupts. Therefore, using `cbStopBackground()` to abort background operations *does not* generate `ON_END_OF_AI_SCAN` or `ON_END_OF_AO_SCAN` events. However, the event handlers can be called immediately after calling `cbStopBackground()`.
- `cbEnableEvent()` is intended for use with Windows applications. Use with console applications can produce unpredictable results.

## User Callback function

The User Callback function is called as an argument of the [cbEnableEvent\(\)](#) function. You create the function using the prototype shown below. You call the function by passing either its address or a pointer to the function to the `CallbackFunc` argument of the `cbEnableEvent()` function.

### Callback function prototype:

C/C++:	<code>void __stdcall CallbackFunc(int BoardNum, unsigned EventType, unsigned EventData, void* UserData);</code>
Visual Basic:	<code>Sub CallbackFunc(ByVal BoardNum&amp;, ByVal EventType&amp;, ByVal EventData&amp;, ByRef UserData as UserDataType)</code>  where <code>UserData</code> is the data type of the <code>UserData</code> argument passed in to <a href="#">cbEnableEvent()</a> (refer to page 171).
Delphi:	<code>procedure CallbackFunc(BoardNum:Integer; EventType:Integer; EventData:Integer; UserData:Pointer);</code>

### Arguments:

<code>BoardNum</code>	Indicates which board caused the event.
<code>EventType</code>	Indicates which event occurred.
<code>EventData</code>	Board-specific data associated with this event. Returns the value of the <code>EventType</code> as listed in the "EventData argument values" section below.
<code>UserData</code>	The pointer or reference to data supplied by the <code>UserData</code> parameter in <a href="#">cbEnableEvent()</a> (refer to page 171). Note that before use, this parameter must be cast to the same data type as passed in to <code>cbEnableEvent()</code> .

### EventData argument values:

EventType	Value of EventData
<code>ON_SCAN_ERROR</code>	The <a href="#">Error code</a> of the scan error.
<code>ON_EXTERNAL_INTERRUPT</code>	The number of interrupts generated since enabling the <code>ON_EXTERNAL_INTERRUPT</code> event.
<code>ON_PRETRIGGER</code>	The number of pretrigger samples available at time of pretrigger.  This value is invalid for some boards when a <code>TOOFEW</code> error occurs. See board details.
<code>ON_DATA_AVAILABLE</code>	The number of samples acquired since the start of scan.
<code>ON_END_OF_AI_SCAN</code>	The total number of samples acquired upon scan completion or end.
<code>ON_END_OF_AO_SCAN</code>	The total number of samples output upon scan completion or end.

## **cbFlashLED()**

Causes the LED on a USB device to flash.

### **Function prototype:**

C/C++:	<code>int cbFlashLED(int BoardNum);</code>
Visual Basic:	<code>Function cbFlashLED(ByVal BoardNum&amp;) as Long</code>
Delphi:	<code>function cbFlashLED(BoardNum:Integer):Integer;</code>

### **Arguments:**

BoardNum	The board number of the USB device whose LED will flash.
----------	--

## cbFromEngUnits()

Converts a single precision voltage (or current) value in engineering units to an integer count value. This function is typically used to obtain a data value from a voltage value for output to a D/A with functions such as `cbAOut()`.

### Function prototype:

**C/C++:** `int cbFromEngUnits(int BoardNum, int Range, float EngUnits, unsigned short *DataVal)`

**Visual Basic:** `Function cbFromEngUnits(ByVal BoardNum&, ByVal Range&, ByVal EngUnits!, DataVal%) As Long`

**Delphi:** `function cbFromEngUnits(BoardNum:Integer; Range:Integer; EngUnits:Single; var DataVal:Word):Integer;`

### Arguments:

BoardNum	The board number associated with the board when it was installed with the <i>InstaCal</i> configuration program. This function uses the board number to determine the resolution and polarity values to use in the conversion. BoardNum may be 0 to 99.
Range	The voltage (or current) range to use for the conversion to counts. When using this function to obtain a value to send to a D/A board, keep in mind that some D/A boards have programmable voltage ranges, and others set the voltage range via switches on the board. In either case, the desired range must be passed to this function. Refer to Table 3 on page 28 for a list of valid range settings.
EngUnits	The single precision voltage (or current) value to use for the conversion to counts. Set the value to be within the range specified by the Range argument.
DataVal	The function returns an integer count to this variable that is equivalent to the EngUnits argument using the resolution of the D/A on the board referenced by BoardNum (if any).

### Returns:

[Error code](#) or 0 if no errors.

DataVal – the integer count equivalent to EngUnits is returned here.

### Notes:

This function is not supported for hardware with resolution greater than 16 bits.

The default resolution of this function is 12 bits, so if the device referenced by BoardNum has neither analog input nor analog output, the result will be a 12 bit conversion.

If the device referenced by BoardNum has both analog input and analog output, the resolution and transfer function of the D/A converter on the device is used.

## cbGetBoardName()

Returns the board name of a specified board.

### Function prototype:

C/C++:	<code>int cbGetBoardName(int BoardNum, char *BoardName)</code>
Visual Basic:	<code>Function cbGetBoardName(ByVal BoardNum&amp;, ByVal BoardName\$) As Long</code>
Delphi:	<code>function cbGetBoardName(BoardNum:Integer; BoardName:PChar):Integer;</code>

### Arguments:

BoardNum	Refers either to the board number associated with a board when it was installed, or GETFIRST or GETNEXT. BoardNum may be 0 to 99, GETFIRST or GETNEXT
BoardName	A null-terminated string variable that the board name will be returned to. This string variable must be pre-allocated to be at least as large as BOARDNAMELEN. This size is guaranteed to be large enough to hold the longest board name string. The "Appendix" in the <i>Universal Library User Guide</i> (available on our web site at <a href="http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf">www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf</a> ) lists the board names and associated device ID codes.

### Returns:

[Error code](#) or 0 if no errors.

BoardName - return string containing the board name.

### Notes:

There are two distinct ways of using this function:

- Pass a board number as the BoardNum argument. The string that is returned describes the board type of the installed board.
- Set BoardNum to GETFIRST or GETNEXT to get a list of all board types that are supported by the library. Set BoardNum to GETFIRST to get the first board type in the list of supported boards. Subsequent calls with Board=GETNEXT returns each of the other board types supported by the library. When you reach the end of the list, BoardName is set to an empty string. Refer to the ulgt04 example program in the installation directory for more details.



## cbGetStatus()

Returns the status about the background operation currently running.

### Function prototype:

C/C++:	<code>int cbGetStatus(int BoardNum, int *Status, long *CurCount, long *CurIndex, int FunctionType)</code>
Visual Basic:	<code>Function cbGetStatus(ByVal BoardNum&amp;, Status%, CurCount&amp;, CurIndex&amp;, FunctionType&amp;) As Long</code>
Delphi:	<code>function cbGetStatus(BoardNum:Integer; var Status:SmallInt; var CurCount:Longint; var CurIndex:Longint; FunctionType:Integer):Integer;</code>

### Arguments:

BoardNum	The board number associated with the board when it was installed with the <i>InstaCal</i> configuration program. BoardNum may be 0 to 99.
Status	The Status argument indicates whether or not a background process is currently executing.
CurCount	The CurCount argument specifies how many points have been input or output since the Background process started. Use it to gauge how far along the operation is towards completion. Generally, CurCount returns the total number of samples transferred between the DAQ board and the Windows data buffer at the time <code>cbGetStatus()</code> was called.

When you set both the CONTINUOUS and BACKGROUND options, CurCount's behavior depends on the board model. Refer to the board-specific information in the *Universal Library User's Guide* for the behavior of your board.

With recent MCC DAQ designs, the CurCount argument continually increases in increments of the packet size as Windows' circular data buffer recycles, until it reaches  $2^{31}$ . Since the Count argument is a signed integer, at 2,147,483,647 + 1, the Count argument rolls back to a negative number (-2,147,483,647). The Count argument resumes incrementing, eventually reaching 0 and increasing back up to 2,147,483,647.

The CurIndex argument is usually more useful than the CurCount argument in managing data collected when you set both the CONTINUOUS and BACKGROUND options.

CurIndex	The CurIndex argument is an index into the Windows data buffer. This index points to the start of the last completed channel scan that was transferred between the DAQ board and the Windows data buffer. If a scan is running but no points in the buffer have been transferred, CurIndex equals -1 in most cases.
----------	---

For CONTINUOUS operations, CurIndex rolls over when the Windows data buffer is full. This rollover indicates that "new" data is now overwriting "old" data. Your goal is to process the old data before it gets overwritten. You can keep ahead of the data flow by copying the old data out of the buffer before new data overwrites it.

The CurIndex argument can help you access the most recently transferred data. Your application does not have to process the data exactly when it becomes available in the buffer – in fact, you should avoid doing so unless absolutely necessary. The CurIndex argument generally increments by the packet size, but in some cases the CurIndex increment can vary within the same scan. One instance of a variable increment is when the packet size is not evenly divisible by the number of channels.

You should determine the best size of the "chunks" of data that your application can most efficiently process, and then periodically check on the `CurIndex` argument value to determine when that amount of additional data has been transferred.

Refer to the *Universal Library User's Guide* for information on your board, particularly when using Pre-Trigger.

`FunctionType` Specifies which scan to retrieve status information about. Set it to one of the constants in the "FunctionType argument values" section below.

#### Returns:

[Error code](#) or 0 if no errors

`Status` - Returns the status of the operation:

`IDLE` - No background operation is running.

`RUNNING` - Background operation is still executing.

`CurCount` - Returns the current number of samples collected.

`CurIndex` - Returns the Current sample index.

#### FunctionType argument values:

<code>AIFUNCTION</code>	Specifies analog input scans started with <a href="#">cbAInScan()</a> or <a href="#">cbAPretrig()</a> .
<code>AOFUNCTION</code>	Specifies analog output scans started with <a href="#">cbAOutScan()</a> .
<code>DIFUNCTION</code>	Specifies digital input scans started with <a href="#">cbDInScan()</a> .
<code>DOFUNCTION</code>	Specifies digital output scans started with <a href="#">cbDOutScan()</a> .
<code>CTRFUNCTION</code>	Specifies counter background operations started with <a href="#">cbCStoreOnInt()</a> or <a href="#">cbCInScan()</a> .
<code>DAQIFUNCTION</code>	Specifies a synchronous input scan started with <a href="#">cbDagInScan()</a> .
<code>DAQOFUNCTION</code>	Specifies a synchronous output scan started with <a href="#">cbDagOutScan()</a> .

## cbGetTCValues()

Converts raw thermocouple data collected using the [cbDagInScan\(\)](#) function to data on a temperature scale (Celsius, Fahrenheit or Kelvin).

### Function prototype:

C/C++:	<pre>int cbGetTCValues(int BoardNum, short *ChanArray, short *ChanTypeArray, int ChanCount, int MemHandle, int FirstPoint, long Count, int Scale, float *TempValArray)</pre>
Visual Basic:	<pre>Function cbGetTCValues(ByVal BoardNum&amp;, ChanArray%, ChanTypeArray%, ByVal ChanCount&amp;, ByVal MemHandle&amp;, ByVal FirstPoint&amp;, ByVal Count&amp;, ByVal CBScale&amp;, TempValArray!) As Long</pre>
Delphi:	<pre>function cbGetTCValues(BoardNum:Integer; var ChanArray:SmallInt; var ChanTypeArray:SmallInt; ChanCount:Integer; MemHandle:Integer; FirstPoint:Integer; CBCount:LongInt; Scale:Integer; var TempValArray:Single):Integer;</pre>

### Arguments:

BoardNum	The board number used to collect the data. BoardNum may be 0 to 99. Refers to the number associated with the board used to collect the data when it was installed with the <i>InstaCal</i> configuration program. The specified board must support synchronous input.
ChanArray	Array containing channel values. Valid channel values are analog and temperature input channels and digital ports. ChanArray must match the channel array used with the <a href="#">cbDagInScan()</a> function.
ChanTypeArray	Array containing channel types. Each element of this array defines the type of the corresponding element in the ChanArray. ChanTypeArray must match the channel type settings used with the <a href="#">cbDagInScan()</a> function.
ChanCount	Number of elements in ChanArray.
MemHandle	This must be a memory handle that was returned by <a href="#">cbWinBufAlloc()</a> when the buffer was allocated. The buffer should contain the data that you want to convert.
FirstPoint	The index into the raw data memory buffer that holds the first sample of the first channel to be converted. The index into the raw memory is (FirstPoint x ChanCount) so that converted data always starts with the first channel specified in the scan. For example, if firstPoint is 14 and the number of channels is 8, the index of the first converted sample is 112.
Count	The number of samples per channel to convert to engineering units. Count should not exceed Windows buffer size / ChanCount – FirstPoint.
Scale	Specifies the temperature scale that the input will be converted to. Choices are CELSIUS, FAHRENHEIT and KELVIN.
TempValArray	The array to hold the converted data. This array must be allocated by the user, and must be large enough to hold count samples x the number of temperature channels.

### Returns:

[Error code](#) or 0 if no errors

TempValArray – Converted data.

## cbInByte()

Reads a byte from a hardware register on a board.

### Function prototype:

C/C++:	<code>int cbInByte(int BoardNum, int PortNum)</code>
Visual Basic:	<code>Function cbInByte(ByVal BoardNum&amp;, ByVal PortNum&amp;) As Long</code>
Delphi:	<code>function cbInByte(BoardNum:Integer; PortNum:Integer):Integer;</code>

### Arguments:

BoardNum	The number associated with the board when it was installed with the <i>InstaCal</i> configuration program. BoardNum may be 0 to 99.
PortNum	Register within the board. Boards are set to a particular base address. The registers on the boards are at addresses that are offsets from the base address of the board (BaseAdr + 0, BaseAdr + 2, etc). Set this argument to the offset for the desired register. This function takes care of adding the base address to the offset, so that the board's address can be changed without changing the code.

### Returns:

The current value of the specified register

### Notes:

`cbInByte()` is used to read 8 bit ports. [cbInWord\(\)](#) is used to read 16-bit ports.

This function was designed for use with ISA bus boards. Use with PCI bus boards is not recommended.

## cbInWord()

Reads a word from a hardware register on a board.

### Function prototype:

C/C++:	<code>int cbInWord(int BoardNum, int PortNum)</code>
Visual Basic:	<code>Function cbInWord(ByVal BoardNum&amp;, ByVal PortNum&amp;) As Long</code>
Delphi:	<code>function cbInWord(BoardNum:Integer; PortNum:Integer):Integer;</code>

### Arguments:

BoardNum	The number associated with the board when it was installed with the <i>InstaCal</i> configuration program. BoardNum may be 0 to 99.
PortNum	Register within the board. Boards are set to a particular base address. The registers on the boards are at addresses that are offsets from the base address of the board (BaseAdr + 0, BaseAdr + 2, etc). Set this argument to the offset for the desired register. This function takes care of adding the base address to the offset, so that the board's address can be changed without changing the code.

### Returns:

The current value of the specified register.

### Notes:

[cbInByte\(\)](#) is used to read 8-bit ports. `cbInWord()` is used to read 16-bit ports.

This function was designed for use with ISA bus boards. Use with PCI bus boards is not recommended.

## cbOutByte()

Writes a byte to a hardware register on a board.

### Function prototype:

C/C++:	<code>int cbOutByte(int BoardNum, int PortNum, int PortVal)</code>
Visual Basic:	<code>Function cbOutByte(ByVal BoardNum&amp;, ByVal PortNum&amp;, ByVal PortVal%) As Long</code>
Delphi:	<code>function cbOutByte(BoardNum:Integer; PortNum:Integer; PortVal:Integer):Integer;</code>

### Arguments:

BoardNum	The number associated with the board when it was installed with the <i>InstaCal</i> configuration program. BoardNum may be 0 to 99.
PortNum	A register within the board. Boards are set to a particular base address. The registers on the boards are at addresses that are offsets from the base address of the board (BaseAdr + 0, BaseAdr + 2, etc).  Set this argument to the offset for the desired register. This function takes care of adding the base address to the offset, so that the board's address can be changed without changing the code.
PortVal	The value that is written to the register.

### Returns:

[Error code](#) or 0 if no errors

### Notes:

`cbOutByte()` is used to write to 8-bit ports. [cbOutWord\(\)](#) is used to write to 16-bit ports.

This function was designed for use with ISA bus boards, and is not recommended for use with PCI-bus boards.

## cbOutWord()

Writes a word to a hardware register on a board.

### Function prototype:

C/C++:	<code>int cbOutWord(int BoardNum, int PortNum, int PortVal)</code>
Visual Basic:	<code>Function cbOutByte(ByVal BoardNum&amp;, ByVal PortNum&amp;, ByVal PortVal%) As Long</code>
Delphi:	<code>function cbOutWord(BoardNum:Integer; PortNum:Integer; PortVal:Integer):Integer;</code>

### Arguments:

BoardNum	The number associated with the board when it was installed with the <i>InstaCal</i> configuration program. BoardNum may be 0 to 99.
PortNum	A register within the board. Boards are set to a particular base address. The registers on the boards are at addresses that are offsets from the base address of the board (BaseAdr + 0, BaseAdr + 2, etc).  Set this argument to the offset for the desired register. This function takes care of adding the base address to the offset, so that the board's address can be changed without changing the code.
PortVal	The value that is written to the register.

### Returns:

[Error code](#) or 0 if no errors

### Notes:

[cbOutByte\(\)](#) is used to write to 8-bit ports. `cbOutWord()` is used to write to 16-bit ports.

This function was designed for use with ISA bus boards, and is not recommended for use with PCI bus boards.

## cbRS485()

Sets the direction of RS-485 communications port buffers.

### Function prototype:

C/C++:	<code>int cbRS485(int BoardNum, int Transmit, int Receive)</code>
Visual Basic:	<code>Function cbRS485(ByVal BoardNum&amp;, ByVal Transmit&amp;, ByVal Receive&amp;) As Long</code>
Delphi:	<code>function cbRS485(BoardNum:Integer; Transmit:Integer; Receive:Integer):Integer;</code>

### Arguments:

BoardNum	The number associated with the board when it was installed with the <i>InstaCal</i> configuration program. BoardNum may be 0 to 99.
Transmit	Set to ENABLED or DISABLED (CBENABLED or CBDISABLED in Visual Basic or Delphi). The transmit RS-485 line driver is turned on. Data written to the RS-485 UART chip is transmitted to the cable connected to that port.
Receive	Set to ENABLED or DISABLED (CBENABLED or CBDISABLED in Visual Basic or Delphi). The receive RS-485 buffer is turned on. Data present on the cable connected to the RS-485 port is received by the UART chip.

### Returns:

[Error code](#) or 0 if no errors

### Notes:

You can simultaneously enable or disable the transmit and receive buffers. If both are enabled, data written to the port is also received by the port. For a complete discussion of RS485 network construction and communication, refer to the CIO-COM485 or PCM-COM485 hardware manual.



## cbStopBackground()

Stops one or more subsystem background operations that are in progress for the specified board. Use this function to stop any function that is running in the background. This includes any function that was started with the BACKGROUND option, as well as [cbCStoreOnInt\(\)](#) (which always runs in the background).

Execute `cbStopBackground()` after normal termination of all background functions to clear variables and flags.

### Function prototype:

C/C++:	<code>int cbStopBackground(int BoardNum, int FunctionType)</code>
Visual Basic:	<code>Function cbStopBackground(ByVal BoardNum&amp;, ByVal FunctionType&amp;) As Long</code>
Delphi:	<code>function cbStopBackground(BoardNum:Integer, FunctionType:Integer):Integer;</code>

### Arguments:

BoardNum	The number associated with the board when it was installed with the <i>InstaCal</i> configuration program. BoardNum may be 0 to 99.
FunctionType	Specifies which background operation to stop. Set it to one of the constants in the "FunctionType argument values" section below.

### Returns:

[Error code](#) or 0 if no errors

### FunctionType argument values:

AIFUNCTION	Specifies analog input scans started with <a href="#">cbAInScan()</a> or <a href="#">cbAPretrig()</a>
AOFUNCTION	Specifies analog output scans started with <a href="#">cbAOutScan()</a> .
DIFUNCTION	Specifies digital input scans started with <a href="#">cbDInScan()</a> .
DOFUNCTION	Specifies digital output scans started with <a href="#">cbDOutScan()</a> .
CTRFUNCTION	Specifies counter background operations started with <a href="#">cbCStoreOnInt()</a> or <a href="#">cbCInScan()</a> .
DAQIFUNCTION	Specifies a synchronous input scan started with <a href="#">cbDagInScan()</a> .
DAQOFUNCTION	Specifies a synchronous output scan started with <a href="#">cbDagOutScan()</a> .

## cbToEngUnits()

Converts an integer count value to an equivalent single precision voltage (or current) value. This function is typically used to obtain a voltage value from data received from an A/D with functions such as `cbAIn()`.

### Function prototype:

C/C++:	<code>int cbToEngUnits(int BoardNum, int Range, unsigned short DataVal, float *EngUnits)</code>
Visual Basic:	<code>Function cbToEngUnits(ByVal BoardNum&amp;, ByVal Range&amp;, ByVal DataVal%, EngUnits!) As Long</code>
Delphi:	<code>function cbToEngUnits(BoardNum:Integer; Range:Integer; DataVal:Word; var EngUnits:Single):Integer;</code>

### Arguments:

BoardNum	The board number associated with the board when it was installed with the <i>InstaCal</i> configuration program. This function uses the board number to determine the resolution and polarity values to use for the conversion. BoardNum may be 0 to 99.
Range	Voltage (or current) range to use for the conversion to engineering units. When using this function to obtain engineering units from a value received from an A/D board, keep in mind that some A/D boards have programmable voltage ranges, and others set the voltage range via switches on the board. In either case, the desired range must be passed to this function. Refer to Table 3 on page 28 for a list of valid range settings.
DataVal	An integer count value (typically, one returned from an A/D board).
EngUnits	The single precision voltage (or current) value that is equivalent to DataVal is returned to this variable. The value will be within the range specified by the Range argument using the resolution of the A/D on the board referenced by BoardNum (if any).

### Returns:

[Error code](#) or 0 if no errors.

EngUnits – the engineering units value equivalent to DataVal is returned to this variable.

### Notes:

This function is not supported for hardware with resolution greater than 16 bits.

The default resolution of this function is 12 bits, so if the device referenced by BoardNum has neither analog input nor analog output, the result will be a 12 bit conversion.

If the device referenced by BoardNum has both analog input and analog output, the resolution and transfer function of the A/D converter on the device is used.

# **Universal Library for .NET**

## **Classes, Methods, and Properties**

---

## UL for .NET Class Library Overview

The Microsoft .NET platform provides a framework that allows for the development of Windows applications using a wide range of new programming languages. These languages include VB .NET, C#, managed C++, JScript, and any other language that is compliant with the .NET Common Language Runtime (CLR). The CLR is a multi-language execution environment.

The interface to the Universal Library consists of standard "C" functions. These functions are not CLR-compliant. Therefore, the Universal Library for .NET was developed. This library enables the various .NET programming languages to call into the Universal Library.

The Universal Library for .NET consists of a set of classes. For the most part, the methods within each class have a corresponding function in the standard UL. Each UL for .NET method has virtually the same parameter set as their UL counterparts.

### MccDaq namespace

The MccDaq namespace contains the classes and enumerated constants by which your UL for .NET applications can access the Universal Library data types and functions.

### MccDaq classes

The MccDaq namespace contains five main classes:

- MccBoard class
- ErrorInfo class
- MccService class
- GlobalConfig class
- DataLogger class

The MccDaq namespace also contains the following four secondary classes:

cBoardConfig	Contains all of the members for setting and getting board-level configuration.
cCtrConfig	Contains all of the members for setting and getting the counter-level configuration of a board.
cDioConfig	Contains all of the members for getting the digital configuration of a board.
cExpansionConfig	Contains all of the members for setting and getting expansion board configuration.

These classes include methods that are accessible from properties of the [MccBoard class](#) (explained below).

### MccBoard class

The MccBoard class provides access to all of the methods for data acquisition and properties providing board information and configuration for a particular board.

The MccBoard class is a member of the MccDaq namespace. Refer to the "[MccDaq namespace](#)" above for an explanation of the MccDaq namespace.

**Class constructors:**

The MccBoard class provides two constructors; one which accepts a board number argument and one with no arguments.

The following code examples demonstrate how to create a new instance of the MccBoard class using the latter version with a default board number of 0.

```
VB .NET:      Private DaqBoard As MccDaq.MccBoard
               DaqBoard = New MccDaq.MccBoard()

C# .NET:      private MccDaq.MccBoard DaqBoard;
               DaqBoard = new MccDaq.MccBoard();
```

The following code examples demonstrate how to create a new instance of the MccBoard class with the board number passed to it.

```
VB .NET:      Private DaqBoard As MccDaq.MccBoard
               DaqBoard = New MccDaq.MccBoard(BoardNumber)

C# .NET:      private MccDaq.MccBoard DaqBoard;
               DaqBoard = new MccDaq.MccBoard(BoardNumber);
```

**Properties and methods**

The MccBoard class includes close to 100 methods for data acquisition. The MccBoard class methods are equivalents of the function calls used in the standard Universal Library. The MccBoard class methods have virtually the same parameter set as their UL counterparts.

The MccBoard class also includes six properties that you can use to examine or change the configuration of your board. The configuration information for all boards is stored in the CB.CFG file, and is loaded from CB.CFG by all programs that use the library.

Each MccBoard property and method is explained briefly later in this chapter, and in detail in the remaining chapters of the reference manual.

**ErrorInfo class**

Contains all of the members for storing and reporting error codes and messages. This class also includes error code enumerated constants, which define the error number and associated message which can be returned when you call a method.

Most UL for .NET methods return ErrorInfo objects. Error information is stored internally on the return from calling the low-level UL function. The error is reported when the user calls the class library methods.

The ErrorInfo class is a member of the MccDaq namespace. Refer to the "[MccDaq namespace](#)" section on page 188 for an explanation of the MccDaq namespace.

**Enumerated constants**

ErrorCode	Lists the named constants for all error codes. For a full explanation of the error associated with each error code and error constant, refer to the " <a href="#">Error Codes</a> " appendix on page 401.
-----------	---

**Properties and methods**

The ErrorInfo class also includes two properties that you can use to examine error information. Each property is explained briefly later in this chapter, and in detail in the remaining chapters of the reference manual.

## MccService class

Contains all of the members for calling utility UL functions.

The MccService class is a member of the MccDaq namespace. Refer to the "[MccDaq namespace](#)" on page 188 for an explanation of the MccDaq namespace.

### Methods

The MccService class contains nine static methods. You do not need to create an instance of the MccService class to call these methods.

## GlobalConfig class

Contains all of the members for getting global board configuration information.

The GlobalConfig class is a member of the MccDaq namespace. Refer to the "[MccDaq namespace](#)" on page 188 for an explanation of the MccDaq namespace.

### Properties and methods

The GlobalConfig class includes three properties that you can use to examine global board configuration information. Each property is explained briefly later in this chapter, and in detail in the remaining chapters of the reference manual.

## DataLogger class

Contains all of the members for reading and converting binary log files.

The DataLogger class is a member of the MccDaq namespace. Refer to the "[MccDaq namespace](#)" on page 188 for an explanation of the MccDaq namespace.

### Property and methods

The DataLogger class provides one property to get a reference to the file name associated with the current instance of the DataLogger.

The DataLogger class includes 14 methods used to read and convert the data contained in a binary log file. These methods are equivalents of the function calls used in the standard Universal Library. The methods have virtually the same parameter set as their UL counterparts.

The DataLogger property and methods are explained briefly later in this chapter, and in detail in the remaining chapters of the reference manual.

## Analog I/O methods

The analog I/O methods available from the [MccBoard class](#) are explained below. These methods perform analog input and output and convert analog data.

- **MccBoard.AIn()** - Takes a single reading from an analog input channel (A/D).
- **MccBoard.AInScan()** - Repeatedly scans a range of analog input (A/D) channels. You can specify the channel range, the number of iterations, the sampling rate, and the A/D range. The data that is collected is stored in an array.
- **MccBoard.ALoadQueue()** - Loads a series of channel/gain pairs into an A/D board's queue. These channel/gains are used with all subsequent analog input methods.
- **MccBoard.AOut()** - Outputs a single value to an analog output (D/A).

- **MccBoard.AOutScan()** - Repeatedly scans a range of analog output (D/A) channels. You can specify the channel range, the number of iterations, and the rate. The data from consecutive elements of an array are sent to each D/A channel in the scan.
- **MccBoard.APretrig()** - Repeatedly scans a range of analog input (A/D) channels waiting for a trigger signal. When a trigger occurs, it returns the specified number of samples and points before the trigger occurred. You can specify the channel range, the sampling rate, and the A/D range. All of the data that is collected is stored in an array.
- **MccBoard.ATrig()** - Reads analog input and waits until it goes above or below a specified threshold. When the trigger condition is met, the current sample is returned.
- **MccBoard.AConvertData()** - Converts analog data from data plus channel tags to separate data and channel tags.

Each raw sample from analog input is a 16-bit value. On some 12-bit A/D boards it consists of a 12-bit A/D value along with a four bit channel number. This method is not intended for use with 16-bit A/D boards.

This conversion is done automatically by the [MccBoard.AIn\(\)](#) method. It can also be done automatically by the [MccBoard.AInScan\(\)](#) method with the `ConvertData` option. In some cases though, it may be useful or necessary to collect the data and then do the conversion sometime later. The [MccBoard.AConvertData\(\)](#) method takes a buffer full of unconverted data and converts it.

- **MccBoard.ACalibrateData()** - Calibrates analog data.

Each raw sample from a board with software calibration factors that must be applied to the sample may be acquired and calibrated, then passed to an array. Alternatively, they can be acquired then passed to the array without calibration. This technique applies the calibration factors to an array of data after the acquisition is complete. When this second technique is used, `ACalibrateData()` may be used to apply the calibration factors to an array of data after the acquisition is complete. The only case where you would withhold calibration until after the acquisition run was complete is on slower CPUs, or when the processing time is at a premium. Applying calibration factors in real time on a per sample basis does eat up machine cycles.

To disable the automatic calibration so that you may apply the calibration later, specify the `NoCalibrateData` option when collecting data with the [MccBoard.AInScan\(\)](#) method.

- **MccBoard.AConvertPretrigData()** - Converts and re-orders pre-trigger data from data plus channel tags to separate data and channel tags.

For devices with a hardware implementation of pretrigger, when data is collected with the [MccBoard.APretrig\(\)](#) method the same data conversion needs to be done as is performed by the [MccBoard.AConvertData\(\)](#) method. There is a further complication because `MccBoard.APretrig()` collects analog data into an array. It treats the array like a circular buffer. While it is waiting for the trigger to occur, it fills the array. When it gets to the end it resets to the start and begins again. When the trigger signal occurs it continues collecting data into the circular buffer until the requested number of samples have been collected.

When the data acquisition is complete, all of the data is in the array but it is in the wrong order. The first element of the array does not contain the first data point. The data has to be rotated in the correct order.

This conversion can be done automatically by the [MccBoard.APretrig\(\)](#) method with the `ConvertData` option. In some cases though, it may be useful or necessary to collect the data and then do the conversion sometime later. The [MccBoard.AConvertPretrigData\(\)](#) method takes a buffer full of unconverted data and converts it.

- **VIn()** - Reads an A/D input channel, and returns a voltage value.
- **VOut()** - Sets the value of a D/A output.

## Configuration methods and properties

The configuration methods and properties available from the [MccBoard](#) class, [cBoardConfig](#) class, [cCtrConfig](#) class, [cDioConfig](#) class, and the [cExpansionConfig](#) class are explained below.

The configuration information for all boards is stored in the configuration file CB.CFG. This information is loaded from CB.CFG by all programs that use the library. The library includes the following classes and methods that retrieve or change configuration options.

- **MccBoard.BoardNum** property – Retrieves the number of the board associated with an instance of the MccBoard class.
- **MccBoard.GetSignal()** - Retrieves the configured auxiliary or DAQ Sync connection and polarity for the specified timing and control signal. This method is intended for advanced users.
- **MccBoard.SelectSignal()** - Configures timing and control signals to use specific auxiliary or DAQ Sync connections as a source or destination. This method is intended for advanced users.
- **MccBoard.SetTrigger()** - Sets up trigger parameters used with the ExtTrigger option for [MccBoard.AInScan\(\)](#).
- **MccBoard.BoardConfig** property - Gets an instance of a cBoardConfig object.
- **MccBoard.BoardConfig.DACUpdate()** - Updates the voltage values on analog output channels.
- **MccBoard.BoardConfig.GetAdRetrigCount()** - Gets the number of samples to acquire during a trigger event when ScanOptions.RetrigMode is set.
- **MccBoard.BoardConfig.GetBaseAdr()** - Gets the base address of a board.
- **MccBoard.BoardConfig.GetBoardType()** - Gets the unique number (device ID) assigned to the board (between 0 and 8000h) indicating the type of board installed.
- **MccBoard.BoardConfig.GetCiNumDevs()** - Gets the number of counter devices on the board.
- **MccBoard.BoardConfig.GetClock()** - Gets the clock frequency in MHz (40, 10, 8, 6, 5, 4, 3, 2, 1), or 0 for not supported.
- **MccBoard.BoardConfig.GetDACStartup()** - Gets the board's configuration register STARTUP bit setting.
- **MccBoard.BoardConfig.GetDACUpdateMode()** - Gets the update mode for a digital-to-analog converter (DAC).
- **MccBoard.BoardConfig.GetDeviceId()** – Gets the name that identifies the instance of a device.
- **MccBoard.BoardConfig.GetDeviceNotes()** – Gets the device notes that are stored in the device's memory.
- **MccBoard.BoardConfig.GetDInMask()** - Determines the bits on a specified port that are configured for input.
- **MccBoard.BoardConfig.GetDiNumDevs()** - Gets the number of digital devices on the board.
- **MccBoard.BoardConfig.GetDmaChan()** - Gets the DMA channel (0, 1 or 3) set for the board.
- **MccBoard.BoardConfig.GetDOutMask()** - Determines the bits on a specified port that are configured for output.
- **MccBoard.BoardConfig.GetDtBoard()** - Gets the number of the board with the DT connector used to connect to external memory boards.
- **MccBoard.BoardConfig.GetIntLevel()** - Gets the interrupt level set for the board (0 for none, or 1 to 15).
- **MccBoard.BoardConfig.GetNumAdChans()** - Gets the number of A/D channels
- **MccBoard.BoardConfig.GetNumDaChans()** - Gets the number of D/A channels.
- **MccBoard.BoardConfig.GetNumExps()** - Gets the number of expansion boards.
- **MccBoard.BoardConfig.GetNumIoPorts()** - Gets the number of I/O ports used by the board.



- **MccBoard.BoardConfig.GetPANID()** - Gets the Personal Area Network (PAN) identifier for wireless communication.
- **MccBoard.BoardConfig.GetRange()** - Gets the selected voltage range.
- **MccBoard.BoardConfig.GetRFChannel()** - Gets the RF channel number that a wireless device uses to communicate.
- **MccBoard.BoardConfig.GetRSS()** - Gets the signal strength in dBm of a signal received by a remote device.
- **MccBoard.BoardConfig.GetUsesExps()** - Gets the True/False value indicating support of expansion boards.
- **MccBoard.BoardConfig.GetWaitState()** - Gets the value of the Wait State jumper (1-enabled, 0-disabled).
- **MccBoard.BoardConfig.SetAdRetrigCount()** - Sets the number of samples to acquire during a trigger event when `ScanOptions.RetrigMode` is set.
- **MccBoard.BoardConfig.SetBaseAdr()** - Sets the base address of a board
- **MccBoard.BoardConfig.SetClock()** - Sets the clock source by the frequency (40, 10, 8, 6, 5, 4, 3, 2, 1), or 0 for not supported.
- **MccBoard.BoardConfig.SetDACStartup()** - Sets the board's configuration register STARTUP bit to 0 or 1 to enable/disable the storing of digital-to-analog converter (DAC) startup values.
- **MccBoard.BoardConfig.SetDACUpdateMode()** - Sets the update mode for a digital-to-analog converter (DAC).
- **MccBoard.BoardConfig.SetDeviceId()** - Sets the name that identifies the instance of a device.
- **MccBoard.BoardConfig.SetDeviceNotes()** - Sets the device notes to store in a device's memory.
- **MccBoard.BoardConfig.SetDmaChan()** - Sets the DMA channel (0, 1 or 3).
- **MccBoard.BoardConfig.SetIntLevel()** - Sets the interrupt level: 0 for none, or 1 to 15.
- **MccBoard.BoardConfig.SetNumAdChans()** - Sets the number of A/D channels available on the board.
- **MccBoard.BoardConfig.SetPanID()** - Sets the Personal Area Network (PAN) identifier used for wireless communication.
- **MccBoard.BoardConfig.SetRange()** - Sets the selected voltage range.
- **MccBoard.BoardConfig.SetRFChannel()** - Sets the RF channel number used for wireless communications.
- **MccBoard.BoardConfig.SetWaitState()** - Sets the value of the Wait State jumper (1 = enabled, 0 = disabled).
- **MccBoard.CtrConfig** property - Gets an instance of a `cCtrConfig` object.
- **MccBoard.CtrConfig.GetCtrType()** - Gets the counter device number of counter type specified with the `configVal` parameter.
- **MccBoard.DioConfig** property - Gets an instance of a `cDioConfig` object.
- **MccBoard.DioConfig.GetConfig()** - Gets the configuration of a digital device (digital input or digital output).
- **MccBoard.DioConfig.GetCurVal()** - Gets the current value of digital outputs.
- **MccBoard.DioConfig.GetDevType()** - Gets the device type of the digital port (AUXPORT, FIRSTPORTA, etc.).
- **MccBoard.DioConfig.GetDInMask()** - Determines the bits on a specified port that are configured for input.
- **MccBoard.DioConfig.GetDOutMask()** - Determines the bits on a specified port that are configured for output.
- **MccBoard.DioConfig.GetNumBits()** - Gets the number of bits in the digital port value.
- **MccBoard.ExpansionConfig** property - Gets an instance of a `cExpansionConfig` object.

- **MccBoard.ExpansionConfig.GetBoardType()** - Gets the expansion board type.
- **MccBoard.ExpansionConfig.GetCjcChan()** - Gets the channel that the CJC is connected to.
- **MccBoard.ExpansionConfig.GetMuxAdChan1()** - Gets the first A/D channel that the board is connected to.
- **MccBoard.ExpansionConfig.GetMuxAdChan2()** - Gets the second A/D channel that the board is connected to.
- **MccBoard.ExpansionConfig.GetNumExpChans()** - Gets the number of expansion board channels.
- **MccBoard.ExpansionConfig.GetRange1()** - Gets the range/gain of the low 16 channels.
- **MccBoard.ExpansionConfig.GetRange2()** - Gets the range/gain of the high 16 channels.
- **MccBoard.ExpansionConfig.GetThermType()** - Gets the type of thermocouple configuration for the board (J, K, E, T, R, S, and B types).
- **MccBoard.ExpansionConfig.SetCjcChan()** - Sets the channel that the CJC is connected to.
- **MccBoard.ExpansionConfig.SetMuxAdChan1()** - Sets the first A/D channel the board is connected to.
- **MccBoard.ExpansionConfig.SetMuxAdChan2()** - Sets the second A/D channel that the board is connected to.
- **MccBoard.ExpansionConfig.SetRange1()** - Sets the range/gain of the low 16 channels.
- **MccBoard.ExpansionConfig.SetRange2()** - Sets the range/gain of the high 16 channels.
- **MccBoard.ExpansionConfig.SetThermType()** - Sets the type of thermocouple configuration for the board (J, K, E, T, R, S, and B types).
- **GlobalConfig.NumBoards** property - Returns the maximum number of boards you can install at one time.
- **GlobalConfig.NumExpBoards** property - Returns the maximum number of expansion boards you can install on a board.
- **GlobalConfig.Version** property - Information used by the library to determine compatibility.

## Counter methods

The counter methods available from the [MccBoard class](#) are explained below. These methods load, read, and configure counters. There are five types of counter chips used in MCC counter boards: 8254s, 8536s, 7266s, 9513s, and generic event counters. Some of the counter commands only apply to one type of counter.

- **MccBoard.C7266Config()** - Selects the basic operating mode of an LS7266 counter.
- **MccBoard.C8254Config()** - Selects the basic operating mode of an 8254 counter.
- **MccBoard.C8536Config()** - Selects the basic operating mode of an 8536 counter chip.
- **MccBoard.C8536Init()** - Initializes and selects all of the chip level features for a 8536 counter board. The options that are set by this command are associated with each counter chip, not the individual counters within it.
- **MccBoard.C9513Config()** - Sets the basic operating mode of a 9513 counter. This method sets all of the programmable options that are associated with a 9513 counter. It is similar in purpose to C8254Config() except that it is used with a 9513 counter.
- **MccBoard.C9513Init()** - Initializes and selects all of the chip level features for a 9513 counter board. The options that are set by this command are associated with each counter chip, not the individual counters within it.
- **MccBoard.CClear()** - Clears a scan counter value (sets it to zero).
- **MccBoard.CConfigScan()** - Configures a scan counter channel. This method only works with counter boards that have counter scan capability.
- **MccBoard.CFreqIn()** - Measures the frequency of a signal by counting it for a specified period of time (GatingInterval), and then converting the count to count/sec (Hz). This method only works with 9513 counters.

- **MccBoard.CIn()** - Reads a counter's current value as a 16-bit integer. (CIn32 () is the preferred counter read method.)
- **MccBoard.CIn32()** - Reads a counter's current value as a 32-bit integer.
- **MccBoard.CInScan()** - Scans a range of scan counter channels, and stores the samples in an array.
- **MccBoard.CLoad()** - Loads a counter with an initial count value as a 16-bit integer. (CLoad32 () is the preferred counter loading method.)
- **MccBoard.CLoad32()** - Loads a counter with a 32-bit integer initial value.
- **MccBoard.CStatus()** - Read the counter status of a counter. Returns various bits that indicate the current state of a counter (currently only applies to LS7266 counters).
- **MccBoard.CStoreOnInt()** - Installs an interrupt handler that stores the current count whenever an interrupt occurs. This method only works with 9513 counters.
- **MccBoard.TimerOutStart()** - Starts a timer square wave output. This method only works with counter boards that have a timer-type counter.
- **MccBoard.TimerOutStop()** - Stops a timer square wave output. This method only works with counter boards that have a timer-type counter.

## Data Logger methods and properties

The methods and property available from the [DataLogger class](#) are explained below. These class members read and convert binary log files.

- **DataLogger.ConvertFile()** - Converts a binary log file to a comma-separated values (.CSV) text file or another text file format that you specify.
- **DataLogger.GetAIChannelCount()** - Retrieves the total number of analog channels that were logged in a binary file.
- **DataLogger.GetAllInfo()** - Retrieves the channel number and unit value of each analog input channel logged in a binary file.
- **DataLogger.GetCJCInfo()** - Retrieves the number of CJC temperature channels logged in a binary file.
- **DataLogger.GetDIOInfo()** - Retrieves the number of digital I/O channels logged in a binary file.
- **DataLogger.GetFileInfo()** - Retrieves the version level and byte size of a binary file.
- **DataLogger.GetFileName()** - Retrieves the name of the n<sup>th</sup> file in the directory containing binary log files.
- **DataLogger.GetPreferences()** - Retrieves API preference settings for time stamp data, analog temperature data, and CJC temperature data. Returns the default values unless changed using `SetPreferences ()`.
- **DataLogger.GetSampleInfo()** - Retrieves the sample interval, sample count, and the date and time of the first data point in a binary file.
- **DataLogger.ReadAIChannels()** - Retrieves analog input data from a binary file, and stores the values in an array.
- **DataLogger.ReadCJCChannels()** - Retrieves CJC temperature data from a binary file, and stores the values in an array.
- **DataLogger.ReadDIOChannels()** - Retrieves digital I/O channel data from a binary file, and stores the values in an array.
- **DataLogger.ReadTimeTags()** - Retrieves date and time values logged in a binary file. This method stores date values in the `dateTags` array, and time values in the `timeTags` array.
- **DataLogger.SetPreferences()** - Sets preferences for returned time stamp data, analog temperature data, and CJC temperature data.
- **DataLogger.FileName** property - Returns the file name associated with an instance of the `DataLogger` class.

## Digital I/O methods

The digital methods available from the [MccBoard class](#) are explained below. These methods perform digital input and output on various types of digital I/O ports.

- **MccBoard.DBitIn()** - Reads a single bit from a digital input port.
- **MccBoard.DBitOut()** - Sets a single bit on a digital output port.
- **MccBoard.DConfigBit()** - Configures a specific digital bit as input or output.
- **MccBoard.DConfigPort()** - Selects whether a digital port is an input or an output.
- **MccBoard.DIn()** - Reads a specified digital input port.
- **MccBoard.DInScan()** - Reads a set number of bytes or words from a digital input port at a specific rate.
- **MccBoard.DOut()** - Writes a byte to a digital output port.
- **MccBoard.DOutScan()** - Writes a series of bytes or words to a digital output port at a specified rate.

## Error Handling method and properties

Most UL for .NET methods return [ErrorInfo](#) objects. The [MccService](#) class includes one method that determines how errors are handled internally by the library. The [ErrorInfo](#) class includes two properties that provide information returned by the method called.

- **MccService.ErrHandling()** - Sets the manner of reporting and handling errors for all method calls.
- **ErrorInfo.Message** property - Gets the text of the error message associated with a specific error code.
- **ErrorInfo.Value** property - Gets the error constant associated with an [ErrorInfo](#) object.

## Memory board methods

The memory board methods available from the [MccBoard class](#) read and write data to and from a memory board, and also set modes that control memory boards (MEGA-FIFO).

The most common use for memory boards is to store large amounts of data from an A/D board via a DT-Connect cable between the two boards. To do this, use the `ExtMemory` option with the [MccBoard.AInScan\(\)](#) or [APretrig\(\)](#) methods.

Once the data has been transferred to the memory board you can use the memory methods to retrieve the data.

- **MccBoard.MemSetDTMode()** - Set DT-Connect mode on a memory board. Memory boards have a DT-Connect interface which can be used to transfer data through a cable between two boards rather than through the PC's system memory. The DT-Connect port on the memory board can be configured as either an input (from an A/D) or as an output (to a D/A). This method configures the port.
- **MccBoard.MemReset()** - Resets the memory board address. The memory board is organized as a sequential device. When data is transferred to the memory board it is automatically put in the next address location. This method resets the current address to the location 0.
- **MccBoard.MemRead()** - Reads a specified number of points from a memory board starting at a specified address.
- **MccBoard.MemWrite()** - Writes a specified number of points to a memory board starting at a specified address.
- **MccBoard.MemReadPretrig()** - Reads data collected with [MccBoard.APretrig\(\)](#). The [MccBoard.APretrig\(\)](#) method writes the pre-triggered data to the memory board in a scrambled order. This method unscrambles the data and returns it in the correct order.

## Revision control methods and properties

The revision control methods and property explained below are available from the [MccBoard class](#). As new revisions of the library are released, bugs from previous revisions are fixed, and occasionally new properties and methods are added. It is Measurement Computing's goal to preserve the programs you have written so that you never change the order or number of arguments in a method. However, sometimes it is not possible to achieve this goal.

The revision control methods initialize the DLL so that the functions are interpreted according to the format of the revision you wrote and compiled your program in.

- **MccBoard.DeclareRevision()** - Declares the revision number of the Universal Library for .NET that your program was written with.
- **MccBoard.GetRevision()** - Returns the version number of the installed Universal Library for .NET.

## Streamer file methods

The streamer file methods available from the [MccBoard class](#) create, fill, and read streamer files.

- **MccBoard.FileAInScan()** - Transfer analog input data directly to file. Very similar to AInScan() except that the data is stored in a file instead of an array.
- **MccBoard.FilePretrig()** - Pre-triggered analog input to a file. Very similar to APretrig() except that the data is stored in a file instead of an array.
- **MccBoard.FileGetInfo()** - Reads streamer file information on how much data is in the file, and the conditions under which it was collected (sampling rate, channels, etc.).
- **MccBoard.FileRead()** - Reads a selected number of data points from a streamer file into a one-dimensional or two-dimensional array.

## Synchronous I/O methods

The synchronous I/O methods available from the [MccBoard class](#) synchronously read, set, or write data from analog channels, counter channels, thermocouple channels, and digital ports.

- **MccBoard.DaqInScan()** – Scans analog, digital, temperature, and counter inputs synchronously, and stores the values in an array.
- **MccBoard.DaqOutScan()** – Outputs values synchronously to analog output channels and digital output ports.
- **MccBoard.DaqSetSetpoints()** – Configures up to 16 detection setpoints associated with the input channels within a scan group.
- **MccBoard.DaqSetTrigger()** – Selects a trigger source and sets up its parameters. This method starts or stops a synchronous data acquisition operation using [MccBoard.DaqInScan\(\)](#) with the ExtTrigger option.

## Temperature input methods

The temperature input methods available from the [MccBoard class](#) convert a raw analog input from an EXP or other temperature sensor board to temperature.

- **MccBoard.TIn()** - Reads a channel from a digital input board, filters it (if specified), does the cold junction compensation, linearizes and converts it to temperature.
- **MccBoard.TInScan()** - Scans a range of temperature inputs. Reads temperatures from a range of channels and returns the temperature values to an array.

## Windows memory management methods

The Windows memory management methods available from the [MccService class](#) take care of allocating, freeing, and copying to/from Windows global memory buffers.

- **MccService.WinBufAlloc()** - Allocates a Windows memory buffer.
- **MccService.WinBufAlloc32()** - Allocates a Windows global memory buffer for use with 32-bit scan methods, and returns a memory handle for the buffer.
- **MccService.WinBufFree()** - Frees a Windows buffer.
- **MccService.WinArrayToBuf()** - Copies data from a one-dimensional or two-dimensional array into a Windows buffer.
- **MccService.WinBufToArray()** - Copies data from a Windows memory buffer into a one-dimensional or two-dimensional array.
- **MccService.WinBufToArray32()** - Copies 32-bit data from a Windows global memory buffer into an array. This method is typically used to retrieve data from the buffer after executing an input scan method.

## Miscellaneous methods, properties, and delegates

The methods explained below are available from the [MccBoard class](#). These methods do not as a group fit into a single category. They get and set board information, convert units, manage events and background operations, copy two-dimensional arrays to/from Windows global memory buffers, and perform serial communication operations.

- **MccBoard.DeviceLogin()** - Opens a device session with a shared device.
- **MccBoard.DeviceLogout()** - Releases the device session with a shared device.
- **MccBoard.DisableEvent()** - Disables one or more events set up with EnableEvent() and disconnects their user-defined handlers.
- **MccBoard.EnableEvent()** - Binds one or more event conditions to a user-defined callback function.
- **EventCallback** delegate - Defines the prototype for the user function for EnableEvent(). This defines the format for the user-defined handlers to be called when the events set up using EnableEvent() occurs.
- **MccBoard.EngArrayToWinBuf()** - Transfers a 2D array of engineering unit values to a Windows buffer as integer values.
- **MccBoard.FlashLED()** - Causes the LED on a USB device to flash.
- **MccBoard.FromEngUnits()** - Converts a single precision voltage (or current) value in engineering units to an integer D/A count value for output to a D/A.
- **MccBoard.GetBoardName()** - Returns the name of a specified board.
- **MccBoard.GetStatus()** - Returns the status of a background operation. Once a background operation starts, your program must periodically check on its progress. This method returns the current status of the operation.
- **MccBoard.GetTCValues()** - Converts raw thermocouple data gathered with DaqInScan() to Celsius, Fahrenheit, or Kelvin.
- **MccBoard.HideLoginDialog()** - Prevents the default login dialog from being shown when a protected function is called while not logged in.
- **MccBoard.InByte()** - Reads a byte from a hardware register on a board.
- **MccBoard.InWord()** - Reads a word from a hardware register on a board.
- **MccBoard.OutByte()** - Writes a byte to a hardware register on a board.
- **MccBoard.OutWord()** - Writes a byte or word to a hardware register on a board.
- **MccBoard.RS485()** - Sets the transmit and receive buffers on an RS485 port.

- **MccBoard.StopBackground()** - Stop a background process. It is sometimes necessary to stop a background process even though the process has been set up to run continuously. This method stops a background process that is running. [StopBackground\(\)](#) should be executed after normal termination of all background functions in order to clear variables and flags.
- **MccBoard.ToEngUnits()** - Converts an integer A/D count value to an equivalent single precision voltage (or current) value.
- **MccBoard.WinBufToEngArray()** - Transfers integer values from a Windows buffer to a 2D array as engineering unit values.
- **MccBoard.BoardName** property - Name of the board associated with an instance of the [MccBoard class](#).

## Universal Library for .NET example programs

The Universal Library for .NET contains many example programs to help you learn and apply UL for .NET methods. We strongly recommend running appropriate example programs before attempting to use the methods.

Table 5 lists the UL for .NET example programs sorted by program name. It includes their featured method calls, special aspects, and other method calls included in the program. All example programs include the `DeclareRevision()` and `ErrHandling()` methods. Table 6 on page 203 lists the UL for .NET example programs sorted by the method name.

Table 5. UL for .NET Example Programs – sorted by program name

Program name	Featured UL for .NET method call	Notes	Other UL for .NET method calls
CInScan01	cbCInScan ()	Scans a range of counter input channels, and writes the data to an array. Board 0 must support counter scans.	WinBuffAlloc32 () WinBufToArray32 () WinBufFree ()
CInScan02	CInScan () CConfigScan ()	Scans a counter input channel in decrement mode, and writes the data to an array. Board 0 must support counter scans.	WinBuffAlloc32 () WinBufToArray32 () WinBufFree ()
DaqInScan01	DaqInScan ()	Synchronously scans analog input channels, digital input ports and counter input channels in the foreground. Board 0 must support synchronous input.	DConfigPort () CConfigScan ()
DaqInScan02	DaqInScan ()	Synchronously scans analog input channels, digital input ports, and counter input channels in the background. Board 0 must support synchronous input.	DConfigPort () CConfigScan () GetStatus () StopBackground ()
DaqInScan03	DaqInScan () GetTCValues ()	Synchronously scans analog input channels, digital input ports and thermocouple input channels in the foreground. Board 0 must support synchronous input.	DConfigPort ()
DaqOutScan01	DaqOutScan ()	Synchronously writes to an analog output channel and a digital output port in the background. Board 0 must support synchronous output.	DConfigPort ()
DaqSetSetpoints 01	DaqSetSetpoints ()	Configures setpoints, adds the setpoint status to the scanlist, and performs asynchronous reads of the setpoint status. Board 0 must support DaqInScan ().	DaqInScan () DConfigPort () GetStatus () StopBackground ()



Program name	Featured UL for .NET method call	Notes	Other UL for .NET method calls
DaqSetTrigger01	DaqSetTrigger()	Configures start and stop triggers. These triggers are used to initiate and terminate A/D conversion using DaqInScan() with the ExtTrigger option selected. Board 0 must support synchronous output.	DConfigPort() GetStatus() StopBackground()
TimerOut01	TimerOutStart() TimerOutStop()	Sends a frequency output to an output timer channel. Board 0 must have a timer output.	
ULAI01	AIn()		ToEngUnits()
ULAI02	AInScan()	Default mode	WinBufToArray() WinBufFree() WinBufAlloc()
ULAI03	AInScan()	Background mode	GetStatus() StopBackground() WinBufToArray() WinBufFree() WinBufAlloc()
ULAI04	AConvertData()		AInScan() GetStatus() StopBackground() WinBufToArray() WinBufFree() WinBufAlloc()
ULAI05	AInScan()	With manual data conversion	GetStatus() StopBackground() WinBufToArray() WinBufFree() WinBufAlloc()
ULAI06	AInScan()	Continuous Background mode	AConvertData() GetStatus() StopBackground() WinBufToArray() WinBufFree() WinBufAlloc()
ULAI07	ATrig()		FromEngUnits()
ULAI08	APretrig()		WinBufToArray() WinBufFree() WinBufAlloc()
ULAI09	ConvertPretrigData()	Background	APretrig() GetStatus() StopBackground() WinBufToArray() WinBufFree() WinBufAlloc()
ULAI10	cbALoadQueue()		AInScan() WinBufToArray() WinBufFree() WinBufAlloc()
ULAI11	cbToEngUnits()		AIn()



Program name	Featured UL for .NET method call	Notes	Other UL for .NET method calls
ULAI12	cbAInScan ()	ExtClock mode	WinBufToArray () WinBufFree () WinBufAlloc ()
ULAI13	cbAInScan ()	Various sampling mode options	WinBufToArray () WinBufFree () WinBufAlloc ()
ULAI14	SetTrigger ()	With ExtTrigger selected	AInScan () FromEngUnits () WinBufToArray () WinBufFree () WinBufAlloc ()
ULAI001	AInScan () AOutScan ()	Concurrent analog input and analog output scans	GetStatus () StopBackground () WinArrayToBuf () WinBufAlloc () WinBufFree () WinBufToArray ()
ULAO01	AOut ()		FromEngUnits () AOut ()
ULAO02	AOutScan ()		WinBufToArray () WinBufFree () WinBufAlloc ()
ULAO03	AOut () DACUpdate () SetDACUpdateMode ( )	Demonstrates the difference between BoardConfig.DACUpdate.Immediate and BoardConfig.DACUpdate.OnCommand D/A update modes. Board 0 must support DAC update mode settings, such as the PCI-DAC6700 Series boards.	FromEngUnits ()
ULCT01	C8254Config ()		CLoad (), CIn ()
ULCT02	C9513Init () C9513Config ()		CLoad (), CIn ()
ULCT03	CStoreOnInt ()		C9513Init (), CLoad () C9513Config (), CIn ()
ULCT04	CFreqIn ()		C9513Init ()
ULCT05	C8536Init () C8536Config ()		CLoad () CIn ()
ULCT06	C7266Config ()		CLoad32 (), CIn32 () CStatus ()
ULDI01	DIn ()		DConfigPort ()
ULDI02	DBitIn ()		DConfigPort ()
ULDI03	DInScan ()		DConfigPort () GetStatus () StopBackground () WinBufToArray () WinBufFree () WinBufAlloc ()
ULDI04	DIn ()	Using the AuxPort	DioConfig() DConfigPort ()

Program name	Featured UL for .NET method call	Notes	Other UL for .NET method calls
ULDI05	DBitIn ()	Using the AuxPort	DioConfig() DConfigPort ()
ULDI06	DConfigBit ()		DBitIn () DioConfig () DConfigPort ()
ULDO01	DOut ()		DConfigPort ()
ULDO02	DBitOut ()		DOut () , DConfigPort ()
ULDO04	DOut ()	Using the AuxPort	DioConfig () DConfigPort ()
ULDO05	DBitOut ()	Using the AuxPort	DOut () DioConfig () DConfigPort ()
ULEV01	EnableEvent () DisableEvent ()	Using OnExternalInterrupt	DConfigPort () DIn ()
ULEV02	EnableEvent () DisableEvent ()	Using OnDataAvailable and OnEndOfAiScan	AInScan () StopBackground () ToEngUnits () WinBufAlloc () WinBufFree () WinBufToArray ()
ULEV03	EnableEvent () DisableEvent ()	Using OnPretrig and OnEndOfAiScan	APretrig () AConvertPretrigData () DConfigPort () DOut () StopBackground () ToEngUnits () WinBufAlloc () WinBufFree () WinBufToArray ()
ULEV04	EnableEvent () DisableEvent ()	Using OnEndOfAoScan	AOutScan () DConfigPort () DOut () FromEngUnits () StopBackground () WinArrayToBuf () WinBufAlloc () WinBufFree ()
ULFI01	FileAInScan ()		FileGetInfo ()
ULFI02	FileRead ()		FileAInScan () FileGetInfo ()
ULFI03	FilePretrig ()		FileGetInfo () FileRead ()
ULGT01	GetErrMsg ()		AIn ()
ULGT03	MccDaq () .MccBoard () class() properties: BoardConfig, DioConfig and ExpansionConfig	Use the MccBoard class properties to get configuration information for a board.	GetBoardName ()

Program name	Featured UL for .NET method call	Notes	Other UL for .NET method calls
ULGT04	GetBoardName ()		MccDaq.MccBoard.BoardName property MccDaq.GlobalConfig.NumBoards property
ULLOG01	GetFileName ()	Retrieves the name of a binary log file.	
ULLOG02	GetFileInfo () GetSampleInfo () GetAIChannelCount () GetCJCInfo () GetDIOInfo ()	Retrieves information about the analog data, CJC data, and digital I/O channel data contained in a binary log file.	GetFileName ()
ULLOG03	ReadAIChannels () ReadCJCChannels () ReadDIOChannels () ReadTimeTags ()	Retrieves the analog input data, CJC temperature data, digital I/O channel data, date values, and time values logged in a binary file, and writes the data to separate arrays.	GetFileName () GetSampleInfo () GetAIChannelCount () GetCJCInfo () GetDIOInfo () SetPreferences ()
ULLOG04	ConvertFile ()	Converts a binary log file to a comma-separated values (.csv) text file or another text file format that you specify.	GetSampleInfo ()
ULMM01	MemReadPretrig ()		APretrig ()
ULMM02	MemRead () MemWrite ()		
ULMM03	AINScan ()	With ExtMemory option	MemReset () MemRead ()
ULTI01	TIn ()		
ULTI02	TInScan ()		
VIn01	VIn ()	Reads an A/D input channel.	
VOut01	VOut ()	Writes to a D/A output channel.	

Table 6. UL for .NET Example Programs – sorted by method name

UL for .NET method call	UL for .NET example program Name	UL for .NET special features/notes
AConvertData ()	ULAI04 ULAI06	
AConvertPretrigData ()	ULAI09 ULEV03	
ACalibrateData ()	None	No example programs at this time
AIN ()	ULAI01 ULGT01 ULAI11	
AINScan ()	ULAI02 ULAI10 ULAI03 ULAI12 ULAI04 ULAI13 ULAI05 ULAI14 ULAI06 ULMM03 ULEV02	Default, Background mode with manual data conversion Continuous Background mode ExtClock mode Various sampling mode options
ALoadQueue ()	ULAI10	
AOut ()	ULAO01 ULAO03	Demonstrates the difference between BoardConfig.DACUpdate.Immediate and BoardConfig.DACUpdate.OnCommand D/A update modes. Board 0 must support DAC update mode settings, such as the PCI-DAC6700 Series boards.

UL for .NET method call	UL for .NET example program Name	UL for .NET special features/notes
AOutScan ()	ULAO02 ULAIO01 ULEV04	Concurrent AInScan () and AOutScan ()
APretrig ()	ULAI08 ULFI03 ULAI09 ULMM01 ULEV03	
ATrig ()	ULAI07 ULMM01	
C7266Config ()	ULCT06	
C8254Config ()	ULCT01	
C8536Config ()	ULCT05	
C8536Init ()	ULCT05	
C9513Config ()	ULCT02 ULCT03	
C9513Init ()	ULCT02 ULCT04 ULCT03	
CClear ()		
CConfigScan ()	CInScan02	Demonstrates how to scan a counter input channel in decrement mode, and then write the data to an array. Board 0 must support counter scans.
CFreqIn ()	ULCT04	
CIn ()	ULCT01 ULCT05 ULCT02	
CIn32 ()	ULCT06	
CInScan ()	CInScan01 CInScan02	Demonstrates how to scan a range of counter channels and then write the data to an array. Board 0 must support counter scans.
CLoad ()	ULCT01 ULCT03 ULCT02 ULCT05	
CLoad32 ()	ULCT06	
ConvertFile ()	ULLOG04	Demonstrates how to convert a binary log file to a different format. You enter the extension of the file type that you want to create.
CStoreOnInt ()	ULCT03	
CStatus ()	ULCT06	
DaqInScan ()	DaqInScan01 DaqInScan02 DaqInScan03	Demonstrates how to synchronously scan analog, counter, and thermocouple input channels, and digital input ports. Board 0 must support synchronous input.
DaqOutScan ()	DaqOutScan01	Demonstrates how to synchronously write to an analog output channel and digital output port in the background. Board 0 must support synchronous output.
DaqSetSetpoints ()	DaqSetSetpoints01	Demonstrates how to configure and use setpoints, including how to add the setpoint status to the scanlist and perform asynchronous reads of the setpoint status. Board 0 must support DaqInScan ().
DaqSetTrigger ()	DaqSettrigger01	Demonstrates how to set up start and stop trigger events and display input channel data.
DBitIn ()	ULDI02 ULDI06 ULDI05	
DBitOut ()	ULDO02 ULDO05	
DConfigBit ()	ULDI06	

UL for .NET method call	UL for .NET example program Name	UL for .NET special features/notes
DConfigPort ()	ULDI01 ULDO01 ULDI02 ULDO02 ULDI03 ULDO05 ULEV01 ULEV03 ULEV04	
DIn ()	ULDI01 ULDI04 ULDI03 ULEV01	
DInScan ()	ULDI03	
DOut ()	ULDO01 ULDO05 ULDO02 ULDO04 ULEV03 ULEV04	
DOutScan ()	None	No example programs at this time
EnableEvents () DisableEvents ()	ULEV01 ULEV03 ULEV02 ULEV04	OnExternalInterrupt OnDataAvailable OnPretrigger OnEndOfAoScan OnScanError OnEndOfAiScan
ErrHandling ()	All samples	All example programs use this method
FileAInScan ()	ULFI01 ULFI02	
FilePretrig ()	ULFI03	ULFI01 ULFI02
FileRead ()	ULFI02 ULFI03	
FlashLED ()	ULFI01	Flashes the onboard LED for visual identification (board 0 must have an external LED, such as the miniLAB 1008 or the USB-1208LS).
FromEngUnits	ULAO01 ULAO03 ULAI07 ULAI14 ULEV04	
GetAIChannelCount ()	ULLOG02 ULLOG03 ULLOG04	Demonstrates how to retrieve the number of analog channels contained in a binary log file.
GetAIInfo ()	ULLOG02	Demonstrates how to retrieve information about the analog input data contained in a binary log file.
GetBoardName	ULGT03 ULGT04	
GetCJCInfo ()	ULLOG02 ULLOG03	Demonstrates how to retrieve information about CJC temperature data contained in a binary log file.
GetDACStartup ()	None	No example programs at this time
GetDACUpdateMode ()	None	No example programs at this time
GetDIOInfo ()	ULLOG02 ULLOG03	Demonstrates how to retrieve information about digital I/O channel data contained in a binary log file.
GetErrMsg ()	ULGT01	
GetFileInfo ()	ULLOG02	Demonstrates how to retrieve the version level and byte size of a binary log file.
GetFileName ()	ULLOG01 ULLOG02 ULLOG03	Demonstrates how to retrieve the name of a binary log file.
GetPreferences ()	None	No example programs at this time.
GetRevision ()	None	No example programs at this time

UL for .NET method call	UL for .NET example program Name	UL for .NET special features/notes
GetSampleInfo ()	ULLOG02 ULLOG03	Demonstrates how to retrieve the sample interval, sample count, and the date and time of the first data point logged in a binary file.
GetStatus ()	ULAI03   ULAI06 ULAI04   ULAI09 ULAI05   ULCT03 ULAIO01   ULDI0	
GetTCValues ()	DaqInScan03	Demonstrates how to retrieve analog, thermocouple, and digital data from a synchronous scan operation. Board 0 must support synchronous input.
InByte ()	None	No example programs at this time
InWord ()	None	No example programs at this time
MccDaq.MccBoard class properties: BoardConfig, DioConfig, and ExpansionConfig	ULGT03 ULGT04	Use the MccBoard class properties to get configuration information for a board.
MemRead ()	ULMM01 ULMM03 ULMM02	
MemReadPretrig ()	ULMM01	
MemReset ()	ULMM03	
MemSetDTMode ()	None	No example programs at this time
MemWrite ()	ULMM02	
ReadAIChannels () ReadCJCChannels () ReadDIOChannels () ReadTimeTags ()	ULLOG03	Demonstrates how to retrieve analog input data, CJC data, DIO port data, and date/time values contained in a binary log file, and then store the data in separate arrays.
RS485 ()	None	No example programs at this time
SetTrigger ()	ULAI14	
SetPreferences ()	ULLOG03	Demonstrates how to store preference settings for time stamped data, analog data, and CJC temperature data.
StopBackground ()	ULAI03   ULAI06 ULAI04   ULAI09 ULAI05   ULCT03 ULAIO01   ULDI03 ULEV02   ULEV03 ULEV04	Concurrent AInScan () and AOutScan ()
TimerOutStart () TimerOutStop ()	TimerOut01	Demonstrates how to start and stop a timer square wave output. The board must have a timer output.
TIn ()	ULTI01	
TInScan ()	ULTI02	
ToEngUnits ()	ULAI01   ULAI11 ULAI07   ULEV02 ULEV03	
VIn ()	VIn01	Demonstrates how to read an A/D input channel.
Vout ()	VOut01	Demonstrates how to write to a D/A output channel.
WinArrayToBuf ()	ULAIO01   ULAIO02 ULEV04	
WinBufAlloc32 ()	CInScan01   CInScan02	Demonstrates how to allocate a Windows global memory buffer for use with 32-bit scan methods.
WinBufToArray32 ()	CInScan01   CInScan02	Demonstrates how to copy 32-bit data from a Windows global memory buffer into an array.

UL for .NET method call	UL for .NET example program Name	UL for .NET special features/notes
WinBufAlloc() WinBufFree() WinBufToArray()	ULAI01 – ULAI06 ULAI08 – ULAI10 ULAI12 – ULAI14 ULAO02 ULCT03 ULDI03 ULEV02 – ULEV04 (ULEV04:WinBufAlloc and WinBufFree only)	

## Analog I/O Methods

### Introduction

This chapter covers Universal Library for .NET methods that handle analog input, analog output and analog data manipulation. These methods are available from the [MccBoard class](#).

Most analog I/O methods include options that may not be compatible with your hardware. To determine which of these methods are compatible with your hardware, refer to the *Universal Library User's Guide* (available in PDF format on our website at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf)).

Table 7 lists the MccDaq.Range enumerated constants you can use in the `range` parameter found in most of the methods explained in this chapter. These values are also used in the [ALoadQueue\(\)](#) method's `gainArray` parameter. Valid ranges for your hardware are listed in the *Universal Library User's Guide*.

Table 7. MccDaq.Range enumerated constants

UL .NET settings	Value	UL .NET settings	Value
MccDaq.Bip20Volts	±20 volts(V)	MccDaq.Uni10Volts	0 to 10 V
MccDaq.Bip10Volts	±10 V	MccDaq.Uni5Volts	0 to 5 V
MccDaq.Bip5Volts	±5 V	MccDaq.Uni2Pt5Volts	0 to 2.5 V
MccDaq.Bip4Volts	±4 V	MccDaq.Uni2Volts	0 to 2 V
MccDaq.Bip2Pt5Volts	±2.5 V	MccDaq.Uni1Pt25Volts	0 to 1.25 V
MccDaq.Bip2Volts	±2 V	MccDaq.Uni1Pt67Volts	0 to 1.67 V
MccDaq.Bip1Pt25Volts	±1.25 V	MccDaq.Uni1Volts	0 to 1 V
MccDaq.Bip1Volts	±1 V	MccDaq.UniPt5Volts	0 to 0.5 V
MccDaq.Bip1Pt67Volts	±1.67 V	MccDaq.UniPt25Volts	0 to 0.25 V
MccDaq.BipPt625Volts	±0.625 V	MccDaq.UniPt2Volts	0 to 0.2 V
MccDaq.BipPt5Volts	±0.5 V	MccDaq.UniPt1Volts	0 to 0.1 V
MccDaq.BipPt25Volts	±0.25 V	MccDaq.UniPt01Volts	0 to 0.01 V
MccDaq.BipPt2Volts	±0.2 V	MccDaq.UniPt02Volts	0 to 0.02 V
MccDaq.BipPt1Volts	±0.1 V	MccDaq.UniPt05Volts	0 to 0.05 V
MccDaq.BipPt05Volts	±0.05 V	MccDaq.Ma0To20	0 to 20 milliamperes (mA)
MccDaq.BipPt01Volts	±0.01 V	MccDaq.Ma4To20	4 to 20 mA
MccDaq.BipPt005Volts	±0.005 V	MccDaq.Ma2To10	2 to 10 mA
		MccDaq.Ma1To5	1 to 5 mA
NotUsed	-1	MccDaq.MaPt5To2Pt5	0.5 to 2.5 mA



## AConvertData()

Converts the raw data collected by [AInScan\(\)](#) into 12-bit A/D values. The `AInScan()` method can return either raw A/D data or converted data, depending on whether or not the `ConvertData()` option is used. For many 12-bit A/D boards, the raw data is a 16-bit value that contains a 12-bit A/D value and a 4-bit channel tag (refer to board-specific information in the *Universal Library User's Guide*). The data returned to `adData` consists of just the 12-bit A/D value. The data returned to `chanTags` consists of just the channel numbers.

Member of the [MccBoard class](#).

### Function prototype:

```
VB .NET:      Public Function AConvertData(ByVal numPoints As Integer, ByRef
               adData As Short, ByRef chanTags As Short) As MccDaq.ErrorInfo

               Public Function AConvertData(ByVal numPoints As Integer, ByRef
               adData As System.UInt16, ByRef chanTags As System.UInt16) As
               MccDaq.ErrorInfo

C# .NET:      public MccDaq.ErrorInfo AConvertData(int numPoints, ref ushort
               adData, out ushort chanTags)

               public MccDaq.ErrorInfo AConvertData(int numPoints, ref short
               adData, out short chanTags)
```

### Parameters:

<code>numPoints</code>	Number of samples to convert
<code>adData</code>	Reference to start of data array
<code>chanTags</code>	Reference to start of channel tag array

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

`adData` - converted data

`chanTags` - channel tags if available.

When collecting data using [AInScan\(\)](#) without the `ConvertData` option, use this method to convert the data after it has been collected. There are cases where the `ConvertData` option is not allowed. For example - if you are using both the `DmaIo` and `Background` option with `AInScan()` on some devices, the `ConvertData` option is not allowed. In those cases this function should be used to convert the data after the data collection is complete.

In those cases, use `AConvertData()` to convert the data after the data collection is complete.

For some boards, each raw data point consists of a 12-bit A/D value with a 4-bit channel number. This method pulls each data point apart and puts the A/D value into the `adData` array and the channel number into the `chanTags` array.

### Notes:

**12-bit A/D boards:** The name of the array must match that used in [AInScan\(\)](#) or [WinBufToArray\(\)](#). Upon returning from `AConvertData()`, `adData` array contains only 12-bit A/D data.

## AConvertPretrigData()

For products with pretrigger implemented in hardware (most products), this function converts the raw data collected by [APretrig\(\)](#). The `APretrig()` method can return either raw A/D data or converted data, depending on whether or not the `ConvertData` option was used. The raw data is not in the correct order as it is collected. After the data collection is completed, it must be rearranged into the correct order. This method also orders the data, starting with the first pretrigger data point and ending with the last post-trigger point.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**

```
Public Function AConvertPretrigData(ByVal preTrigCount As Integer,
    ByVal totalCount As Integer, ByRef adData As Short, ByRef chanTags
    As Short) As MccDaq.ErrorInfo

    Public Function AConvertPretrigData(ByVal preTrigCount As Integer,
    ByVal totalCount As Integer, ByRef adData As System.UInt16, ByRef
    chanTags As System.UInt16) As MccDaq.ErrorInfo
```

**C# .NET:**

```
public MccDaq.ErrorInfo AConvertPretrigData(int preTrigCount, int
    totalCount, ref ushort adData, out ushort chanTags)

public MccDaq.ErrorInfo AConvertPretrigData(int preTrigCount, int
    totalCount, ref short adData, out short chanTags)
```

### Parameters:

<code>preTrigCount</code>	Number of pre-trigger samples (this value must match the value returned by the <code>PretrigCount</code> parameter in the <a href="#">APretrig()</a> method)
<code>totalCount</code>	Total number of samples that were collected
<code>adData</code>	Reference to data array (must match array name used in <a href="#">APretrig()</a> method)
<code>chanTags</code>	Reference to channel tag array, or a NULL reference may be passed if using 16-bit boards or if channel tags are not desired (see the note regarding 16-bit boards below).

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

`adData` - converted data

When you collect data with [APretrig\(\)](#) and you don't use the `ConvertData` option, you must use this method to convert the data after it is collected. There are cases where the `ConvertData` option is not allowed: for example, if you use the `Background` option with `APretrig()` on some devices, the `ConvertData` option is not allowed. In those cases, this method should be used to convert the data after the data collection is complete.

### Notes:

**12-Bit A/D Boards:** On some 12-bit boards, each raw data point consists of a 12-bit A/D value with a 4-bit channel number. This method pulls each data point apart and puts the A/D value into the `adData` and the channel number into the `chanTags` array.

Upon returning from `AConvertPretrigData()`, `adData` array contains only 12-bit A/D data.

**16-Bit A/D Boards:** This method is for use with 16-bit A/D boards only insofar as ordering the data. No channel tags are returned.

Name of the `ADData` array must match that used in [AInScan\(\)](#) or [WinBufToArray\(\)](#).

**Visual Basic programmers:**

After the data is collected with [APretrig\(\)](#), it must be copied to a BASIC array with [WinBufToArray\(\)](#).

**Important**

The entire array must be copied, which includes the extra 512 samples needed by `APretrig()`. Example code is provided below.

```
SampleCount& = 10000
Dim A_D_Data%(SampleCount& + 512)
Dim Chan_Tags%(SampleCount& + 512)
APretrig%(LowChan, HighChan, PretrigCount&, SampleCount&...)
WinBufToArray%(MemHandle%, A_D_Data%, SampleCount& + 512)
AConvertPretrigData%(Pretrig_Count&, SampleCount&, A_D_Data%, Chan_Tags%)
```

## ACalibrateData()

Calibrates the raw data collected by [AInScan\(\)](#) from boards with real time software calibration when the real time calibration has been turned off. The `AInScan()` method can return either raw A/D data or calibrated data, depending on whether or not the `NoCalibrateData` option was used.

Member of the [MccBoard class](#).

### Function prototype:

**VB .NET:**

```
Public Function ACalibrateData(ByVal numPoints As Integer, ByVal
range As MccDaq.Range, ByRef adData As Short) As MccDaq.ErrorInfo

Public Function ACalibrateData(ByVal numPoints As Integer, ByVal
range As MccDaq.Range, ByRef adData As System.UInt16) As
MccDaq.ErrorInfo
```

**C# .NET:**

```
public MccDaq.ErrorInfo ACalibrateData(int numPoints, MccDaq.Range,
ref ushort adData)

public MccDaq.ErrorInfo ACalibrateData(int numPoints, MccDaq.Range
range, ref short adData)
```

### Parameters:

<code>numPoints</code>	Number of samples to convert
<code>range</code>	The programmable gain/range used when the data was collected. Refer to Table 7 on page 208 for a list of valid range settings.
<code>adData</code>	Reference to data array

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

`adData` - converted data

### Notes:

When collecting data using [AInScan\(\)](#) with the `NoCalibrateData` option, use this method to calibrate the data after it is collected.

- The name of the array must match that used in [AInScan\(\)](#) or [WinBufToArray\(\)](#).
- Applying software calibration factors in real time on a per sample basis eats up machine cycles. If your CPU is slow, or if processing time is at a premium, withhold calibration until after the acquisition run is complete. Turning off real time software calibration saves CPU time during a high speed acquisition run.

Processor speed is a factor for DMA transfers and for real time software calibration. Processors of less than 150 MHz Pentium class may impose speed limits below the capability of the board (refer to specific board information.) If your processor is less than a 150 MHz Pentium, and you need an acquisition speed in excess of 200 kHz, use the `NoCalibrateData` option to turn off real-time software calibration and save CPU time. After the acquisition is run, calibrate the data with `ACalibrateData()`.

## AIn()

Reads an A/D input channel. This method reads the specified A/D channel from the specified board. If the specified A/D board has programmable gain then it sets the gain to the specified range. The raw A/D value is converted to an A/D value and returned to `DataValue`.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**

```
Public Function AIn(ByVal channel As Integer, ByVal range As
MccDaq.Range , ByRef dataValue As Short) As MccDaq.ErrorInfo

Public Function AIn(ByVal channel As Integer, ByVal range As
MccDaq.Range, ByRef dataValue As System.UInt16) As MccDaq.ErrorInfo
```

**C# .NET:**

```
public MccDaq.ErrorInfo AIn(int channel, MccDaq.Range range, out
ushort DataValue)

public MccDaq.ErrorInfo AIn(int channel, MccDaq.Range range, out
short DataValue)
```

### Parameters:

channel	A/D channel number. The maximum allowable channel depends on which type of A/D board is being used. For boards with both single ended and differential inputs, the maximum allowable channel number also depends on how the board is configured. For example, a CIO-DAS1600 has 8 channels for differential, 16 for single ended. Expansion boards are also supported by this method, so this parameter can contain values up to 272. See board specific information for EXP boards if you are using an expansion board.
range	A/D Range code. If the selected A/D board does not have a programmable gain feature, this parameter is ignored. If the A/D board does have programmable gain, set the <code>range</code> parameter to the desired A/D range. Refer to board specific information for a list of the supported A/D ranges of each board. Refer to Table 7 on page 208 for a list of valid range settings.
dataValue	Reference to data value.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

`dataValue` - The value of the A/D sample.

## AInScan()

Scans a range of A/D channels and stores the samples in an array. `AInScan()` reads the specified number of A/D samples at the specified sampling rate from the specified range of A/D channels from the specified board. If the A/D board has programmable gain, then it sets the gain to the specified range. The collected data is returned to the data array.

Member of the [MccBoard class](#).

### Function prototype:

**VB .NET:**

```
Public Function AInScan(ByVal lowChan As Integer, ByVal highChan As Integer, ByVal numPoints As Integer, ByRef rate As Integer, ByVal range As MccDaq.Range, ByVal memHandle As Integer, ByVal options As MccDaq.ScanOptions) As MccDaq.ErrorInfo
```

**C# .NET:**

```
public MccDaq.ErrorInfo AInScan(int lowChan , int highChan, int numPoints, ref int rate, MccDaq.Range range, int memHandle, MccDaq.ScanOptions options)
```

### Parameters:

<code>lowChan</code>	First A/D channel of the scan. When <code>ALoadQueue()</code> is used, the channel count is determined by the total number of entries in the channel gain queue. <code>lowChan</code> is ignored.
<code>highChan</code>	Last A/D channel of the scan. When <code>ALoadQueue()</code> is used, the channel count is determined by the total number of entries in the channel gain queue. <code>highChan</code> is ignored.  <b>low / high Channel #</b> - The maximum allowable channel depends on which type of A/D board is being used. For boards that have both single ended and differential inputs the maximum allowable channel number also depends on how the board is configured. For example, a CIO-DAS1600 has 8 channels for differential, 16 for single ended.
<code>numPoints</code>	Number of A/D samples to collect. Specifies the total number of A/D samples that will be collected. If more than one channel is being sampled then the number of samples collected per channel is equal to $\text{count} / (\text{highChan} - \text{lowChan} + 1)$ .
<code>rate</code>	The rate at which samples are acquired, in samples per second per channel. For example, sampling four channels, 0-3, at a rate of 10,000 scans per second (10 kHz) results in an A/D converter rate of 40 kHz: four channels at 10,000 samples per channel per second. With other software, you specify the total A/D chip rate. In those systems, the per channel rate is equal to the A/D rate divided by the number of channels in a scan.  The channel count is determined by the <code>lowChan</code> and <code>highChan</code> parameters. $\text{Channel Count} = (\text{highChan} - \text{lowChan} + 1)$ .  When <code>ALoadQueue()</code> is used, the channel count is determined by the total number of entries in the channel gain queue. <code>lowChan</code> and <code>highChan</code> are ignored.  <code>rate</code> also returns the value of the actual rate set, which may be different from the requested rate because of pacer limitations.
<code>range</code>	A/D range code. If the selected A/D board does not have a programmable range feature, this parameter is ignored. Otherwise, set the <code>range</code> parameter to any range that is supported by the selected A/D board. Refer to board-specific information for a list of the supported A/D ranges of each board. Refer to Table 7 on page 208 for a list of valid range settings.

memHandle	Handle for Windows buffer to store data. This buffer must have been previously allocated with the <a href="#">WinBufAlloc()</a> method.
options	Bit fields that control various options. Refer to the constants in the "options parameter values" section below.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

rate - actual sampling rate used.

memHandle - collected A/D data returned via the Windows buffer.

**options parameter values:**

All of the options settings are `MccDaq.ScanOptions` enumerated constants. To set a variable to one of these constants, you must refer to the `MccDaq` object and the `ScanOptions` enumeration (*variable* = `MccDaq.ScanOptions.SingleIo`, *variable* = `MccDaq.ScanOptions.DmaIo`, etc.).

**Transfer method options:** The following four options determine how data is transferred from the board to PC memory. If none of these options are specified (recommended), the optimum sampling mode is automatically chosen based on board type and sampling speed.

SingleIo	A/D transfers to memory are initiated by an interrupt. One interrupt per conversion. Rates attainable using <code>SingleIo</code> are PC-dependent and generally less than 10 kHz. Use the default method unless you have a reason to select a specific transfer method.
DmaIo	A/D transfers are initiated by a DMA request.
BlockIo	A/D transfers are handled in blocks (by REP-INSW for example).  <b>BlockIo is not recommended for slow acquisition rates:</b> If the rate of acquisition is very slow (for example less than 200 Hz) <code>BlockIo</code> is probably not the best choice for transfer mode. The reason for this is that status for the operation is not available until one packet of data has been collected (typically 512 samples). The implication is that if acquiring 100 samples at 100 Hz using <code>BlockIo</code> , the operation will not complete until 5.12 seconds has elapsed.
BurstIo	Allows higher sampling rates for sample counts up to full FIFO. Data is collected into the local FIFO. Data transfers to the PC are held off until after the scan is complete. For <code>Background</code> scans, the count and index returned by <a href="#">GetStatus</a> remain 0 and the status equals <code>Running</code> until the scan finishes. When the scan is complete and the data is retrieved, the count and index are updated and the status equals <code>Idle</code> . <code>BurstIo</code> is the default mode for non-Continuous fast scans (aggregate sample rates above 1000 Hz) with sample counts up to full-FIFO. To avoid the <code>BurstIo</code> default, specify <code>BlockIo</code> . <code>BurstIo</code> is not a valid option for most boards. It is used mainly for USB products.
BurstMode	Enables burst mode sampling. Scans from <code>lowChan</code> to <code>highChan</code> are clocked at the maximum A/D rate between samples in order to minimize channel to channel skew. Scans are initiated at the rate specified by rate.  <code>BurstMode</code> is not recommended for use with the <code>SingleIo</code> option. If this combination is used, the count value should be set as low as possible, preferably to the number of channels in the scan. Otherwise, overruns may occur.

ConvertData	<p>If the <code>ConvertData</code> option is used for 12 bit boards then the data that is returned to the buffer will automatically be converted to 12 bit A/D values. If <code>ConvertData</code> is not used then the data from 12 bit A/D boards will be return unmodified (16-bit values that contain both a 12 bit A/D value and a 4 bit channel number). After the data collection is complete you can call <a href="#">AConvertData()</a> to convert the data after the fact. On some devices, <code>ConvertData</code> may not be specified if you are using the <code>Background</code> option and DMA transfers. This option is ignored for the 16-bit boards.</p>
Background	<p>If the <code>Background</code> option is not used, the <a href="#">AInScan()</a> method will not return to your program until all of the requested data has been collected and returned to the buffer. When the <code>Background</code> option is used, control will return immediately to the next line in your program and the data collection from the A/D into the buffer will continue in the background. Use <a href="#">GetStatus()</a> with <code>AiFunction</code> to check on the status of the background operation. Alternatively, some boards support <a href="#">EnableEvent()</a> for event notification of changes in status of <code>Background</code> scans. Use <a href="#">StopBackground()</a> with <code>AiFunction</code> to stop the background process before it has completed. <code>StopBackground()</code> should be executed after normal termination of all background functions in order to clear variables and flags.</p>
Continuous	<p>This option puts the method in an endless loop. Once it collects the required number of samples, it resets to the start of the buffer and begins again. The only way to stop this operation is with <a href="#">StopBackground()</a>. Normally this option should be used in combination with <code>Background</code> so that your program will regain control.</p> <p><b>numPoints parameter settings in Continuous mode:</b> For some DAQ hardware, <code>numPoints</code> must be an integer multiple of the <i>packet size</i>. Packet size is the amount of data that a DAQ device transmits back to the PC's memory buffer during each data transfer. Packet size can differ among DAQ hardware, and can even differ on the same DAQ product depending on the transfer method.</p> <p>In some cases, the minimum value for the <code>numPoints</code> parameter may change when the <code>Continuous</code> option is used. This can occur for several reasons; the most common is that in order to trigger an interrupt on boards with FIFOs, the circular buffer must occupy at least half the FIFO. Typical half-FIFO sizes are 256, 512 and 1024.</p> <p>Another reason for a minimum <code>numPoints</code> value is that the buffer in memory must be periodically transferred to the user buffer. If the buffer is too small, data will be overwritten during the transfer resulting in garbled data.</p> <p>Refer to the board-specific section in the <i>Universal Library User's Guide</i> (available on our web site at <a href="http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf">www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf</a>) for packet size information for your particular DAQ hardware.</p>
ExtClock	<p>If this option is used then conversions will be controlled by the signal on the external clock input rather than by the internal pacer clock. Each conversion will be triggered on the appropriate edge of the clock input signal (refer to board-specific information contained in the <i>Universal Library Users Guide</i>). In most cases, when this option is used the <code>rate</code> parameter is ignored. The sampling rate is dependent on the clock signal. Options for the board will default to a transfer mode that will allow the maximum conversion rate to be attained unless otherwise specified.</p> <p>In some cases, such as with the PCI-DAS4020/12, an approximation of the rate is used to determine the size of the packets to transfer from the board. Set the <code>rate</code> parameter to an approximate maximum value.</p>



	<p><b>SingleIo is recommended for slow external clock rates:</b> If the rate of the external clock is very slow (say less than 200 Hz) and the board you are using supports <code>BlockIo</code>, you may want to include the <code>SingleIo</code> option. This is because that the status for the operation is not available until one packet of data has been collected (typically 512 samples). The implication is that, if acquiring 100 samples at 100 Hz using <code>BlockIo</code> (the default for boards that support it if <code>ExtClock</code> is used), the operation will not complete until 5.12 seconds has elapsed.</p>
<code>ExtMemory</code>	<p>Causes the command to send the data to a connected memory board via the DT-Connect interface rather than returning the data to the buffer. Data for each call to this method will be appended unless <a href="#">MemReset()</a> is called. The data should be unloaded with the <a href="#">MemRead()</a> method before collecting new data. When <code>ExtMemory</code> option is used, the reference to the buffer(<code>memHandle</code>) may be set to null or 0. <code>Continuous</code> option cannot be used with <code>ExtMemory</code>. Do not use <code>ExtMemory</code> and <code>DtConnect</code> together. The transfer modes <code>DmaIo</code>, <code>SingleIo</code> and <code>BlockIo</code> have no meaning when used with this option.</p>
<code>ExtTrigger</code>	<p>If this option is specified, the sampling will not begin until the trigger condition is met. On many boards, this trigger condition is programmable (refer to <a href="#">SetTrigger()</a> and to board-specific info for details). On other boards, only 'polled gate' triggering is supported. In this case assuming active high operation, data acquisition will commence immediately if the trigger input is high. If the trigger input is low, acquisition will be held off until it goes high. Acquisition will then continue until <code>numPoints</code> samples have been taken regardless of the state of the trigger input. This option is most useful if the signal is a pulse with a very low duty cycle (trigger signal in TTL low state most of the time) so that triggering will be held off until the occurrence of the pulse.</p>
<code>NoTodInts</code>	<p>Disables the system's time-of-day interrupts for the duration of the scan. These interrupts are used to update the systems real time clock and are also used by various other programs.</p> <p>These interrupts can limit the maximum sampling speed of some boards - particularly the PCM-DAS08. If the interrupts are turned off using this option, the real-time clock will fall behind by the length of time that the scan takes.</p>
<code>NoCalibrateData</code>	<p>Turns off real-time software calibration for boards which are software calibrated, by applying calibration factors to the data on a sample by sample basis as it is acquired. Examples are the PCM-DAS16/330 and PCM-DAS16x/12.</p> <p>Turning off software calibration saves CPU time during a high speed acquisition run. This may be required if your processor is less than a 150 MHz Pentium and you desire an acquisition speed in excess of 200 kHz. These numbers may not apply to your system. Only trial will tell for sure. DO NOT use this option if you do not have to. If this option is used, the data must be calibrated after the acquisition run with the <a href="#">ACalibrateData()</a> method.</p>
<code>DtConnect</code>	<p>All A/D values will be sent to the A/D board's DT-Connect port. This option is incorporated into the <code>ExtMemory</code> option. Use <code>DtConnect</code> only if the external board is not supported by the Universal Library.</p>
<code>RetrigMode</code>	<p>Re-arms the trigger after a trigger event is performed. With this mode, the scan begins when a trigger event occurs. When the scan completes, the trigger is re-armed to acquire the next the batch of data. You can specify the number of samples in the scan for each trigger event (described below). The <code>RetrigMode</code> option can be used with the <code>Continuous</code> option to continue arming the trigger until <code>StopBackground()</code> is called.</p>

You can specify the number of samples to acquire with each trigger event. This is the trigger count (`retrigCount`). Use `SetAdRetrigCount()` to set the trigger count. If you specify a trigger count that is either zero or greater than the value of the `AInScan numPoints` argument, the trigger count is set to the value of `numPoints`.

Specify the `Continuous` option with the trigger count set to zero to fill the buffer with `numPoints` samples, re-arm the trigger, and refill the buffer upon the next trigger.

**Caution!** You will generate an error if you specify a total A/D rate beyond the capability of the board. For example, if you specify `LowChan = 0`, `HighChan = 7` (8 channels total), and `Rate = 20,000`, and you are using a CIO-DAS16/JR, you will get an error — you have specified a total rate of  $8 * 20,000 = 160,000$ , but the CIO-DAS16/JR is capable of converting only 120,000 samples per second. The maximum sampling rate depends on the A/D board that is being used. It is also dependent on the sampling mode options.

**Important**

In order to understand the functions, read the board-specific information contained in the *Universal Library User's Guide* (available on our web site at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf)). Review and run the example programs before attempting any programming of your own. Following this advice will save you hours of frustration, and possibly time wasted holding for technical support.

This note, which appears elsewhere, is especially applicable to this method. Now is the time to read board-specific information for your board (see the *Universal Library User's Guide*). We suggest that you make a copy of that page to refer to as you read this manual and examine the example programs.

## ALoadQueue()

Loads the A/D board's channel/gain queue. This method only works with A/D boards that have channel/gain queue hardware.

Some products do not support channel / gain queue, and some that do support it are limited on the order of elements, number of elements, and gain values that can be included, etc. Please refer to the device-specific information in the *Universal Library User's Guide* (available on our web site at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf)) to find details for your particular product.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**                   Public Function ALoadQueue(ByVal chanArray As Short( ), ByVal gainArray As MccDaq.Range(), ByVal count As Integer) As MccDaq.ErrorInfo

**C# .NET:**                   public MccDaq.ErrorInfo ALoadQueue(short[] chanArray, MccDaq.Range[] gainArray, int count)

### Parameters:

chanArray	Array containing channel values. This array should contain all of the channels that will be loaded into the channel gain queue.
gainArray	Array containing A/D range values. This array should contain each of the A/D ranges that will be loaded into the channel gain queue. Refer to Table 7 on page 208 for a list of valid A/D range settings.
count	Number of elements in chanArray and gainArray or 0 to disable channel/gain queue. Specifies the total number of channel/gain pairs that will be loaded into the queue.  chanArray and gainArray should contain at least count elements. Set count = 0 to disable the board's channel/gain queue. The maximum value is specific to the queue size of the A/D boards channel gain queue.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

### Notes:

Normally, the [AInScan\(\)](#) method scans a fixed range of channels (from lowChan to highChan) at a fixed A/D range. If you load the channel gain queue with this method then all subsequent calls to AInScan() will cycle through the channel/range pairs that you have loaded into the queue.

## AOut()

Sets the value of a D/A output.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**

```
Public Function AOut(ByVal channel As Integer, ByVal range As
MccDaq.Range, ByVal dataValue As Short) As MccDaq.ErrorInfo

Public Function AOut(ByVal channel As Integer, ByVal range As
MccDaq.Range, ByVal dataValue As System.UInt16) As MccDaq.ErrorInfo
```

**C# .NET:**

```
public MccDaq.ErrorInfo AOut(int channel, MccDaq.Range range, ushort
dataValue)

public MccDaq.ErrorInfo AOut(int channel, MccDaq.Range range, short
dataValue)
```

### Parameters:

channel	D/A channel number. The maximum allowable channel depends on which type of D/A board is being used.
range	D/A range code. The output range of the D/A channel can be set to any of those supported by the board. If the D/A board does not have programmable ranges then this parameter will be ignored. Refer to Table 7 on page 208 for a list of valid range settings.
dataValue	Value to set D/A to. Must be in the range 0 - N where N is the value $2^{\text{Resolution}} - 1$ of the converter  <b>Exception:</b> using 16-bit boards with Basic range is -32768 to 32767. Refer to the discussion on Basic signed integers for more information.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

### Notes:

**Simultaneous Update Boards:** If you set the simultaneous update jumper for simultaneous operation, use [AOutScan\(\)](#) for simultaneous update of multiple channels. `AOut()` always writes the D/A data then reads the D/A, which causes the D/A output to be updated.

## AOutScan()

Outputs values to a range of D/A channels. This method can be used for paced analog output on hardware that supports paced output. It can also be used to update all analog outputs at the same time when the `Simultaneous` option is used.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**                   Public Function AOutScan(ByVal lowChan As Integer, ByVal highChan As Integer, ByVal numPoints As Integer, ByRef rate As Integer, ByVal range As MccDaq.Range , ByVal memHandle As Integer, ByVal options As MccDaq.ScanOptions) As MccDaq.ErrorInfo

**C# .NET:**                   public MccDaq.ErrorInfo AOutScan(int lowChan, int highChan, int numPoints, ref int rate, MccDaq.Range range, int memHandle, MccDaq.ScanOptions options)

### Parameters:

<code>lowChan</code>	First D/A channel of scan.
<code>highChan</code>	Last D/A channel of scan.  <code>lowChan/highChan</code> - The maximum allowable channel depends on which type of D/A board is being used.
<code>numPoints</code>	Number of D/A values to output. Specifies the total number of D/A values that will be output. Most D/A boards do not support timed outputs. For these boards, set the count to the number of channels in the scan.
<code>rate</code>	Sample rate in scans per second. For many D/A boards the <code>rate</code> is ignored and can be set to <code>NotUsed</code> . For D/A boards with trigger and transfer methods which allow fast output rates, such as the CIO-DAC04/12-HS, <code>rate</code> should be set to the D/A output rate (in scans/sec). This parameter also returns the value of the actual rate set. This value may be different from the user specified rate because of pacer limitations.  If supported, this is the rate at which scans are triggered. If you are updating 4 channels, 0-3, then specifying a rate of 10,000 scans per second (10 kHz) will result in the D/A converter rates of 10 kHz — (one D/A per channel). The data transfer rate will be 40,000 words per second — 4 channels * 10,000 updates per scan.  The maximum update rate depends on the D/A board that is being used, and the sampling mode options.
<code>range</code>	D/A range code. The output range of the D/A channel can be set to any of those supported by the board. If the D/A board does not have a programmable gain this parameter is ignored. Refer to Table 7 on page 208 for a list of valid range settings.
<code>memHandle</code>	Handle for Windows buffer from which data will be output. This buffer must have been previously allocated with the <a href="#">WinBufAlloc()</a> method and data values loaded (perhaps using <a href="#">WinArrayToBuf()</a> ).
<code>options</code>	Bit fields that control various options. Refer to the constants in the "options parameter values" section on page 222.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

`rate` - actual sampling rate used.

**options parameter values:**

All of the options settings are `MccDaq.ScanOptions` enumerated constants. To set a variable to one of these constants, you must refer to the `MccDaq` object and the `ScanOptions` enumeration (*variable* = `MccDaq.ScanOptions.Continuous`, *variable* = `MccDaq.ScanOptions.Background`, etc.).

Continuous	This option may only be used with boards which support interrupt, DMA or REP-INSW transfer methods. This option puts the method in an endless loop. Once it outputs the specified number ( <code>numPoints</code> ) of D/A values, it resets to the start of the buffer and begins again. The only way to stop this operation is by calling <a href="#">StopBackground()</a> with <code>AoFunction</code> . This option should only be used in combination with <code>Background</code> so that your program can regain control.
Background	This option may only be used with boards which support interrupt, DMA or REP-INSW transfer methods. When this option is used the D/A operations will begin running in the background and control will immediately return to the next line of your program. Use <a href="#">GetStatus()</a> with <code>AoFunction</code> to check the status of background operation. Alternatively, some boards support <a href="#">EnableEvent()</a> for event notification of changes in status of <code>Background</code> scans. Use <a href="#">StopBackground()</a> with <code>AoFunction</code> to terminate background operations before they are completed. <code>StopBackground()</code> should be executed after normal termination of all background functions in order to clear variables and flags.
Simultaneous	When this option is used (if the board supports it and the appropriate switches are set on the board) all of the D/A voltages will be updated simultaneously when the last D/A in the scan is updated. This generally means that all the D/A values will be written to the board, then a read of a D/A address causes all D/As to be updated with new values simultaneously.
ExtClock	If this option is used then conversions will be paced by the signal on the external clock input rather than by the internal pacer clock. Each conversion will be triggered on the appropriate edge of the clock input signal (refer to board-specific information contained in the <i>Universal Library Users Guide</i> ). When this option is used the <code>Rate</code> parameter is ignored. The sampling rate is dependent on the clock signal. Options for the board will default to transfer types that allow the maximum conversion rate to be attained unless otherwise specified.
ExtTrigger	If this option is specified the sampling will not begin until the trigger condition is met. On many boards, this trigger condition is programmable. Refer to the <a href="#">SetTrigger()</a> method and to board-specific information contained in the <i>Universal Library Users Guide</i> for details.
ADCClockTrig	Triggers a data output operation when the ADC clock starts.
ADCClock	Paces the data output operation using the ADC clock.

**Caution!** You will generate an error if you specify a total D/A rate beyond the capability of the board. For example: If you specify `LowChan = 0` and `HighChan = 3` (4 channels total) and `Rate = 100,000`, and you are using a cSBX-DDA04, you will get an error. You have specified a total rate of  $4 \times 100,000 = 400,000$ . The cSBX-DDA04 is rated to 330,000 updates per second. The maximum update rate depends on the D/A board that is being used. It is also dependent on the sampling mode options.

## APretrig()

Waits for a trigger to occur and then returns a specified number of analog samples before and after the trigger occurred. If only 'polled gate' triggering is supported, the trigger input line (refer to the user's manual for the board) must be at TTL low before this method is called, or a `TrigState` error will occur. The trigger occurs when the trigger condition is met. Refer to the [SetTrigger\(\)](#) method for more details.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**

```
Public Function APretrig(ByVal lowChan As Integer, ByVal highChan As Integer, ByRef pretrigCount As Integer, ByRef totalCount As Integer, ByRef rate As Integer, ByVal range As MccDaq.Range, ByVal memHandle As Integer, ByVal options As MccDaq.ScanOptions) As MccDaq.ErrorInfo
```

**C# .NET:**

```
public MccDaq.ErrorInfo APretrig(int lowChan, int highChan, ref int pretrigCount, ref int totalCount, ref int rate, MccDaq.Range range, int memHandle, MccDaq.ScanOptions options)
```

### Parameters:

`lowChan` First A/D channel of scan.

`highChan` Last A/D channel of scan.

**lowChan/highChan** - The maximum allowable channel depends on which type of A/D board is being used. For boards with both single ended and differential inputs, the maximum allowable channel number also depends on how the board is configured (e.g., 8 channels for differential inputs, 16 for single ended inputs).

`pretrigCount` Number of pre-trigger A/D samples to collect. Specifies the number of samples to collect before the trigger occurs.

For products using a hardware implementation of pretrigger (most products), `pretrigCount` must be less than the `(totalCount - 512)`. For these devices, if the trigger occurs too early, fewer than the requested number of pre-trigger samples will be collected, and a `TooFew` error will occur. The `pretrigCount` will be set to indicate how many samples were actually collected. The post trigger samples will still be collected.

For software implementations of pretrigger, `pretrigCount` must be less than `totalCount`. For these devices, triggers that occur before the requested number of pre-trigger samples are collected are ignored. Refer to board-specific information contained in the *Universal Library Users Guide*.

`totalCount` Total number of A/D samples to collect. Specifies the total number of samples that will be collected and stored in the buffer.

For products using a hardware implementation of pretrigger (most products), `totalCount` must be greater than or equal to the `pretrigCount + 512`. If the trigger occurs too early, fewer than the requested number of samples will be collected, and a `TooFew` error will occur. The `totalCount` will be set to indicate how many samples were actually collected.

For software implementations of pretrigger, `totalCount` must be greater than `pretrigCount`. For these devices, triggers that occur before the requested number of pre-trigger samples are collected are ignored. Refer to board-specific information contained in the *Universal Library Users Guide*.

`totalCount` must be evenly divisible by the number of channels being scanned. If it is not, this method will adjust the number (down) to the next valid value and return that value to the `totalCount` parameter.

	<p>pretrigCount must also be evenly divisible by the number of channels being scanned. If it is not, this function will adjust the number (up) to the next valid value and return that value to the pretrigCount parameter.</p>
rate	Sample rate in scans per second.
range	A/D Range code. If the selected A/D board does not have a programmable gain feature, this parameter is ignored. Otherwise, set to any range that is supported by the selected A/D board. Refer to board specific information for a list of the supported A/D ranges of each board. Refer to Table 7 on page 208 for a list of valid range settings.
memHandle	<p>Handle for Windows buffer to store data. This buffer must have been previously allocated with the <a href="#">WinBufAlloc()</a> method.</p> <p>For hardware trigger types, the buffer referenced by memHandle must be big enough to hold at least totalCount + 512 integers.</p>
options	Bit fields that control various options . Refer to the constants in the "options parameter values" section below.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

pretrigCount - Number of pre-trigger samples

totalCount - Total number of samples collected

rate - Actual sampling rate

memHandle - Collected A/D data returned via the Windows buffer

**options parameter values:**

All of the options settings are MccDaq.ScanOptions enumerated constants. To set a variable to one of these constants, you must refer to the MccDaq object and the ScanOptions enumeration (*variable* = MccDaq.ScanOptions.DtConnect, *variable* = MccDaq.ScanOptions.ExtMemory, etc.).

ConvertData	<p>For hardware trigger types, the data is collected into a "circular" buffer. When the data collection is complete, the data is in the wrong order. If you use the ConvertData option, the data is automatically rotated into the correct order (and converted to 12 bit values if required) when the data acquisition is complete. Otherwise, call <a href="#">AConvertPretrigData()</a> to rotate the data. You cannot use the ConvertData option in combination with the Background option for this method. The ConvertData option is not required for software triggered types.</p>
Background	<p>If the Background option is not used, the APretrig() method will not return to your program until all of the requested data has been collected and returned to the buffer. When the Background option is used, control returns immediately to the next line in your program, and the data collection from the A/D into the buffer will continue in the background. Use <a href="#">GetStatus()</a> with AiFunction to check on the status of the background operation. Alternatively, some boards support <a href="#">EnableEvent()</a> for event notification of changes in status of Background scans. Use <a href="#">StopBackground()</a> with AiFunction to terminate the background process before it has completed.</p> <p>Call StopBackground() after normal termination of all background functions to clear variables and flags.</p> <p>For hardware trigger types, you cannot use the ConvertData option in combination with the Background option for this method. To correctly order and parse the data, use <a href="#">AConvertPretrigData()</a> after the function completes.</p>



ExtClock	This option is available only for boards that have separate inputs for external pacer and external trigger. Refer to your hardware manual or board-specific information.
ExtMemory	<p>Causes this method to send the data to a connected memory board via the DT-Connect interface rather than returning the data to the buffer. If you use this option to send the data to a MEGA-FIFO memory board, then you must use <a href="#">MemReadPretrig()</a> to later read the pre-trigger data from the memory board. If you use <a href="#">MemRead()</a>, the data will NOT be in the correct order.</p> <p>Every time this option is used, it overwrites any data already stored in the memory board. All data should be read from the board (with <a href="#">MemReadPretrig()</a>) before collecting any new data. When this option is used, the memHandle parameter is ignored. The MEGA-FIFO memory must be fully populated in order to use the <code>APretrig()</code> method with the <code>ExtMemory</code> option.</p>
DTConnect	When the <code>DtConnect</code> option is used with this method the data from ALL A/D conversions is sent out the DT-Connect interface. While this method is waiting for a trigger to occur, it will send data out the DT-Connect interface continuously. If you have a Measurement Computing memory board plugged into the DT-Connect interface then you should use the <code>ExtMemory</code> option rather than this option.

**Important:**

For hardware trigger types, the buffer referenced by `memHandle` must be big enough to hold at least `totalCount + 512` integers.

## ATrig()

Waits for a specified analog input channel to go above or below a specified value. ATrig continuously reads the specified channel and compares its value to trigValue. Depending on whether trigType is set to TrigAbove or TrigBelow, it waits for the first A/D sample that is above or below trigValue. The first sample that meets the trigger criteria is returned to dataValue.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**

```
Public Function ATrig(ByVal chan As Integer, ByVal trigType As MccDaq.TriggerType, ByVal trigValue As Short, ByVal range As MccDaq.Range, ByRef dataValue As Short) As MccDaq.ErrorInfo
Public Function ATrig(ByVal chan As Integer, ByVal trigType As MccDaq.TriggerType, ByVal trigValue As System.UInt16, ByVal range As MccDaq.Range, ByRef dataValue As System.UInt16) As MccDaq.ErrorInfo
```

**C# .NET:**

```
public MccDaq.ErrorInfo ATrig(int chan, MccDaq.TriggerType trigType, short trigValue, MccDaq.Range range, out short dataValue) public MccDaq.ErrorInfo ATrig(int chan, MccDaq.TriggerType trigType, ushort trigValue, MccDaq.Range range, out ushort dataValue)
```

### Parameters:

chan	A/D channel number. The maximum allowable channel depends on which type of A/D board is being used. For boards with both single ended and differential inputs, the maximum allowable channel number also depends on how the board is configured. For example a CIO-DAS1600 has eight channels for differential inputs and 16 channels for single-ended inputs.
trigType	MccDaq.TriggerType.TrigAbove or MccDaq.TriggerType.TrigBelow. Specifies whether to wait for the analog input to be above or below the specified trigger value.
trigValue	The threshold value that all A/D values are compared to. Must be in the range 0 - 4095 for 12 bit A/D boards, or 0-65,535 for 16-bit A/D boards. Refer to your BASIC manual for information on signed BASIC integer data types.
range	Gain code. If the selected A/D board does not have a programmable gain feature, this parameter is ignored. Otherwise, set to any range that is supported by the selected A/D board. Refer to Table 7 on page 208 for a list of valid range settings. Refer to board-specific information for a list of the supported A/D ranges of each board.
dataValue	Returns the value of the first A/D sample to meet the trigger criteria.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

dataValue - value of the first A/D sample to match the trigger criteria.

### Notes:

Ctrl-C will not terminate the wait for an analog trigger that meets the specified condition. There are only two ways to terminate this call: satisfy the trigger condition or reset the computer.

**Caution!** Use caution when using this method in Windows programs. All active windows will lock on the screen until the trigger condition is satisfied. All keyboard and mouse activity will also lock until the trigger condition is satisfied.

## VIn()

Reads an A/D input channel, and returns a voltage value. If the specified A/D board has programmable gain, then this function sets the gain to the specified range. The voltage value is returned to `dataValue`.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**                   Public Function VIn(ByVal channel As Integer, ByVal range As MccDaq.Range, ByRef dataValue As Single, ByVal options As MccDaq.VInOptions) As MccDaq.ErrorInfo

**C# .NET:**                   public MccDaq.ErrorInfo VIn (System.Int32 channel, MccDaq.Range range, System.Single dataValue, MccDaq.VInOptions options)

### Parameters:

<code>channel</code>	A/D channel number. The maximum allowable channel depends on which type of A/D board is being used. For boards with both single ended and differential inputs, the maximum allowable channel number also depends on how the board is configured.
<code>range</code>	A/D range code. If the board has a programmable gain, it will be set according to this parameter value. Keep in mind that some A/D boards have a programmable gain feature, and others set the gain via switches on the board. In either case, the range that the board is configured for must be passed to this method. Refer to Table 7 on page 208 for a list of valid range settings.
<code>dataValue</code>	Reference to the data value.
<code>options</code>	Reserved for future use.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

`dataValue` - The value in volts of the A/D sample.

### options parameter values:

<code>Default</code>	Reserved for future use.
----------------------	--------------------------

## VOut()

Sets the value of a D/A output.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**                   Public Function VOut(ByVal channel As Integer, ByVal range As MccDaq.Range, ByVal dataValue As Single, ByVal options As MccDaq.VOutOptions) As MccDaq.ErrorInfo

**C# .NET:**                   public MccDaq.ErrorInfo VOut(System.Int32 channel, MccDaq.Range range, System.Single dataValue, MccDaq.VOutOptions options)

### Parameters:

channel	The D/A channel number. The maximum allowable channel depends on which type of D/A board is being used.
range	The D/A range code. If the board has a programmable gain, it will be set according to this parameter value. The output range of the D/A channel can be set to any of those supported by the board. Keep in mind that some D/A boards have programmable voltage ranges, and others set the voltage range via switches on the board. In either case, the range that the board is configured for must be passed to this method.
dataValue	The voltage value to be written.
options	Reserved for future use.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

### options parameter values:

Default	Reserved for future use.
---------	--------------------------

---

## Configuration Methods and Properties

### Introduction

This section covers Universal Library for .NET methods and properties that retrieve or change configuration options on a board. The configuration information for all boards is stored in the configuration file CB.CFG. This information is loaded from CB.CFG by all programs that use the library.

To determine which of these methods are compatible with your hardware, refer to the board-specific information contained in the *Universal Library User's Guide* (available in PDF format on our website at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf)).

## BoardConfig property

Represents an instance of the [cBoardConfig](#) class. Use this property to call the board configuration methods.

Member of the [MccBoard](#) class.

### Property prototype:

```
VB .NET:      Public ReadOnly Property BoardConfig As MccDaq.cBoardConfig
C# .NET      public MccDaq.cBoardConfig BoardConfig [get]
```

### Methods:

Over 20 UL for .NET configuration methods are accessible only from the BoardConfig property. Before you call any of these methods, you need to create an instance of an MccBoard object.

```
Dim MyBoard As MccDaq.MccBoard
MyBoard = New MccDaq.MccBoard (MyBoardNum)
```

To call a method from the BoardConfig property, use the notation shown in the example below.

```
MyErrorInfo = MyBoard.BoardConfig.GetBoardType (MyBoardType)
```

Each method available from the BoardConfig property is explained below.

## BoardConfig.DACUpdate()

Updates the voltage values on analog output channels. This method is usually called after a [SetDACUpdateMode\(\)](#) method call with its configVal parameter set to 1 (on command).

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

### Function prototype:

```
VB .NET:      Public Function DACUpdate() As MccDaq.ErrorInfo
C# .NET:      public MccDaq.ErrorInfo DACUpdate()
```

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

## BoardConfig.GetAdRetrigCount()

Gets the number of samples to acquire during each trigger event when ScanOptions.RetrigMode is enabled.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

### Function prototype:

```
VB .NET:      Public Function GetAdRetrigCount (ByRef retrigCount As Integer) As
MccDaq.ErrorInfo
C# .NET:      public MccDaq.ErrorInfo GetAdRetrigCount (System.Int32 retrigCount)
```

### Parameters:

retrigCount	Specifies the number of samples to acquire per trigger event when RetrigMode is set.
-------------	--

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

## BoardConfig.GetBaseAdr()

Gets the base address used by the Universal Library to communicate with a board. This is recommended for use only with ISA bus boards.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

### Function prototype:

VB .NET:                   Public Function GetBaseAdr(ByVal devNum As Integer, ByRef configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo GetBaseAdr(int devNum, out int configVal)

### Parameters:

devNum	Number of the base address to return (PCI boards may have several address ranges).
configVal	The board's base address.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

## BoardConfig.GetBoardType()

Gets the unique number (device ID) assigned to the board (between 0 and 8000h) indicating the type of board installed.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

### Function prototype:

VB .NET:                   Public Function GetBoardType(ByRef configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo GetBoardType(out int configVal)

### Parameters:

configVal	Returns a number indicating the board type.
-----------	---

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

## BoardConfig.GetCiNumDevs()

Gets the number of counter devices on the board.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

### Function prototype:

VB .NET:                   Public Function GetCiNumDevs(ByRef configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo GetCiNumDevs(out int configVal)

### Parameters:

configVal	Returns the number of counter devices.
-----------	--

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

## BoardConfig.GetClock()

Gets the counter's clock frequency in MHz (40, 10, 8, 6, 5, 4, 3, 2, 1), or 0 for not supported.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

### Function prototype:

VB .NET:                   Public Function GetClock(ByRef configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo GetClock(out int configVal)

### Parameters:

configVal                   Clock frequency in MHz.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

## BoardConfig.GetDACStartup()

Returns the board's configuration register STARTUP bit setting. Refer to the "[Notes](#)" section for the [SetDACStartup\(\)](#) method on page 240 for more information.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

### Function prototype:

VB .NET:                   Public Function GetDACStartup(ByVal devNum As Integer, ByRef configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo GetDACStartup(int devNum, out int configVal)

### Parameters:

devNum                    The number of the DAC channel whose startup bit setting you want to get.

configVal                   Returns the setting of the startup bit (0 or 1).

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

configVal                   Returns 0 if startup bit is disabled, or 1 if startup bit is enabled.

## BoardConfig.GetDACUpdateMode()

Returns the update mode for a digital-to-analog converter (DAC).

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

### Function prototype:

VB .NET:                   Public Function GetDACUpdateMode(ByVal configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo GetDACUpdateMode(out int configVal)

### Parameters:

configVal                   Returns a number indicating the DAC update mode (0 = *immediate*, 1 = *on command*).



**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

configVal                      If ConfigVal returns 0, the DAC update mode is immediate. Values written with [AOut \(\)](#) or [AOutScan \(\)](#) are automatically output by the DAC channels. If ConfigVal returns 1, the DAC update mode is set to *on command*. Values written with [AOut \(\)](#) or [AOutScan \(\)](#) are not output by the DAC channels until a [DACUpdate \(\)](#) method call is made.

**BoardConfig.GetDeviceID()**

Returns the name that identifies the instance of a device.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:                      Public Function GetDeviceId(ByRef configVal As String, ByRef maxLen As Integer) As MccDaq.ErrorInfo

C# .NET:                      public MccDaq.ErrorInfo GetDeviceId(System.String configVal, System.Int32 maxLen)

**Parameters:**

configVal                      Returns a string containing the name that identifies the device.

maxLen                         Specifies the maximum number of bytes to read, and returns the number of bytes that were actually read.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**BoardConfig.GetDeviceNotes()**

Returns the device notes that are stored in the device's memory.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:                      Public Function GetDeviceNotes(ByVal start As Integer, ByRef configVal As String, ByRef maxLen As Integer) As MccDaq.ErrorInfo

C# .NET:                      public MccDaq.ErrorInfo GetDeviceNotes(System.Int32 start , System.String configVal, System.Int32 maxLen)

**Parameters:**

start                            The start address of the device's memory to begin reading.

maxLen                         The maximum number of bytes to read from the device's memory. Returns the number of bytes actually read.

configVal                      Returns the text stored in the device's memory.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**BoardConfig.GetDiNumDevs()**

Gets the number of digital devices on the board.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:                   Public Function GetDiNumDevs(ByRef configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo GetDiNumDevs(out int configVal)

**Parameters:**

configVal                   Returns the number of digital devices.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**BoardConfig.GetDmaChan()**

Gets the DMA channel (0, 1 or 3) set for the board.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:                   Public Function GetDmaChan(ByRef configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo GetDmaChan(out int configVal)

**Parameters:**

configVal                   Returns DMA channel. 0, 1, or 3

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**BoardConfig.GetDtBoard()**

Gets the number of the board with the DT-Connect interface used to connect to external memory boards.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:                   Public Function GetDtBoard(ByRef configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo GetDtBoard(out int configVal)

**Parameters:**

configVal                   Returns the board number of the board that the external memory board is connected to.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**BoardConfig.GetIntLevel()**

Gets the interrupt level set for the board (0 for none, or 1 to 15).

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:           Public Function GetIntLevel (ByRef configVal As Integer) As MccDaq.ErrorInfo

C# .NET:           public MccDaq.ErrorInfo GetIntLevel(out int configVal)

**Parameters:**

configVal                 Returns the interrupt level (0 for none, or 1 – 15).

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**BoardConfig.GetNumAdChans()**

Gets the number of A/D channels.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:           Public Function GetNumAdChans (ByRef configVal As Integer) As MccDaq.ErrorInfo

C# .NET:           public MccDaq.ErrorInfo GetNumAdChans(out int configVal)

**Parameters:**

configVal                 Returns the number of A/D channels.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**BoardConfig.GetNumDaChans()**

Gets the number of D/A channels.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:           Public Function GetNumDaChans (ByRef configVal As Integer) As MccDaq.ErrorInfo

C# .NET:           public MccDaq.ErrorInfo GetNumDaChans(out int configVal)

**Parameters:**

configVal                 Returns the number of D/A channels.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**BoardConfig.GetNumExps()**

Gets the number of expansion boards.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:           Public Function GetNumExps (ByRef configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                    public MccDaq.ErrorInfo GetNumExps(out int configVal)

**Parameters:**

configVal                    Returns the number of expansion boards attached to the board.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**BoardConfig.GetNumIoPorts()**

Gets the number of I/O ports used by the board.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:                    Public Function GetNumIoPorts(ByRef configVal As Integer) As  
                              MccDaq.ErrorInfo

C# .NET:                    public MccDaq.ErrorInfo GetNumIoPorts(out int configVal)

**Parameters:**

configVal                    Returns the number of I/O ports used by the board.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**BoardConfig.GetPANID()**

Returns the Personal Area Network (PAN) identifier for wireless communication.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:                    Public Function GetPANID(ByRef configVal As Integer) As  
                              MccDaq.ErrorInfo

C# .NET:                    public MccDaq.ErrorInfo GetPANID(System.Int32 configVal)

**Parameters:**

configVal                    Returns a number from 0 to 65534 that identifies the Personal Area Network used  
                              for wireless communication.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**BoardConfig.GetRange()**

Gets the selected voltage range. For switch-selectable gains only.

If the selected A/D board does not have a programmable gain feature, this method returns the range as defined by the installed *InstaCal* settings. If *InstaCal* and the board are installed correctly, the range returned corresponds to the input range set by switches on the board. Refer to board-specific information for a list of the A/D ranges supported by each board.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:           Public Function GetRange(ByRef configVal As MccDaq.Range) As MccDaq.ErrorInfo

C# .NET:           public MccDaq.ErrorInfo GetRange(out MccDaq.Range configVal)

**Parameters:**

configVal               Returns the selected voltage range. Refer to Table 7 on page 208 for a list of valid configVal settings.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**BoardConfig.GetRFChannel()**

Returns the RF channel number that a wireless device uses to communicate.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:           Public Function GetRFChannel(ByRef configVal As Integer) As MccDaq.ErrorInfo

C# .NET:           public MccDaq.ErrorInfo GetRFChannel(System.Int32 configVal)

**Parameters:**

configVal               Returns the number of the RF channel selected for wireless communication. Valid channel numbers are 12 to 23.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**BoardConfig.GetRSS()**

Returns the signal strength in dBm of a signal received by a remote device.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:           Public Function GetRSS(ByRef configVal As Integer) As MccDaq.ErrorInfo

C# .NET:           public MccDaq.ErrorInfo GetRSS(System.Int32 configVal)

**Parameters:**

configVal               Returns the received signal strength in dBm of the remote device. In general, values are negative.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**BoardConfig.GetUsesExps()**

Gets the *True/False* value indicating support of expansion boards.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:           Public Function GetUsesExps(ByRef configVal As Integer) As MccDaq.ErrorInfo

C# .NET:           public MccDaq.ErrorInfo GetUsesExps(out int configVal)

**Parameters:**

configVal               Returns *True* if the board supports expansion boards, or *False* if the board does not support expansion boards.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**BoardConfig.GetWaitState()**

Gets the value of the Wait State jumper (1-enabled, 0-disabled).

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:           Public Function GetWaitState(ByRef configVal As Integer) As MccDaq.ErrorInfo

C# .NET:           public MccDaq.ErrorInfo GetWaitState(out int configVal)

**Parameters:**

configVal               Returns the wait state of the board.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**BoardConfig.SetAdRetrigCount()**

Sets the number of samples to acquire during each trigger event when ScanOptions.RetrigMode is enabled.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:           Public Function SetAdRetrigCount(ByRef retrigCount As Integer) As MccDaq.ErrorInfo

C# .NET:           public MccDaq.ErrorInfo SetAdRetrigCount(System.Int32 retrigCount)

**Parameters:**

retrigCount           Specifies the number of samples to acquire per trigger event when RetrigMode is set. Set to zero to use the value of the numPoints argument of the scan function.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**BoardConfig.SetBaseAdr()**

Sets the base address used by the Universal Library to communicate with a board. This is recommended for use only with ISA bus boards.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:                   Public Function SetBaseAdr(ByVal devNum As Integer, ByVal configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo SetBaseAdr(int devNum, int configVal)

**Parameters:**

devNum                   Number of the base address to configure (should always be 0 – can't configure PCI base addresses).

configVal               Sets the base address of the board.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**BoardConfig.SetClock()**

Sets the counter's clock source by the frequency (40, 10, 8, 6, 5, 4, 3, 2, 1), or 0 for not supported.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:                   Public Function SetClock(ByVal configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo SetClock(int configVal)

**Parameters:**

configVal               Sets the clock frequency in MHz.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**BoardConfig.SetDmaChan()**

Sets the DMA channel (0, 1 or 3).

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:                   Public Function SetDmaChan(ByVal configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo SetDmaChan(int configVal)

**Parameters:**

configVal               Sets the DMA channel to 0, 1 or 3.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**BoardConfig.SetDACStartup()**

Sets the board's configuration register STARTUP bit to 0 or 1 to enable/disable the storing of digital-to-analog converter (DAC) startup values. Each time the DAC board is powered up, the stored values are written to the DACs. New DAC start-up values are stored in memory by calling `AOut()` or `AOutScan()` after calling `SetDACStartup()` with the argument set to 1. Refer to the "[Notes](#)" section below for more information.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:                   Public Function SetDACStartup(ByVal configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo SetDACStartup(int configVal)

**Parameters:**

configVal                   Set to 0 to disable, or 1 to enable the storing of startup values for the channel.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**Notes:**

Use the SetDACStartup() method to store the DAC values you would like each DAC channel to be set to each time the board is powered up.

To store the current DAC values as start-up values, call SetDACStartup() with a configVal value of 1. Then, each time you call [AOut\(\)](#) or [AOutScan\(\)](#), the value written for each channel is stored in NV RAM. The last value written to a particular channel while SetDACStartup() is set to 1 is the value that that channel will be set to at power up. Call SetDACStartup() again with a configVal value of 0 to stop storing values in NV RAM.

**Example:**

```
DacBoard.BoardConfig.SetDACStartup(1);
for (int i =1; i <8; i++)
{
    DacBoard.AOut(i, BIP5VOLTS, DACValue[i]);
}
DacBoard.BoardConfig.SetDACStartup(0);
```

## BoardConfig.SetDACUpdateMode()

Sets the update mode for a digital-to-analog converter (DAC).

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:                   Public Function SetDACUpdateMode(ByVal devNum as Integer, ByVal configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo SetDACUpdateMode(int devNum, int configVal)

**Parameters:**

devNum                   Number of the channel whose update mode you want set.

configVal                   When set to 0, the DAC update mode is *immediate*. Values written with [AOut\(\)](#) or [AOutScan\(\)](#) are automatically output by the DAC channels.

When set to 1, the DAC update mode is *on command*. Values written with AOut() or AOutScan() are not output by the DAC channel(s) until a [DACUpdate\(\)](#) method call is made.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.



## BoardConfig.SetDeviceID()

Sets the name that identifies the instance of a device.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

### Function prototype:

VB .NET:                   Public Function SetDeviceId(ByVal configVal As String, ByRef maxLen As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo SetDeviceId(System.String configVal, System.Int32 maxLen)

### Parameters:

configVal                   Sets the string that contains the name identifying a device.

maxLen                     Specifies the maximum number of bytes to write, and returns the number of bytes that were actually written. For WLS Series devices, the string can contain up to 20 characters.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

## BoardConfig.SetDeviceNotes()

Sets the device notes to store in the device's memory.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

### Function prototype:

VB .NET:                   Public Function SetDeviceNotes(ByVal start As Integer, ByVal configVal As String, ByRef maxLen As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo SetDeviceNotes(System.Int32 start, System.String configVal, System.Int32 maxLen)

### Parameters:

start                       The start address of the device's memory to begin writing.

maxLen                     The maximum number of bytes to write to the device's memory. Returns the number of bytes actually written.

configVal                   The text to store in the device's memory.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

## BoardConfig.SetIntLevel()

Sets the interrupt level: 0 for none, or 1 to 15. Recommended for use only with ISA bus boards.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

### Function prototype:

VB .NET:                   Public Function SetIntLevel(ByVal configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo SetIntLevel(int configVal)

**Parameters:**

configVal                      Sets the interrupt level. Valid settings are 0 for none, or 1 – 15.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**BoardConfig.SetNumAdChans()**

Sets the number of A/D channels available on the board.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:                      Public Function SetNumAdChans (ByVal configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                      public MccDaq.ErrorInfo SetNumAdChans(int configVal)

**Parameters:**

configVal                      Sets the number of A/D channels on the board. Check board specific info for valid numbers. Note that this setting affects the single-ended/differential input mode of boards for which this setting is programmable.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**BoardConfig.SetPANID()**

Sets the Personal Area Network (PAN) identifier used for wireless communication.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:                      Public Function SetPANID (ByVal configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                      public MccDaq.ErrorInfo SetPANID(System.Int32 configVal)

**Parameters:**

configVal                      Sets the number (from 0 to 65534) that identifies the Personal Area Network used for wireless communication.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**BoardConfig.SetRange()**

Sets the selected voltage range. For use with boards for which the range is manually selected.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:                      Public Function SetRange (ByVal configVal As MccDaq.Range ) As MccDaq.ErrorInfo

C# .NET:                      public MccDaq.ErrorInfo SetRange (MccDaq.Range configVal)

**Parameters:**

configVal                      Range code.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**BoardConfig.SetRFChannel()**

Sets the RF channel number used for wireless communications.

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:                      Public Function SetRFChannel(ByVal configVal As Integer) As  
                                 MccDaq.ErrorInfo

C# .NET:                      public MccDaq.ErrorInfo SetRFChannel(System.Int32 configVal)

**Parameters:**

configVal                      Sets the number of the RF channel to use for wireless communications. Valid  
                                 channel numbers are 12 to 23.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**BoardConfig.SetWaitState()**

Sets the value of the Wait State jumper (1 = enabled, 0 = disabled).

Member of the [cBoardConfig](#) class. Accessible from the [MccBoard.BoardConfig](#) property.

**Function prototype:**

VB .NET:                      Public Function SetWaitState(ByVal configVal As Integer) As  
                                 MccDaq.ErrorInfo

C# .NET:                      public MccDaq.ErrorInfo SetWaitState(int configVal)

**Parameters:**

configVal                      Sets the wait state on the board.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

## BoardNum property

Number of the board associated with an instance of the `MccBoard` class.

Member of the [MccBoard](#) class.

### Property prototype:

VB .NET:	<code>Public ReadOnly Property BoardNum As Integer</code>
C# .NET:	<code>public int BoardNum [get]</code>

## CtrConfig property

Represents an instance of the [cCtrConfig](#) class. Use this property to call counter chip configuration methods.

Member of the [MccBoard](#) class.

### Property prototype:

```
VB .NET:      Public ReadOnly Property CtrConfig As MccDaq.cCtrConfig
C# .NET      public MccDaq.cCtrConfig CtrConfig [get]
```

### Methods:

The `GetCtrType()` configuration method is accessible only from the `CtrConfig` property. Before you call this method, you need to create an instance of an `MccBoard` object.

```
Dim MyBoard As MccDaq.MccBoard
MyBoard = New MccDaq.MccBoard(MyBoardNum)
```

To call this method from the `CtrConfig` property, use the notation shown in the example below:

```
MyErrorInfo = MyBoard.CtrConfig.GetCtrType(MyCtrNum, MyCtrType)
```

This method is explained below.

## CtrConfig.GetCtrType()

Gets the value that indicates the counter type.

Member of the [cCtrConfig](#) class. Accessible from the [MccBoard.CtrConfig](#) property.

### Function prototype:

```
VB .NET:      Public Function GetCtrType(ByVal devNum As Integer, ByRef configVal
As Integer) As MccDaq.ErrorInfo
C# .NET:      public MccDaq.ErrorInfo GetCtrType(int devNum, out int configVal )
```

### Parameters:

<code>devNum</code>	Number of the counter device.
<code>configVal</code>	Returns the type of counter where: 1 = 8254, 2 = 9513, 3 = 8536, 4 = 7266, 5 = event counter, 6 = scan counter, and 7 = timer counter.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

## DioConfig property

Represents an instance of the [cDioConfig](#) class. Use this property to call various digital I/O configuration methods.

Member of the [MccBoard](#) class.

### Property prototype:

```
VB .NET:      Public ReadOnly Property DioConfig As MccDaq.cDioConfig
C# .NET      public MccDaq.cDioConfig DioConfig [get]
```

### Methods:

Six configuration methods are accessible only from the DioConfig property. Before you call any of these methods, you need to create an instance of an MccBoard object.

```
Dim MyBoard As MccDaq.MccBoard
MyBoard = New MccDaq.MccBoard (MyBoardNum)
```

To call these methods from the DioConfig property, use the notation shown in the example below.

```
MyErrorInfo = MyBoard.DioConfig.GetNumBits (MyDevNum, MyNumBits)
```

These methods are explained below.

## DioConfig.GetDInMask()

Determines the bits on a specified port that are configured for input.

Member of the [cDioConfig](#) class. Accessible from the [MccBoard.DioConfig](#) property.

### Function prototype:

```
VB .NET:      Public Function GetDInMask(ByVal devNum As Integer, ByRef configVal
As Integer) As MccDaq.ErrorInfo
C# .NET:      public MccDaq.ErrorInfo GetDInMask(int devNum, out int configVal)
```

### Parameters:

devNum	Number of the port whose input bit configuration you want to determine.
configVal	Returns a bit mask showing the bit configuration of the specified port. Any of the lower eight bits that return a value of 1 are configured for input. Each of the upper eight bits always return 0.

### Returns:

An [ErrorInfo object](#) that indicates the status of the operation.

### Notes:

Use GetDInMask() with the GetDOutMask() method to determine if an AuxPort is configurable. If you apply both methods to the same port, and both configVal parameters returned have input and output bits that overlap, the port is not configurable. You can determine overlapping bits by *Anding* both parameters.

For example, the PCI-DAS08 has seven bits of digital I/O (four outputs and three inputs). For this board, the configVal parameter returned by GetDInMask() is always 7 (0000 0111), while the configVal parameter returned by GetDOutMask() is always 15 (0000 1111). When you *And* both configVal parameters together, you get a non-zero number (7). Any non-zero number indicates that input and output bits overlap for the specified port, and that port is a non-configurable AuxPort.

## DioConfig.GetDOutMask()

Determines the bits on a specified port that are configured for output.

Member of the [cDioConfig](#) class. Accessible from the [MccBoard.DioConfig](#) property.

### Function prototype:

VB .NET:                   Public Function GetDOutMask(ByVal devNum As Integer, ByRef configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo GetDOutMask(int devNum, out int configVal)

### Parameters:

devNum                   Number of the port whose output bit configuration you want to determine.

configVal               Returns a bit mask showing the bit configuration of the specified port. Any of the lower eight bits that return a value of 1 are configured for output. Each of the upper eight bits always return 0.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

### Notes:

Use [GetDInMask\(\)](#) with the [GetDOutMask\(\)](#) method to determine if an AuxPort is configurable. If you apply both methods to the same port, and both configVal parameters returned have input and output bits that overlap, the port is not configurable. You can determine overlapping bits by *Anding* both parameters.

For example, the PCI-DAS08 has seven bits of digital I/O (four outputs and three inputs). For this board, the configVal parameter returned by [GetDInMask\(\)](#) is always 7 (0000 0111), while the configVal parameter returned by [GetDOutMask\(\)](#) is always 15 (0000 1111). When you *And* both configVal parameters together, you get a non-zero number (7). Any non-zero number indicates that input and output bits overlap for the specified port, and that port is a non-configurable AuxPort.

## DioConfig.GetConfig()

Gets the configuration of a digital device (digital input or digital output).

Member of the [cDioConfig](#) class. Accessible from the [MccBoard.DioConfig](#) property.

### Function prototype:

VB .NET:                   Public Function GetConfig(ByVal devNum As Integer, ByRef configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo GetConfig(int devNum, out int configVal)

### Parameters:

devNum                   Number of the digital device.

configVal               Current configuration (1 = DigitalOut, 2 = DigitalIn).

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

## DioConfig.GetCurVal()

Gets the current value of digital outputs.

Member of the [cDioConfig](#) class. Accessible from the [MccBoard.DioConfig](#) property.

**Function prototype:**

VB .NET:                   Public Function GetCurVal(ByVal devNum As Integer, ByRef configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo GetCurVal(int devNum, out int configVal)

**Parameters:**

devNum                    Number of the digital device.  
configVal                 Current value of the digital output.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**DioConfig.GetDevType()**

Gets the device type of the digital port (AuxPort, FirstPortA, etc.).

Member of the [cDioConfig](#) class. Accessible from the [MccBoard.DioConfig](#) property.

**Function prototype:**

VB .NET:                   Public Function GetDevType(ByVal devNum As Integer, ByRef configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo GetDevType(int devNum, out int configVal)

**Parameters:**

devNum                    Number of the digital device.  
configVal                 Constant that indicates the type of device (AuxPort, FirstPortA, etc.).

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**DioConfig.GetNumBits()**

Gets the number of bits in the digital port.

Member of the [cDioConfig](#) class. Accessible from the [MccBoard.DioConfig](#) property.

**Function prototype:**

VB .NET:                   Public Function GetNumBits(ByVal devNum As Integer, ByRef configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo GetNumBits(int devNum, out int configVal)

**Parameters:**

devNum                    Number of the digital device.  
configVal                 Number of bits in the digital port.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.



## ExpansionConfig property

Represents an instance of the [cExpansionConfig](#) class. Use this property to call various expansion board configuration methods.

Member of the [MccBoard](#) class.

### Property prototype:

```
VB .NET:      Public ReadOnly Property ExpansionConfig As MccDaq.cExpansionConfig
C# .NET      public MccDaq.cExpansionConfig ExpansionConfig {get}
```

### Methods:

Over a dozen configuration methods are accessible only from the ExpansionConfig property. Before you call any of these methods, you need to create an instance of an MccBoard object.

```
Dim MyBoard As MccDaq.MccBoard
MyBoard = New MccDaq.MccBoard (MyBoardNum)
```

To call these methods from the ExpansionConfig property, use the notation shown in the example below.

```
MyErrorInfo = MyBoard.ExpansionConfig.GetBoardType (MyExpNum, MyExpType)
```

These methods are explained below.

## ExpansionConfig.GetBoardType()

Gets the expansion board type.

Member of the [cExpansionConfig](#) class. Accessible from the [MccBoard.ExpansionConfig](#) property.

### Function prototype:

```
VB .NET:      Public Function GetBoardType (ByVal devNum As Integer, ByRef
              configVal As Integer) As MccDaq.ErrorInfo
C# .NET:      public MccDaq.ErrorInfo GetBoardType (int devNum, out int configVal)
```

### Parameters:

devNum	Number of the expansion board.
configVal	Returns a number indicating the expansion board type (refer to the "BoardType Codes" topic in the <i>Universal Library User's Guide</i> ).

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

## ExpansionConfig.GetCjcChan()

Gets the channel that the CJC is connected to.

Member of the [cExpansionConfig](#) class. Accessible from the [MccBoard.ExpansionConfig](#) property.

### Function prototype:

```
VB .NET:      Public Function GetCjcChan (ByVal devNum As Integer, ByRef configVal
              As Integer) As MccDaq.ErrorInfo
C# .NET:      public MccDaq.ErrorInfo GetCjcChan (int devNum, out int configVal)
```

**Parameters:**

devNum	Number of the expansion board.
configVal	Returns a number indicating the channel on the A/D board that the CJC is connected to.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**ExpansionConfig.GetMuxAdChan1()**

Gets the first A/D channel that the EXP board is connected to.

Member of the [cExpansionConfig](#) class. Accessible from the [MccBoard.ExpansionConfig](#) property.

**Function prototype:**

VB .NET:	Public Function GetMuxAdChan1(ByVal devNum As Integer, ByRef configVal As Integer) As MccDaq.ErrorInfo
C# .NET:	public MccDaq.ErrorInfo GetMuxAdChan1(int devNum, out int configVal)

**Parameters:**

devNum	Number of the expansion board.
configVal	Number indicating the first A/D channel that the EXP board is connected to.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**ExpansionConfig.GetMuxAdChan2()**

Gets the second A/D channel that the EXP board is connected to.

Member of the [cExpansionConfig](#) class. Accessible from the [MccBoard.ExpansionConfig](#) property.

**Function prototype:**

VB .NET:	Public Function GetMuxAdChan2(ByVal devNum As Integer, ByRef configVal As Integer) As MccDaq.ErrorInfo
C# .NET:	public MccDaq.ErrorInfo GetMuxAdChan2(int devNum, out int configVal)

**Parameters:**

devNum	Number of the expansion board.
configVal	Number indicating the second A/D channel that the EXP board is connected to.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**ExpansionConfig.GetNumExpChans()**

Gets the number of expansion board channels.

Member of the [cExpansionConfig](#) class. Accessible from the [MccBoard.ExpansionConfig](#) property.

**Function prototype:**

VB .NET:	Public Function GetNumExpChans(ByVal devNum As Integer, ByRef configVal As Integer) As MccDaq.ErrorInfo
----------	---

C# .NET:                    `public MccDaq.ErrorInfo GetNumExpChans(int devNum, out int configVal)`

**Parameters:**

devNum                    Number of the expansion board.  
 configVal                Number of channels on the expansion board.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**ExpansionConfig.GetRange1()**

Gets the range/gain of the low 16 channels.

Member of the [cExpansionConfig](#) class. Accessible from the [MccBoard.ExpansionConfig](#) property.

**Function prototype:**

VB .NET:                    `Public Function GetRange1(ByVal devNum As Integer, ByRef configVal As Integer) As MccDaq.ErrorInfo`

C# .NET:                    `public MccDaq.ErrorInfo GetRange1(int devNum, out int configVal)`

**Parameters:**

devNum                    Number of the expansion board.  
 configVal                Returns the range (gain) of the low 16 channels.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**ExpansionConfig.GetRange2()**

Gets the range/gain of the high 16 channels.

Member of the [cExpansionConfig](#) class. Accessible from the [MccBoard.ExpansionConfig](#) property.

**Function prototype:**

VB .NET:                    `Public Function GetRange2(ByVal devNum As Integer, ByRef configVal As Integer) As MccDaq.ErrorInfo`

C# .NET:                    `public MccDaq.ErrorInfo GetRange2(int devNum, out int configVal)`

**Parameters:**

devNum                    Number of the expansion board.  
 configVal                Returns the range (gain) of the high 16 channels.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**ExpansionConfig.GetThermType()**

Gets the type of thermocouple or RTD configuration for the board (J, K, E, T, R, S, and B types).

Member of the [cExpansionConfig](#) class. Accessible from the [MccBoard.ExpansionConfig](#) property.



## ExpansionConfig.SetCjcChan()

Sets the channel that the CJC is connected to.

Member of the [cExpansionConfig](#) class. Accessible from the [MccBoard.ExpansionConfig](#) property.

### Function prototype:

VB .NET:                   Public Function SetCjcChan(ByVal devNum As Integer, ByVal configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo SetCjcChan(int devNum, int configVal)

### Parameters:

devNum                    Number of the expansion board.  
configVal                 Sets the A/D channel to connect to the CJC.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

## ExpansionConfig.SetMuxAdChan1()

Sets the first A/D channel that the EXP board is connected to.

Member of the [cExpansionConfig](#) class. Accessible from the [MccBoard.ExpansionConfig](#) property.

### Function prototype:

VB .NET:                   Public Function SetMuxAdChan1(ByVal devNum As Integer, ByVal configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo SetMuxAdChan1(int devNum, int configVal)

### Parameters:

devNum                    Number of the expansion board.  
configVal                 Number indicating the first A/D channel that the EXP board is connected to.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

## ExpansionConfig.SetMuxAdChan2()

Sets the second A/D channel that the EXP board is connected to.

Member of the [cExpansionConfig](#) class. Accessible from the [MccBoard.ExpansionConfig](#) property.

### Function prototype:

VB .NET:                   Public Function SetMuxAdChan2 (ByVal devNum As Integer, ByVal configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo SetMuxAdChan2(int devNum, int configVal)

### Parameters:

devNum                    Number of the expansion board.

configVal                 Number indicating the second A/D channel that the EXP board is connected to.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

## ExpansionConfig.SetRange1()

Sets the range/gain of the low 16 channels.

Member of the [cExpansionConfig](#) class. Accessible from the [MccBoard.ExpansionConfig](#) property.

### Function prototype:

VB .NET:                   Public Function SetRange1 (ByVal devNum As Integer, ByVal configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo SetRange1(int devNum, int configVal)

### Parameters:

devNum                    Number of the expansion board.

configVal                 Sets the range (gain) of the low 16 channels.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

## ExpansionConfig.SetRange2()

Sets the range/gain of the high 16 channels.

Member of the [cExpansionConfig](#) class. Accessible from the [MccBoard.ExpansionConfig](#) property.

### Function prototype:

VB .NET:                   Public Function SetRange2 (ByVal devNum As Integer, ByVal configVal As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo SetRange2(int devNum, int configVal)

### Parameters:

devNum                    Number of the expansion board.

configVal                 Sets the range (gain) of the high 16 channels.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**ExpansionConfig.SetThermType()**

Sets the type of thermocouple or RTD configuration for the board (J, K, E, T, R, S, and B types).

Member of the [cExpansionConfig](#) class. Accessible from the [MccBoard.ExpansionConfig](#) property.

**Function prototype:**

VB .NET:	Public Function SetThermType(ByVal devNum As Integer, ByVal configVal As Integer) As MccDaq.ErrorInfo
C# .NET:	public MccDaq.ErrorInfo SetThermType(int devNum, int configVal)

**Parameters:**

devNum	Number of the expansion board.
configVal	Number that sets the type of thermocouple configured for the board. (J = 1, K = 2, T = 3, E = 4, R = 5, S = 6, B = 7, Platinum .00392 = 257, Platinum .00391 = 258, Platinum .00385 = 259, Copper .00427 = 260, Nickel/Iron .00581 = 261, Nickel/Iron .00527 = 262)

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

## GetSignal()

Retrieves the configured Auxiliary or DAQ Sync connection and polarity for the specified timing and control signal.

This method is intended for advanced users. Except for the SYNC\_CLK input, you can easily view the settings for the timing and control signals using *InstaCal*.

Member of the [MccBoard](#) class.

**Note:** This method is not supported by all board types. Refer to the board-specific information contained in the *Universal Library User's Guide* (available in PDF format on our website at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf)).

### Function prototype:

VB .NET:      `Public Function GetSignal(ByVal direction As MccDaq.SignalDirection, ByVal signalType As MccDaq.SignalType, ByVal index As Integer, ByRef connectionPin As MccDaq.ConnectionPin, ByRef signalPolarity As MccDaq.SignalPolarity ) As MccDaq.ErrorInfo`

C# .NET:      `public MccDaq.ErrorInfo GetSignal(MccDaq.SignalDirection direction , MccDaq.SignalType signalType, int index, out MccDaq.ConnectionPin connectionPin, out MccDaq.SignalPolarity signalPolarity )`

### Parameters:

direction	Specifies whether retrieving the source (MccDaq.SignalDirection.SignalIn) or destination (MccDaq.SignalDirection.SignalOut).
signalType	Signal type whose connection is to be retrieved. Refer to "signalType parameter values" under the <a href="#">SelectSignal()</a> method section on page 260 for valid signal types.
index	Used to indicate which connection to reference when there is more than one connection associated with the output Signal type. When querying output signals, increment this value until BadIndex is returned or 0 is returned via the connection parameter to determine all the output connectionPins for the specified output Signal. The first connectionPin is indexed by 0.  For input signals (direction= MccDaq.SignalDirection.SignalIn), always set index to 0.
connectionPin	The specified connection is returned through this variable. Note that this is set to 0 if no connection is associated with the signalType, or if the index is set to an invalid value. Refer to "direction, connectionPin, and polarity parameter values" under the <a href="#">SelectSignal()</a> method section on page 260 for expected return values.
signalPolarity	Holds the polarity for the associated signalType and connectionPin.  For output signals assigned an AuxOut connectionPin, the return value is either MccDaq.SignalPolarity.Inverted or MccDaq.SignalPolarity.NonInverted.  For AdcConvert, DacUpdate, AdcTbSrc and DacTbSrc, input signals, either MccDaq.SignalPolarity.PositiveEdge or MccDaq.SignalPolarity.NegativeEdge are returned.  All other signals return 0.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.



**Notes:**

The above timing and control configuration information can also be viewed and edited inside *InstaCal*: Open *InstaCal*, click on the board, and press the **Configure...** button or menu item. If the board supports DAQ Sync and Auxiliary Input/Output signal connections, a button labeled **Advanced Timing & Control Configuration** displays. Press this button to open a display for viewing and modifying the above timing and control signals.

## NumBoards property

Returns the maximum number of boards you can install at one time.

Member of the [GlobalConfig](#) class.

### Property prototype:

VB .NET:                   Public Shared ReadOnly Property NumBoards As Integer  
C# .NET:                   public int NumBoards [get]

## NumExpBoards property

Returns the maximum total number of expansion boards you can install.

Member of the [GlobalConfig](#) class.

### Property prototype:

VB .NET:                   Public Shared ReadOnly Property NumExpBoards As Integer  
C# .NET:                   public static int NumExpBoards [get]

## SelectSignal()

Configures timing and control signals to use specific Auxiliary or DAQ Sync connections as a source or destination.

This method is intended for advanced users. Except for the SyncClk input, you can easily configure all the timing and control signals using *InstaCal*.

Member of the [MccBoard](#) class.

### SelectSignal is not supported by all boards

This method is not supported by all board types. Refer to the board-specific information contained in the *Universal Library User's Guide* (available in PDF format on our website at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf)).

### Function prototype:

**VB .NET:**

```
Public Function SelectSignal(ByVal direction As
MccDaq.SignalDirection , ByVal signalType As MccDaq.SignalType,
ByVal connectionPin As MccDaq.ConnectionPin, ByVal polarity As
MccDaq.SignalPolarity) As MccDaq.ErrorInfo
```

**C# .NET:**

```
public MccDaq.ErrorInfo SelectSignal(MccDaq.SignalDirection
direction, MccDaq.SignalType signal, MccDaq.ConnectionPin
connectionPin, MccDaq.SignalPolarity polarity)
```

### Parameters:

direction	<p>Direction of the specified signal type to be assigned a connector pin. For most signal types, this should be either <code>MccDaq.SignalDirection.SignalIn</code> or <code>MccDaq.SignalDirection.SignalOut</code>.</p> <p>For the <code>SyncClk</code>, <code>AdcTbSrc</code> and <code>DacTbSrc</code> signals, the external source can also be disabled by specifying <code>Disabled (=0)</code>, such that it is neither input nor output. Set it in conjunction with the <code>signalType</code>, <code>connectionPin</code>, and <code>polarity</code> arguments. Refer to the "direction, connectionPin, and polarity parameter values" section starting on page 260.</p>
signalType	<p>Signal type to be associated with a connector pin. Set it to one of the constants in the "signalType parameter values" section on page 260.</p>
connectionPin	<p>Designates the connector pin to associate the signal type and direction. Since individual pin selection is not allowed for the DAQ-Sync connectors, all DAQ-Sync pin connections are referred to as <code>DsConnector</code>. The <code>MccDaq.ConnectionPin.AuxIn</code> and <code>MccDaq.ConnectionPin.AuxOut</code> settings match their corresponding hardware pin names.</p>
polarity	<p><code>AdcTbSrc</code> and <code>DacTbSrc</code> input signals (<code>direction = MccDaq.SignalDirection.SignalIn</code>) can be set for either rising edge (<code>MccDaq.SignalPolarity.PositiveEdge</code>) or falling edge (<code>MccDaq.SignalPolarity.NegativeEdge</code>) signals. The <code>AuxOut</code> connections can be set to <code>MccDaq.SignalPolarity.Inverted</code> or <code>MccDaq.SignalPolarity.NonInverted</code> from their internal polarity.</p>

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

**signalType parameter values:**

All of the signalType settings are MccDaq.SignalType enumerated constants. To set a variable to one of these constants, you must refer to the MccDaq object and the SignalType enumeration (*variable* = MccDaq.SignalType.AdcConvert, *variable* = MccDaq.SignalType.AdcGate, etc.).

AdcConvert	A/D conversion pulse or clock.
AdcGate	External gate for A/D conversions.
AdcScanClk	A/D channel scan signal.
AdcScanStop	A/D scan completion signal.
ADC_SSH	A/D simultaneous sample and hold signal.
AdcStartScan	Start of A/D channel-scan sequence signal.
AdcStartTrig	A/D scan start trigger.
AdcStopTrig	A/D stop- or pre- trigger.
AdcTbSrc	A/D pacer timebase source.
Ctr1Clk	CTR1 clock source.
Ctr2Clk	CTR2 clock source.
DacStartTrig	D/A start trigger.
DacTbSrc	D/A pacer timebase source.
DacUpdate	D/A update signal.
DGnd	Digital ground.
SyncClk	STC timebase signal.

**direction, connectionPin, and polarity parameter values:**

- All of the direction settings are MccDaq.SignalDirection enumerated constants. To set a variable to one of these constants, you must refer to the MccDaq object and the SignalDirection enumeration (*variable* = MccDaq.SignalDirection.SignalIn, *variable* = MccDaq.SignalDirection.SignalOut, etc.).
- All of the connectionPin settings are MccDaq.ConnectionPin enumerated constants. To set a variable to one of these constants, you must refer to the MccDaq object and the ConnectionPin enumeration (*variable* = MccDaq.ConnectionPin.AuxIn0, *variable* = MccDaq.ConnectionPin.DsConnector, etc.).
- All of the polarity settings are MccDaq.SignalPolarity enumerated constants. To set a variable to one of these constants, you must refer to the MccDaq object and the SignalPolarity enumeration (*variable* = MccDaq.SignalPolarity.PositiveEdge, *variable* = MccDaq.ConnectionPin.Negative, etc.).

Valid input (direction= MccDaq.SignalDirection.SignalIn) settings include:		
signalType	connectionPin	polarity
AdcConvert	AuxIn0 to AuxIn5 DsConnector	PositiveEdge or NegativeEdge
AdcGate	AuxIn0 to AuxIn5 DsConnector	See <a href="#">SetTrigger</a> .
AdcStartTrig	AuxIn0 to AuxIn5 DsConnector	
AdcStopTrig	AuxIn0 to AuxIn5 DsConnector	

AdcTbSrc	AuxIn0 to AuxIn5	PositiveEdge or NegativeEdge
DacStartTrig	AuxIn0 to AuxIn5 DsConnector	Not assigned here.
DscTbSrc	AuxIn0 to AuxIn5	PositiveEdge or NegativeEdge
DacUpdate	AuxIn0 to AuxIn5 DsConnector	PositiveEdge or NegativeEdge
SyncClk	DsConnector	Not assigned here.

Valid output (direction= MccDaq.SignalDirection.SignalOut) settings include:		
signalType	connectionPin	polarity
AdcConvert	AuxIn0 to AuxIn5 DsConnector	Inverted* or NonInverted
AdcScanClk	AuxOut0..AuxOut2	
AdcScanStop	AuxOut0..AuxOut2	
AdcSsh	AuxOut0..AuxOut2 DsConnector	
AdcStartScan	AuxOut0..AuxOut2 DsConnector	
AdcStartTrig	AuxOut0..AuxOut2 DsConnector	
AdcStopTrig	AuxOut0..AuxOut2 DsConnector	
CtrlClk	AuxOut0..AuxOut2	
Ctr2Clk	AuxOut0..AuxOut2	
DacStartTrig	AuxOut0..AuxOut2 DsConnector	
DacUpdate	AuxOut0..AuxOut2 DsConnector	
DGND	AuxOut0..AuxOut2	Not assigned here.
SyncClk	DsConnector	
* Inverted is only valid for Auxiliary Output (AuxOut) connections.		

Valid disabled settings (direction = MccDaq.SignalDirection.Disabled):		
signalType	connectionPin	polarity
AdcTbSrc	Not assigned here.	Not assigned here.
DacTbSrc		
SyncClk		

**Notes:**

- You can view and edit the above timing and control configuration information from *InstaCal*. Open *InstaCal*, click on the board, and press the **Configure...** button or menu item. If the board supports DAQ Sync and Auxiliary Input/Output signal connections, an **Advanced Timing & Control Configuration** button displays. Press that button to open a display for viewing and modifying the above timing and control signals.
- Except for the AdcTbSrc, DacTbSrc and SyncClk signals, selecting an input signal connection does not necessarily activate it. Alternately, assigning an output signal to a connection does activate the signal upon performing the respective operation. For instance, when running an ExtClock AInScan(), AdcConvert SignalIn selects the connection to use as an external clock to pace the A/D conversions; if AInScan() is run without setting the ExtClock option, however, the selected connection is not activated and the signal at that connection is ignored. In both cases, the AdcConvert signal is output the connection(s) selected for the AdcConvert SignalOut. Since there are no scan options for enabling the Timebase Source and the SyncClk, selecting an input for the A/D or D/A Timebase Source, or SyncClk does activate the input source for the next respective operations.

- Multiple input signals can be mapped to the same `AuxIn` connection by successive calls to `SelectSignal()`; however, only one connection can be mapped to each input signal. If another connection had already been assigned to an input signal, the former selection is de-assigned and the new connection is assigned.
- Only one output signal can be mapped to the same `AuxOut` connection; however, multiple connections can be mapped to the same output signal by successive calls to `SelectSignal()`. If an output signal had already been assigned to a connection, then the former output signal is de-assigned and the new output signal is assigned to the connection.
- When selecting `DsConnector` for a signal, only one direction per signal type can be defined at a given time. Attempting to assign both Directions of a signal to the `DsConnector` results in only the latest selection being applied. If the signal type had formerly been assigned an input direction from the `DsConnector`, assigning the output direction for that signal type results in the input signal being reassigned to its default connection.
- `Adc_Tb_Src` and `Dac_Tb_Src` are intended to synchronize the timebase of the analog input and output pacers across two or more boards. Internal calculations of sampling and update rates assume that the external timebase has the same frequency as its internal clock. Adjust sample rates to compensate for differences in clock frequencies.

For instance, if the external timebase has a frequency of 10 MHz on a board that has a internal clock frequency of 40 MHz, the scan function samples or updates at a rate of about 1/4 the rate entered. However, while compensating for differences in external timebase and internal clock frequency, if the rate entered results in an invalid pacer count, the method returns a `BADRATE` error.

## SetTrigger()

Selects the trigger source and sets up its parameters. This trigger is used to initiate analog to digital conversions using the following Universal Library for .NET methods:

- [AInScan\(\)](#), if the `ExTrigger` option is selected.
- [CInScan\(\)](#), if the `ExTrigger` option is selected.
- [DInScan\(\)](#), if the `ExTrigger` option is selected.
- [APretrig\(\)](#)
- [FilePretrig\(\)](#)

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**

```
Public Function SetTrigger(ByVal trigType As MccDaq.TriggerType ,
    ByVal lowThreshold As Short, ByVal highThreshold As Short) As
    MccDaq.ErrorInfo

    Public Function SetTrigger(ByVal trigType As MccDaq.TriggerType,
    ByVal lowThreshold As System.UInt16, ByVal highThreshold As
    System.UInt16) As MccDaq.ErrorInfo
```

**C# .NET:**

```
public MccDaq.ErrorInfo SetTrigger(MccDaq.TriggerType trigType,
    short lowThreshold, short highThreshold)
public MccDaq.ErrorInfo SetTrigger(MccDaq.TriggerType trigType,
    ushort lowThreshold, ushort highThreshold)
```

### Parameters:

<code>trigType</code>	Specifies the type of triggering based on the external trigger source. Set it to one of the constants in the "trigType parameter values" section below.
<code>lowThreshold</code>	Selects the low threshold used when the trigger input is analog. The range depends upon the resolution of the trigger circuitry. Must be 0 to 255 for 8-bit trigger circuits, 0 to 4095 for 12-bit trigger circuits, and 0 to 65535 for 16-bit trigger circuits. Refer to the " <a href="#">Notes</a> " section on page 264.
<code>highThreshold</code>	Selects the high threshold used when the trigger input is analog. The range depends upon the resolution of the trigger circuitry. Must be 0 to 255 for 8-bit trigger circuits, 0 to 4095 for 12-bit trigger circuits, and 0 to 65535 for 16-bit trigger circuits. Refer to the " <a href="#">Notes</a> " section on page 264.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

### trigType parameter values:

All of the `trigType` settings are `MccDaq.TriggerType` enumerated constants. To set a variable to one of these constants, you must refer to the `MccDaq` object and the `TriggerType` enumeration (*variable* = `MccDaq.TriggerType.GateNegHys`, *variable* = `MccDaq.TriggerType.GatePosHys`, etc.).

Trigger Source	trigType	Explanation
Analog	GateNegHys	AD conversions are enabled when the external analog trigger input is more positive than <code>highThreshold</code> . AD conversions are disabled when the external analog trigger input more negative than <code>lowThreshold</code> . Hysteresis is the level between <code>lowThreshold</code> and <code>highThreshold</code> .
	GatePosHys	AD conversions are enabled when the external analog trigger input is more negative than <code>lowThreshold</code> . AD conversions are disabled when the external analog trigger input is more positive than <code>highThreshold</code> . Hysteresis is the level between <code>lowThreshold</code> and <code>highThreshold</code> .
	GateAbove	AD conversions are enabled as long as the external analog trigger input is more positive than <code>highThreshold</code> .
	GateBelow	AD conversions are enabled as long as the external analog trigger input is more negative than <code>lowThreshold</code> .
Analog	TrigAbove	AD conversions are enabled when the external analog trigger makes a transition from below <code>highThreshold</code> to above. Once conversions are enabled, the external trigger is ignored.
	TrigBelow	AD conversions are enabled when the external analog trigger input makes a transition from above <code>lowThreshold</code> to below. Once conversions are enabled, the external trigger is ignored.
	GateInWindow	AD conversions are enabled as long as the external analog trigger is inside the region defined by <code>lowThreshold</code> and <code>highThreshold</code> .
	GateOutWindow	AD conversions are enabled as long as the external analog trigger is outside the region defined by <code>lowThreshold</code> and <code>highThreshold</code> .
Digital	GateHigh	AD conversions are enabled as long as the external digital trigger input is 5 V (logic HIGH or 1).
	GateLow	AD conversions are enabled as long as the external digital trigger input is 0 V (logic LOW or 0).
	TrigHigh	AD conversions are enabled when the external digital trigger is 5 V (logic HIGH or '1'). Once conversions are enabled, the external trigger is ignored.
	TrigLow	AD conversions are enabled when the external digital trigger is 0 V (logic LOW or '0'). Once conversions are enabled, the external trigger is ignored.
	TrigPosEdge	AD conversions are enabled when the external digital trigger makes a transition from 0 V to 5 V (logic LOW to HIGH). Once conversions are enabled, the external trigger is ignored.
	TrigNegEdge	AD conversions are enabled when the external digital trigger makes a transition from 5 V to 0 V (logic HIGH to LOW). Once conversions are enabled, the external trigger is ignored.

**Notes:**

The value of the threshold must be within the range of the analog trigger circuit associated with the board. Refer to the board-specific information in the *Universal Library User's Guide*. For example, on the PCI-DAS1602/16, the analog trigger circuit handles  $\pm 10$  V. A value of 0 corresponds to -10 V, whereas a value of 65535 corresponds to +10 V.

If you are using signed integer types, the thresholds range from -32768 to 32767 for 16-bit boards, instead of from 0 to 65535. In this case, the unsigned value of 65535 corresponds to a value of -1, 65534 corresponds to -2, ..., 32768 corresponds to -32768.

For most boards that support analog triggering, you can pass the required trigger voltage level and the appropriate `Range` to `FromEngUnits` to calculate the `highThreshold` and `lowThreshold` values.



For some boards (refer to the "Analog Input Boards" chapter in the *Universal Library User's Guide*), you must manually calculate the threshold by first calculating the least significant bit (LSB) for a particular range for the trigger resolution of your hardware. You then use the LSB to find the threshold in counts based on an analog voltage trigger threshold.

To calculate the threshold, do the following:

1. Calculate the LSB by dividing the full scale range (FSR) by  $2^{\text{resolution}}$ . FSR is the entire span from  $-FS$  to  $+FS$  of your hardware for a particular range. For example, the full scale range of  $\pm 10$  V is 20 V.
2. Calculate how many times you need to add the LSB calculated in step 1 to the negative full scale ( $-FS$ ) to reach the trigger threshold value.

The maximum threshold value is  $2^{\text{resolution}} - 1$ . The formula is shown here:

$$\text{Abs}(-FS - \text{threshold in volts}) \div (\text{LSB}) = \text{threshold in counts}$$

Here are two examples that use this formula—one for 8-bit trigger resolution and one for 12-bit trigger resolution.

- 8-bit example using the  $\pm 10$  V range with a  $-5$  V threshold:

**Calculate LSB:**  $\text{LSB} = 20 \div 2^8 = 20 \div 256 = 0.078125$

**Calculate threshold:**  $\text{Abs}(-10 - (-5)) \div .078125 = 5 \div 0.078125 = 64$  (round this result if it is not an integer). A count of 64 translates to a voltage threshold of  $-5.0$  V.

- 12-bit example using the  $\pm 10$  V range with a  $+1$  V threshold:

**Calculate LSB:**  $\text{LSB} = 20 \div 2^{12} = 20 \div 4096 = 0.00488$

**Calculate threshold:**  $\text{Abs}(-10 - 1) \div .00488 = 11 \div 0.00488 = 2254$  (rounded from 2254.1). A count of 2254 translates to a voltage threshold of  $0.99952$  V.

## Version property

This information is used by the library to determine compatibility.

Member of the [GlobalConfig class](#).

### Property prototype:

VB .NET:	<code>Public Shared ReadOnly Property Version As Integer</code>
C# .NET:	<code>public int Version [get]</code>

---

## **Counter Methods**

### **Introduction**

This section covers Universal Library for .NET methods that load, read, and configure counters. There are five types of counter chips used in MCC counter boards: 8254's, 8536's, 7266's, 9513's and generic event counters. Some of the counter methods apply to only one type of counter.

## C7266Config()

Configures a 7266 counter for desired operation. This method can only be used with boards that contain a 7266 counter chip (Quadrature Encoder boards). For more information, refer to the LS7266R1 data sheet (ls7266r1.pdf) located in the "Documents" subdirectory of the installation.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**

```
Public Function C7266Config(ByVal counterNum As Integer, ByVal
    quadrature As MccDaq.Quadrature , ByVal countingMode As
    MccDaq.CountingMode , ByVal dataEncoding As MccDaq.DataEncoding,
    ByVal indexMode As MccDaq.IndexMode , ByVal invertIndex As
    MccDaq.OptionState , ByVal flagPins As MccDaq.FlagPins , ByVal
    gateState As MccDaq.OptionState ) As MccDaq.ErrorInfo
```

**C# .NET:**

```
public MccDaq.ErrorInfo C7266Config(int counterNum,
    MccDaq.Quadrature quadrature, MccDaq.CountingMode countingMode,
    MccDaq.DataEncoding dataEncoding, MccDaq.IndexMode indexMode,
    MccDaq.OptionState invertIndex, MccDaq.FlagPins flagPins,
    MccDaq.OptionState gateState)
```

### Parameters:

counterNum	Counter Number (1 - n), where n is the number of counters on the board.
quadrature	Selects the resolution multiplier for quadrature input, or disables quadrature input (NoQuad) so that the counters can be used as standard TTL counters. NoQuad, X1Quad, X2Quad or X4Quad.
countingMode	Selects operating mode for the counter. NormalMode, RangeLimit, NoRecycle, ModuloN. Set it to one of the constants in the "countingMode parameter values" section on page 269.
dataEncoding	Selects the format of the data that is returned by the counter - either Binary or BCD format. BinaryCount or BCDCount.
indexMode	Selects which action will be taken when the Index signal is received. The IndexMode must be set to IndexDisabled whenever a Quadrature is set to NOQuad or when GateState is set to Enabled. Set it to one of the constants in the "indexMode parameter values" section on page 269.
invertIndex	Selects the polarity of the Index signal. If set to Disabled, the Index signal is assumed to be positive polarity. If set to Enabled, the Index signal is assumed to be negative polarity.
flagPins	Selects which signals will be routed to the FLG1 and FLG2 pins. Set it to one of the constants in the "flagPins parameter values" section on page 269.
gateState	If gateState is set to ENABLED, then the channel INDEX input is routed to the RCNTR pin of the LS7266 chip, and is used as a gating signal for the counter. Whenever gateState = ENABLED the indexMode must be set to "IndexDisabled".

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

**countingMode parameter values:**

All of the `countingMode` settings are `MccDaq.CountingMode` enumerated constants. To set a variable to one of these constants, you must refer to the `MccDaq` object and the `CountingMode` enumeration (*variable* = `MccDaq.CountingMode.NormalMode`, *variable* = `MccDaq.CountingMode.NormalMode`, `CountingMode.RangeLimit`, etc.).

<code>NormalMode</code>	Each counter operates as a 24-bit counter that rolls over to 0 when the maximum count is reached.
<code>RangeLimit</code>	In range limit count mode, an upper and lower limit is set, mimicking limit switches in the mechanical counterpart. The upper limit is set by loading the PRESET register with the <a href="#">CLoad()</a> method after the counter has been configured. The lower limit is always 0. When counting up, the counter freezes whenever the count reaches the value that was loaded into the PRESET register. When counting down, the counter freezes at 0. In either case the counting is resumed only when the count direction is reversed.
<code>NoRecycle</code>	In non-recycle mode, the counter is disabled whenever a count overflow or underflow takes place. The counter is re-enabled when a reset or load operation is performed on the counter.
<code>ModuloN</code>	In <code>ModuloN</code> mode, an upper limit is set by loading the PRESET register with a maximum count. Whenever counting up, when the maximum count is reached, the counter will roll-over to 0 and continue counting up. Likewise when counting down, whenever the count reaches 0, it will roll over to the maximum count (in the PRESET register) and continue counting down.

**indexMode parameter values:**

All of the `indexMode` settings are `MccDaq.IndexMode` enumerated constants. To set a variable to one of these constants, you must refer to the `MccDaq` object and the `IndexMode` enumeration (*variable* = `MccDaq.IndexMode.IndexDisabled`, *variable* = `MccDaq.IndexMode.LoadCtr`, etc.).

<code>IndexDisabled</code>	The Index signal is ignored.
<code>LoadCtr</code>	The channel <code>INDEX</code> input is routed to the <code>LCNTR</code> pin of the LS7266 counter chip. The counter is loaded whenever the signal occurs.
<code>LoadOutLatch</code>	The channel <code>INDEX</code> input is routed to the <code>LCNTR</code> pin of the LS7266 counter chip. The current count is latched whenever the signal occurs. When this mode is selected, the <a href="#">CIn()</a> method will return the same count value each time it is called until the Index signal occurs.
<code>ResetCtr</code>	The channel <code>INDEX</code> input is routed to the <code>RCNTR</code> pin of the LS7266 counter chip. The counter is reset whenever the signal occurs.

**flagPins parameter values:**

All of the `flagPins` settings are `MccDaq.FlagPins` enumerated constants. To set a variable to one of these constants, you must refer to the `MccDaq` object and the `FlagPins` enumeration (*variable* = `MccDaq.FlagPins.CarryBorrow`, *variable* = `MccDaq.FlagPins.CompareBorrow`, etc.).

<code>CarryBorrow</code>	FLG1 pin is Carry output, FLG2 is Borrow output.
<code>CompareBorrow</code>	FLG1 pin is Compare output, FLG2 is Borrow output.
<code>CarryBorrowUpDown</code>	FLG1 pin is Carry/Borrow output, FLG2 is Up/Down signal.
<code>IndexError</code>	FLG1 pin is Index output, FLG2 is Error output.

## C8254Config()

Configures 8254 counter for desired operation. This method can only be used with 8254 counters. For more information, see the 82C54 data sheet in accompanying 82C54.pdf file located in the "Documents" subdirectory of the installation.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**                   Public Function C8254Config(ByVal counterNum As Integer, ByVal config As MccDaq.C8254Mode ) As MccDaq.ErrorInfo

**C# .NET:**                   public MccDaq.ErrorInfo C8254Config(int counterNum, MccDaq.C8254Mode config)

### Parameters:

counterNum	Selects one of the counter channels. An 8254 has 3 counters. The value may be 1 - n, where n is the number of 8254 counters on the board (refer to board-specific info in the ).
config	Refer to the 8254 data sheet for a detailed description of each of the configurations. Set it to one of the constants in the "config parameter values" section below.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

### config parameter values:

All of the config settings are MccDaq.C8254Mode enumerated constants. To set a variable to one of these constants, you must refer to the MccDaq object and the C8254Mode enumeration (*variable* = MccDaq.C8254Mode.HighOnLastCount, *variable* = MccDaq.C8254Mode.LastShot, etc.).

HighOnLastCount	Output of counter (OUT N) transitions from low to high on terminal count and remains high until reset. See Mode 0 in the 8254 data sheet in accompanying 82C54.pdf file located in the <i>Documents</i> subdirectory of the installation.
OneShot	Output of counter (OUT N) transitions from high to low on rising edge of GATE N, then back to high on terminal count. See mode 1 in the 8254 data sheet in accompanying 82C54.pdf file located in the <i>Documents</i> subdirectory of the installation.
RateGenerator	Output of counter (OUT N) pulses low for one clock cycle on terminal count, reloads counter and recycles. See mode 2 in the 8254 data sheet in accompanying 82C54.pdf file located in the <i>Documents</i> subdirectory of the installation.
SquareWave	Output of counter (OUT N) is high for count < 1/2 terminal count then low until terminal count, whereupon it recycles. This mode generates a square wave. See mode 3 in the 8254 data sheet in the accompanying 82C54.pdf file located in the <i>Documents</i> subdirectory of the installation.
SoftWareStrobe	Output of counter (OUT N) pulses low for one clock cycle on terminal count. Count starts after counter is loaded. See mode 4 in the 8254 data sheet in the accompanying 82C54.pdf file located in the <i>Documents</i> subdirectory of the installation.
HardwareStrobe	Output of counter (OUT N) pulses low for one clock cycle on terminal count. Count starts on rising edge at GATE N input. See mode 5 in the 8254 data sheet in accompanying 82C54.pdf file located in the <i>Documents</i> subdirectory of the installation.

## C8536Config()

Configures 8536 counter for desired operation. This method can only be used with 8536 counters. For more information, refer to the *Zilog 8536* product specification. The document is available on our web site at [www.mccdaq.com/PDFmanuals/Z8536.pdf](http://www.mccdaq.com/PDFmanuals/Z8536.pdf).

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**                      Configure for software triggering:

```
Public Shared Function C8536Config(ByVal counterNum As Integer,
    ByVal outputControl As Mccdaq.C8536OutputControl, ByVal recycleMode As
    MccDaq.RecycleMode, ByVal trigType As MccDaq.C8536TriggerType) As
    MccDaq.ErrorInfo
```

Configure for hardware triggering; use when existing code includes  
MccDaq.OptionState.

```
Public Function C8536Config(ByVal counterNum As Integer, ByVal
    outputControl As MccDaq.C8536OutputControl, ByVal recycleMode As
    MccDaq.RecycleMode, ByVal retrigger As MccDaq.OptionState) As
    MccDaq.ErrorInfo
```

**C# .NET:**                      Configure for software triggering:

```
public MccDaq.ErrorInfo C8536Config(int counterNum,
    MccDaq.C8536OutputControl outputControl, MccDaq.RecycleMode
    recycleMode, MccDaq.C8536TriggerType trigType)
```

Configure for hardware triggering; use when existing code includes  
MccDaq.OptionState.

```
public MccDaq.ErrorInfo C8536Config(int counterNum,
    MccDaq.C8536OutputControl outputControl, MccDaq.RecycleMode
    recycleMode, MccDaq.OptionState retrigger)
```

### Parameters:

counterNum	Selects one of the counter channels. An 8536 has three counters. The value may be 1, 2 or 3. INT32 Series boards have two chips installed, so the counterNum value may be 1 to 6.
outputControl	Specifies the action of the output signal. Set it to one of the constants in the "outputControl parameter values" section on page 272.
recycleMode	If set to Recycle (as opposed to OneTime), the counter automatically reloads to the starting count every time it reaches 0, and then counting continues.
retrigger	If set to Enabled, every trigger on the counter's trigger input initiates loading of the initial count. Counting proceeds from the initial count.
trigType	Specifies the trigger type. Set it to one of the constants in the "trigType parameter values" section on page 272.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

**outputControl parameter values:**

All of the `outputControl` settings are `MccDaq.C8536OutputControl` enumerated constants. To set a variable to one of these constants, refer to the `MccDaq` object and the `C8536OutputControl` enumeration (*variable* = `MccDaq.C8536OutputControl.HighPulseOnTc`, *variable* = `MccDaq.C8536OutputControl.ToggleOnTc`, etc.).

<code>HighPulseOnTc</code>	Output transitions from low to high for one clock pulse on terminal count.
<code>ToggleOnTc</code>	Output changes state on terminal count.
<code>HighUntilTc</code>	Output transitions to high at the start of counting, and then goes low on terminal count.

**trigType parameter values:**

All of the `trigType` settings are `MccDaq.C8536TriggerType` enumerated constants. To set a variable to one of these constants, refer to the `MccDaq` object and the `C8536TriggerType` enumeration (*variable* = `MccDaq.C8536OutputControl.HighPulseOnTc`, *variable* = `MccDaq.C8536OutputControl.ToggleOnTc`, etc.).

<code>HWStartTrig</code>	The first trigger on the counter's trigger input initiates loading of the initial count. Counting proceeds from the initial count.
<code>HWRetrig</code>	Every trigger on the counter's trigger input initiates loading of the initial count. Counting proceeds from the initial count.
<code>SWStartTrig</code>	The <a href="#">CLoad()</a> method initiates loading of the initial count. Counting proceeds from the initial count.



## C8536Init()

Initializes the counter linking features of an 8536 counter chip. The linking of counters 1 and 2 must be accomplished prior to enabling the counters.

Refer to the *Zilog 8536 product specification* for a description of the hardware affected by this mode. The document is available on our web site at [www.mccdaq.com/PDFmanuals/Z8536.pdf](http://www.mccdaq.com/PDFmanuals/Z8536.pdf).

Member of the [MccBoard](#) class.

### Function prototype:

VB .NET:                   Public Function C8536Init(ByVal chipNum As Integer, ByVal ctrlOutput As MccDaq.CtrlOutput) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo C8536Init(int chipNum, MccDaq.CtrlOutput ctrlOutput)

### Parameters:

chipNum	Selects one of the 8536 chips on the board, 1 to <i>n</i> .
ctrlOutput	Specifies how the counter 1 is to be linked to counter 2, if at all. Set it to one of the constants in the "ctrlOutput parameter values" section below.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

### ctrlOutput parameter values:

All of the ctrlOutput settings are MccDaq.CtrlOutput enumerated constants. To set a variable to one of these constants, you must refer to the MccDaq object and the CtrlOutput enumeration (*variable* = MccDaq.CtrlOutput.NotLinked, *variable* = MccDaq.CtrlOutput.GateCtr2, etc.).

NotLinked	Counter 1 is not connected to any other counter's inputs.
GateCtr2	Output of counter 1 is connected to the GATE of counter #2.
TrigCtr2	Output of counter 1 is connected to the trigger of counter #2.
InCtr2	Output of counter 1 is connected to counter #2 clock input.

## C9513Config()

Sets all of the configurable options of a 9513 counter. For more information, see the AM9513A data sheet in accompanying 9513A.pdf file located in the *Documents* subdirectory of the installation.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**

```
Public Function C9513Config(ByVal counterNum As Integer, ByVal
gateControl As MccDaq.GateControl , ByVal counterEdge As
MccDaq.CountEdge , ByVal counterSource As MccDaq.CounterSource ,
ByVal specialGate As MccDaq.OptionState , ByVal reload As
MccDaq.Reload , ByVal recycleMode As MccDaq.RecycleMode , ByVal
bcdMode As MccDaq.BCDMode , ByVal countDirection As
MccDaq.CountDirection , ByVal outputControl As
MccDaq.C9513OutputControl ) As MccDaq.ErrorInfo
```

**C# .NET:**

```
public MccDaq.ErrorInfo C9513Config(int counterNum,
MccDaq.GateControl gateControl, MccDaq.CountEdge counterEdge,
MccDaq.CounterSource counterSource, MccDaq.OptionState specialGate,
MccDaq.Reload reload, MccDaq.RecycleMode recycleMode, MccDaq.BCDMode
bcdMode, MccDaq.CountDirection countDirection,
MccDaq.C9513OutputControl outputControl)
```

### Parameters:

counterNum	Counter number (1 - n) where n is the number of counters on the board. (For example, a CIO-CTR5 has 5, a CIO-CTR10 has 10, etc. See board specific info).
gateControl	Sets the gating response for level, edge, etc. Set it to one of the constants in the "gateControl parameter values" section on page 275.
counterEdge	Which edge to count. Referred to as "Source Edge" in 9513 data book. Can be set to POSITIVEEDGE (count on rising edge) or NEGATIVEEDGE (count on falling edge).
counterSource	Each counter may be set to count from one of 16 internal or external sources. Set it to one of the constants in the "counterSource parameter values" section on page 275.
specialGate	Special gate may be enabled (MccDaq.OptionState.Enabled) or disabled (MccDaq.OptionState.Disabled).
reload	Reload the counter from the load register (reload = MccDaq.Reload.LoadReg) or alternately load from the load register, then the hold register (reload = MccDaq.Reload.LoadAndHoldReg).
recycleMode	Execute once (MccDaq.RecycleMode.OneTime) or reload and recycle (MccDaq.RecycleMode.Recycle).
bcdMode	Counter may operate in <i>binary coded decimal count</i> (MccDaq.BCDMode.BCDCount) or <i>binary count</i> (MccDaq.BCDMode.BinaryCount).
countDirection	AM9513 may count up (MccDaq.CountDirection.CountUp) or down (MccDaq.CountDirection.CountDown).
outputControl	The type of output desired. Set it to one of the constants in the "outputControl parameter values" on page 275.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

**gateControl parameter values :**

All of the `gateControl` settings are `MccDaq.GateControl` enumerated constants. To set a variable to one of these constants, you must refer to the `MccDaq` object and the `GateControl` enumeration (*variable* = `MccDaq.GateControl.NoGate`, *variable* = `MccDaq.GateControl.AhlTcPrevCtr`, etc.).

<code>NoGate</code>	No gating
<code>AhlTcPrevCtr</code>	Active high TCN -1
<code>AhlNextGate</code>	Active High Level GATE N + 1
<code>AhlPrevGate</code>	Active High Level GATE N - 1
<code>AhlGate</code>	Active High Level GATE N
<code>AllGate</code>	Active Low Level GATE N
<code>AheGate</code>	Active High Edge GATE N
<code>Alegate</code>	Active Low Edge GATE N

**counterSource parameter values:**

All of the `counterSource` settings are `MccDaq.CounterSource` enumerated constants. To set a variable to one of these constants, you must refer to the `MccDaq` object and the `CounterSource` enumeration (*variable* = `MccDaq.CounterSource.TcPrevCtr`, *variable* = `MccDaq.CounterSource.CtrInput1`, etc.).

<code>TcPrevCtr</code>	TCN - 1 (Terminal count of previous counter)
<code>CtrInput1</code>	SRC 1 (Counter Input 1)
<code>CtrInput2</code>	SRC 2 (Counter Input 2)
<code>CtrInput3</code>	SRC 3 (Counter Input 3)
<code>CtrInput4</code>	SRC 4 (Counter Input 4)
<code>CtrInput5</code>	SRC 5 (Counter Input 5)
<code>Gate1</code>	GATE 1
<code>Gate2</code>	GATE 2
<code>Gate3</code>	GATE 3
<code>Gate4</code>	GATE 4
<code>Gate5</code>	GATE 5
<code>Freq1</code>	F1
<code>Freq2</code>	F2
<code>Freq3</code>	F3
<code>Freq4</code>	F4
<code>Freq5</code>	F5

**outputControl parameter values :**

All of the `outputControl` settings are `MccDaq.9513OutputControl` enumerated constants. To set a variable to one of these constants, you must refer to the `MccDaq` object and the `9513OutputControl` enumeration (*variable* = `MccDaq.9513OutputControl.AlwaysLow`, *variable* = `MccDaq.9513OutputControl.HighPulseOnTc`, etc.).

<code>AlwaysLow</code>	Inactive, Output Low
------------------------	----------------------

HighPulseOnTc	High pulse on Terminal Count
ToggleOnTc	TC Toggled
Disconnected	Inactive, Output High Impedance
LowPulseOnTc	Active Low Terminal Count Pulse
3, 6, 7	(numeric values) Illegal

**Notes:**

The information provided here and in [C9513Init\(\)](#) will only help you understand how Universal Library syntax corresponds to the 9513 data sheet (refer to the accompanying 9513A.pdf file located in the *Documents* subdirectory of the installation). It is not a substitute for the data sheet. You cannot program and use a 9513 counter/timer without the data sheet.

## C9513Init()

Initializes all of the chip-level features of a 9513 counter chip. This method can only be used with 9513 counters. For more information, refer to the AM9513A data sheet in accompanying 9513A.pdf file located in the *Documents* subdirectory of the installation.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**

```
Public Function C9513Init(ByVal chipNum As Integer, ByVal
    foutDivider As Integer, ByVal foutSource As MccDaq.CounterSource,
    ByVal compare1 As MccDaq.CompareValue , ByVal compare2 As
    MccDaq.CompareValue , ByVal timeOfDay As MccDaq.TimeOfDay ) As
    MccDaq.ErrorInfo
```

**C# .NET:**

```
public MccDaq.ErrorInfo C9513Init(int chipNum, int foutDivider,
    MccDaq.CounterSource foutSource, MccDaq.CompareValue compare1,
    MccDaq.CompareValue compare2, MccDaq.TimeOfDay timeOfDay)
```

### Parameters:

chipNum	Specifies which 9513 chip is to be initialized. For a CTR05 board, set to 1. For a CTR10 board, set to either 1 or 2, and for a CTR20 set to 1-4.
foutDivider	F-Out divider (0-15). If set to 0, foutDivider is the rate of foutSource divided by 16. If set to a number between 1 and 15, foutDivider is the rate of foutSource divided by foutDivider.
foutSource	Specifies source of the signal for F-Out signal. Set it to one of the constants in the "foutSource parameter values" section on page 278.
compare1	MccDaq.CompareValue.Enabled or MccDaq.CompareValue.Disabled
compare2	MccDaq.CompareValue.Enabled or MccDaq.CompareValue.Disabled.
timeOfDay	MccDaq.TimeOfDay.Disabled, or three different enabled settings. Set it to one of the constants in the "timeOfDay parameter values" section on page 278.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

**foutSource parameter values:**

All of the `foutSource` settings are `MccDaq.CounterSource` enumerated constants. To set a variable to one of these constants, you must refer to the `MccDaq` object and the `CounterSource` enumeration (*variable* = `MccDaq.CounterSource.CtrInout1`, *variable* = `MccDaq.CounterSource.CtrInput2`, etc.).

<b>foutSource</b>	<b>9513 Data Sheet Equivalent</b>	<b>foutSource</b>	<b>9513 Data Sheet Equivalent</b>
<code>CtrInput1</code>	SRC 1 (Counter Input 1)	<code>Gate3</code>	GATE3
<code>CtrInput2</code>	SRC 2 (Counter Input 2)	<code>Gate4</code>	GATE4
<code>CtrInput3</code>	SRC 3 (Counter Input 3)	<code>Gate5</code>	GATE5
<code>CtrInput4</code>	SRC 4 (Counter Input 4)	<code>Freq1</code>	F1
<code>CtrInput5</code>	SRC 5 (Counter Input 5)	<code>Freq2</code>	F2
<code>Gate1</code>	GATE1	<code>Freq3</code>	F3
<code>Gate2</code>	GATE2	<code>Freq4</code>	F4

**timeOfDay parameter values:**

All of the `timeOfDay` settings are `MccDaq.TimeOfDay` enumerated constants. To set a variable to one of these constants, you must refer to the `MccDaq` object and the `TimeOfDay` enumeration (*variable* = `MccDaq.TimeOfDay.Disable`, *variable* = `MccDaq.TimeOfDay.One`, etc.).

<b>timeOfDay</b>	<b>9513 Data Sheet Equivalent</b>
<code>Disabled</code>	TOD Disabled
<code>One</code>	TOD Enabled/5 Input
<code>Two</code>	TOD Enabled/6 Input
<code>Three</code>	TOD Enabled/10 Input
<b>No parameters for</b>	<b>9513 Data Sheet Equivalent</b>
0 (FOUT on)	FOUT Gate
0 (Data bus matches board)	Data Bus Width
1 (Disable Increment)	Data Pointer Control
1 (BCD Scaling)	Scalar Control

**Notes:**

The information provided here and in [C9513Config\(\)](#) will only help you understand how Universal Library for .NET syntax corresponds to the 9513 data sheet (refer to the accompanying 9513A.pdf file located in the *Documents* subdirectory of the installation). It is not a substitute for the data sheet. You cannot program and use a 9513 counter/timer without the data sheet.

## CClear()

Clears a scan counter value (sets it to zero). This method only works with counter boards that have counter scan capability.

Member of the [MccBoard](#) class.

### Function prototype:

VB .NET:           Public Function CClear(ByVal counterNum As Integer) As  
                          MccDaq.ErrorInfo

C# .NET:           public MccDaq.ErrorInfo CClear(int counterNum)

### Parameters:

counterNum           The counter to clear.

**Note:** This parameter is zero-based (the first counter number to clear is "0").

### Returns:

[Error code](#) or 0 if no errors

## CConfigScan()

Configures a counter channel. This method only works with counter boards that have counter scan capability.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**

```
Public Function CConfigScan(ByVal counterNum As Integer, ByVal mode As MccDaq.CounterModes, ByVal debounceTime As MccDaq.DebounceTimes,
    ByVal debounceMode As MccDaq.DebounceModes, ByVal edgeDetection As MccDaq.EdgeDetection, ByVal tickSize As Integer, ByVal mapCounter As Integer) As MccDaq.ErrorInfo
```

**C# .NET:**

```
public MccDaq.ErrorInfo CConfigScan(int counterNum,
    MccDaq.CounterModes mode, MccDaq.DebounceTimes debounceTime,
    MccDaq.DebounceModes debounceMode, MccDaq.EdgeDetection
    edgeDetection, int tickSize, int mapCounter)
```

### Parameters:

counterNum	The counter to set up.  <b>Note:</b> This parameter is zero-based (the first counter number to set up is "0").
mode	Bit fields that control various options. All of the mode settings are MccDaq.CounterModes enumerated constants. Set it to one of the constants in the "mode parameter values" section below.
debounceTime	Used to bypass the debounce mode, or to set a channel's comparator output to one of 16 debounce times. Debounce is used to eliminate switch-induced transients typically associated with electro-mechanical devices including relays, proximity switches, and encoders.  All of the debounceTime settings are MccDaq.DebounceTimes enumerated constants. Set it to one of the constants in the "debounceTimes parameter values" section on page 281.
debounceMode	Sets the mode of the debounce module. The debounceTrigger settings are MccDaq.DebounceModes enumerated constants. Set it to one of the constants in the "debounceMode parameter values" section on page 282.
edgeDetection	Determines whether the rising edge or falling edge is to be detected. The edgeDetection settings are MccDaq.EdgeDetection enumerated constants. The choices are RisingEdge and FallingEdge.
tickSize	Reserved for future use.
mapCounter	Used to select the mapped counter. A mapped channel is one of the counter input channels other than counterNum that can participate with the input signal of the counter defined by counterNum by gating the counter or decrementing the counter.

### Returns:

[Error code](#) or 0 if no errors

### mode parameter values:

ClearOnRead	The counter is cleared at the beginning of every sample.
StopAtMax	The counter will stop at the top of its count. For the CIn32 () method, the top of the count depends on whether the Bit32 option is used. If it is, the top of the count is FFFFFFFF hex. If not, the top of the count is FFFF hex. By default, the counter counts upward and rolls over on the 32-bit boundary.



DecrementOn	<p>Allows the mapped channel to decrement the counter. With this option, the main counter will increment the counter, and the mapped counter can be used to decrement the counter. By default, the counter decrement option is set to "off".</p> <p>This mode is not compatible with <code>CIn()</code> or <code>CIn32()</code>. If a counter is configured for <code>DecrementOn</code>, calling <code>CIn()</code> or <code>CIn32()</code> for that counter will result in a <code>BADCOUNTERMODE</code> error.</p>
GatingOn	<p>Selects gating "on." When "on", the counter is enabled when the mapped channel to gate the counter is high. When the mapped channel is low, the counter is disabled but holds the count value.</p> <p>This mode is not compatible with <code>CIn()</code> or <code>CIn32()</code>. If a counter is configured for <code>GatingOn</code>, calling <code>CIn()</code> or <code>CIn32()</code> for that counter will result in a <code>BADCOUNTERMODE</code> error.</p>
LatchOnMap	<p>Causes the count to be latched by the signal on the mapped counter. By default, the count is latched by the internal "start of scan" signal, so the count is updated each time it's read.</p> <p>This mode is not compatible with <code>CIn()</code> or <code>CIn32()</code>. If a counter is configured for <code>LatchOnMap</code>, calling <code>CIn()</code> or <code>CIn32()</code> for that counter will result in a <code>BADCOUNTERMODE</code> error.</p>
Bit32	<p>Selects a 32-bit counter. This mode affects only <code>CIn32()</code> and <code>CIn()</code> and only when the counter is configured for <code>StopAtMax</code>. Recommended for use only with <code>CIn32()</code>. (Using the <code>Bit32</code> option with <code>CIn()</code> is not very useful, since the value returned by <code>CIn()</code> is only 16 bits. The effect is that the value returned by <code>CIn()</code> rolls over at 64k 65,535 times before stopping.)</p>
Encoder	Sets the specified counter to encoder mode.
EncoderModeX1	Sets the encoder measurement mode to X1.
EncoderModeX2	Sets the encoder measurement mode to X2.
EncoderModeX4	Sets the encoder measurement mode to X4.
LatchOnZ	Selects the Encoder Z mapped signal to latch the counter outputs. This allows the user to know the exact counter value when an edge is present on another counter.
ClearOnZOn	Selects "clear on Z" on. The counter is cleared on the rising edge of the mapped (Z) counter. By default, the "ClearOnZ" option is off, and the counter is not cleared.

**debounceTimes parameter values:**

Debounce500ns	Sets the counter channel's comparator output to 500 ns.
Debounce 1500ns	Sets the counter channel's comparator output to 1500 ns.
Debounce 3500ns	Sets the counter channel's comparator output to 3500 ns.
Debounce 7500ns	Sets the counter channel's comparator output to 7500 ns.
Debounce 15500ns	Sets the counter channel's comparator output to 15500 ns.
Debounce 31500ns	Sets the counter channel's comparator output to 31500 ns.
Debounce 63500ns	Sets the counter channel's comparator output to 63500 ns.
Debounce 127500ns	Sets the counter channel's comparator output to 127500 ns.
Debounce 100us	Sets the counter channel's comparator output to 100 us.
Debounce 300us	Sets the counter channel's comparator output to 300 us.
Debounce 700us	Sets the counter channel's comparator output to 700 us.

Debounce 1500us	Sets the counter channel's comparator output to 1500 us.
Debounce 3100us	Sets the counter channel's comparator output to 3100 us.
Debounce 6300us	Sets the counter channel's comparator output to 6300 us.
Debounce 12700us	Sets the counter channel's comparator output to 12700 us.
Debounce 25500us	Sets the counter channel's comparator output to 25500 us.

**debounceMode parameter values:**

- TriggerAfterStable: This mode rejects glitches and only passes state transitions after a specified period of stability (the debounce time). This mode is used with electro-mechanical devices like encoders and mechanical switches to reject switch bounce and disturbances due to a vibrating encoder that is not otherwise moving. The debounce time should be set short enough to accept the desired input pulse but longer than the period of the undesired disturbance.
- CtrTriggerBeforeStable: Use this mode when the input signal has groups of glitches and each group is to be counted as one. The trigger before stable mode will recognize and count the first glitch within a group but reject the subsequent glitches within the group if the debounce time is set accordingly. In this case the debounce time should be set to encompass one entire group of glitches.

## CFreqIn()

Measures the frequency of a signal. This method can only be used with 9513 counters. This method uses internal counters #5 and #4.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**

```
Public Function CFreqIn(ByVal signalSource As MccDaq.SignalSource ,
    ByVal gateInterval As Integer, ByRef count As Short, ByRef freq As
    Integer) As MccDaq.ErrorInfo

Public Function CFreqIn(ByVal signalSource As MccDaq.SignalSource,
    ByVal gateInterval As Integer, ByRef count As System.UInt16, ByRef
    freq As Integer) As MccDaq.ErrorInfo
```

**C# .NET:**

```
public MccDaq.ErrorInfo CFreqIn(MccDaq.SignalSource signalSource,
    int gateInterval, out short count, out int freq)

public MccDaq.ErrorInfo CFreqIn(MccDaq.SignalSource signalSource,
    int gateInterval, out ushort count, out int freq)
```

### Parameters:

signalSource	<p>Specifies the source of the signal to calculate the frequency from.</p> <p>The signal to be measured is routed internally from the source specified by signalSource to the clock input of counter 5. On boards with more than one 9513 chip, there is more than one counter 5. Which counter 5 is used is also determined by SigSource. Set it to one of the constants in the "signalSource parameter values" section on page 284.</p> <p>The value of signalSource determines which chip will be used. CtrInput6 through CtrInput10, Freq6 through Freq10 and Gate6 through Gate9 indicate chip two will be used. The signal to be measured must be present at the chip two input specified by SigSource.</p> <p>Note: The gating connection from counter 4 output to counter 5 gate must be made between counters 4 and 5 of <i>this chip</i> (see below). Refer to board-specific information to determine valid values for your board.</p>
gateInterval	<p>Gating interval in milliseconds (must be &gt; 0). Specifies the time, in milliseconds, that the counter will count. The optimum gateInterval depends on the frequency of the measured signal. The counter can count up to 65535. If the gating interval is too low, then the count will be too low and the resolution of the frequency measurement will be poor. For example, if the count changes from 1 to 2 the measured frequency doubles.</p> <p>If the gating interval is too long, the counter will overflow and a FreqOverflow error will occur.</p> <p>This method will not return until the gateInterval has expired. There is no background option. Under Windows, this means that window activity will stop for the duration of the call. Adjust the gateInterval so this does not pose a problem to your user interface.</p>
count	The raw count.
freq	The measured frequency in Hz.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

Count - Count that the frequency calculation is based on.

Freq - Measured frequency in Hz

**signalSource parameter values:**

All of the `signalSource` settings are `MccDaq.SignalSource` enumerated constants. To set a variable to one of these constants, you must refer to the `MccDaq` object and the `SignalSource` enumeration (*variable* = `MccDaq.SignalSource.CtrInput1`, *variable* = `MccDaq.SignalSource.Gate1`, etc.).

**One 9513 chip (Chip 1 used):**

- `CtrInput1` through `CtrInput5`
- `Gate1` through `Gate4`
- `Freq1` through `Freq5`

**Two 9513 chips (Chip 1 or Chip 2 used):**

- `CtrInput1` through `CtrInput10`
- `Gate1` through `Gate9` (excluding `Gate5`)
- `Freq1` through `Freq10`

**Four 9513 chips (Chips 1- 4 may be used):**

- `CtrInput1` through `CtrInput20`
- `Gate1` through `Gate19` (excluding gates 5, 10 & 15)
- `Freq1` through `Freq20`

**Notes:**

- This method requires an electrical connection between counter 4 output and counter 5 gate. This connection must be made between counters 4 and 5 on the chip specified by `signalSource`.
- [C9513Init\(\)](#) must be called for each `chipNum` that will be used by this method. The values of `foutDivider`, `foutSource`, `compare1`, `compare2`, and `timeOfDay` are irrelevant to this method and may be any value shown in the `C9513Init()` method description.
- If you select an external clock source for the counters, the `gateInterval`, `count`, and `freq` settings are only valid if the external source is 1 MHz. Otherwise, you need to scale the values according to the frequency of the external clock source.

For example, for an external clock source of 2 MHz, increase your `gateInterval` setting by a factor of 2, and also double the `count` and `freq` values returned when analyzing your results.

## CIn()

Reads the current count from a counter.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**

```
Public Function CIn(ByVal counterNum As Integer, ByRef count As Short) As MccDaq.ErrorInfo

Public Function CIn(ByVal counterNum As Integer, ByRef count As System.UInt16) As MccDaq.ErrorInfo
```

**C# .NET:**

```
public MccDaq.ErrorInfo CIn(int counterNum, out ushort count)

public MccDaq.ErrorInfo CIn(int counterNum, out short count)
```

### Parameters:

counterNum	The counter to read the current count from. Valid values are 1 to 20, up to the number of counters on the board.
count	Counter value returned here.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

### Notes:

count: Refer to your BASIC manual for information on BASIC integer data types. -32,768 to 32,767 for BASIC languages. BASIC reads counters as:

- -1 reads as 65535
- -21768 reads as 32768
- 32767 reads as 32767
- 2 reads as 2
- 0 reads as 0

**CIn () vs. CIn32 ():** Although the CIn () and [CIn32 \(\)](#) methods perform the same operation, CIn32 () is the preferred method to use.

The only difference between the two is that CIn () returns a 16-bit count value and CIn32 () returns a 32-bit value. Both CIn () and CIn32 () can be used, but CIn32 () is required whenever you need to read count values greater than 16-bits (counts > 65535).

## CIn32()

Reads the current count from a counter, and returns it as a 32 bit integer.

Member of the [MccBoard](#) class.

### Function prototype:

```
VB .NET:      Public Function CIn32(ByVal counterNum As Integer, ByRef count As
               Integer) As MccDaq.ErrorInfo

               Public Function CIn32(ByVal counterNum As Integer, ByRef count As
               System.UInt32) As MccDaq.ErrorInfo

C# .NET:      public MccDaq.ErrorInfo CIn32(int counterNum, out uint count)
               public MccDaq.ErrorInfo CIn32(int counterNum, out int count)
```

### Parameters:

counterNum	The counter to read the current count from. Valid values are 1 to $n$ , where $n$ is the number of counters on the board.
count	Current count value from selected counter.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

### Notes:

**CIn() vs. CIn32():** Although the [CIn\(\)](#) and CIn32 () methods perform the same operation, CIn32 () is the preferred method to use.

The only difference between the two is that CIn () returns a 16-bit count value and CIn32 () returns a 32-bit value. Both CIn () and CIn32 () can be used, but CIn32 () is required whenever you need to read count values greater than 16-bits (counts > 65535).

## CInScan()

Scans a range of scan counter channels, and stores the samples in an array. This method only works with counter boards that have counter scan capability.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**                   Public Function CInScan(ByVal firstCtr As Integer, ByVal lastCtr As Integer, ByVal numPoints As Integer, ByRef rate As Integer, ByVal memHandle As Integer, ByVal options As MccDaq.ScanOptions) As MccDaq.ErrorInfo

**C# .NET:**                   public MccDaq.ErrorInfo CInScan(int firstCtr, int lastCtr, int numPoints, int rate, int memHandle, MccDaq.ScanOptions options)

### Parameters:

firstCtr	First counter channel of the scan.  This parameter is zero-based, so the first counter number is "0".
lastCtr	Last counter channel of the scan.  This parameter is zero-based, so the first counter number is "0".  The maximum allowable channel for both firstCtr and lastCtr depends on how many scan counters are available on the Measurement Computing device in use.
numPoints	Number of counter samples to collect. Specifies the total number of counter samples that will be collected. If more than one channel is being sampled then the number of samples collected per channel is equal to Count / (firstCtr - lastCtr + 1).
rate	The rate at which samples are taken – the counts are latched and saved in board memory, in samples per second.  Rate also returns the value of the actual rate set, which may be different from the requested rate because of pacer limitations.
memHandle	The handle for the Windows buffer to store data (Windows). This buffer must have been previously allocated with the <a href="#">WinBufAlloc32()</a> method.
options	Bit fields that control various options. Set it to one of the constants in the "options parameter values" section on page 288.

### Returns:

[Error code](#) or 0 if no errors

rate – the actual sampling rate used.

memHandle – the collected counter data returned via the Windows buffer.

**options parameter values:**

All of the options settings are `MccDaq.ScanOptions` enumerated constants. To set a variable to one of these constants, you must refer to the `MccDaq` object and the `ScanOptions` enumeration (*variable* = `MccDaq.ScanOptions.Continuous`, *variable* = `MccDaq.ScanOptions.Background`, etc.).

Background	<p>When the <code>Background</code> option is used, control returns immediately to the next line in your program, and the data collection from the counters into the buffer continues in the background. If the <code>Background</code> option is not used, the <code>CInScan()</code> method does not return to your program until all of the requested data has been collected and returned to the buffer.</p> <p>Use <a href="#">GetStatus()</a> to check on the status of the background operation. Use <a href="#">StopBackground()</a> with <code>CtrFunction</code> to terminate the background process before it has completed. Execute <code>StopBackground()</code> after normal termination of all background functions in order to clear variables and flags.</p>
Continuous	<p>This option puts the function in an endless loop. Once it collects the required number of samples, it resets to the start of the buffer and begins again. The only way to stop this operation is with <code>StopBackground()</code> with <code>CtrFunction</code>. Normally, you should use this option with <code>Background</code> so that your program regains control.</p>
ExtTrigger	<p>If this option is specified, sampling does not begin until the trigger condition is met. You can set the trigger condition to rising edge, falling edge, or the level of the digital trigger input with the <a href="#">SetTrigger()</a> method. Refer to board-specific information in the <i>UL User's Guide</i>.</p>
ExtClock	<p>If this option is specified, conversions will be controlled by the signal on the external clock input rather than by the internal pacer clock. Each conversion will be triggered on the appropriate edge of the clock input signal (refer to board-specific information in the <i>UL User's Guide</i>). When this option is used the <code>rate</code> parameter is ignored. The sampling rate is dependent on the clock signal. Options for the board will default to a transfer mode that will allow the maximum conversion rate to be attained unless otherwise specified.</p>



## CLoad()

Loads the specified counter's Load, Hold, Alarm, QuadCount, QuadPreset or PreScaler register with a count. When loading a counter with a starting value, it is never loaded directly into the counter's count register. Rather, it is loaded into the load or hold register. From there, the counter, after being enabled, loads the count from the appropriate register, generally on the first valid pulse.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**

```
Public Function CLoad(ByVal regNum As MccDaq.CounterRegister, ByVal
loadValue As Integer) As MccDaq.ErrorInfo

Public Function CLoad(ByVal regNum As MccDaq.CounterRegister, ByVal
loadValue As System.UInt32) As MccDaq.ErrorInfo
```

**C# .NET:**

```
public MccDaq.ErrorInfo CLoad(MccDaq.CounterRegister regNum, uint
loadValue)

public MccDaq.ErrorInfo CLoad(MccDaq.CounterRegister regNum, int
loadValue)
```

### Parameters:

regNum	The register to load the count to. Set it to one of the constants in the "regNum parameter values" section below.
loadValue	The value to be loaded. This value must be between 0 and $2^{\text{resolution}} - 1$ of the counter. Refer to the discussion of Basic signed integers in the "16-bit values using a signed integer data type" section in the "Universal Library Description & Use" chapter of the <i>Universal Library User's Guide</i> (available on our web site at <a href="http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf">www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf</a> ).

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

### regNum parameter values:

All of the regNum settings are MccDaq.CounterRegister enumerated constants. To set a variable to one of these constants, you must refer to the MccDaq object and the CounterRegister enumeration (*variable* = MccDaq.CounterRegister.LoadReg1, *variable* = MccDaq.CounterRegister.HoldReg1, etc.).

LoadReg1 ... 20	Load registers 1 to 20. Can span many chips.
HoldReg1 ... 20	Hold registers 1 to 20. Can span several chips. (9513 only)
Alarm1Chip1	Alarm register 1 of the first counter chip. (9513 only)
Alarm2Chip1	Alarm register 2 of the first counter chip. (9513 only)
Alarm1Chip2	Alarm register 1 of the 2nd counter chip. (9513 only)
Alarm2Chip2	Alarm register 2 of the 2nd counter chip. (9513 only)
Alarm1Chip3	Alarm register 1 of the third counter chip. (9513 only)
Alarm2Chip3	Alarm register 2 of the third counter chip. (9513 only)
Alarm1Chip4	Alarm register 1 of the four counter chip. (9513 only)
Alarm2Chip4	Alarm register 2 of the four counter chip. (9513 only)
QuadCount1 to QuadCount4	Current Count (LS7266 only)
QuadPreset1 to QuadPreset4	Preset register (LS7266 only)

QuadPrescaler1 to QuadPrescaler4      Prescaler register (LS7266 only)

**Notes:**

You cannot load a count-down-only counter with less than 2.

**Counter types:** There are several counter types supported. Please refer to the data sheet for the registers available for a counter type.

**CLoad() vs. CLoad32():** The CLoad() and [CLoad32\(\)](#) perform the same operation. These methods differ in that CLoad() loads a 16-bit count value, while CLoad32() loads a 32-bit value. The only time you need to use CLoad32() is to load counts that are larger than 32 bits (counts > 65535).

## CLoad32()

Loads the specified counter's COUNT, PRESET or PRESCALER register with a count.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**

```
Public Function CLoad32(ByVal regNum As MccDaq.CounterRegister ,
    ByVal loadValue As Integer) As MccDaq.ErrorInfo
Public Function CLoad32(ByVal regNum As MccDaq.CounterRegister,
    ByVal loadValue As System.UInt32) As MccDaq.ErrorInfo
```

**C# .NET:**

```
public MccDaq.ErrorInfo CLoad32(MccDaq.CounterRegister regNum, uint
    loadValue)
public MccDaq.ErrorInfo CLoad32(MccDaq.CounterRegister regNum, int
    loadValue)
```

### Parameters:

regNum	The register to load the value into. Set it to one of the constants in the "regNum parameter values" section below.
loadValue	The value to be loaded into regNum.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

### regNum parameter values:

All of the regNum settings are MccDaq.CounterRegister enumerated constants. To set a variable to one of these constants, you must refer to the MccDaq object and the CounterRegister enumeration (variable = MccDaq.CounterRegister.LoadReg1, variable = MccDaq.CounterRegister.HoldReg1, etc.).

LoadReg1 ... 20	Load registers 1 to 20. Can span many chips.
HoldReg1 ... 20	Hold registers 1 to 20. Can span several chips. (9513 only)
Alarm1Chip1	Alarm register 1 of the first counter chip. (9513 only)
Alarm2Chip1	Alarm register 2 of the first counter chip. (9513 only)
Alarm1Chip2	Alarm register 1 of the 2nd counter chip. (9513 only)
Alarm2Chip2	Alarm register 2 of the 2nd counter chip. (9513 only)
Alarm1Chip3	Alarm register 1 of the third counter chip. (9513 only)
Alarm2Chip3	Alarm register 2 of the third counter chip. (9513 only)
Alarm1Chip4	Alarm register 1 of the four counter chip. (9513 only)
Alarm2Chip4	Alarm register 2 of the four counter chip. (9513 only)
QuadCount1 to QuadCount4	Used to initialize the counter
QuadPreset1 to QuadPreset4	Used to set upper limit of counter in some modes.
QuadPrescaler1 to QuadPrescaler4	Used for clock filtering (valid values: 0 to 255).

### Notes:

CLoad() vs. CLoad32(): Although the CLoad() and CLoad32() methods perform the same operation, CLoad32() is the preferred method to use.

The only difference between the two is that CLoad() loads a 16-bit count value, and CLoad32() loads a 32-bit value. The only time you need to use CLoad32() is to load counts that are larger than 32 bits (counts > 65535).

## CStatus()

Returns status information about the specified counter (7266 counters only)

### Function prototype:

**VB .NET:**                   Public Function CStatus(ByVal counterNum As Integer, ByRef statusBits As MccDaq.StatusBits ) As MccDaq.ErrorInfo

**C# .NET:**                   public MccDaq.ErrorInfo CStatus(int counterNum, out MccDaq.StatusBits statusBits)

### Parameters:

counterNum	The number of the counter whose status bits you want to read. Valid values are 1 to <i>n</i> , where <i>n</i> is the number of counters on the board.
statusBits	Current status from selected counter is returned here. The status consists of individual bits that indicate various conditions within the counter. Set it to one of the constants in the "statusBits parameter values" section below.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

All of the statusBits settings are MccDaq.StatusBits enumerated constants. To set a variable to one of these constants, you must refer to the MccDaq object and the StatusBits enumeration (*variable* = MccDaq.StatusBits.UnderFlow, *variable* = MccDaq.StatusBits.Overflow, etc.).

### statusBits parameter values:

Underflow	set to 1 whenever the count decrements past 0. Is cleared to 0 whenever CStatus() is called.
Overflow	Set to 1 whenever the count increments past it's upper limit. Is cleared to 0 whenever CStatus() is called.
Compare	Set to 1 whenever the count matches the preset register. Is cleared to 0 whenever CStatus() is called.
Sign	Set to 1 when the MSB of the count is 1. Is cleared to 0 whenever the MSB of the count is set to 0.
Error	Set to 1 whenever an error occurs due to excessive noise on the input. Is cleared to 0 by calling <a href="#">C7266Config()</a> .
UpDown	Set to 1 when counting up. Is cleared to 0 when counting down
Index	Set to 1 when index is valid. Is cleared to 0 when index is not valid.

## CStoreOnInt()

Installs an interrupt handler that will store the current count whenever an interrupt occurs. This method can only be used with 9513 counters. This method will continue to operate in the background until either `intCount` is satisfied or [StopBackground\(\)](#) with `CtrFunction` is called.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**                    `Public Function CStoreOnInt(ByVal intCount As Integer, ByRef  
cnterControl As MccDaq.CounterControl, ByVal memHandle As Integer) As  
MccDaq.ErrorInfo`

**C# .NET:**                    `public MccDaq.ErrorInfo CStoreOnInt(int intCount, ref  
MccDaq.CounterControl cnterControl, int memHandle)`

### Parameters:

<code>intCount</code>	The counters will be read every time an interrupt occurs, until <code>IntCount</code> number of interrupts have occurred. If <code>intCount = 0</code> , the method will run until <a href="#">StopBackground()</a> is called. (refer to <code>memHandle</code> below).
<code>cnterControl</code>	The array should have an element for each counter on the board. (5 elements for CTR-05 board, 10 elements for a CTR-10, etc.). Each element corresponds to a counter channel. Each element should be set to either <code>MccDaq.CounterControl.Disabled</code> or <code>MccDaq.CounterControl.Enabled</code> .  All channels set to <code>MccDaq.CounterControl.Enabled</code> will be read when an interrupt occurs.
<code>memHandle</code>	Handle for Windows buffer. If <code>intCount</code> is non-zero, the buffer referenced by <code>memHandle</code> must be of sufficient size to hold ( <code>intCount * Number of Counters</code> ) points.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

### Notes:

If the library revision is set to 4.0 or greater, the following code changes are required:

- If `intCount` is non-zero, the buffer referenced by `memHandle` must be able to hold (`intCount * Number of Counters`) points.  
  
For example, if you set `intCount` to 100 for a CTR-05 board, you must allocate the size of the buffer to be  $(100 * 5) = 500$ . This new functionality keeps the user application from having to move the data out of the buffer for every interrupt, before it is overwritten. Now, for each interrupt, the counter values will be stored in adjacent memory locations in the buffer.

#### Allocate the proper buffer size for non-zero `IntCount` settings

Specifying `intCount` as a non-zero value and failing to allocate the proper sized buffer results in a runtime error. There is no way for the Universal Library to determine if the buffer has been allocated with the proper size.

- If `intCount = 0`, the functionality is unchanged.

## TimerOutStart()

Starts a timer square wave output. Use [TimerOutStop\(\)](#) to stop the output. This method only works with counter boards that have a timer-type counter.

Member of the [MccBoard](#) class.

### Function prototype:

```
VB .NET:      Public Function TimerOutStart(ByVal timerNum As Integer, ByRef
frequency As Double) As MccDag.ErrorInfo
```

```
C# .NET:      public MccDag.ErrorInfo TimerOutStart(int timerNum, double
               frequency)
```

### Parameters:

timerNum	The timer to output the square wave from. Valid values are zero up to the number of timers on the board - 1.
----------	--

frequency	The desired square wave frequency. The timers clock will be divided down by integer values to produce the frequency. The actual frequency output will be returned. Valid values are dependent on the timer's clock and the timer resolution.
-----------	--

## Returns:

Error code or 0 if no errors

frequency – the actual frequency set.

## TimerOutStop()

Stops a timer square wave output. Use [TimerOutStart\(\)](#) to start the output. This method only works with counter boards that have a timer-type counter.

Member of the [MccBoard](#) class.

### Function prototype:

VB .NET:                   Public Function TimerOutStop(ByVal timerNum As Integer) As  
                              MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo TimerOutStop(int timerNum)

### Parameters:

timerNum                   The timer to stop. Valid values are zero up to the number of timers on the  
                              board – 1.

### Returns:

[Error code](#) or 0 if no errors

---

## Data Logger Methods and Property

### Introduction

This section covers Universal Library for .NET methods and a file name property used to read and convert data logged to a binary file from a data acquisition product equipped with data logging functionality. The data is typically logged to a CompactFlash® memory card, which may then be inserted into a media reader for reading and conversion using these methods.

Data is stored in a binary file. The data may consist of analog data, CJC temperature data, digital I/O data, time stamped data, and information about the device configuration. You can use the data logger methods and property to read this information, apply conversions to the data, and convert the files to a comma separated values (.CSV) text file or another specified text file format.



## FileName property

Returns the file name associated with the current instance of the DataLogger class.

Member of the [DataLogger](#) class.

### Property prototype:

VB .NET:	Public Shared ReadOnly Property DataLogger As String
C# .NET:	public string FileName [get]

## ConvertFile()

Converts a binary log file to a comma-separated values (.CSV) text file or another text file format that you specify.

Member of the [DataLogger](#) class.

### Function prototype:

**VB .NET:**                   Public Function ConvertFile(ByRef destFileName As String, ByVal startSample As Integer, ByVal count As Integer, ByVal delimiter As MccDaq.FieldDelimiter) As MccDaq.ErrorInfo

**C# .NET:**                   public MccDaq.ErrorInfo ConvertFile(string destFileName, int startSample, int count, MccDaq.FieldDelimiter delimiter)

### Parameters:

destFileName	The name and destination path of the converted file. Use the file extension of the file type that you want to create.
startSample	The first sample to read.
count	The number of samples to read.
delimiter	Specifies the character to use between fields in the converted file.  All of the delimiter settings are MccDaq.FieldDelimiter enumerated constants. Choices are MccDaq.FieldDelimiter.Comma, MccDaq.FieldDelimiter.Semicolon, MccDaq.FieldDelimiter.Space, and MccDaq.FieldDelimiter.Tab.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

### Notes:

- Time stamp data is stored according to the `timeZone` preference and `timeFormat` preference. Refer to [SetPreferences\(\)](#) on page 314.
- Time stamps in the converted file may be in either 12-hour or 24-hour format based on the value of the `timeFormat` preference. Time stamps can optionally be converted to local time based on the value of the `timeZone` preference.
- AI temperature data is returned according to the `Units` preference. Refer to [SetPreferences\(\)](#) on page 314.
- The `units` preference is only applied to the AI data if the data was logged as temperature data. Refer to [GetAllInfo\(\)](#) on page 300. This value is ignored if the AI data was logged as raw data.
- The `units` preference is always applied to CJC data, since it is always logged as temperature data.
- If the `destFileName` argument ends with a .CSV extension, the `delimiter` parameter must be set to `MccDaq.FieldDelimiter.Comma`. Otherwise, an `INVALIDDELIMITER` error is returned.
- You can open a comma-separated values text file (.CSV) directly in Microsoft Excel. Text files with extensions other than .CSV can only be imported into Excel.

## GetAIChannelCount()

Retrieves the total number of analog input channels logged in a binary file.

Member of the [DataLogger](#) class.

### Function prototype:

VB .NET:                   Public Function GetAIChannelCount (ByRef aiCount As Integer) As  
                              MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo GetAIChannelCount (ref int aiCount)

### Parameters:

aiCount                   The number of analog input channels logged in the file.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

aiCount – Returns the number of analog input channels logged in the binary file.

## GetAllInfo()

Retrieves the channel number and unit value of each analog input channel logged in a binary file.

Member of the [DataLogger](#) class.

### Function prototype:

**VB .NET:**                   Public Function GetAllInfo(ByRef channelNumbers As Integer, ByRef units As MccDaq.LoggerUnits, As MccDaq.ErrorInfo

**C# .NET:**                   public MccDaq.ErrorInfo GetAllInfo(ref int channelNumbers, ref MccDaq.LoggerUnits units)

### Parameters:

channelNumbers	An array that contains the analog input channel numbers logged in the file.
units	An array that contains the unit values set by the device in InstaCal for each analog input channel logged in the file.  The units settings are MccDaq.LoggerUnits enumerated constants. Choices are MccDaq.LoggerUnits.Temperature and MccDaq.LoggerUnits.Raw.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

channelNumbers – Returns the analog input channel numbers logged in the binary file.

units – Returns the unit value set by the device in *InstaCal* for each analog input channel logged in the binary file (MccDaq.LoggerUnits.Temperature or MccDaq.LoggerUnits.Raw.)

## GetCJCInfo()

Retrieves the number of CJC temperature channels logged in a binary file.

Member of the [DataLogger](#) class.

### Function prototype:

VB .NET:                   Public Function GetCJCInfo (ByRef cjcCount As Integer) As  
                                  MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo GetCJCInfo(ref int cjcCount)

### Parameters:

cjcCount                   The number of CJC temperature channels logged in the file.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

cjcCount – Returns the number of CJC temperature channels logged in the binary file.

## GetDIOInfo()

Retrieves the number of digital I/O channels logged in a binary file.

Member of the [DataLogger](#) class.

### Function prototype:

VB .NET:                   Public Function GetDIOInfo (ByRef dioCount As Integer) As  
                              MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo GetDIOInfo(ref int dioCount)

### Parameters:

dioCount                   The number of digital I/O channels logged in the file.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

dioCount – Returns the number of digital I/O channels logged in the binary file.

## GetFileInfo()

Gets file information from the file associated with the current instance of the DataLogger.

Member of the [DataLogger](#) class.

### Function prototype:

VB .NET:                   Public Function GetFileInfo(ByRef version As Integer, ByRef size As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo GetFileInfo(ref int version, ref int size)

### Parameters:

version	The version level of the file.
size	The size in bytes of the file.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

version – Returns the version level of the binary file.

size – Returns the size in bytes of the binary file.

## GetFileName()

Retrieves the name and path of the  $n^{\text{th}}$  file in the directory containing binary log files.

Member of the [DataLogger](#) class.

### Function prototype:

**VB .NET:**                   Public Shared Function GetFileName(ByVal fileName As Integer,  
ByRef path As String, ByRef fileName As String) As MccDaq.ErrorInfo

**C# .NET:**                   public static MccDaq.ErrorInfo GetFileName(int fileName, ref  
string path, ref string fileName)

### Parameters:

fileName	Index of the file whose name you want to return. Specify one of the following:  The number ( $n$ ) that represents the location of the file in the directory (where $n = 0, 1, 2,$ and so on), or  <code>MccService.GetFirst</code> – get the first file in the directory, or  <code>MccService.GetNext</code> – get the next file in the directory, based on the current index.  This parameter is the index of the file in the directory, and is not part of the filename.
path	The full path of the directory containing the log files.
fileName	The full path and name of the binary file. The path must be NULL terminated and cannot be longer than 256 characters.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

fileName – Returns the file name and path of the binary file.

### Notes:

To access all of the files in a directory, first call `GetFileName()` with the `fileName` set to `MccService.GetFirst`, then again with the `fileName` set to `MccService.GetNext` until the method returns the error code [NOMOREFILES](#).



## GetPreferences()

Retrieves API preference settings for time stamp data, analog temperature data, and CJC temperature data. Returns the default values unless changed using [SetPreferences\(\)](#).

Member of the [DataLogger](#) class.

### Function prototype:

**VB .NET:**                   Public Shared Function GetPreferences(ByRef timeFormat As MccDaq.TimeFormat, ByRef timeZone As MccDaq.TimeZone, ByRef units As MccDaq.TempScale) As MccDaq.ErrorInfo

**C# .NET:**                   public MccDaq.ErrorInfo GetPreferences(ref MccDaq.TimeFormat timeFormat, ref MccDaq.TimeZone timeZone, ref MccDaq.TempScale units)

### Parameters:

timeFormat	Returns the format used to display time stamp data.  All of the timeFormat settings are MccDaq.TimeFormat enumerated constants. Choices are MccDaq.TimeFormat.12Hour (for example 2:32:51PM) and MccDaq.TimeFormat.TwentyFourHour (for example 14:32:51).
timeZone	Returns the time zone to store time stamp data.  All of the timeZone settings are MccDaq.TimeZone enumerated constants. Choices are MccDaq.TimeZone.Local and MccDaq.TimeZone.GMT.
units	Returns the unit to use for analog temperature data. This value is ignored if raw data values are logged.  All of the units settings are MccDaq.TempScale enumerated constants. Choices are MccDaq.TempScale.Celsius, MccDaq.TempScale.Fahrenheit, and MccDaq.TempScale.Kelvin.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

timeFormat – Returns the format to apply to time stamp data from API functions that return time data.

timeZone – Returns the time zone to apply to time stamp data from API functions that return time data.

units – Returns the unit to use when converting temperature data from API functions that return temperature data.

## GetSampleInfo()

Retrieves the sample interval, sample count, and the date and time of the first data point in a binary file.

Member of the [DataLogger](#) class.

### Function prototype:

**VB .NET:**                   Public Function GetSampleInfo (ByRef sampleInterval As Integer, ByRef sampleCount As Integer, ByRef startDate As Integer, ByRef startTime As Integer) As MccDaq.ErrorInfo

**C# .NET:**                   public MccDaq.ErrorInfo GetFileInfo(ref int sampleInterval, ref int sampleCount, ref int startDate, ref int startTime)

### Parameters:

sampleInterval	The time, in seconds, between samples.
sampleCount	The number of samples contained in the file.
startDate	The date of the first data point logged in the file. Date values are packed in the following format:  Byte 0:            day Byte 1:           month Byte 2 - 3:       year
startTime	The time when the first data point was logged in the file. Time values are packed in the following format:  Byte 0:           seconds Byte 1:           minutes Byte 2:           hours Byte 3:           0xff = 24hour format 0x0 = AM 0x1 = PM

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

sampleInterval – Returns the time, in seconds, between samples.

sampleCount – Returns the number of samples in the file.

startDate – Returns the date of the first data point logged in the file.

startTime – Returns the time when the first data point was logged in the file.

### Notes:

Time stamped data is returned according to the `timeZone` and `timeFormat` preferences. Refer to [SetPreferences\(\)](#) on page 314.

## ReadAIChannels()

Retrieves analog input data from a binary file, and stores the values in an array.

Member of the [DataLogger](#) class.

### Function prototype:

**VB .NET:**                   Public Function ReadAIChannels(ByVal startSampleAs Integer, ByVal count Integer, ByRef aiChannels As Single) As MccDaq.ErrorInfo

**C# .NET:**                   public MccDaq.ErrorInfo ReadAIChannels (int startSample, int count, ref float [] aiChannels)

### Parameters:

**startSample**               The first sample to read from the binary file.

**count**                    The number of samples to read from the binary file.

**aiChannels**               Receives the analog input values.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

**aiChannels** - Returns the analog input values logged in the file.

### Notes:

The units of the analog input data that is returned is set by the value of the **Units** preference. Refer to [SetPreferences\(\)](#) on page 314.


The **Units** preference is only applied if the logged data is temperature data. This value is ignored if the data logged is raw.

### Analog array:

The user is responsible for allocating the size of the analog data array, and ensuring that it is large enough to hold the data that will be returned. You can calculate the array allocation using the **SampleCount** value from [GetSampleInfo\(\)](#), and the **AICount** value from [GetAIChannelCount\(\)](#):

```
float* aiChannels = new float[sampleCount * aiCount];
```

The figure below shows the layout of the analog array, and how the elements should be indexed.

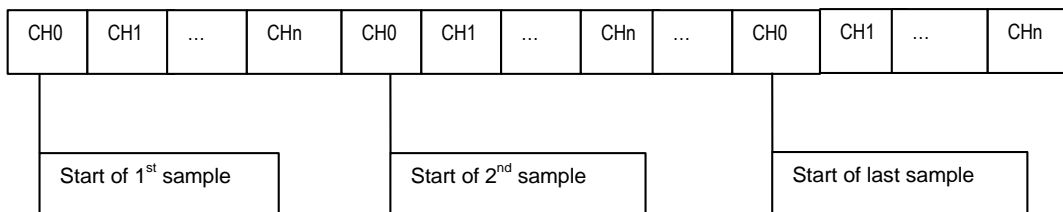

Where n is (numberOfChannels - 1). CH0 - CHn refer to the channels in the array, not the input channels of the device.

For example, assume that all of the even number input channels are logged. The analog array channels are mapped as shown here:

Array Channel	Device Input Channel
0	0
1	2
2	4
3	6

Use the following code fragment to access the elements of the analog array:

```
for (i=0; i<numberOfSamples; i++)
{
    for (j=0; j<numberOfAIChannels; j++)
    {
        a = analogArray[(i *numberOfAIChannels) + j];
    }
}
```

where

the `numberOfSamples` is set by the `sampleCount` value from [GetSampleInfo\(\)](#)

the `numberOfAIChannels` is set by the `aiCount` value from [GetAIChannelCount\(\)](#)

## ReadCJCChannels()

Retrieves CJC temperature data from a binary file, and stores the values in an array.

Member of the [DataLogger](#) class.

### Function prototype:

**VB .NET:**                   Public Function ReadCJCChannels(ByVal startSampleAs Integer, ByVal count Integer, ByRef cjcChannels As Single) As MccDaq.ErrorInfo

**C# .NET:**                   public MccDaq.ErrorInfo ReadCJCChannels(int startSample, int count, ref float [] cjcChannels)

### Parameters:

**startSample**               The first sample to read from the binary file.

**count**                    The number of samples to read from the binary file.

**cjcChannels**              Receives the CJC temperature values.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

**cjcChannels** - Returns the CJC temperature values logged in the file.

### Notes:

The unit of the CJC temperature data that is returned is set by the value of the **Units** preference. Refer to [SetPreferences\(\)](#) on page 314.


The **Units** preference is only valid if the logged data is temperature data. This value is ignored if the data logged is raw.

### CJC array:

The user is responsible for allocating the size of the CJC array, and ensuring that it is large enough to hold the data that will be returned. You can calculate the array allocation using the **SampleCount** value from [GetSampleInfo\(\)](#), and the **cjcCount** value from [GetCJCInfo\(\)](#):

```
float* cjcChannels = new float[SampleCount * CJCCount];
```

The figure below shows the layout of the CJC array, and how the elements should be indexed.

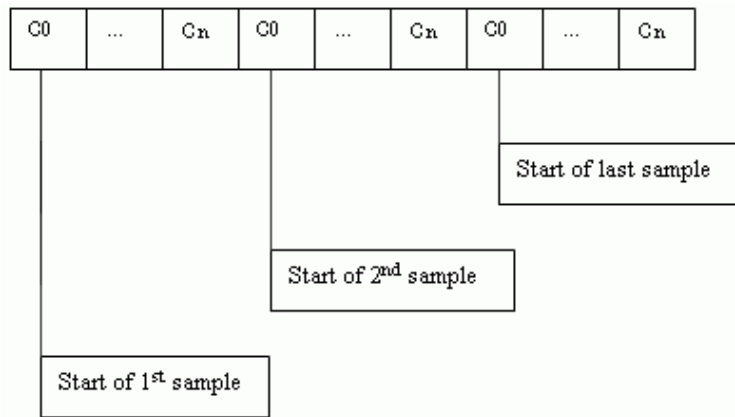

where  $n$  is  $(CJCCount - 1)$

Use the following code fragment to access the elements of the CJC array.

```
for (i=0; i<numberOfSamples; i++)
{
    for (j=0; j<numberOfCJCChannels; j++)
    {
        c = cjcArray[(i * numberOfCJCChannels) + j];
    }
}
```

where

numberOfSamples is set by the `sampleCount` value from `GetSampleInfo()`.

numberOfCJCChannels is set by the `cjcCount` value from `GetCJCInfo()`.



## ReadTimeTags()

Retrieves date and time values logged in a binary file. This method stores date values in the dateTags array, and time values in the timeTags array.

Member of the [DataLogger](#) class.

### Function prototype:

**VB .NET:**                   Public Function ReadTimeTags(ByVal startSample As Integer, ByVal count Integer, ByRef dateTags As Integer, ByRef timeTags As Integer) As MccDaq.ErrorInfo

**C# .NET:**                   public MccDaq.ErrorInfo ReadTimeTags(int startSample, int count, ref int [] dateTags, ref int [] timeTags)

### Parameters:

startSample	The first sample to read from the binary file.
count	The number of samples to read from the binary file.
dateTags	Receives the date tag values. Dates are packed in the following format:  Byte 0:            day Byte 1:           month Byte 2 - 3:       year
timeTags	Receives the time tag values. Times are packed in the following format:  Byte 0:           seconds Byte 1:           minutes Byte 2:           hours Byte 3:           0xff = 24hour format 0x0 = AM 0x1 = PM

### Returns:

- An [ErrorInfo](#) object that indicates the status of the operation.
- dateTags - Returns the date value for each sample logged in the file.
- timeTags - Returns the time value for each sample logged in the file.

### Notes:

Time stamped data is stored according to the `timeZone` preference and the `timeFormat` preference. Refer to [SetPreferences\(\)](#) on page 314.

Time stamped data is logged in the file if *InstaCal* is configured to do so. If time stamps are not logged, the time array is filled with values calculated from the file header information.

### Array size:

The user is responsible for allocating the size of the date and time arrays, and ensuring that they are large enough to hold the data that is returned. You can calculate the array allocation using the `sampleCount` value from [GetSampleInfo\(\)](#) on page 306.

```
int*   dates = new int[sampleCount];
int*   times = new int[sampleCount];
```



**dateTags array**

The figure below shows the layout of the dateTags array, and how the elements should be indexed.

D0	D1	D2	...	Dn
----	----	----	-----	----

where:  $n$  is  $(\text{numberOfSamples} - 1)$

Each sample has only one date. Use the following code fragment to access the elements of the dateTags array:

```
for (i=0; i<numberOfSamples; i++)
{
    d = dateTagsArray[i];
}
```

**timeTags array**

The figure below shows the layout of the timeTags array, and how the elements should be indexed.

T0	T1	T2	...	Tn
----	----	----	-----	----

where:  $n$  is  $(\text{numberOfSamples} - 1)$

Each sample has only one time stamp. Use the following code fragment to access the elements of the timeTags array:

```
for (i=0; i<numberOfSamples; i++)
{
    t = timeTagsArray[i];
}
```

## SetPreferences()

Sets preferences for returned time stamped data, analog temperature data, and CJC temperature data.

Member of the [DataLogger](#) class.

### Function prototype:

**VB .NET:**                   Public Shared Function SetPreferences(ByVal timeFormat As MccDaq.TimeFormat, ByVal timeZone As MccDaq.TimeZone, ByVal units As MccDaq.TempScale) As MccDaq.ErrorInfo

**C# .NET:**                   public MccDaq.ErrorInfo SetPreferences(MccDaq.TimeFormat timeFormat, ref MccDaq.TimeZone timeZone, ref MccDaq.TempScale units)

### Parameters:

timeFormat	<p>Specifies the time format to apply when returning time stamped data (when using <a href="#">ReadTimeTags()</a> for example).</p> <p>All of the timeFormat settings are MccDaq.TimeFormat enumerated constants. Choices are MccDaq.TimeFormat.12Hour (for example 2:32:51) and MccDaq.TimeFormat.TwentyFourHour (for example 14:32:51).</p> <p>timeFormat defaults to MccDaq.TimeFormat.12Hour.</p>
timeZone	<p>Specifies whether to convert time stamped data that is returned (when using <a href="#">ReadTimeTags()</a> for example) to the local time zone or to return the time stamps as they are stored in the file (in the GMT time zone).</p> <p>All of the timeZone settings are MccDaq.TimeZone enumerated constants. Choices are MccDaq.TimeZone.Local and MccDaq.TimeZone.GMT.</p> <p>timeZone defaults to MccDaq.TimeZone.Local.</p>
units	<p>Specifies the unit for analog data. This value is ignored if counts are logged.</p> <p>All of the Units settings are MccDaq.TempScale enumerated constants. Choices are MccDaq.TempScale.Celsius, MccDaq.TempScale.Fahrenheit, and MccDaq.TempScale.Kelvin.</p> <p>units defaults to MccDaq.TempScale.Fahrenheit.</p>

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

### Notes:

- The timeFormat and timeZone preferences are applied to all time data returned using API methods that return time data.
- The units preference specifies the temperature scale that the API applies when reading and converting analog, CJC, and time stamped data.

---

## Digital I/O Methods

### Introduction

Use the methods explained in this chapter to read and set digital values. Most digital ports are configurable, while some others are non-configurable. Some types of hardware allow readback of the values that output ports are set to on configurable port types. Devices using 8255 chips for digital I/O are one example. For these devices, input methods such as `DIn()` are valid for ports configured as output.

## DBitIn()

Reads the state of a single digital input bit. This method treats all of the DIO ports of a particular type on a board as a single port. It lets you read the state of any individual bit within this port. Note that for some port types, such as 8255 ports, if the port is configured for `DigitalOut`, this method provides readback of the last output value.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**                    `Public Function DBitIn(ByVal portType As MccDaq.DigitalPortType, ByVal bitNum As Integer, ByRef bitValue As MccDaq.DigitalLogicState) As MccDaq.ErrorInfo`

**C# .NET:**                    `public MccDaq.ErrorInfo DBitIn(MccDaq.DigitalPortType portType, int bitNum, out MccDaq.DigitalLogicState bitValue)`

### Parameters:

<code>portType</code>	There are three general types of digital ports - ports that are programmable as input or output, ports that are fixed input or output and ports for which each bit may be programmed as input or output. For the first of these types, set <code>portType</code> to <code>FirstPortA</code> . For the latter two types, set <code>portType</code> to <code>AuxPort</code> . Some boards have both types of digital ports (DAS1600). Set <code>portType</code> to either <code>FirstPortA</code> or <code>AuxPort</code> , depending on which digital inputs you wish to read.
<code>bitNum</code>	This specifies the bit number within the single large port.
<code>bitValue</code>	Place holder for return value of bit. Value will be 0 or 1. A 0 indicates a logic low reading, a 1 indicates a logic high reading. Logic high does not necessarily mean 5 V. See the board manual for chip input specifications.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

`BitValue` - value (0 or 1) of specified bit returned here.

## DBitOut()

Sets the state of a single digital output bit. This method treats all of the DIO chips of a particular type on a board as a single very large port. It lets you set the state of any individual bit within this large port. If the port type is not `AuxPort`, you **must** use [DConfigPort\(\)](#) to configure the port for output first. If the port type is `AuxPort`, you **may** need to use [DConfigBit\(\)](#) or [DConfigPort\(\)](#) to configure the bit for output first. Check the board specific information in the *Universal Library User's Guide* (available on our web site at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf)) to determine if `AuxPort` should be configured for your hardware.

Member of the [MccBoard](#) class.

### Function prototype:

VB .NET:           Public Function DBitOut(ByVal portType As MccDaq.DigitalPortType ,  
                                  ByVal bitNum As Integer, ByVal bitValue As MccDaq.DigitalLogicState)  
                                  As MccDaq.ErrorInfo

C# .NET:           public MccDaq.ErrorInfo DBitOut(MccDaq.DigitalPortType portType, int  
                                  bitNum, MccDaq.DigitalLogicState bitValue)

### Parameters:

portType	There are three general types of digital ports - ports that are programmable as input or output, ports that are fixed input or output and ports for which each bit may be programmed as input or output. For the first of these types, set <code>PortType</code> to <code>FirstPortA</code> . For the latter two types, set <code>portType</code> to <code>AuxPort</code> . Some boards have both types of digital ports (DAS1600). Set <code>portType</code> to either <code>FirstPortA</code> or <code>AuxPort</code> depending on which digital port you wish to write to.
bitNum	This specifies the bit number within the single large port. The specified bit must be in a port that is currently configured as an output.
bitValue	The value to set the bit to. Value will be 0 or 1. A 0 indicates a logic low output, a 1 indicates a logic high output. Logic high does not necessarily mean 5V. Refer to the board's user's guide for chip specifications.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

## DConfigBit()

Configures a specific digital bit as Input or Output. This method treats all DIO ports of the `AuxPort` type on a board as a single port. This method is NOT supported by 8255 type DIO ports. Please refer to board specific information for details.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**                    `Public Function DConfigBit(ByVal portNum As MccDaq.DigitalPortType, ByVal bitNum As Integer, ByVal direction As MccDaq.DigitalPortDirection) As MccDaq.ErrorInfo`

**C# .NET:**                    `public MccDaq.ErrorInfo DConfigBit(MccDaq.DigitalPortType portNum, int bitNum, MccDaq.DigitalPortDirection direction)`

### Parameters:

<code>portNum</code>	The port ( <code>AuxPort</code> ) whose bits are to be configured. The port specified must be bitwise configurable. See board specific information for details.
<code>bitNum</code>	The bit number to configure as input or output. See board specific information for details.
<code>direction</code>	<code>MccDaq.DigitalPortDirection DigitalOut</code> or <code>DigitalIn</code> configures the specified bit for output or input, respectively.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

## DConfigPort()

Configures a digital port as input or output. This method is for use with ports that may be programmed as input or output, such as those on the 82C55 chips and 8536 chips. Refer to the board's hardware User Guide for details of chip operation.

Member of the [MccBoard](#) class.

### Function prototype:

VB .NET:                   Public Function DConfigPort(ByVal portNum As MccDaq.DigitalPortType, ByVal direction As MccDaq.DigitalPortDirection) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo DConfigPort(MccDaq.DigitalPortType portNum, MccDaq.DigitalPortDirection direction)

### Parameters:

portNum	The specified port must be configurable. For most boards, AuxPort is not configurable; consult your board-specific documentation.
direction	MccDaq.DigitalPortDirection.DigitalOut or MccDaq.DigitalPortDirection.DigitalIn configures the entire eight-bit or four-bit port for output or input.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

### Notes:

When used on ports within an 8255 chip, this method will reset all ports on that chip configured for output to a zero state. This means that if you set an output value on FirstPortA and then change the configuration on FirstPortB from Output to Input, the output value at FirstPortA will be all zeros. You can, however, set the configuration on SecondPortX without affecting the value at FirstPortA. For this reason, this method is usually called at the beginning of the program for each port requiring configuration.

## DIn()

Reads a digital input port. Note that for some port types, such as 8255 ports, if the port is configured for DigitalOut, this method will provide readback of the last output value.

Member of the [MccBoard](#) class.

### Function prototype:

```
VB .NET:      Public Function DIn(ByVal portNum As MccDaq.DigitalPortType , ByRef
               dataValue As Short) As MccDaq.ErrorInfo

               Public Function DIn(ByVal portNum As MccDaq.DigitalPortType, ByRef
               dataValue As System.UInt16) As MccDaq.ErrorInfo

C# .NET:      public MccDaq.ErrorInfo DIn(MccDaq.DigitalPortType portNum, out
               ushort dataValue)

               public MccDaq.ErrorInfo DIn(MccDaq.DigitalPortType portNum, out
               short dataValue)
```

### Parameters:

portNum	Specifies which digital I/O port to read. Some hardware does allow readback of the state of the output using this method. Check the board-specific information in the <i>Universal Library User's Guide</i> .
dataValue	Digital input value returned here.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

dataValue - Digital input value returned here

### Notes:

The size of the ports vary. If it is an eight bit port, the returned value is in the 0 - 255 range. If it is a four bit port, the value is in the 0 - 15 range.

Refer to the board-specific information contained in the *Universal Library User's Guide* for clarification of valid portNum values (available in PDF format on our website at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf))



## DInScan()

Performs multiple reads of a digital input port of a high speed digital port on a board with a pacer clock - such as the CIO-PDMA16.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**

```
Public Function DInScan(ByVal portNum As MccDaq.DigitalPortType ,
    ByVal numPoints As Integer, ByRef rate As Integer, ByVal memHandle
    As Integer, ByVal options As MccDaq.ScanOptions ) As
    MccDaq.ErrorInfo
```

**C# .NET:**

```
public MccDaq.ErrorInfo DInScan(MccDaq.DigitalPortType portNum, int
    numPoints, ref int rate, int memHandle, MccDaq.ScanOptions options)
```

### Parameters:

portNum	Specifies which digital I/O port to read (usually FirstPortA or FirstPortB). The specified port must be configured as an input.
numPoints	The number of times to read digital input.
rate	Number of times per second (Hz) to read the port. The actual sampling rate in some cases will vary a small amount from the requested rate. The actual rate will be returned to the rate parameter.
memHandle	Handle for Windows buffer to store data. This buffer must have been previously allocated with the <a href="#">WinBufAlloc()</a> or <a href="#">WinBufAlloc32()</a> method.
options	Bit fields that control various options. Set it to one of the constants in the "options parameter values" section below.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

rate - actual sampling rate returned.

memHandle - digital input value returned via allocated Windows buffer.

### options parameter values:

All of the options settings are MccDaq.ScanOptions enumerated constants. To set a variable to one of these constants, you must refer to the MccDaq object and the ScanOptions enumeration (*variable* = MccDaq.ScanOptions.Background, *variable* = MccDaq.ScanOptions.Continuous, etc.).

Background	<p>If the Background option is not used, the <a href="#">DInScan()</a> method will not return to your program until all of the requested data has been collected and returned to memHandle.</p> <p>When the Background option is used, control will return immediately to the next line in your program and the transfer from the digital input port to memHandle will continue in the background. Use <a href="#">GetStatus()</a> with DiFunction to check on the status of the background operation. Use <a href="#">StopBackground()</a> with DiFunction to terminate the background process before it has completed.</p>
Continuous	<p>This option puts the method in an endless loop. Once it transfers the required number of bytes it resets to the start of the buffer and begins again. The only way to stop this operation is by calling StopBackground() with DiFunction. Normally this option should be used in combination with Background so that your program will regain control.</p>

ExtClock	If this option is used then transfers will be controlled by the signal on the trigger input line rather than by the internal pacer clock. Each transfer will be triggered on the appropriate edge of the trigger input signal (refer to board-specific information in the <i>Universal Library User's Guide</i> ). When this option is used, the <code>rate</code> parameter is ignored. The transfer rate is dependent on the trigger signal.
ExtTrigger	If this option is used then the scan will not begin until the signal on the trigger input line meets the trigger criteria.
WordXfer	Normally this method reads a single (byte) port. If <code>WordXfer</code> is specified, it will read two adjacent ports on each read, and store the value of both ports together as the low and high byte of a single array element in the buffer.  When <code>WordXfer</code> is used, it is generally required to set <code>portNum</code> to <code>FirstPortA</code> .

**Notes:**

**Transfer method** - May not be specified. DMA is used.

## DOut()

Writes a byte to a digital output port. If the port type is not `AuxPort`, you **must** use [DConfigPort\(\)](#) to configure the port for output first. If the port type is `AuxPort`, you **may** need to use [DConfigPort\(\)](#) to configure the port for output first. Check the board specific information in the *Universal Library User's Guide* (available on our web site at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf)) to determine if `AuxPort` should be configured for your hardware.

Member of the [MccBoard](#) class.

### Function prototype:

```
VB .NET:      Public Function DOut(ByVal portNum As MccDaq.DigitalPortType, ByVal
               dataValue As Short) As MccDaq.ErrorInfo
               Public Function DOut(ByVal portNum As MccDaq.DigitalPortType, ByVal
               dataValue As System.UInt16) As MccDaq.ErrorInfo

C# .NET:      public MccDaq.ErrorInfo DOut(MccDaq.DigitalPortType portNum, ushort
               dataValue)
               public MccDaq.ErrorInfo DOut(MccDaq.DigitalPortType portNum, short
               dataValue)
```

### Parameters:

portNum	There are three general types of digital ports - ports that are programmable as input or output, ports that are fixed input or output, and ports for which each bit may be programmed as input or output. For the first of these types, set portNum to FirstPortA. For the latter two types, set portNum to AuxPort. Some boards have both types of digital ports (DAS1600). Set portNum to either FirstPortA or AuxPort depending on which digital port you wish to write to.
dataValue	Digital input value to be written.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

### Notes:

The size of the ports vary. If it is an eight bit port, the output value is in the 0 - 255 range. If it is a four bit port, the value is in the 0 - 15 range. Refer to the board-specific information in the *Universal Library User's Guide* for valid portNum values (available in PDF format on our website at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf))

## DOutScan()

Performs multiple writes to a digital output port of a high speed digital port on a board with a pacer clock, such as the CIO-PDMA16 or CIO-PMA32.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**                   Public Function DOutScan(ByVal portNum As MccDaq.DigitalPortType, ByVal count As Integer, ByRef rate As Integer, ByVal memHandle As Integer, ByVal options As MccDaq.ScanOptions ) As MccDaq.ErrorInfo

**C# .NET:**                   public MccDaq.ErrorInfo DOutScan(MccDaq.DigitalPortType portNum, int count, ref int rate, int memHandle, MccDaq.ScanOptions options)

### Parameters:

portNum	Specifies which digital I/O port to write (usually FirstPortA or FirstPortB). The specified port must be configured as an output.
count	The number of times to write digital output.
rate	Number of times per second (Hz) to write to the port. The actual update rate in some cases will vary a small amount from the requested rate. The actual rate will be returned to the rate parameter.
memHandle	Handle for Windows buffer to store data in (Windows). This buffer must have been previously allocated with the <a href="#">WinBufAlloc()</a> method.
options	Bit fields that control various options. Set it to one of the constants in the "options parameter values" section below.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

rate - actual sampling rate returned.

### options parameter values:

All of the options settings are MccDaq.ScanOptions enumerated constants. To set a variable to one of these constants, you must refer to the MccDaq object and the ScanOptions enumeration (*variable* = MccDaq.ScanOptions.Background, *variable* = MccDaq.ScanOptions.Continuous, etc.).

Background	<p>If the Background option is not used, the DOutScan() method will not return control to your program until all of the requested data has been output.</p> <p>When the Background option is used, control will return immediately to the next line in your program and the transfer to the digital output port from memHandle will continue in the background. Use <a href="#">GetStatus()</a> with DoFunction to check on the status of the background operation. Use <a href="#">StopBackground()</a> with DoFunction to terminate the background process before it has completed.</p>
Continuous	<p>This option puts the method in an endless loop. Once it transfers the required number of bytes it resets to the start of the buffer and begins again. The only way to stop this operation is with <a href="#">StopBackground()</a> with DoFunction. Normally this option should be used in combination with Background so that your program will regain control.</p>

ExtClock	If this option is used then transfers will be controlled by the signal on the trigger input line rather than by the internal pacer clock. Each transfer will be triggered on the appropriate edge of the trigger input signal (refer to board-specific information contained in the <i>Universal Library Users Guide</i> ). When this option is used the rate parameter is ignored. The transfer rate is dependent on the trigger signal.
WordXfer	<p>Normally this method writes a single (byte) port. If WordXfer is specified, it will write two adjacent ports as the low and high byte of a single array element in <code>dataBuffer</code>.</p> <p>When WordXfer is used, it is generally required to set <code>portNum</code> to <code>FirstPortA</code>.</p>
NonStreamedIO	<p>When this option is used, you can output non-streamed data to a specific DAC output channel.</p> <p>To load the data output buffer into the device's internal output FIFO, the aggregate size of the data output buffer must be <math>\leq</math> the size of the internal data output FIFO in the device. Once the sample data are transferred or downloaded to the device, the device is responsible for outputting the data. You can't make any changes to the output buffer once the output begins.</p> <p>With NonStreamedIO mode, you do not have to periodically feed output data through the program to the device for the data output to continue. However, the size of the buffer is limited.</p>
ADCClockTrig	Triggers a data output operation when the ADC clock starts.
ADCClock	Paces the data output operation using the ADC clock.

**Notes:**

- `MccDaq.ScanOptions.ByteXfer` is the default option. Make sure you are using an array when your data is arranged in bytes. Use the `MccDaq.ScanOptions.WordXfer` option for word array transfers.
- `NonStreamedIO` can only be used with the number of samples (`count`) set equal to the size of the FIFO or less.
- Transfer method may not be specified. DMA is used.

---

## Error Handling Methods and Properties

### Introduction

Use the methods and properties explained in this chapter to get information from error codes returned by other UL for .NET methods. Most library methods return `ErrorInfo` objects. These objects contain properties that provide information on the status of the method called. The different routines built into the methods for handling errors include stopping the program when an error occurs, and printing error messages versus error codes.

## ErrHandling()

Sets the error handling for all subsequent method calls. Most methods return error codes after each call. In addition, other error handling features are built into the library. This method controls those features. If the Universal Library cannot find the configuration file CB.CFG, it always terminates the program, regardless of the `ErrHandling()` setting.

Member of the [MccService](#) class.

### Function prototype:

**VB .NET:**                   Public Shared Function ErrHandling(ByVal errorReporting As MccDaq.ErrorReporting, ByVal errorHandling As MccDaq.ErrorHandling) As MccDaq.ErrorInfo

**C# .NET:**                   public static MccDaq.ErrorInfo ErrHandling(MccDaq.ErrorReporting errorReporting, MccDaq.ErrorHandling errorHandling)

### Parameters:

<code>errorReporting</code>	This parameter controls when the library will print error messages on the screen. The default is <code>DontPrint</code> . Set it to one of the constants in the "errorReporting parameter values" section below.
<code>errorHandling</code>	This parameter specifies what class of error will cause the program to halt. Set it to one of the constants in the "errorHandling parameter values" section below.

### Returns:

Returns an [ErrorInfo](#) object that always has [ErrorInfo.Value](#) = `NoErrors`.

### errorReporting parameter values:

All of the `errorReporting` settings are `MccDaq.ErrorReporting` enumerated constants. To set a variable to one of these constants, you must refer to the `MccDaq` object and the `ErrorReporting` enumeration (*variable* = `MccDaq.ErrorReporting.DontPrint`, *variable* = `MccDaq.ErrorReporting.PrintWarnings`, etc.).

<code>DontPrint</code>	Errors will not generate a message to the screen. In that case your program must always check the returned error code after each library call to determine if an error occurred.
<code>PrintWarnings</code>	Only warning errors will generate a message to the screen. Your program will have to check for fatal errors.
<code>PrintFatal</code>	Only fatal errors will generate a message to the screen. Your program must check for warning errors.
<code>PrintAll</code>	All errors will generate a message to the screen.

### errorHandling parameter values:

All of the `errorReporting` settings are `MccDaq.ErrorHandling` enumerated constants. To set a variable to one of these constants, you must refer to the `MccDaq` object and the `ErrorHandling` enumeration (*variable* = `MccDaq.ErrorHandling.DontStop`, *variable* = `MccDaq.ErrorHandling.StopFatal`, etc.).

<code>DontStop</code>	The program will always continue executing when an error occurs.
<code>StopFatal</code>	The program will halt if a "fatal" error occurs.
<code>StopAll</code>	Will stop whenever any error occurs. You can check error codes to determine the cause of the error.

**Notes:**

Warnings vs. fatal errors: All errors that can occur are classified as either "warnings" or "fatal."

- Errors that can occur in normal operation in a bug free program (disk is full, too few samples before trigger occurred) are classified as "warnings."
- All other errors indicate a more serious problem and are classified as "fatal."



## Message property

Use the `ErrorInfo.Message` property to get the error message associated with an `ErrorInfo` object. Most UL for .NET methods return an `ErrorInfo` object. If no error occurred, an `ErrorInfo` object is returned with the `Message` property set to "No error has occurred".

Member of the [ErrorInfo](#) class.

### Property prototype:

```
VB .NET:      Public ReadOnly Property Message As String
C# .NET:      public string Message [get]
```

### Notes:

Refer to the [ErrHandling\(\)](#) method for an alternate method of handling errors.

### Value property

Use the `ErrorInfo.Value` property to get the error constant associated with an `ErrorInfo` object. Most UL for .NET methods return an `ErrorInfo` object. If an error occurs, an `ErrorInfo` object is returned with a non-zero value in the `Value` property.

Member of the [ErrorInfo](#) class.

### Property prototype:

```
VB .NET:      Public ReadOnly Property Value As MccDaq.ErrorInfo.ErrorCode
C# .NET:      public MccDaq.ErrorInfo.ErrorCode Value [get]
```

### Notes:

Refer to the [ErrHandling\(\)](#) method for an alternate method of handling errors.

---

## Memory Board Methods

Use the methods explained in this chapter to read and write data to and from a memory board, and also set modes that control memory boards (MEGA-FIFO).

The most common use for the memory boards is to store large amounts of data from an A/D board via a DT-Connect cable to a memory board. To do this, use the `ExtMemory` option with [AInScan\(\)](#) or [APretrig\(\)](#).

Once the data is transferred to the memory board, you can use the memory methods to retrieve it.

## MemRead()

Reads data from a memory board into an array. Member of the [MccBoard](#) class.

### Function prototype:

```
VB .NET:      Public Function MemRead(ByRef dataBuffer As Short, ByVal firstPoint
               As Integer, ByVal numPoints As Integer) As MccDaq.ErrorInfo

               Public Function MemRead(ByRef dataBuffer As System.UInt16, ByVal
               firstPoint As Integer, ByVal numPoints As Integer) As
               MccDaq.ErrorInfo

C# .NET:      public MccDaq.ErrorInfo MemRead(out short dataBuffer, int
               firstPoint, int numPoints)

               public MccDaq.ErrorInfo MemRead(out ushort dataBuffer, int
               firstPoint, int numPoints)
```

### Parameters:

dataBuffer	Reference to the data array.
firstPoint	Index of first point to read, or FromHere. Use the firstPoint parameter to specify the first point to be read. For example, to read data sample numbers 200 through 250, set firstPoint = 200 and numPoints = 50.
numPoints	Number of data points (words) to read.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

dataBuffer - data read from the memory board.

### Notes:

If you are going to read a large amount of data from the board in small chunks, set firstPoint to FromHere to read each successive chunk. Using FromHere speeds up the operation of MemRead() when working with large amounts of data.

For example, to read 300,000 points in 100,000 point chunks, the calls would look like this:

```
DaqBoard0.MemRead(DataBuffer, 0, 100000)
DaqBoard0.MemRead(DataBuffer, FROMHERE, 1000000)
DaqBoard0.MemRead(DataBuffer, FROMHERE, 1000000)
```

**DT-Connect Conflicts** - The [MemRead\(\)](#) method can not be called while a DT-Connect transfer is in progress. For example, if you start collecting A/D data to the memory board in the background (by calling [AInScan\(\)](#) with the DTConnect + Background options) you cannot call MemRead() until the AInScan() has completed. If you do you will get a DtActive error.

## MemReadPretrig()

Reads pre-trigger data from a memory board that has been collected with the [APretrig\(\)](#) method and re-arranges the data in the correct order (pre-trigger data first, then post-trigger data). This method can only be used to retrieve data that has been collected with the [APretrig\(\)](#) method with `ExtMemory` set in the options parameter. After each [APretrig\(\)](#) call, all data must be unloaded from the memory board with this method. If any more data is sent to the memory board then the pre-trigger data will be lost.

Member of the [MccBoard](#) class.

### Function prototype:

```
VB .NET:      Public Function MemReadPretrig(ByRef dataBuffer As Short, ByVal
               firstPoint As Integer, ByVal numPoints As Integer) As
               MccDaq.ErrorInfo

               Public Function MemReadPretrig(ByRef dataBuffer As System.UInt16,
               ByVal firstPoint As Integer, ByVal numPoints As Integer) As
               MccDaq.ErrorInfo

C# .NET:      public MccDaq.ErrorInfo MemReadPretrig(out short dataBuffer, int
               firstPoint, int numPoints)

               public MccDaq.ErrorInfo MemReadPretrig(out ushort dataBuffer, int
               firstPoint, int numPoints)
```

### Parameters:

<code>dataBuffer</code>	Reference to the data array
<code>firstPoint</code>	Index of first point to read or <code>FromHere</code> . Use the <code>FirstPoint</code> parameter to specify the first point to be read. For example, to read data sample numbers 200 through 250, set <code>firstPoint = 200</code> and <code>numPoints = 50</code> .
<code>numPoints</code>	Number of data samples (words) to read.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

`dataBuffer` - data read from memory board

### Notes:

If you are going to read a large amount of data from the board in small chunks, set `FirstPoint` to `FromHere` to read each successive chunk. Using `FromHere` speeds up the operation of `MemRead()` when working with large amounts of data.

For example, to read 300,000 points in 100,000 chunks, the calls would look like this:

```
DaqBoard0.MemReadPretrig(DataBuffer, 0, 100000)
DaqBoard0.MemReadPretrig(DataBuffer, FROMHERE, 1000000)
DaqBoard0.MemReadPretrig(DataBuffer, FROMHERE, 1000000)
```

**DT-Connect Conflicts** - The [MemReadPretrig\(\)](#) method cannot be called while a DT-Connect transfer is in progress. For example, if you start collecting A/D data to the memory board in the background (by calling [AInScan\(\)](#) with the `DTConnect + Background` options) you cannot call `MemReadPretrig()` until the `AInScan()` has completed. If you do you will get a `DTACTIVE` error.

## MemReset()

Resets the memory board reference to the start of the data. The memory boards are sequential devices. They contain a counter which points to the 'current' word in memory. Every time a word is read or written this counter increments to the next word.

Member of the [MccBoard](#) class.

### Function prototype:

```
VB .NET:      Public Function MemReset() As MccDaq.ErrorInfo
C# .NET:      public MccDaq.ErrorInfo MemReset()
```

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

### Notes:

This method is used to reset the counter back to the start of the memory. Between successive calls to [AInScan\(\)](#), you should call this method so that the second `AInScan()` overwrites the data from the first call. Otherwise, the data from the first `AInScan()` will be followed by the data from the second `AInScan()` in the memory on the card.

Likewise, anytime you call [MemRead\(\)](#) or [MemWrite\(\)](#), it will leave the counter pointing to the next memory location after the data that you read or wrote. Call `MemReset()` to reset back to the start of the memory buffer before the next call to `AInScan()`.

## MemSetDTMode()

Sets the DT-Connect Mode of a memory board.

Member of the [MccBoard](#) class.

### Function prototype:

VB .NET:                   Public Function MemSetDTMode(ByVal mode As MccDaq.DTMode ) As  
                                  MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo MemSetDTMode(MccDaq.DTMode mode)

### Parameters:

mode                       Must be set to either DTIn or DTOut. Set the mode on the memory board to DTIn to transfer data from an A/D board to the memory board. Set mode = DTOut to transfer data from a memory board to a D/A board.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

### Notes:

This command only controls the direction of data transfer between the memory board and its parent board that is connected to it via a DT-Connect cable.

If using the ExtMemory option for AInScan(), etc., *this method should not be used*. The memory board mode is already set through the ExtMemory option.

Use this method only if the parent board is not supported by the Universal Library.

## MemWrite()

Writes data from an array to the memory card.

Member of the [MccBoard](#) class.

### Function prototype:

```
VB .NET:      Public Function MemWrite(ByRef dataBuffer As Short, ByVal firstPoint
               As Integer, ByVal numPoints As Integer) As MccDaq.ErrorInfo

               Public Function MemWrite(ByRef dataBuffer As System.UInt16, ByVal
               firstPoint As Integer, ByVal numPoints As Integer) As
               MccDaq.ErrorInfo

C# .NET:      public MccDaq.ErrorInfo MemWrite(ref short dataBuffer, int
               firstPoint, int numPoints)

               public MccDaq.ErrorInfo MemWrite(ref ushort dataBuffer, int
               firstPoint, int numPoints)
```

### Parameters:

dataBuffer	Reference to the data array.
firstPoint	Index of first point to write or FromHere. Use the firstPoint parameter to specify where in the board's memory to write the first point. For example, to write to location numbers 200 through 250, set firstPoint = 200 and numPoints = 50.
numPoints	Number of data points (words) to write.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

### Notes:

To write large amounts of data to the board in small chunks, set firstPoint to FromHere to write each successive chunk. Using FromHere speeds up the operation of MemWrite() when working with large amounts of data.

For example, to write 300,000 points in 100,000 point chunks, the calls would look like this:

```
DaqBoard1.MemWrite(DataBuffer, 0, 100000)
DaqBoard1.MemWrite(DataBuffer, FROMHERE, 100000)
DaqBoard1.MemWrite(DataBuffer, FROMHERE, 100000)
```

**DT-Connect Conflicts** - The MemWrite() method cannot be called while a DT-Connect transfer is in progress. For example, if you start collecting A/D data to the memory board in the background (by calling [AInScan\(\)](#) with the DTCONNECT + BACKGROUND options). You cannot call MemWrite() until the AInScan() has completed. If you do, you will get a DTACTIVE error.

---

## Revision Control Methods

### Introduction

Use the methods and properties explained in this chapter to initialize the Universal Library DLL so that the underlying functions are interpreted according to the format of the revision you wrote and compiled your program in.

As new revisions of the library are released, bugs from previous revisions are fixed and occasionally new methods are added. It is our goal to preserve existing programs you have written and therefore to never change the order or number of arguments in a method. However, sometimes it is not possible to achieve this goal.



## DeclareRevision()

Initializes the Universal Library with the revision number of the library used to write your program. Must be the first Universal Library for .NET method to be called by your program.

Member of the [MccService class](#).

### Function prototype:

```
VB .NET:      Public Shared Function DeclareRevision(ByRef revNum As Single) As
                MccDag.ErrorInfo
```

C# .NET: `public static MccDaq.ErrorInfo DeclareRevision(ref float revNum)`

### Parameters:

revNum	Revision number of the Universal Library used to interpret method parameters.
--------	---

## Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

**Notes:**

**Default:** Any program using the 32-bit library and not containing this line of code will be defaulted to revision 5.4 parameter assignments.

As new revisions of the library are released, bugs from previous revisions are fixed and occasionally new properties and methods are added. It is Measurement Computing's goal to preserve existing programs you have written and therefore to never change the order or number of parameters in a method.

With the `DeclareRevision()` method, programs do not have to be rewritten in each line where new functions are used, and the program then recompiled. The revision control method initializes the DLL so that the functions are interpreted according to the format of the revision that you wrote and compiled your program in. The method works by interpreting the UL function call from your program and filling in any arguments needed to run with the new revision.

If your program has declared you are running code written for an earlier revision and you call a new method, you must rewrite your program to include the new parameter, and declare the current revision in `DeclareRevision()`.

## GetRevision()

Gets the revision level of Universal Library DLL and the VXD.

Member of the [MccService class](#).

**Function prototype:**

```
VB.NET:      Public Shared Function GetRevision(ByRef revNum As Single, ByRef
vxdRevNum As Single) As MccDag.ErrorInfo
```

```
C# .NET:      public static MccDag.ErrorInfo GetRevision(out float revNum, out
               float vxdRevNum)
```

### Parameters:

revNum Place holder for the revision number of Library DLL.

vxdRevNum Place holder for the revision number of Library VXD.

## Returns:

revNum - Revision number of the Library DLL

**vxdRevNum** - Revision number of the Library VXD

An [ErrorInfo](#) object that indicates if the revision levels of VXD and DLL are incompatible.

---

## Streamer File Methods

### Introduction

Use the streamer file methods explained in the chapter to create, fill, and read streamer files.

## FileAInScan()

Scans a range of A/D channels and stores the samples in a disk file. `FileAInScan()` reads the specified number of A/D samples at the specified sampling rate from the specified range of A/D channels from the specified board. If the A/D board has programmable gain, it sets the gain to the specified range. The collected data is returned to a file in binary format. Use [FileRead\(\)](#) to load data from that file into an array. See board specific information to determine if this method is supported on your board.

Member of the [MccBoard class](#).

### Function prototype:

```
VB .NET:      Public Function FileAInScan(ByVal lowChan As Integer, ByVal highChan
               As Integer, ByVal numPoints As Integer, ByRef rate As Integer, ByVal
               range As MccDaq.Range , ByVal fileName As String, ByVal options As
               MccDaq.ScanOptions) As MccDaq.ErrorInfo

C# .NET:      public MccDaq.ErrorInfo FileAInScan(int lowChan, int highChan, int
               numPoints, ref int rate, MccDaq.Range range, string fileName,
               MccDaq.ScanOptions options)
```

### Parameters:

<code>lowChan</code>	First A/D channel of scan.
<code>highChan</code>	Last A/D channel of scan.
<code>numPoints</code>	The maximum allowable channel depends on which type of A/D board is being used. For boards with both single ended and differential inputs, the maximum allowable channel number also depends on how the board is configured (for example, eight channels for differential, 16 for single ended).
<code>rate</code>	Specifies the total number of A/D samples that will be collected. If more than one channel is being sampled, the number of samples collected per channel is equal to $\text{Count} / (\text{HighChan} - \text{LowChan} + 1)$ .
<code>range</code>	Sample rate in samples per second (Hz) per channel. The maximum sampling rate depends on the A/D board that is being used (refer to the <code>rate</code> description in <a href="#">AInScan()</a> ).
<code>filename</code>	If the selected A/D board does not have a programmable range feature, this parameter is ignored. Otherwise set the <code>range</code> parameter to any range that is supported by the selected A/D board. Refer to Table 7 on page 208 for a list of valid range settings. Refer to board-specific information in the <i>Universal Library User's Guide</i> for a list of the supported A/D ranges of each board.
<code>options</code>	The name of the file in which to store the data. If the file doesn't exist, it will be created.
<code>options</code>	Bit fields that control various options. Set it to one of the constants in the "options parameter values" section below.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

`rate` = actual sampling rate

### options parameter values:

All of the `options` settings are `MccDaq.ScanOptions` enumerated constants. To set a variable to one of these constants, you must refer to the `MccDaq` object and the `ScanOptions` enumeration (*variable* = `MccDaq.ScanOptions.ExtClock`, *variable* = `MccDaq.ScanOptions.ExtTrigger`, etc.).

ExtClock	If this option is used, conversions are controlled by the signal on the trigger input line rather than by the internal pacer clock. Each conversion is triggered on the appropriate edge of the trigger input signal (see board specific info). Additionally, the <code>rate</code> parameter is ignored. The sampling rate is dependent on the trigger signal.
ExtTrigger	<p>If this option is specified, the sampling does not begin until the trigger condition is met.</p> <p>On many boards, this trigger condition is programmable (see <a href="#">SetTrigger()</a> method and board specific info for details) and can be programmed for rising or falling edge or an analog level.</p> <p>On other boards, only "polled gate" triggering is supported. Assuming active high operation, data acquisition commences immediately if the trigger input is high. If the trigger input is low, acquisition is held off until it goes high. Acquisition continues until <code>numPoints</code> samples are taken, regardless of the state of the trigger input. For 'polled gate' triggering, this option is most useful if the signal is a pulse with a very low duty cycle (trigger signal in TTL low state most of the time) to hold off triggering until the pulse occurs.</p>
DtConnect	Samples are sent to the DT-Connect port if the board is equipped with one.

#### Notes:

##### Important

In order to understand the functions, you must read the board-specific information contained in the *Universal Library User's Guide* (available on our web site at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf)). Review and run the example programs before attempting to program yourself. Following this advice will save you hours of frustration, and possibly time wasted holding for technical support. This note, which appears elsewhere, is especially applicable to this method. Read the board-specific information for your board. We suggest that you make a copy of that page to refer to as you read this manual and examine the example programs.

**OVERRUN Error** - (Error code 29): This error indicates that the data was not written to the file as fast as the data was sampled. Consequently some data was lost. The value returned from [FileGetInfo\(\)](#) in `*TotalCount` is the number of points that were successfully collected.

## FileGetInfo()

This method returns information about a streamer file. When [FileAInScan\(\)](#) or [FilePretrig\(\)](#) fills the streamer file, information is stored about how the data was collected (sample rate, channels sampled etc.). This method returns that information. Refer to board-specific information in the *Universal Library User's Guide* to determine if your board supports `FileAInScan()` and/or `FilePretrig()`.

Member of the [MccService](#) class.

### Function prototype:

**VB .NET:**

```
Public Shared Function FileGetInfo(ByVal fileName As String, ByRef
lowChan As Short, ByRef highChan As Short, ByRef pretrigCount As
Integer, ByRef totalCount As Integer, ByRef rate As Integer, ByRef
range As MccDaq.Range ) As MccDaq.ErrorInfo
```

**C# .NET:**

```
public static MccDaq.ErrorInfo FileGetInfo(string fileName, out
short lowChan, out short highChan, out int pretrigCount, out int
totalCount, out int rate, out MccDaq.Range range)
```

### Parameters:

fileName	Name of streamer file.
lowChan	Variable to return lowChan to.
highChan	Variable to return highChan to.
pretrigCount	Variable to return pretrigCount to.
totalCount	Variable to return totalCount to.
rate	Variable to return sampling rate to.
range	Variable to return A/D range code to. Refer to Table 7 on page 208 for a list of valid range settings.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

lowChan - low A/D channel of scan

highChan - high A/D channel of scan

totalCount - total number of points collected

pretrigCount - number of pre-trigger points collected

rate - sampling rate when data was collected

range - Range of A/D when data was collected

## FilePretrig()

Scan a range of channels continuously while waiting for a trigger.

Once the trigger occurs, `FilePretrig()` returns the specified number of samples, including the specified number of pre-trigger samples to a disk file. This method waits for a trigger signal to occur on the Trigger Input. Once the trigger occurs, it returns the specified number (`TotalCount`) of A/D samples, including the specified number of pre-trigger points. It collects the data at the specified sampling rate (`rate`) from the specified range (`lowChan-highChan`) of A/D channels from the specified board. If the A/D board has programmable gain then it sets the gain to the specified range. The collected data is returned to a file. See board specific info to determine if this method is supported by your board.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**                    `Public Function FilePretrig(ByVal lowChan As Integer, ByVal highChan As Integer, ByRef pretrigCount As Integer, ByRef totalCount As Integer, ByRef rate As Integer, ByVal range As MccDaq.Range , ByVal fileName As String, ByVal options As MccDaq.ScanOptions ) As MccDaq.ErrorInfo`

**C# .NET:**                    `public MccDaq.ErrorInfo FilePretrig(int lowChan, int highChan, ref int pretrigCount, ref int totalCount, ref int rate, MccDaq.Range range, string fileName, MccDaq.ScanOptions options )`

### Parameters:

<code>lowChan</code>	First A/D channel of scan
<code>highChan</code>	Last A/D channel of scan
	The maximum allowable channel depends on which type of A/D board is being used. For boards that have both single ended and differential inputs the maximum allowable channel number also depends on how the board is configured. Refer to board-specific information for the maximum number of channels allowed in differential and single ended modes.
<code>pretrigCount</code>	Specifies the number of samples before the trigger that will be returned. <code>PretrigCount</code> must be less than 16000, and <code>PretrigCount</code> must also be less than <code>TotalCount - 512</code> . If the trigger occurs too early, then fewer than the requested number of pre-trigger samples will be collected. In that case a <code>TooFew</code> error will occur. The <code>PretrigCount</code> will be set to indicate how many samples were collected and the post trigger samples will still be collected.
<code>totalCount</code>	Sets the total number of samples to be collected and stored in the file. <code>TotalCount</code> must be greater than or equal to <code>PretrigCount + 512</code> . If the trigger occurs too early, fewer than the requested number of samples will be collected and a <code>TooFew</code> error will occur. The <code>TotalCount</code> will be set to indicate how many samples were actually collected.
<code>rate</code>	Sample rate in samples per second (Hz) per channel. The maximum sampling rate depends on the A/D board that is being used. This is the rate at which scans are triggered.

	<p>If you are sampling 4 channels, 0 - 3, then specifying a rate of 10,000 scans per second (10 kHz) will result in the A/D converter rate of 40 kHz: 4 channels at 10,000 samples per channel per second. This is different from some software, where you specify the total A/D chip rate. In those systems, the per channel rate is equal to the A/D rate divided by the number of channels in a scan. This parameter also returns the value of the actual set. This may be different from the requested rate because of pacer limitations.</p>
range	<p>If the selected A/D board does not have a programmable range feature, this parameter is ignored. Otherwise, set the <code>range</code> parameter to any range that is supported by the selected A/D board. Refer to Table 7 on page 208 for a list of valid range settings. Refer to board-specific information in the <i>Universal Library User's Guide</i> for a list of the supported A/D ranges of each board.</p>
filename	<p>The name of the file in which to store the data. If the file doesn't exist, it will be created.</p>
options	<p>Bit fields that control various options. Set it to one of the constants in the "options parameter values" section below.</p>

#### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

`preTrigCount` - actual number of pre-trigger samples collected

`totalCount` - actual number of samples collected

`rate` = actual sampling rate

#### options parameter values:

All of the `options` settings are `MccDaq.ScanOptions` enumerated constants. To set a variable to one of these constants, you must refer to the `MccDaq` object and the `ScanOptions` enumeration (`variable = MccDaq.ScanOptions.ExtClock` or `variable = MccDaq.ScanOptions.DtConnect`).

ExtClock	<p>If this option is used then conversions will be controlled by the signal on the trigger input line rather than by the internal pacer clock. Each conversion will be triggered on the appropriate edge of the trigger input signal (see board specific info). When this option is used the <code>rate</code> parameter is ignored. The sampling rate is dependent on the trigger signal.</p>
DtConnect	<p>Samples are sent to the DT-Connect port if the board is equipped with one.</p>

#### Notes:

[OVERRUN Error](#) - (Error code 29): This error indicates that the data was not written to the file as fast as the data was sampled. Consequently some data was lost. The value in `TotalCount` will be the number of points that were successfully collected.



## FileRead()

This method reads data from a streamer file, and returns the data in a one-dimensional or two-dimensional array. When [FileAInScan\(\)](#) or [FilePreTrig\(\)](#) fills the streamer file, this method returns the content of that file. Refer to information on your board in the *Universal Library User's Guide* to determine if your board supports [FileAInScan\(\)](#) and/or [FilePreTrig\(\)](#).

Member of the [MccService](#) class.

### Function prototype:

VB .NET:

Returns a one-dimensional array of short values:

```
Public Shared Function FileRead(ByVal fileName As String, ByVal
firstPoint As Integer, ByRef numPoints As Integer, ByRef dataBuffer
As Short) As MccDaq.ErrorInfo
```

Returns a one-dimensional array of System.UInt16 values:

```
Public Shared Function FileRead(ByVal fileName As String, ByVal
firstPoint As Integer, ByRef numPoints As Integer, ByRef dataBuffer
As System.UInt16) As MccDaq.ErrorInfo
```

Returns a two-dimensional array of double values:

```
Public Shared Function FileRead(ByVal fileName As String, ByVal
firstPoint As Integer, ByRef numPoints As Integer, ByRef dataBuffer
As Double(), ByVal numChannels As Integer) As MccDaq.ErrorInfo
```

C# .NET:

Returns a one-dimensional array of short values:

```
public static MccDaq.ErrorInfo FileRead(string fileName, int
firstPoint, ref int numPoints, out short dataBuffer)
```

Returns a one-dimensional array of System.UInt16 values:

```
public static MccDaq.ErrorInfo FileRead(string fileName, int
firstPoint, ref int numPoints, out ushort dataBuffer)
```

Returns a two-dimensional array of double values:

```
public static MccDaq.ErrorInfo FileRead(string fileName, int
firstPoint, ref int numPoints, out double dataBuffer[,], int
numChannels)
```

### Parameters:

filename	Name of the streamer file.
firstPoint	Index of the first point to read.
numPoints	The number of points to read from the file.
dataBuffer	Reference to the array in the data buffer that data is read into.
numChannels	The number of channels to read into dataBuffer.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

dataBuffer - data read from a file.

numPoints - number of points actually read.

numPoints may be less than the requested number of points if an error occurs.

**Notes:**

**Data format:** The data is returned as 16-bits. The 16-bits may represent 12 bits of analog, 12-bits of analog plus 4 bits of channel, or 16-bits of analog.

**Loading portions of files:** The file may contain much more data than can fit in `dataBuffer`. In those cases, use `numPoints` and `firstPoint` to read a selected piece of the file into `dataBuffer`. Call [FileGetInfo\(\)](#) first to find out how many points are in the file.

---

## Synchronous I/O Methods

### Introduction

Use the methods discussed in this chapter to synchronously acquire and output analog, digital, counter, and temperature data. These functions can be used with hardware equipped with synchronous input and output capability.

## DaqInScan()

Scans analog, digital, counter, and temperature input channels synchronously, and stores the samples in an array. This method only works with boards that support synchronous input.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**

```
Public Function DaqInScan(ByVal chanArray As Short( ), ByVal
chanTypeArray As MccDaq.ChannelType(), ByVal gainArray As
MccDaq.Range(), ByVal chanCount As Integer, ByRef rate As Integer,
ByRef pretrigCount As Integer, ByRef totalCount As Integer, ByVal
memHandle As Integer, ByVal options MccDaq.ScanOptions) As
MccDaq.ErrorInfo
```

**Visual Basic:**

```
public ErrorInfo DaqInScan(short[] chanArray, MccDaq.ChannelType[]
chanTypeArray, MccDaq.Range[] gainArray, int chanCount, ref int
rate, ref int pretrigCount, ref int totalCount, int memHandle,
MccDaq.ScanOptions options)
```

### Parameters:

chanArray	Array containing channel values. Valid channel values are analog input channels, digital ports, counter input channels, and temperature input channels of the device.
chanTypeArray	<p>Array containing channel types. Each element of this array defines the type of the corresponding element in the chanArray.</p> <p>All of the chanTypeArray settings are Mccdaq.ChannelType enumerated constants. Set it to one of the constants in the "chanTypeArray parameter values" section on page 349.</p>
gainArray	<p>Array containing A/D range codes. If the corresponding element in the chanArray is not an analog input channel, the range code for this channel is ignored.</p> <p>All of the gainArray settings are MccDaq.Range enumerated constants. Set to any range that is supported by the selected A/D board. Refer to board-specific information in the <i>Universal Library User's Guide</i> for a list of the supported A/D ranges of each board. Refer to Table 7 on page 208 for a list of valid range settings.</p>
chanCount	Number of elements in each of the three arrays - chanArray, chanTypeArray and gainArray.
rate	<p>The sample rate at which samples are acquired, in samples per second per channel.</p> <p>rate also returns the value of the actual rate set, which may be different from the requested rate because of pacer limitations.</p>
pretrigCount	<p>Sets the number of pre-trigger samples to collect. Specifies the number of samples to collect before the trigger occurs. This method won't run in pre-trigger mode if preTrigCount is set to zero. preTrigCount is ignored if the ExtTrigger option is not specified.</p> <p>preTrigCount also returns the value of the actual pre-trigger count set, which may be different from the set pre-trigger count because pre-trigger count must be a multiple of the channel count (chanCount).</p> <p>preTrigCount must be evenly divisible by the number of channels being scanned (chanCount). If it is not, this method adjusts the number (down) to the next valid value and returns that value to the preTrigCount parameter.</p>

totalCount	<p>Total number of samples to collect. Specifies the total number of samples to collect and store in the buffer. <code>totalCount</code> must be greater than <code>pretrigCount</code>.</p> <p><code>totalCount</code> also returns the value of the actual total count set, which may be different from the requested total count because total count must be a multiple of the channel count (<code>chanCount</code>).</p> <p><code>totalCount</code> must be evenly divisible by the number of channels being scanned (<code>chanCount</code>). If it is not, this method adjusts the number (down) to the next valid value and returns that value to the <code>totalCount</code> parameter.</p>
memHandle	Handle for the Windows buffer to store data. This buffer must have been previously allocated with the <a href="#">WinBufAlloc()</a> method.
options	Bit fields that control various options. All of the <code>options</code> settings are <code>Mccdaq.ScanOptions</code> enumerated constants. This field may contain any combination of non-contradictory choices in the " <a href="#">options parameter values</a> " section below.

**chanTypeArray parameter values :**

Analog	Analog input channel.
Digital8	8-bit digital input port.
Digital16	16-bit digital input port. (FIRSTPORTA only)
Ctr16	16-bit counter.
Ctr32Low	Lower 16-bits of a 32-bit counter.
Ctr32High	Upper 16-bits of a 32-bit counter.
CJC	CJC channel.
TC	<p>Thermocouple channel.</p> <p>The <a href="#">GetTCValues()</a> method can be used to convert raw thermocouple data to data on a temperature scale (Celsius, Fahrenheit or Kelvin). <b>Note:</b> If at least one TC channel is listed in the channel array, and averaging is enabled for that channel, the averaging will be applied to all of the channels listed in the channel array.</p>
SetpointStatus	The setpoint status register. This is a bitfield indicating the state of each of the setpoints. A "1" indicates that the setpoint criteria has been met.

**chanTypeArray flag values:**

SetpointEnable	<p>Enables a setpoint. When this option is specified, it must be OR'ed with the <code>ChanTypeArray</code> parameter values.</p> <p>You set the setpoint criteria with the <a href="#">DaqSetSetpoints()</a> method. The number of channels set with the <code>SetpointEnable</code> flag must match the number of setpoints set by the <a href="#">DaqSetSetpoints()</a> method's <code>setpointCount</code> parameter.</p>
----------------	--

**options parameter values :**

Background	When the <code>Background</code> option is used, control returns immediately to the next line in your program and the data collection from the counters into the buffer continues in the background. If the <code>Background</code> option is not used, the <code>DaqInScan()</code> method does not return to your program until all of the requested data has been collected and returned to the buffer.
------------	--

Use [GetStatus\(\)](#) with `DaqiFunction` to check on the status of the background operation. Use [StopBackground\(\)](#) with `DaqiFunction` to terminate the background process before it has completed. Execute `StopBackground()` after normal termination of all background functions in order to clear variables and flags.

If the `Background` option is not used, the `DaqInScan()` method will not return to your program until all of the requested data has been collected and returned to the buffer. When the `Background` option is used, control will return immediately to the next line in your program and the data collection from the A/D into the buffer will continue in the background. Use [GetStatus\(\)](#) with `DaqiFunction` to check on the status of the background operation. Use [StopBackground\(\)](#) with `DaqiFunction` to terminate the background process before it has completed. Execute `StopBackground()` after normal termination of all background functions in order to clear variables and flags.

Continuous	This option puts the function in an endless loop. Once it collects the required number of samples, it resets to the start of the buffer and begins again. The only way to stop this operation is to use <a href="#">StopBackground()</a> with the <code>DaqiFunction</code> . Normally this option should be used in combination with <code>Background</code> so that your program will regain control.
ExtClock	If this option is used, conversions will be controlled by the signal on the external clock input rather than by the internal pacer clock. Each conversion will be triggered on the appropriate edge of the clock input signal. When this option is used the <code>rate</code> argument is ignored. The sampling rate is dependent on the clock signal. Options for the board will default to a transfer mode that will allow the maximum conversion rate to be attained unless otherwise specified.
ExtTrigger	If this option is specified, the sampling will not begin until the trigger condition is met (refer to the <a href="#">DaqSetTrigger()</a> method).

#### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

`rate` – Actual sampling rate used.

`preTrigCount` – Actual pre-trigger count used.

`totalCount` – Actual total count used.

`memHandle` - Collected data returned via the Windows buffer.

## DaqOutScan()

Outputs values synchronously to analog output channels and digital output ports. This method only works with boards that support synchronous output.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**

```
Public Function DaqOutScan(ByVal chanArray As Short(), ByVal
chanTypeArray As MccDaq.ChannelType(), ByVal gainArray As
MccDaq.Range(), ByVal chanCount As Integer, ByRef rate As Integer,
ByVal count As Integer, ByVal memHandle As Integer, ByVal options As
MccDaq.ScanOptions) As MccDaq.ErrorInfo
```

**Visual Basic:**

```
public ErrorInfo DaqOutScan(short[] chanArray, MccDaq.ChannelType[]
chanTypeArray, MccDaq.Range[] gainArray, int chanCount, ref int rate, int count,
int memHandle, MccDaq.ScanOptions options)
```

### Parameters:

chanArray	Array containing channel values. Valid channel values are analog output channels and digital ports.
chanTypeArray	Array containing channel types. Each element of this array defines the type of the corresponding element in the chanArray. The chanTypeArray settings are Mccdaq.ChannelType enumerated constants. Choices are:  <div style="display: flex; justify-content: space-between;"> <span>Analog</span> <span>Analog output channel.</span> </div> <div style="display: flex; justify-content: space-between;"> <span>Digital16</span> <span>16-bit digital output port. (FirstPortA only)</span> </div>
gainArray	Array containing D/A range codes. If the corresponding element in the ChanArray is not an analog output channel, the range code for this channel is ignored. If the board does not have programmable gain, this parameter is ignored, and therefore can be set to null.
chanCount	Number of elements in each of the three arrays - chanArray, chanTypeArray and gainArray.
rate	Sample rate in scans per second. rate also returns the value of the actual rate set, which may be different from the requested rate because of pacer limitations.
count	Sets the total number of values to output. count must be a multiple of chanCount.
memHandle	Handle for the Windows buffer from which data is output. This buffer must have been previously allocated with the <a href="#">WinBufAlloc()</a> method and data values loaded (for example using <a href="#">WinArrayToBuf()</a> .
options	Bit fields that control various options. All of the options settings are Mccdaq.ScanOptions enumerated constants. This field may contain any combination of non-contradictory choices in the "options parameter values" section below.

### options parameter values:

Background	When this option is used, the output operations begin running in the background, and control immediately returns to the next line of your program.  Use <a href="#">GetStatus()</a> with DaqoFunction to check the status of background operation. Use <a href="#">StopBackground()</a> with DaqoFunction to terminate background operations before they are completed. Execute StopBackground() with DaqoFunction after normal termination of all background functions in order to clear variables and flags.
------------	--

Continuous	<p>This option puts the method in an endless loop. Once it outputs the specified number (<code>Count</code>) of output values, it resets to the start of the buffer and begins again. The only way to stop this operation is by calling <a href="#">StopBackground()</a> with <code>DaqFunction</code>. This option should only be used in combination with <code>Background</code> so that your program regains control.</p>
ExtClock	<p>If this option is used, conversions are paced by the signal on the external clock input rather than by the internal pacer clock. Each conversion is triggered on the appropriate edge of the clock input signal.</p> <p>When this option is used, the <code>rate</code> parameter is ignored. The sampling rate is dependent on the clock signal. Options for the board default to transfer types that allow the maximum conversion rate to be attained unless otherwise specified.</p>
ADCClockTrig	<p>If this option is used, the data output operation is triggered upon the start of the ADC clock.</p>
ADCClock	<p>When this option is used, the data output operation is paced by the ADC clock.</p>
NonStreamedIO	<p>This option allows non-streamed data output to be generated to a specified output channel.</p> <p>In this mode, the aggregate size of data output buffer must be less than or equal to the size of the internal data output FIFO on the Measurement Computing device. This allows the data output buffer to be loaded into the device's internal output FIFO.</p> <p>Once the sample updates are transferred (or downloaded) to the device, the device is responsible for outputting the data. While the size is limited, and the output buffer cannot be changed once the output is started, this mode has the advantage being able to continue data output without having to periodically feed output data through the program to the device.</p>

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

`rate` – Actual sampling rate used.



## DaqSetSetpoints()

Configures up to 16 detection setpoints associated with the input channels within a scan group. This method only works with boards that support synchronous input.

Member of the [MccBoard](#) class.

### Function Prototype:

**VB .NET:**

```
Public Function DaqSetSetpoints(ByVal limitAArray As Single(), ByVal
limitBArray As Single(), ByVal reserved As Single(), ByVal
setpointFlagsArray As MccDaq.SetpointFlag(), ByVal
setpointOutputArray As MccDaq.SetpointOutput(), ByVal output1Array
As Single(), ByVal output2Array As Single(), outputMask1Array As
Single(), outputMask2Array As Single(), ByVal setpointCount As
Integer) As MccDaq.ErrorInfo
```

**C# .NET:**

```
public MccDaq.ErrorInfo DaqSetSetpoints(float[] limitAArray, float[]
limitBArray, float[] reserved, MccDaq.SetpointFlag[]
setpointFlagsArray, MccDaq.SetpointOutput[] setpointOutputArray,
float[] output1Array, float[] output2Array, float[]
outputMask1Array, float[] outputMask2Array, int setpointCount)
```

### Parameters:

limitAArray	Array containing the limit A values for the input channels used for the setpoint. Limit A specifies a value used to determine if the setpoint criteria are met.
limitBArray	Array containing the limit B values for the input channels used for the setpoint. Limit B specifies a value used to determine if the setpoint criteria are met.
reserved	Reserved for future use.
setpointFlagsArray	Array containing the setpoint flags. All of the setpointFlagsArray settings are MccDaq.SetpointFlag enumerated constants. Set it to one of the constants in the "setpointFlagsArray parameter values" section below.
setpointOutputArray	Array containing output sources. All of the setpointOutputArray settings are MccDaq.SetPointOutput enumerated constants. Set it to one of the constants in the "setpointOutputArray parameter values" section on page 354.
output1Array	Array containing the values for the output channels used for the setpoint.
output2Array	Array containing the values for the output channels used for the setpoint.
outputMask1Array	Array containing the output masks for output value 1 – for FIRSTPORTC only.
outputMask2Array	Array containing the output masks for output value 2 – for FIRSTPORTC only.
setpointCount	Number of setpoints to configure (0 -16). Set to 0 to disable the setpoints.

### setpointFlagsArray parameter values:

Flag	Description
EqualLimitA	Setpoint criteria: The input channel = limit A.
LessThanLimitA	Setpoint criteria: The input channel < limit A.
GreaterThanLimitB	Setpoint criteria: The input channel > limit B.
OutsideLimits	Setpoint criteria: The input channel < limit A and > limit B.
InsideLimits	Setpoint criteria: The input channel > limit A and < limit B.
Hysteresis	Setpoint criteria: If the input channel > limit A then output value 1. If the input channel < limit B then output value 2.

UpdateOnTrueOnly      If the criteria is met then output value 1.

UpdateOnTrueAndFalse      If the criteria is met then output value 1, else output value 2.

**setpointOutputArray parameter values:**

Output Source	Description
None	Perform no outputs.
FirstPortC	Output to FIRSTPORTC when the criteria is met.
DAC0	Output to DAC0 when the criteria is met. You must have a device with DAC0.
DAC1	Output to DAC1 when the criteria is met. You must have a device with DAC1.
DAC2	Output to DAC2 when the criteria is met. You must have a device with DAC2.
DAC3	Output to DAC3 when the criteria is met. You must have a device with DAC3.
TMR0	Output to timer 0 when the criteria is met.
TMR1	Output to timer 1 when the criteria is met.

**Returns:**

Error code or 0 if no errors

## DaqSetTrigger()

Selects the trigger source and sets up its parameters. This trigger is used to initiate or terminate an acquisition using the [DaqInScan\(\)](#) method if the `ExtTrigger` option is selected. This method only works with boards that support synchronous output.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**

```
Public Function DaqSetTrigger(ByVal trigSource As MccDaq.TriggerSource, ByVal trigSense MccDaq.TriggerSensitivity, ByVal trigChan As Integer, ByVal chanType As MccDaq.ChannelType, ByVal gain As MccDaq.Range, ByVal level As Single, ByVal variance As Single, ByVal trigEvent As MccDaq.TriggerEvent) As MccDaq.ErrorInfo
```

**Visual Basic:**

```
public ErrorInfo DaqSetTrigger(MccDaq.TriggerSource trigSource, MccDaq.TriggerSensitivity trigSense, int trigChan, MccDaq.ChannelType chanType, MccDaq.Range gain, float level, float variance, MccDaq.TriggerEvent trigEvent)
```

### Parameters:

trigSource	Specifies the type of triggering based on the external trigger source. All of the <code>trigSource</code> settings are <code>MccDaq.TriggerSource</code> enumerated constants. Set it to one of the constants in the " <a href="#">trigSource parameter values</a> " section on page 356.
trigSense	Specifies the trigger sensitivity. The trigger sensitivity normally defines the way in which a trigger event is detected based upon the characteristics of the trigger input signal. Often, it defines the way in which the trigger input signal(s) should be compared to the trigger level parameter value.  All of the <code>trigSense</code> settings are <code>MccDaq.TriggerSensitivity</code> enumerated constants. Set it to one of the constants in the " <a href="#">trigSense parameter values</a> " section on page 356.
trigChan	The trigger channel. This channel must be a configured channel in the channel array (refer to <a href="#">DaqInScan()</a> ).
chanType	The channel type. All of the <code>chanType</code> settings are <code>MccDaq.ChannelType</code> enumerated constants. <code>chanType</code> should match the channel type setting for the trigger channel configured using the <a href="#">DaqInScan()</a> method.
gain	The trigger channel gain code. If the device has programmable gain, this parameter should match the gain code setting when the channel is configured using the <a href="#">DaqInScan()</a> method. The <code>gain</code> parameter is ignored if <code>trigChan</code> is not an analog channel.
level	A single precision floating point value which represents, in engineering units, the level at or around which the trigger event should be detected. This option is used for trigger types that depend on an input channel comparison to detect the start trigger or stop trigger event.  The actual level at which the trigger event is detected depends upon trigger sensing and variability. Refer to <a href="#">Trigger levels</a> on page 357 for more information.
variance	A single-precision floating point value which represents, in engineering units, the amount that the trigger event can vary from the <code>level</code> parameter.  While the <code>TrigSense</code> parameter indicates the direction of the input signal relative to the <code>level</code> parameter, the <code>variance</code> parameter specifies the degree to which the input signal can vary relative to the <code>level</code> parameter.

trigEvent	<p>Specifies the trigger event type. Valid values indicate either a start trigger event (<code>MccDaq.TriggerEvent.Start</code>) or a stop trigger event (<code>MccDaq.TriggerEvent.Stop</code>).</p> <p><b>Start:</b> The start trigger event defines the conditions under which post-trigger acquisition data collection should be initiated or triggered. The start trigger event can vary in complexity from starting immediately, to starting on complex channel value definitions.</p> <p><b>Stop:</b> The stop trigger event signals the current data acquisition process to terminate. The stop trigger event can be as simple as that of a scan count, or as complex as involving a channel value level condition.</p>
-----------	---

**trigSource parameter values:**

TrigImmediate	Start trigger event only. Acquisition begins immediately upon invocation the <a href="#">DaqInScan()</a> method. No pre-trigger data acquisition is possible with this trigger type.
TrigExtTTL	Start trigger event only. Acquisition begins on the selectable edge of an external TTL signal. No pre-trigger data acquisition is possible with this trigger type.
TrigAnalogHW	Start trigger event only. Acquisition begins upon a selectable criteria of the input signal (above level, below level, rising edge, etc.) <code>trigChan</code> must be defined as the first channel in the channel scan group. No pre-trigger data acquisition is possible with this trigger type.
TrigAnalogSW	Post-trigger data acquisition begins upon a selectable criteria of the input signal (above level, below level, rising edge, etc.)
TrigDigPattern	Post-trigger data acquisition begins upon receiving a specified digital pattern on the specified digital port.
TrigCounter	Post-trigger data acquisition begins upon detection of specified counter criteria.
TrigScanCount	Stop trigger event only. Stops collecting post-trigger data when the specified number of post-trigger scans are completed.

**trigSense parameter values:**

RisingEdge	Triggers when the signal goes from low to high (TTL trigger) or rises through a specified level (hardware analog, software analog, and counter).
FallingEdge	Triggers when the signal goes from high to low (TTL trigger) or falls through a specified level (hardware analog, software analog, and counter).
AboveLevel	Triggers when the signal is above a specified level (hardware analog, software analog, counter, and digital pattern).
BelowLevel	Triggers when the signal is below a specified level (hardware analog, software analog, counter, and digital pattern).
EqLevel	Triggers when the signal equals a specified level (hardware analog, software analog, counter, and digital pattern).
NeLevel	Triggers when the signal does not equal a specified level (hardware analog, software analog, counter, and digital pattern).

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**Notes:**

**Trigger levels:** The actual level at which the trigger event is detected depends upon trigger sensing and variability. The various ranges of possible values for the `level` parameter based on the trigger source are:

<code>TrigAnalogHW</code>	The voltage used to define the trigger level. Trigger detection is performed in hardware.
<code>TrigAnalogSW</code>	The voltage used to define the trigger level. Trigger detection is performed in software.
<code>TrigDigPattern</code>	Sets the bit pattern for the digital channel trigger. Choices are: 0.0 (no bits set): 255.0 (all bits set) for 8-bit digital ports. 0.0 (no bits set): 65,535.0 (all bits set) for 16-bit digital ports.
<code>TrigCounter</code>	Selects either Pulse or Totalize counter values (0.0 – 65,535).
<code>TrigImmediate</code>	Ignored
<code>TrigScanCount</code>	Ignored

**Trigger start and stop criteria:** The table below lists the trigger start and stop criteria based on the selected trigger type and sensitivity.

Table 8. Interaction of trigger variance with trigger level and trigger sensitivities

Trigger Start/Stop Source (TrigSource)	Trigger Sensitivity (TrigSense)	Trigger Start/Stop Criteria
TrigAnalogHW (Start trigger event only)	RisingEdge	Triggers when the signal value < (level – variance) Then, the signal value > level
	FallingEdge	Triggers when the signal value > (level + variance) Then, the signal value < level
	AboveLevel	Triggers when the signal value > (level)
	BelowLevel	Triggers when the signal value < (level)
TrigAnalogSW	RisingEdge	Triggers/stops when the signal value < (level – variance) Then, the signal value > level
	FallingEdge	Triggers/stops when the signal value > (level + variance) Then, the signal value < level
	AboveLevel	Triggers/stops when the signal value > (level)
	BelowLevel	Triggers/stops when the signal value < (level)
	EqLevel	Triggers/stops when the (level – Variance) < signal value < (level + variance)
	NeLevel	Triggers/stops when the signal value < (level – variance) OR when the signal value > (level + variance)
TrigDigPattern	AboveLevel	Triggers/stops when (digital port value AND (bitwise) variance) > (level AND (bitwise) variance)
	BelowLevel	Triggers/stops when the (digital port value AND (bitwise) variance) < (level AND (bitwise) variance)
	EqLevel	Triggers/stops when the (digital port value AND (bitwise) variance) = (level AND (bitwise) variance)
	NeLevel	Triggers/stops when the (digital port value AND (bitwise) variance) != (level AND (bitwise) variance)

Trigger Start/Stop Source (TrigSource)	Trigger Sensitivity (TrigSense)	Trigger Start/Stop Criteria
TrigCounter	RisingEdge	Triggers/stops when the counter channel $< (\text{level} - \text{variance})$ Then, the counter channel $> \text{level}$
	FallingEdge	Triggers/stops when the counter channel $> (\text{level} + \text{variance})$ Then, the counter channel $< \text{level}$
	AboveLevel	Triggers/stops when the counter channel $> (\text{level} - \text{variance})$
	BelowLevel	Triggers/stops when the counter channel $< (\text{level} + \text{variance})$
	EqLevel	Triggers/stops when $(\text{level} - \text{variance}) < \text{counter channel} < (\text{level} + \text{variance})$
	NeLevel	Triggers/stops when the counter channel $< (\text{level} - \text{variance})$ OR when the counter channel $> (\text{level} + \text{variance})$

---

# Temperature Input Methods

## Introduction

Use the methods explained in this chapter to convert a raw analog input from an EXP or other temperature sensor board to temperature.

## TIn()

Reads an analog input channel, linearizes it according to the selected temperature sensor type, and returns the temperature in degrees.

The CJC channel, the gain, and sensor type, are read from the *InstaCal* configuration file. Run the *InstaCal* configuration program to set these items.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**                   Public Function TIn(ByVal chan As Integer, ByVal scale As MccDaq.TempScale, ByRef tempValue As Single, ByVal options As MccDaq.ThermocoupleOptions ) As MccDaq.ErrorInfo

**C# .NET:**                   public MccDaq.ErrorInfo TIn(int chan, MccDaq.TempScale scale, out float tempValue, MccDaq.ThermocoupleOptions options)

### Parameters:

chan	Input channel to read.
scale	Specifies the temperature scale that the input is converted to. Choices are MccDaq.TempScale.Celsius, MccDaq.TempScale.Fahrenheit and MccDaq.TempScale.Kelvin.
tempValue	The temperature in degrees is returned here. Thermocouple resolution is approximately 0.25 °C, depending on scale, range and thermocouple type. RTD resolution is 0.1 °C.
options	Bit fields that control various options. Set it to one of the constants in the "options parameter values" section below.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

tempValue - Temperature returned here

### options parameter values:

All of the options settings are MccDaq.ThermocoupleOptions enumerated constants. To set a variable to one of these constants, you must refer to the MccDaq object and the ThermocoupleOptions enumeration (*variable*= MccDaq.ThermocoupleOptions.Filter or *variable* = MccDaq.ThermocoupleOptions.NoFilter).

Filter	When selected, a smoothing function is applied to temperature readings, very much like the electrical smoothing inherent in all hand held temperature sensor instruments. This is the default. Ten samples are read from the specified channel and averaged. The average is the reading returned. Averaging removes normally distributed signal line noise.
NoFilter	When selected, the temperature readings are not smoothed, resulting in a scattering of readings around a mean.



**Notes:**

**Using CIO-EXP boards:** For CIO-EXP boards, the channel number is calculated using the following formula, where:

- ADChan is the A/D channel that is connected to the multiplexer
  - MuxChan is a number ranging from 0 to 15 that specifies the channel number on a particular bank of the multiplexer board
- $$\text{Chan} = (\text{ADChan} * 16) + (16 + \text{MuxChan})$$

For example, you have an EXP16 connected to a CIO-DAS08 via the CIO-DAS08 channel 0. (Remember that DAS08 channels are numbered 0, 1, 2, 3, 4, 5, 6 & 7). If you connect a thermocouple to channel 5 of the EXP16, the value for chan would be  $(0 * 16) + (16 + 5) = 0 + 21 = \mathbf{21}$ .

**Using 6K-EXP boards:** For 6K-EXP boards, the channel number is calculated using one of the following formulas, where:

- ADChan is the A/D channel that is connected to the multiplexer.
  - MuxChan is a number ranging from 0 to 15 that specifies the channel number on a particular bank of the multiplexer board.
  - If the A/D board has 16 or less single-ended channels:
- $$\text{Chan} = (\text{ADChan} * 16) + (16 + \text{MuxChan})$$

For example, you have a 6K-EXP16 connected to a PCI-DAS6052 via the a PCI-DAS6052 channel 0. If you connect a thermocouple to channel 5 of the 6K-EXP16, the value for Chan would be  $(0 * 16) + (16 + 5) = 0 + 21 = \mathbf{21}$ .

- If the A/D board has 64 single-ended channels and the A/D multiplexer channel is less than or equal to 7:
- $$\text{Chan} = (\text{ADChan} * 16) + (64 + \text{MuxChan})$$

For example, you have a 6K-EXP16 connected to a PCI-DAS6031 via the a PCI-DAS6031 channel 7. If you connect a thermocouple to channel 5 of the 6K-EXP16, the value for Chan would be  $(7 * 16) + (64 + 5) = 112 + 69 = \mathbf{181}$ .

- If the A/D board has 64 single-ended channels and the A/D multiplexer channel is greater than or equal to 31:

$$\text{Chan} = (\text{ADChan} * 16 - 320) + \text{MuxChan}$$

For example, you have a 6K-EXP16 connected to a PCI-DAS6031 via the PCI-DAS6031 channel 32. If you connect a thermocouple to channel 5 of the 6K-EXP16, the value for Chan is  $(32 * 16 - 320) + 5 = 192 + 5 = \mathbf{197}$ .

**CJC Channel:** The Cold Junction Compensation (CJC) channel is set in the *InstaCal* install program. If you have multiple EXP boards, Universal Library will apply the CJC reading to the linearization formula in the following manner:

- First, if you have chosen a CJC channel for the EXP board that the channel you are reading is on, it will use the CJC temp reading from that channel.
- Second, if you left the CJC channel for the EXP board that the channel you are reading is on to NOT SET, the library will use the CJC reading from the next lower EXP board with a CJC channel selected.

For example: You have 4 CIO-EXP16 boards connected to a CIO-DAS08 on channel 0, 1, 2 and 3. You choose CIO-EXP16 #1 (connected to CIO-DAS08 channel 0) to have its CJC read on CIO-DAS08 channel 7, AND, you leave the CIO-EXP16's 2, 3 and 4 CJC channels to NOT SET. Result: The CIO-EXP boards all use the CJC reading from CIO-EXP16 #1, connected to channel 7 for linearization. As you can see, it is important to keep the CIO-EXP boards in the same case and out of any breezes to ensure a clean CJC reading.

**Important**

If the EXP board is connected to an A/D that does not have programmable gain (DAS08, DAS16, DAS16F) then the A/D board range is read from the configuration file (cb.cfg). In most cases, hardware selectable ranges should be set to  $\pm 5$  V for thermocouples and 0 to 10 V for RTDs. Refer to the board-specific information in the *Universal Library User's Guide* (available on our web site at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf)) or in the user manual for your board. If the board does have programmable RTDs gains, the TIn () method will set the appropriate A/D range.

**Specific Errors:** If an OutOfRange or OpenConnection error occurs, the value returned is -9999.0. If a NotReady error occurs, the value returned is -9000.

## TInScan()

Reads a range of channels from an analog input board, linearizes them according to temperature sensor type, and returns the temperatures to an array in degrees.

The CJC channel, the gain, and temperature sensor type are read from the configuration file. Use the *InstaCal* configuration program to change any of these options.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**

```
Public Function TInScan(ByVal lowChan As Integer, ByVal highChan As Integer, ByVal scale As MccDaq.TempScale, ByVal dataBuffer As Single(), ByVal options As MccDaq.ThermocoupleOptions) As MccDaq.ErrorInfo
```

**C# .NET:**

```
public MccDaq.ErrorInfo TInScan(int lowChan, int highChan, MccDaq.TempScale scale, out float dataBuffer, MccDaq.ThermocoupleOptions options)
```

### Parameters:

lowChan	Low mux channel of scan.
highChan	High mux channel of scan.
scale	Specifies the temperature scale that the input is converted to. Choices are MccDaq.TempScale.Celsius, MccDaq.TempScale.Fahrenheit and MccDaq.TempScale.Kelvin.
dataBuffer	The temperature is returned in degrees. Each element in the array corresponds to a channel in the scan. dataBuffer must be at least large enough to hold (highChan - lowChan + 1) temperature values. Thermocouple resolution is approximately 0.25 °C, depending on scale, range and thermocouple type. RTD resolution is 0.1 °C.
options	Bit fields that control various options. Set it to one of the constants in the "options parameter values" section below.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

dataBuffer[] - Temperature values in degrees are returned here for each channel in scan.

### options parameter values:

All of the options settings are MccDaq.ThermocoupleOptions enumerated constants. To set a variable to one of these constants, you must refer to the MccDaq object and the ThermocoupleOptions enumeration (*variable* = MccDaq.ThermocoupleOptions.Filter or *variable* = MccDaq.ThermocoupleOptions.NoFilter).

Filter	When selected, a smoothing function is applied to temperature readings, very much like the electrical smoothing inherent in all hand held temperature sensor instruments. This is the default. Ten samples are read from the specified channel and averaged. The average is the reading returned. Averaging removes normally distributed signal line noise.
NoFilter	When selected, the temperature readings are not smoothed, resulting in a scattering of readings around a mean.

**Notes:**

**Using EXP boards:** For EXP boards, these channel numbers (Chan) are calculated using the following formula:

- ADChan = A/D channel that is connected to the multiplexer
  - MuxChan is a number ranging from 0 to 15 that specifies the channel number on a particular bank of the multiplexer board
- $$\text{Chan} = (\text{ADChan} * 16) + (16 + \text{MuxChan})$$

For example, you have an EXP16 connected to a CIO-DAS08 via the CIO-DAS08 channel 0. (Remember, DAS08 channels are numbered 0, 1, 2, 3, 4, 5, 6 & 7). If you connect thermocouples to channels 5, 6, and 7 of the EXP16, the value for lowChan would be  $(0+1) * 16 + 5 = \mathbf{21}$ , and the value for highChan would be  $(0+1) * 16 + 7 = \mathbf{23}$ .

**Important**

For an EXP board connected to an A/D board that does not have programmable gain (DAS08, DAS16, DAS16F), the A/D board range is read from the configuration file (cb.cfg). In most cases, set hardware-selectable ranges to  $\pm 5$  V for thermocouples, and to 0 to 10 V for RTDs. Refer to the board-specific information in the *Universal Library User's Guide* (available on our web site at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf)) or in the user manual for your board. If the board has programmable RTDs gains, the TIn() method sets the appropriate A/D range.

**Using 6K-EXP boards:** For 6K-EXP boards, the channel number (Chan) is calculated using one of the following formulas, where:

- ADChan is the A/D channel that is connected to the multiplexer.
  - MuxChan is a number ranging from 0 to 15 that specifies the channel number on a particular bank of the multiplexer board.
  - If the A/D board has 16 or less single-ended channels:
- $$\text{Chan} = (\text{ADChan} * 16) + (16 + \text{MuxChan})$$

For example, you have a 6K-EXP16 connected to a PCI-DAS6052 via the a PCI-DAS6052 channel 0. If you connect a thermocouple to channels 5, 6, and 7 of the 6K-EXP16, the value for lowChan would be  $(0 * 16) + (16 + 5) = 0 + 21 = \mathbf{21}$ , and the value for highChan would be  $(0 * 16) + (16 + 5) = 0 + 231 = \mathbf{23}$ .

- If the A/D board has 64 single-ended channels and the A/D multiplexer channel is less than or equal to 7:
- $$\text{Chan} = (\text{ADChan} * 16) + (64 + \text{MuxChan})$$

For example, you have a 6K-EXP16 connected to a PCI-DAS6031 via the PCI-DAS6031 channel 7. Connect a thermocouple to channels 5, 6, and 7 of the 6K-EXP16: the value for lowChan is  $(7 * 16) + (64 + 5) = 112 + 69 = \mathbf{181}$ , and the value for highChan is  $(7 * 16) + (64 + 7) = 112 + 71 = \mathbf{183}$ .

- If the A/D board has 64 single-ended channels and the A/D multiplexer channel is greater than or equal to 32:
- $$\text{Chan} = (\text{ADChan} * 16 - 320) + \text{MuxChan}$$

For example, you have a 6K-EXP16 connected to a PCI-DAS6031 via the PCI-DAS6031 channel 32. If you connect a thermocouple to channels 5, 6, and 7 of the 6K-EXP16, the value for lowChan is  $(32 * 16 - 320) + 5 = 192 + 5 = \mathbf{197}$ , and the value for highChan is  $(32 * 16 - 320) + 7 = 192 + 7 = \mathbf{199}$ .

**CJC Channel:** The Cold Junction Compensation (CJC) channel is set in the InstaCal install program. If you have multiple EXP boards, Universal Library will apply the CJC reading to the linearization formula in the following manner:

- First, if you have chosen a CJC channel for the EXP board that the channel you are reading is on, it will use the CJC temp reading from that channel.

- Second, if you have left the CJC channel for the EXP board that the channel you are reading is on to NOT SET, the library will use the CJC reading from the next lower EXP board with a CJC channel selected.

For example: You have 4 CIO-EXP16 boards connected to a CIO-DAS08 on channel 0, 1, 2 and 3. You choose CIO-EXP16 #1 (connected to CIO-DAS08 channel 0) to have its CJC read on CIO-DAS08 channel 7, AND, you leave the CIO-EXP16's 2, 3 and 4 CJC channels to NOT SET. Result: The CIO-EXP boards all use the CJC reading from CIO-EXP16 #1, connected to channel 7 for linearization. As you can see, it is important to keep the CIO-EXP boards in the same case and out of any breezes to ensure a clean CJC reading.

**Important**

In order to understand the functions, you must read the board-specific information contained in the *Universal Library User's Guide* (available on our web site at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf)).

Review and run the example programs before attempting any programming of your own. Following this advice will save you hours of frustration, and possibly time wasted holding for technical support.

This note, which appears elsewhere, is especially applicable to this method. Read the board-specific information for your board (see the *Universal Library User's Guide*). We suggest that you make a copy of that page to refer to as you read this manual and examine the example programs.

**Specific errors:** For most boards, if an `OUTOFRANGE` or `OPENCONNECTION` error occurs, the value in the array element associated with the channel causing the error returned will be -9999.0 (Refer to board-specific information in the *Universal Library User's Guide*).

---

# Windows Memory Management Methods

## Introduction

Use the methods explained in this section to allocate, free, and copy to/from Windows global memory buffers.

## WinBufAlloc()

Allocates a Windows global memory buffer which can be used with the scan methods and returns a memory handle for it.

Member of the [MccService](#) class.

### Function prototype:

VB .NET:                   Public Shared Function WinBufAlloc(ByVal numPoints As Integer) As Integer

C# .NET:                   public static int WinBufAlloc(int numPoints)

### Parameters:

numPoints                   Size of buffer to allocate. Specifies how many data points (16-bit integers, NOT bytes) can be stored in the buffer.

### Returns:

0 if buffer could not be allocated or a non-zero integer handle to the buffer.

### Notes:

Unlike most other methods in the library, this method does not return an [ErrorInfo](#) object. It returns a Windows global memory handle, which can then be passed to the scan methods in the library. If an error occurs, the handle will come back as 0 to indicate that the buffer was not allocated.

## WinBufAlloc32()

Allocates a Windows global memory buffer for use with 32-bit scan methods, and returns a memory handle for the buffer.

Member of the [MccService](#) class.

### Function prototype:

VB .NET:                   Public Shared Function WinBufAlloc32(ByVal numPoints As Integer) As Integer

C# .NET:                   public int WinBufAlloc32(int numPoints)

### Parameters:

numPoints                   The size of buffer to allocate. Specifies how many data points (32-bit integers, NOT bytes) can be stored in the buffer.

### Returns:

0 if buffer could not be allocated, or a non-zero integer handle to the buffer.

### Notes:

Unlike most other methods in the library, this method does not return an error code. It returns a Windows global memory handle which can then be passed to the scan methods in the library. If an error occurs, the handle will come back as 0 to indicate that the buffer was not allocated.



## WinBufFree()

Frees a Windows global memory buffer which was previously allocated with the [WinBufAlloc\(\)](#) or [WinBufAlloc32\(\)](#) method.

Member of the [MccService](#) class.

### Function prototype:

VB .NET:           Public Shared Function WinBufFree(ByVal memHandle As Integer) As  
                          MccDaq.ErrorInfo

C# .NET:           public static MccDaq.ErrorInfo WinBufFree(int memHandle)

### Parameters:

memHandle           A Windows memory handle. This must be a memory handle that was returned by  
WinBufAlloc() or WinBufAlloc32() when the buffer was allocated.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

## WinArrayToBuf()

Copies data from a one-dimensional or two-dimensional array into a Windows memory buffer.

Member of the [MccService](#) class.

### Function prototype:

VB .NET:

Copies data from a one-dimensional array of short values:

```
Public Shared Function WinArrayToBuf(ByRef dataArray As Short, ByVal memHandle As Integer, ByVal firstPoint As Integer, ByVal numPoints As Integer) As MccDaq.ErrorInfo
```

Copies data from a one-dimensional array of System.UInt16 values:

```
Public Shared Function WinArrayToBuf(ByRef dataArray As System.UInt16, ByVal memHandle As Integer, ByVal firstPoint As Integer, ByVal numPoints As Integer) As MccDaq.ErrorInfo
```

Copies data from a two-dimensional array of double values:

```
Public Shared WinArrayToBuf(ByRef dataArray(,) As Double, ByVal memHandle As Integer, ByVal firstPoint As Integer, ByVal numPoints As Integer, ByVal numChannels As Integer) As MccDaq.ErrorInfo
```

C# .NET:

Copies data from a one-dimensional array of short values:

```
public static MccDaq.ErrorInfo WinArrayToBuf(ref short dataArray, int memHandle, int firstPoint, int numPoints)
```

Copies data from a one-dimensional array of System.UInt16 values:

```
public static MccDaq.ErrorInfo WinArrayToBuf(ref ushort dataArray, int memHandle, int firstPoint, int numPoints)
```

Copies data from a two-dimensional array of double values:

```
public static MccDaq.ErrorInfo WinArrayToBuf(ref double[,] dataArray, int memHandle, int firstPoint, int numPoints, int numChannels)
```

### Parameters:

dataArray	The array containing the data to be copied. The first dimension should equal the number of channels. The second dimension should equal the number of points/channel.
memHandle	This must be a memory handle that was returned by <a href="#">WinBufAlloc()</a> when the buffer was allocated. The data will be copied into this buffer.
firstPoint	Index of the first point in the memory buffer where data will be copied to.
numPoints	Number of data points to copy from dataArray.
numChannels	Number of channels to copy from dataArray.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

### Notes:

This method copies data from an array to a Windows global memory buffer. This would typically be used to initialize the buffer with data before doing an output scan. You can use the `firstPoint` and `numPoints` parameters to fill a portion of the buffer. This is useful if you want to send new data to the buffer after a Background + Continuous output scan has been started, for example during circular buffering.

## WinBufToArray()

Copies data from a Windows memory buffer into a one-dimensional or two-dimensional array.

Member of the [MccService](#) class.

### Function prototype:

VB .NET:

Copies data to a one-dimensional array of short values:

```
Public Shared Function WinBufToArray(ByVal memHandle As Integer,
ByRef dataArray As Short, ByVal firstPoint As Integer, ByVal
numPoints As Integer) As MccDaq.ErrorInfo
```

Copies data to a one-dimensional array of System.UInt16 values:

```
Public Shared Function WinBufToArray(ByVal memHandle As Integer,
ByRef dataArray As System.UInt16( ), ByVal firstPoint As Integer,
ByVal numPoints As Integer) As MccDaq.ErrorInfo
```

Copies data to a two-dimensional array of double values:

```
public static ErrorInfo WinBufToArray(ByVal memHandle As Integer,
ByRef dataArray(,) As Double, ByVal firstPoint As Integer, ByVal
numPoints As Integer, ByVal numChannels As Integer) As
MccDaq.ErrorInfo
```

C# .NET:

Copies data to a one-dimensional array of short values:

```
public static MccDaq.ErrorInfo WinBufToArray(int memHandle, out
short dataArray, int firstPoint, int numPoints)
```

Copies data to a one-dimensional array of System.UInt16 values:

```
public static MccDaq.ErrorInfo WinBufToArray(int memHandle, out
ushort dataArray, int firstPoint, int numPoints)
```

Copies data to a two-dimensional array of double values:

```
public static MccDaq.ErrorInfo WinBufToArray(int memHandle, out
double [,] dataArray, int firstPoint, int numPoints, int
numChannels)
```

### Parameters:

memHandle

This must be a memory handle that was returned by WinBufAlloc() when the buffer was allocated. The data will be copied from this buffer.

dataArray

Array that the data will be copied to. The first dimension should equal the number of channels. The second dimension should equal the number of points/channel.

firstPoint

Index of the first point in the memory buffer that the data will be copied from.

numPoints

Number of data points to copy into dataArray.

numChannels

Number of channels to copy into dataArray.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

### Notes:

This method copies data from a Windows global memory buffer to a single value or into an array of doubles. This would typically be used to retrieve data from the buffer after executing an input scan method. You can use the firstPoint and numPoints parameters to copy only a portion of the buffer to the array. This can be useful if you want foreground code to manipulate previously collected data while a Background scan continues to collect new data.

## WinBufToArray32()

Copies 32-bit data from a Windows global memory buffer into a one-dimensional or two-dimensional array. This method is typically used to retrieve data from the buffer after executing an input scan method.

Member of the [MccService](#) class.

### Function prototype:

VB .NET:

Copies data into a two-dimensional array of double values:

```
Public Shared Function WinBufToArray32 (ByVal memHandle As Integer,
ByRef dataArray(,) As Double, ByVal firstPoint As Integer, ByVal
numPoints As Integer, ByVal numChannels As Integer) As
MccDaq.ErrorInfo
```

Copies data into an array of integer values:

```
Public Shared Function WinBufToArray32 (ByVal memHandle As Integer,
ByRef dataArray As Integer, ByVal firstPoint As Integer, ByVal
numPoints As Integer) As MccDaq.ErrorInfo
```

Copies data into an array of System.UInt32 values:

```
Public Shared Function WinBufToArray32 (ByVal memHandle As Integer,
ByRef dataArray As System.UInt32(), ByVal firstPoint As Integer,
ByVal numPoints As Integer) As MccDaq.ErrorInfo
```

C# .NET:

Copies data into a two-dimensional array of double values:

```
public MccDaq.ErrorInfo WinBufToArray32 (int memHandle, out
double[,] dataArray, int firstPoint, int numPoints, int numChannels)
```

Copies data into an array of integer values:

```
public MccDaq.ErrorInfo WinBufToArray32(int memHandle, out int
dataArray, int firstPoint, int numPoints)
```

Copies data into an array of System.UInt32 values:

```
public MccDaq.ErrorInfo WinBufToArray32(int memHandle, out uint
dataArray, int firstPoint, int numPoints)
```

### Parameters:

memHandle

The memory handle that was returned by [WinBufAlloc32\(\)](#) when the buffer was allocated. The buffer should contain the data that you want to copy.

dataArray

The array where that the data is copied.

firstPoint

The index of first point in the memory buffer that data is copied from.

numPoints

The number of data points to copy.

numChannels

The number of channels to copy into dataArray.

### Returns:

[Error code](#) or 0 if no errors

**Notes:**

You can copy only a portion of the buffer to the array using the `firstPoint` and `numPoints` argument. This is useful if you want foreground code to manipulate previously collected data while a Background scan continues to collect new data.

Although this method is available to both Windows C and Delphi programs, it is not necessary, since you can manipulate the memory buffer directly by casting the `MemHandle` returned from `WinBufAlloc32()` to the appropriate type. This method avoids having to copy the data from the memory buffer to an array. Refer to the following example:

```
/*declare and initialize the variables*/
long numPoints = 1000;
unsigned short *dataArray = NULL;
int MemHandle = 0;

/*allocate the buffer and cast it to a pointer to an unsigned long*/
MemHandle = WinBufAlloc32(numPoints);
dataArray = (unsigned long*)MemHandle;

/*scan in the data*/
CInScan(.....,MemHandle,...);

/*print the results*/
for(int i=0; i<numPoints; ++i)
    printf("Data[%d]=%d\n", i, dataArray[i]);

/*free the buffer and NULL the pointer*/
WinBufFree(MemHandle);
dataArray = NULL;
```

---

## Miscellaneous Methods, Properties, and Delegates

### Introduction

The methods and properties explained in this chapter do not as a group fit into a single category. They get and set board information, convert units, manage events and background operations, and perform serial communication operations.

## BoardName property

Name of the board associated with an instance of the MccBoard class.

Member of the [MccBoard](#) class.

### Function prototype:

VB .NET:	<code>Public ReadOnly Property BoardName As String</code>
C# .NET:	<code>public string BoardName [get]</code>

## DeviceLogin()

Opens a device session with a shared device.

Member of the [MccBoard](#) class.

### Function prototype:

VB .NET                   Public Function DeviceLogin(ByVal userName As String, ByVal password As String) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo DeviceLogin(System.String userName, System.String password)

### Parameters:

userName                   A null-terminated string that identifies the user name used to log in to a device session.

password                   A null-terminated string that identifies the password used to log in to a device session.

### Returns:

[Error code](#) or 0 if no errors.

### Notes:

If the user name or password is invalid, the INVALIDLOGIN error is returned.

If the session is already opened by another user, the SESSIONINUSE error is returned.

## DeviceLogout()

Releases the device session with a shared device.

Member of the [MccBoard](#) class.

### Function prototype:

VB .NET                   Public Function DeviceLogout() As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo DeviceLogout()



## HideLoginDialog()

Prevents the default login dialog from being shown when a protected function is called while not logged in.

Member of the [MccBoard](#) class.

### Function prototype:

VB .NET	Public Function HideLoginDialog(ByVal hide As Boolean) As MccDaq.ErrorInfo
C# .NET:	public MccDaq.ErrorInfo HideLoginDialog(System.Boolean hide)

### Parameters:

hide	If true, default dialog will not be shown when a protected function is called while the user is not logged in.
------	--

### Returns:

[Error code](#) or 0 if no errors.

### Notes:

Overrides *InstaCal* Show Login Dialog prompt setting.

## DisableEvent()

Disables one or more event conditions, and disconnects their user-defined handlers.

Member of the [MccBoard](#) class.

### Function prototype:

VB .NET:           Public Function DisableEvent(ByVal eventType As MccDaq.EventType) As MccDaq.ErrorInfo

C# .NET:           public MccDaq.ErrorInfo DisableEvent(MccDaq.EventType eventType)

### Parameters:

eventType           Specifies one or more event conditions that will be disabled. More than one event type can be specified by bitwise OR'ing the event types. Note that specifying an event that has not been enabled is benign and will not cause any errors. Refer to "[eventType parameter values](#)" on page 380 for a list of valid event types.  
To disable all events in a single call, use AllEventTypes.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

### Notes:

For most event types, this method cannot be called while any background operations ([AInScan\(\)](#), [APretrig\(\)](#), or [AOutScan\(\)](#)) are active. Perform a [StopBackground\(\)](#) before calling `DisableEvent()`. However, for `OnExternalInterrupt` events, you can call `DisableEvent()` while the board is actively generating events.

### Important

In order to understand the functions, you must read the board-specific information contained in the *Universal Library User's Guide* (available on our web site at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf)).

Review and run the example programs before attempting any programming of your own. Following this advice will save you hours of frustration, and possibly time wasted holding for technical support.

This note, which appears elsewhere, is especially applicable to this method. Now is the time to read the board-specific information for your board (see the *Universal Library User's Guide*). We suggest that you make a copy of that page to refer to as you read this manual and examine the example programs.

## EnableEvent()

This method binds one or more event conditions to a user-defined callback function. Upon detection of an event condition, the user-defined function is invoked with board- and event-specific data. Detection of event conditions occurs in response to interrupts. Typically, this method is used in conjunction with interrupt driven processes such as [AInScan](#), [APretrig](#), or [AOutScan](#).

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**

```
Public Function EnableEvent (ByVal eventType As MccDaq.EventType,
    ByVal eventParameter As Integer, ByVal callbackFunc As
    MccDaq.EventCallback, ByVal userData As IntPtr) As MccDaq.ErrorInfo
Public Function EnableEvent (ByVal eventType As MccDaq.EventType,
    ByVal eventParameter As System.UInt32, ByVal callbackFunc As
    MccDaq.EventCallback, ByVal userData As IntPtr) As MccDaq.ErrorInfo
Public Function EnableEvent (ByVal eventType As MccDaq.EventType,
    ByVal eventParameter As MccDaq.EventParameter, ByVal callbackFunc As
    MccDaq.CallbackFunction, ByVal userData As IntPtr) As
    MccDaq.ErrorInfo
```

**C# .NET:**

```
public MccDaq.ErrorInfo EnableEvent (MccDaq.EventType eventType, uint
    eventParameter, MccDaq.EventCallback callbackFunc, System.IntPtr
    userData)

public MccDaq.ErrorInfo EnableEvent (MccDaq.EventType eventType, int
    eventParameter, MccDaq.EventCallback callbackFunc, System.IntPtr
    userData)

public MccDaq.ErrorInfo EnableEvent (MccDaq.EventType eventType,
    MccDaq.EventParameter eventParameter, MccDaq.CallbackFunction
    callbackFunc, System.IntPtr userData)
```

### Parameters:

eventType	Specifies one or more event conditions that will be bound to the user-defined callback function. More than one event type can be specified by bitwise OR'ing the event types. Set it to one of the constants in the "eventType parameter values" section on page 380.
eventParameter	<p>Additional data required to specify some event conditions, such as an OnDataAvailable event or OnExternalInterrupt event.</p> <p>For OnDataAvailable events, eventParameter is used to determine the minimum number of samples to acquire during an analog input scan before generating the event. For OnExternalInterrupt events, eventParameter is used to latch digital bits on supported hardware by setting it to one of the constants in the "eventParameter parameter values" section on page 380.</p> <p>Most event conditions ignore this value.</p>
callbackFunc	<p>A delegate type that is the user-defined callback function to handle the above event type(s). A <i>delegate</i> is a data structure that refers either to a static method, or to a class instance and an instance method of that class.</p> <p>The callbackFunc needs the same parameters as the EventCallback delegate declaration. Refer to the "<a href="#">EventCallback delegate</a>" section on page 381 for proper syntax and return values.</p>
userData	Reference to user-defined data that is passed to the EventCallback delegate. This parameter is NOT de-referenced by the library or its drivers; as a consequence, a NULL pointer can be supplied.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**eventType parameter values:**

OnScanError	Generates an event upon detection of a driver error during Background input and output scans. This includes OverRun, UnderRun, and TooFew errors.
OnExternalInterrupt	For some digital and counter boards, generates an event, latches digital input data, or latches digital output data upon detection of a pulse at the External Interrupt pin.
OnPretrigger	For <a href="#">APretrig()</a> , generates an event upon detection of the first trigger.
OnDataAvailable	Generates an event whenever the number of samples acquired during an analog input scan increases by eventParameter samples or more. Note that for BlockIo scans, events will be generated on packet transfers; for example, even if EventParameter is set to 1, events will only be generated every packet-size worth of data (256 samples for the PCI-DAS1602) for aggregate rates greater than 1 kHz for the default <a href="#">AInScan()</a> mode.  For <a href="#">APretrig()</a> , the first event is not generated until a minimum of EventParameter samples after the pretrigger.
OnEndOfAiScan	Generates an event upon completion or fatal error of a <a href="#">AInScan()</a> or <a href="#">APretrig()</a> . This event is NOT generated when scans are aborted using <a href="#">StopBackground()</a> .
OnEndOfAoScan	Generates an event upon completion or fatal error of a <a href="#">AOutScan()</a> . This event is not generated when scans are aborted using <a href="#">StopBackground()</a> .

**eventParameter parameter values:**

LatchDI	Returns the data that was latched in at the most recent interrupt edge.
LatchDO	Latches out the data most recently written to the hardware.

**Notes:**

- [EnableEvent\(\)](#) cannot be called while any background operations ([AInScan\(\)](#), [APretrig\(\)](#), or [AOutScan\(\)](#)) are active. If a background operation is in progress when [EnableEvent\(\)](#) is called, [EnableEvent\(\)](#) will return the `AlreadyActive` error. You should perform a [StopBackground\(\)](#) before calling [EnableEvent\(\)](#).
- Events can be generated no faster than the user callback function can handle them. If an event type becomes multi-signaled before the event handler returns, events will be merged, such that the event handler is called once per event type, and the event handler is supplied with the event data corresponding to the latest event. In addition, if more than one event type becomes signaled, the event handler for each event type is called in the same order in which they are listed above.
- Events are generated while handling board-generated interrupts. As a consequence, using [StopBackground\(\)](#) to abort background operations will not generate `OnEndOfAoScan` or `OnEndOfAiScan` events. However, the event handlers can be called directly immediately after calling [StopBackground\(\)](#).

**Important**

In order to understand the functions, you must read the board -specific information found in the *Universal Library User's Guide* (available on our web site at [www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf](http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf)).

Review and run the example programs prior to attempting any programming of your own. Following this advice will save you hours of frustration, and possibly time wasted holding for technical support.

This note, which appears elsewhere, is especially applicable to this method. Read the board-specific information for your board (see the *Universal Library User's Guide*). We suggest that you make a copy of that page to refer to as you read this manual and examine the example programs.

## EngArrayToWinBuf()

Transfers a 2D array of engineering unit values to a Windows buffer as integer values.

The conversion from engineering unit values to integer values uses the D/A resolution of the board associated with the MccBoard object.

This method is usually used to obtain values compatible with the `AOutScan()` method or the `DaqOutScan()` method from a 2D array of engineering unit values, such as those provided by Measurement Studio signal generation methods. The converted values are transferred to the buffer based on the `gain`, `firstPoint`, `count`, and `numChannels` parameters.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**

```
Public Function EngArrayToWinBuf(ByVal gain As MccDaq.Range, ByVal
engUnits As Double(,), ByVal memHandle As Integer, ByVal firstPoint
As Integer, ByVal count As Integer, ByVal numChannels As Integer) As
MccDaq.ErrorInfo

Public Function EngArrayToWinBuf(ByVal gainArray As MccDaq.Range(),
ByVal gainCount As Integer, ByVal engUnits As Double(,), ByVal
memHandle As Integer, ByVal firstPoint As Integer, ByVal count As
Integer, ByVal numChannels As Integer) As MccDaq.ErrorInfo
```

**C# .NET:**

```
public MccDaq.ErrorInfo EngArrayToWinBuf(MccDaq.Range gain,
double[,] engUnits, int memHandle, int firstPoint, int count, int
numChannels)

public MccDaq.ErrorInfo EngArrayToWinBuf(MccDaq.Range[] gainArray,
int gainCount, double [,] engUnits, int memHandle, int firstPoint,
int count, int numChannels)
```

### Parameters:

<code>gain</code>	The range to use for converting the data. This range should be the same as the range specified for <code>AOutScan()</code> or <code>DaqOutScan()</code> .
<code>gainArray</code>	<p>The array containing the D/A range values used during the analog output scan.</p> <p>If a gain queue was not used for the scan, this array should only contain 1 element whose value matches the gain used during the scan. If a gain queue was used during the scan, this array should match the <code>gainArray</code> value used in <code>DaqOutScan()</code>.</p> <p>If the corresponding range in the <code>gainArray</code> is set to <code>NotUsed</code> (<code>MccDaq.Range.NotUsed</code>), engineering unit values are returned as integer values.</p>
<code>gainCount</code>	The number of array elements in <code>gainArray</code> . Set <code>gainCount</code> to 1 when no gain queue was used for the scan. If a gain queue was used for the scan, this number should match the number of gain queue pairs defined in <code>DaqOutScan()</code> .
<code>engUnits</code>	The array of data to convert to binary units and store in the windows memory buffer. With the <code>engUnits</code> array, the channel numbers are stored in the first dimension, and the number of points/channel is stored in the second dimension.
<code>memHandle</code>	The handle to the windows memory buffer that holds the binary data that is output. This value should be large enough to hold $(\text{count} \times \text{numChannels})$ samples.

<code>firstPoint</code>	The index into the windows memory buffer that will hold the first sample of the converted first channel. The index into the raw memory is ( <code>firstPoint</code> x <code>numChannels</code> ) so that converted data always starts with the first channel specified in the scan. For example, if <code>firstPoint</code> is 14 and the number of channels is 8, the index of the first converted sample is 112.
<code>count</code>	The number of samples per channel to convert from engineering units. Count should not exceed Windows buffer size / ( <code>numChannels</code> - <code>firstPoint</code> ).
<code>numChannels</code>	The number of channels of data stored in the existing array to be transferred.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**Notes:**

This method stores the samples specified by `firstPoint` in the windows memory buffer. Each sample is converted using the ranges set by `gain`.

If the corresponding range in the `gainArray` is set to `NotUsed`, engineering unit values are returned as integer values.

## EventCallback delegate

The EventCallback delegate is called as a parameter of the [EnableEvent\(\)](#) method. A delegate is a data structure that refers either to a static method, or to a class instance and an instance method of that class.

You create the data structure using the prototype shown below. You call the delegate by passing either its address or a pointer to the delegate to the `callbackFunc` parameter of the [EnableEvent\(\)](#) method.

### Delegate prototype:

**VB .NET:**                    `Public Sub MyCallback(ByVal BoardNum As Integer, ByVal EventType As MccDaq.EventType, ByVal EventData As UInt32, ByVal pUserData As System.IntPtr)`

**C# .NET:**                    `public delegate void EventCallback(int BoardNum, MccDaq.EventType EventType, uint EventData, IntPtr pUserData);`

### Parameters:

<code>BoardNum</code>	Indicates which board caused the event.
<code>EventType</code>	Indicates which event occurred.
<code>EventData</code>	Board-specific data associated with this event. Set it to one of the constants in the "EventData parameter values" section below.
<code>pUserData</code>	Pointer to or reference of data supplied by the <code>userData</code> parameter in the <a href="#">EnableEvent()</a> method. Note that before using this parameter value, it must be cast to the same data type as it was passed to <code>EnableEvent()</code> .

### Returns:

`pUserData` – Returns the value specified by the `userData` parameter in `EnableEvent()`.

### EventData parameter values:

<b>EventType</b>	<b>Value of EventData</b>
<code>OnScanError</code>	The <a href="#">Error code</a> of the scan error.
<code>OnExternalInterrupt</code>	The number of interrupts generated since enabling the <code>ON_EXTERNAL_INTERRUPT</code> event.
<code>OnPretrigger</code>	The number of pretrigger samples available at time of pretrigger. Value is invalid for some boards when a <code>TOOFEW</code> error occurs. See board details.
<code>OnDataAvailable</code>	The number of samples acquired since the start of scan.
<code>OnEndOfAiScan</code>	The total number of samples acquired upon scan completion or end.
<code>OnEndOfAoScan</code>	The total number of samples output upon scan completion or end.

## FlashLED()

Causes the LED on a USB device to flash.

Member of the [MccBoard](#) class.

### Function prototype:

VB .NET:	<code>Public Function FlashLED() As MccDaq.ErrorInfo</code>
C# .NET:	<code>public MccDaq.ErrorInfo FlashLED()</code>



## FromEngUnits()

Converts a single precision voltage (or current) value in engineering units to an integer count value. This method is typically used to obtain a data value from a voltage value for output to a D/A with methods such as `AOut()`.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**

```
Public Function FromEngUnits(ByVal range As MccDaq.Range , ByVal
engUnits As Single, ByRef dataVal As Short) As MccDaq.ErrorInfo

Public Function FromEngUnits(ByVal range As MccDaq.Range, ByVal engUnits As Single,
ByRef dataVal As System.UInt16) As MccDaq.ErrorInfo
```

**C# .NET:**

```
public MccDaq.ErrorInfo FromEngUnits(MccDaq.Range range, float
engUnits, out ushort dataVal)

public MccDaq.ErrorInfo FromEngUnits(MccDaq.Range range, float
engUnits, out short dataVal)
```

### Parameters:

<code>range</code>	The voltage (or current) range to use for the conversion to counts. When using this method to obtain a value to send to a D/A board, keep in mind that some D/A boards have programmable voltage ranges, and others set the voltage range via switches on the board. In either case, the desired range must be passed to this method. Refer to Table 7 on page 208 for a list of valid range settings.
<code>engUnits</code>	The single precision voltage (or current) value to use for the conversion to counts. Set the value to be within the range specified by the <code>range</code> parameter.
<code>dataVal</code>	Returns an integer count to this variable that is equivalent to the <code>engUnits</code> parameter using the resolution of the D/A on the board (if any).

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

`dataVal` – the integer count equivalent to `engUnits` is returned here.

### Note:

This method is not supported for hardware with resolution greater than 16 bits.

The default resolution of this method is 12 bits, so if the device has neither analog input nor analog output, the result will be a 12 bit conversion.

If the device has both analog input and analog output, the resolution and transfer function of the D/A converter on the device is used.

## GetBoardName()

Returns the board name of a specified board.

Member of the [MccService](#) class.

### Function prototype:

```
VB.NET:      Public Shared Function GetBoardName(ByVal boardNumber As Integer,
ByRef boardName As String) As MccDag.ErrorInfo
```

```
C# .NET:      public MccDag.ErrorInfo GetBoardName(int boardNumber, ref string
               boardName)
```

### Parameters:

boardNumber	Refers either to the board number associated with a board when it was installed, or GETFIRST or GETNEXT.
-------------	--

boardName	A null-terminated string variable that the board name is returned to. Refer to the Appendix, "Board Type Codes," in the <i>Universal Library User's Guide</i> (available on our web site at <a href="http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf">www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf</a> ).
-----------	---

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

boardName - return string containing the board name.

**Notes:**

There are two distinct ways of using this method:

- Pass a board number as the `BoardNum` argument. The string that is returned describes the board type of the installed board.
- Set `BoardNum` to `GETFIRST` or `GETNEXT` to get a list of all board types that are supported by the library. Set `BoardNum` to `GETFIRST` to get the first board type in the list of supported boards. Subsequent calls with `Board=GETNEXT` returns each of the other board types supported by the library. When you reach the end of the list, `BoardName` is set to an empty string. Refer to the `ULGT04` example program in the installation directory for more details.

## GetStatus()

Returns the status about the background operation currently running.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**                   Public Function GetStatus(ByRef status As Short, ByRef curCount As Integer, ByRef curIndex As Integer, ByVal functionType As MccDaq.FunctionType) As MccDaq.ErrorInfo

**C# .NET:**                   public MccDaq.ErrorInfo GetStatus(out short status, out int curCount, out int curIndex, MccDaq.FunctionType functionType)

### Parameters:

**status**                   status indicates whether or not a background process is currently executing.

**curCount**               The **curCount** parameter specifies how many points have been input or output since the Background process started. Use it to gauge how far along the operation is towards completion. Generally, **curCount** returns the total number of samples transferred between the DAQ board and the Windows data buffer at the time **GetStatus()** was called.

When you set both the **Continuous** and **Background** options, **curCount**'s behavior depends on the board model. Refer to the board-specific information in the *Universal Library User's Guide* for the behavior of your board.

With recent MCC DAQ designs, the **curCount** parameter continually increases in increments of the packet size as Windows' circular data buffer recycles, until it reaches  $2^{31}$ . Since the **count** parameter is a signed integer, at  $2,147,483,647 + 1$ , the **count** parameter rolls back to a negative number ( $-2,147,483,647$ ). The **count** parameter resumes incrementing, eventually reaching 0 and increasing back up to  $2,147,483,647$ .

The **curIndex** parameter is usually more useful than the **curCount** parameter in managing data collected when you set both the **Continuous** and **Background** options.

**curIndex**               The **curIndex** parameter is an index into the Windows data buffer. This index points to the start of the last completed channel scan that was transferred between the DAQ board and the Windows data buffer. If a scan is running but no points in the buffer have been transferred, **curIndex** equals -1 in most cases.

For **Continuous** operations, **curIndex** rolls over when the Windows data buffer is full. This rollover indicates that "new" data is now overwriting "old" data. Your goal is to process the old data before it gets overwritten. You can keep ahead of the data flow by copying the old data out of the buffer before new data overwrites it.

The **curIndex** parameter can help you access the most recently transferred data. Your application does not have to process the data exactly when it becomes available in the buffer – in fact, you should avoid doing so unless absolutely necessary. The **curIndex** parameter generally increments by the packet size, but in some cases the **curIndex** parameter can vary within the same scan. One instance of a variable increment is when the packet size is not evenly divisible by the number of channels.

You should determine the best size of the "chunks" of data that your application can most efficiently process, and then periodically check on the **curIndex** parameter value to determine when that amount of additional data has been transferred.

Refer to the *Universal Library User's Guide* for information on your board, particularly when using Pre-Trigger.

`functionType` Specifies which scan to retrieve status information about. Set it to one of the constants in the "functionType parameter values" section below.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

`Status`                      `Idle` - No background operation is running.  
                                  `Running` - Background operation still underway

`curCount` - current number of samples collected

`curIndex` - Current sample index

**functionType parameter values:**

`AiFunction`                Specifies analog input scans started with [AInScan\(\)](#) or [APretrig\(\)](#).

`AoFunction`                Specifies analog output scans started with [AOutScan\(\)](#).

`DiFunction`                Specifies digital input scans started with [DInScan\(\)](#).

`DoFunction`                Specifies digital output scans started with [DOutScan\(\)](#).

`CtrFunction`               Specifies counter background operations started with [CStoreOnInt\(\)](#) or [CInScan\(\)](#).

`DaqiFunction`             Specifies a synchronous input scan started with [DaqInScan\(\)](#).

`DaqoFunction`             Specifies a synchronous output scan started with [DaqOutScan\(\)](#).

## GetTCValues()

Converts raw thermocouple data from a Windows global memory buffer collected using the [DaqInScan\(\)](#) method to a one-dimensional or two-dimensional array of data on a temperature scale (Celsius, Fahrenheit or Kelvin).

Member of the [MccBoard](#) class.

### Function prototype:

VB .NET:

Copies data to a one-dimensional array of single values:

```
Public Function GetTCValues(ByVal chanArray As Short(), ByVal
chanTypeArray As MccDaq.ChannelType(), ByVal chanCount As Integer,
ByVal memHandle As Integer, ByVal firstPoint As Integer, ByVal count
As Integer, ByVal scale As MccDaq.TempScale, ByRef tempValArray As
Single) As MccDaq.ErrorInfo
```

Copies data to a two-dimensional array of double values:

```
Public Function GetTCValues(ByVal chanArray As Short(),ByVal
chanTypeArray As MccDaq.ChannelType(), ByVal chanCount As Integer,
ByVal memHandle As Integer, ByVal firstPoint As Integer,
ByVal count As Integer, ByVal scale As MccDaq.TempScale, ByRef
tempValArray As Double(,)) As MccDaq.ErrorInfo
```

C# .NET:

Copies data to a one-dimensional array of single values:

```
public ErrorInfo GetTCValues(short[] chanArray, MccDaq.ChannelType()
chanTypeArray, int chanCount, int memHandle, int firstPoint,
int count, MccDaq.TempScale scale, out float tempValArray)
```

Copies data to a two-dimensional array of double values:

```
public ErrorInfo GetTCValues(short[] chanArray, MccDaq.ChannelType
[] chanTypeArray, int chanCount, int memHandle, int firstPoint,
int count, MccDaq.TempScale scale, out double[,] tempValArray)
```

### Parameters:

chanArray	Array containing channel values. Valid channel values are analog and temperature input channels and digital ports. chanArray must match the channel array used with the <a href="#">DaqInScan()</a> method.
chanTypeArray	Array containing channel types. Each element of this array defines the type of the corresponding element in the chanArray. chanTypeArray must match the channel type settings used with the <a href="#">DaqInScan()</a> method.
chanCount	Number of elements in chanArray.
memHandle	The memory handle that was returned by <a href="#">WinBufAlloc()</a> when the buffer was allocated. The buffer should contain the data that you want to convert.
firstPoint	The index into the raw data memory buffer that holds the first sample of the first channel to be converted. The index into the raw memory is (firstPoint x chanCount) so that converted data always starts with the first channel specified in the scan. For example, if firstPoint is 14 and the number of channels is 8, the index of the first converted sample is 112.
count	The number of samples per channel to convert to engineering units. count should not exceed Windows buffer size / chanCount – firstPoint.
scale	Specifies the temperature scale that the input will be converted to. Choices are Celsius, Fahrenheit and Kelvin.

`tempValArray`      The array to hold the converted data. This array must be allocated by the user, and must be large enough to hold `count` samples x the number of temperature channels.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation

`tempValArray` – Converted data.

## InByte()

Reads a byte from a hardware register on a board.

Member of the [MccBoard](#) class.

### Function prototype:

VB .NET:           Public Function InByte(ByVal portNum As Integer) As Integer  
C# .NET:           public int InByte(int portNum)

### Parameters:

portNum           Register within the board. Boards are set to a particular base address. The registers on the boards are at addresses that are offsets from the base address of the board (BaseAdr + 0, BaseAdr + 2, etc).  
  
Set this parameter to the offset for the desired register. This method takes care of adding the base address to the offset, so that the board's address can be changed without changing the code.

### Returns:

The current value of the specified register

### Notes:

[InByte\(\)](#) is used to read 8 bit ports. [InWord\(\)](#) is used to read 16-bit ports.

This method was designed for use with ISA bus boards. Use with PCI bus boards is not recommended.

## InWord()

Reads a word from a hardware register on a board.

Member of the [MccBoard](#) class.

### Function prototype:

VB .NET:           Public Function InWord(ByVal portNum As Integer) As Integer  
C# .NET:           public int InWord(int portNum)

### Parameters:

portNum           Register within the board. Boards are set to a particular base address. The registers on the boards are at addresses that are offsets from the base address of the board (BaseAdr + 0, BaseAdr + 2, etc).  
  
Set this parameter to the offset for the desired register. This method takes care of adding the base address to the offset, so that the board's address can be changed without changing the code.

### Returns:

The current value of the specified register.

### Notes:

[InByte\(\)](#) is used to read 8-bit ports. `InWord()` is used to read 16-bit ports.

This method was designed for use with ISA bus boards. Use with PCI bus boards is not recommended.



## OutByte()

Writes a byte to a hardware register on a board.

Member of the [MccBoard](#) class.

### Function prototype:

VB .NET:                   Public Function OutByte(ByVal portNum As Integer, ByVal portVal As Integer) As MccDaq.ErrorInfo

C# .NET:                   public MccDaq.ErrorInfo OutByte(int portNum, int portVal)

### Parameters:

portNum                   Register within the board. Boards are set to a particular base address. The registers on the boards are at addresses that are offsets from the base address of the board (BaseAdr + 0, BaseAdr + 2, etc).

Set this parameter to the offset for the desired register. This method takes care of adding the base address to the offset, so that the board's address can be changed without changing the code.

portVal                   Value that is written to the register.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

### Notes:

[OutByte\(\)](#) is used to write to 8-bit ports. [OutWord\(\)](#) is used to write to 16-bit ports.

This method was designed for use with ISA bus boards. Use with PCI bus boards is not recommended.

## OutWord()

Writes a word to a hardware register on a board.

Member of the [MccBoard](#) class.

### Function prototype:

VB .NET:           Public Function OutWord(ByVal portNum As Integer, ByVal portVal As Integer) As MccDaq.ErrorInfo

C# .NET:           public MccDaq.ErrorInfo OutWord(int portNum, int portVal)

### Parameters:

portNum           Register within the board. Boards are set to a particular base address. The registers on the boards are at addresses that are offsets from the base address of the board (BaseAdr + 0, BaseAdr + 2, etc).

Set this parameter to the offset for the desired register. This method takes care of adding the base address to the offset, so that the board's address can be changed without changing the code.

PortVal           Value that is written to the register.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

### Notes:

[OutByte\(\)](#) is used to write to 8-bit ports. `OutWord()` is used to write to 16-bit ports.

This method was designed for use with ISA bus boards. Use with PCI bus boards is not recommended.

## RS485()

Sets the direction of RS-485 communications port buffers.

Member of the [MccBoard](#) class.

### Function prototype:

VB .NET:           Public Function RS485(ByVal transmit As MccDaq.OptionState , ByVal receive As MccDaq.OptionState ) As MccDaq.ErrorInfo

C# .NET:           public MccDaq.ErrorInfo RS485(MccDaq.OptionState transmit, MccDaq.OptionState receive)

### Parameters:

transmit	Set to Enabled or Disabled. The transmit RS-485 line driver is turned on. Data written to the RS-485 UART chip is transmitted to the cable connected to that port.
receive	Set to MccDaq.OptionState.Enabled or MccDaq.OptionState.Disabled. The receive RS-485 buffer is turned on. Data present on the cable connected to the RS-485 port is received by the UART chip.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

### Notes:

You can simultaneously enable or disable the transmit and receive buffers. If both are enabled, data written to the port is also received by the port. For a complete discussion of RS485 network construction and communication, refer to the CIO-COM485 or PCM-COM485 hardware manual.

## StopBackground()

Stops one or more subsystem background operations that are in progress for the specified board. Use this method to stop any method that is running in the background. This includes any method that was started with the `Background` option, as well as [CStoreOnInt\(\)](#) (which always runs in the background).

Execute [StopBackground\(\)](#) after normal termination of all background functions to clear variables and flags.

Member of the [MccBoard](#) class.

### Function prototype:

VB .NET:           Public Function StopBackground(ByVal funcType As  
                          MccDaq.FunctionType) As MccDaq.ErrorInfo

C# .NET:           public MccDaq.ErrorInfo StopBackground(MccDaq.FunctionType funcType)

### Parameters:

functionType	Specifies which background operation to stop. Set it to one of the constants in the "functionType parameter values" section below.
--------------	--

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**functionType** parameter values:

AiFunction:	Specifies analog input scans started with <a href="#">AInScan()</a> or <a href="#">APretrig()</a> .
AoFunction	Specifies analog output scans started with <a href="#">AOutScan()</a> .
DiFunction	Specifies digital input scans started with <a href="#">DInScan()</a> .
DoFunction	Specifies digital output scans started with <a href="#">DOutScan()</a> .
CtrFunction	Specifies counter background operations started with <a href="#">CStoreOnInt()</a> or <a href="#">CInScan()</a> .
DaqiFunction	Specifies a synchronous input scan started with <a href="#">DaqInScan()</a> .
DaqoFunction	Specifies a synchronous output scan started with <a href="#">DaqOutScan()</a> .

## ToEngUnits()

Converts an integer count value to an equivalent single precision voltage (or current) value. This method is typically used to obtain a voltage value from data received from an A/D with methods such as `AIN()`.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**

```
Public Function ToEngUnits(ByVal range As MccDaq.Range , ByVal
dataVal As Short, ByRef engUnits As Single) As MccDaq.ErrorInfo
Public Function ToEngUnits(ByVal range As MccDaq.Range, ByVal
dataVal As System.UInt16, ByRef engUnits As Single) As
MccDaq.ErrorInfo
```

**C# .NET:**

```
Public MccDaq.ErrorInfo ToEngUnits(MccDaq.Range range, ushort
dataVal, out float engUnits)
Public MccDaq.ErrorInfo ToEngUnits(MccDaq.Range range, short
dataVal, out float engUnits)
```

### Parameters:

<code>range</code>	Voltage (or current) range to use for the conversion to engineering units. When using this method to obtain engineering units from a value received from an A/D board, keep in mind that some A/D boards have programmable voltage ranges, and others set the voltage range via switches on the board. In either case, the desired range must be passed to this method. Refer to Table 7 on page 208 for a list of valid range settings.
<code>dataVal</code>	An integer count value (typically, one returned from an A/D board).
<code>engUnits</code>	The single precision voltage (or current) value that is equivalent to <code>dataVal</code> is returned to this variable. The value will be within the range specified by the <code>range</code> parameter.

### Returns:

An [ErrorInfo](#) object that indicates the status of the operation.

`engUnits` – the engineering units value equivalent to `dataVal` is returned to this variable.

### Note:

This method is not supported for hardware with resolution greater than 16 bits.

The default resolution of this method is 12 bits, so if the device has neither analog input nor analog output, the result will be a 12 bit conversion.

If the device has both analog input and analog output, the resolution and transfer function of the D/A converter on the device is used.

## WinBufToEngArray()

Transfers integer values from a Windows buffer to a 2D array as engineering unit values.

The conversion from integer values to engineering unit values uses the A/D resolution of the board associated with the MccBoard object.

This method is usually used to obtain values compatible with those required by Measurement Studio waveform display controls from a Windows buffer containing data from a method such as AInScan() or DaqInScan().

The converted values are transferred to the 2D array based on the gain, firstPoint, count, and numChannels parameters.

Member of the [MccBoard](#) class.

### Function prototype:

**VB .NET:**

```
Public Function WinBufToEngArray(ByVal gain As MccDaq.Range, ByVal memHandle As Integer, ByVal engUnits As Double(), ByVal firstPoint As Integer, ByVal count As Integer, ByVal numChannels As Integer) As MccDaq.ErrorInfo
```

```
Public Function WinBufToEngArray(ByVal gainArray As MccDaq.Range(), ByVal gainCount As Integer, ByVal memHandle As Integer, ByVal engUnits As Double(), ByVal firstPoint As Integer, ByVal count As Integer, ByVal numChannels As Integer) As MccDaq.ErrorInfo
```

**C# .NET:**

```
public MccDaq.ErrorInfo WinBufToEngArray(MccDaq.Range gain, int memHandle, double[,] engUnits, int firstPoint, int count, int numChannels)

public MccDaq.ErrorInfo WinBufToEngArray(MccDaq.Range gainArray, int gainCount, int memHandle, double[,] engUnits, int firstPoint, int count, int numChannels)
```

### Parameters:

gain	The range to use for converting scan data. This range should be the same as the range specified for AInScan() or DaqInScan().
gainArray	The array containing the A/D range values used during the analog input scan.  If a gain queue was not used for the scan, this array should only contain 1 element whose value matches the gain used during the scan. If a gain queue was used during the scan, this array should match the gainArray value used in ALoadQueue() or DaqInScan().  If the corresponding range in the gainArray is set to NotUsed (MccDaq.Range.NotUsed), raw data is returned in engineering unit values.
gainCount	The number of array elements in gainArray. Set gainCount to 1 when no gain queue was used for the scan. If a gain queue was used for the scan, this number should match the number of gain queue pairs defined in ALoadQueue() or DaqInScan().
memHandle	The handle to the memory buffer holding the raw data to be converted to engineering units. This value should be large enough to hold (count x numChannels) samples.

<code>engUnits</code>	The array to hold the converted data. This array must be allocated by the user, and must be large enough to hold <code>count</code> samples. The first dimension should be the number of channels. The second dimension should equal the number of points/channel.
<code>firstPoint</code>	The index into the raw data memory buffer that holds the first sample of the first channel to be converted. The index into the raw memory is ( <code>firstPoint</code> x <code>numChannels</code> ) so that converted data always starts with the first channel specified in the scan. For example, if <code>firstPoint</code> is 14 and the number of channels is 8, the index of the first converted sample is 112.
<code>count</code>	The number of samples per channel to convert to engineering units. <code>count</code> should not exceed Windows buffer size / ( <code>numChannels</code> - <code>firstPoint</code> ).
<code>numChannels</code>	The number of channels of data stored in the existing array to be transferred.

**Returns:**

An [ErrorInfo](#) object that indicates the status of the operation.

**Notes:**

If `gainCount` is greater than one, the conversions cycle through the array until `count` samples have been converted. When only one gain is specified, that gain is applied to all conversions. Data is returned in engineering unit values as a two-dimensional array.

If the corresponding range in the `gainArray` is set to `NotUsed`, raw data is returned in engineering unit values.

# **Appendix**



---

# Error Codes

The following table lists error codes that are returned when running Universal Library or Universal Library for .NET.

Universal Library .NET errors can be referenced from the [MccDaq.ErrorInfo.Message](#) property.

Each entry in the list has four parts: the error code number, its symbolic name, its error message, and an explanation. Both the Universal Library function and its Universal Library .NET equivalent method are referred to when appropriate. Error code and error messages are identical for both programming libraries. The only difference in the error names used by each library is the case—the Universal Library error names are all uppercase (for example `NOERRORS`), while the Universal Library for .NET error names are mixed case (for example `NOERRORS`).

Error number	Error name	Error message
0	<b>NOERRORS</b>	No error has occurred The function executed successfully.
1	<b>BADBOARD</b>	Invalid board number The <code>BoardNum</code> argument that was specified does not match any of the boards that are listed in the configuration file. Run the configuration program to check which board numbers are configured.
2	<b>DEADDIGITALDEV</b>	Digital device is not responding - is base address correct? The digital device on the specified board is not responding. Either the board was installed incorrectly or the board is defective. Run the configuration program and make sure that the correct board was installed.
3	<b>DEADCOUNTERDEV</b>	Counter device is not responding - is base address correct? The counter device on the specified board is not responding. Either the board was installed incorrectly or the board is defective. Run the configuration program and make sure that the correct board was installed.
4	<b>DEADDADEV</b>	D/A is not responding - is base address correct? The D/A device on the specified board is not responding. Either the board was installed incorrectly or the board is defective. Run the configuration program and make sure that the correct board was installed.
5	<b>DEADADDEV</b>	A/D is not responding - is base address correct? The A/D device on the specified board is not responding. Either the board was installed incorrectly or the board is defective. Run the configuration program and make sure that the correct board was installed.
6	<b>NOTDIGITALCONF</b>	Selected board does not have digital I/O A digital I/O function or method was called with a board number that referred to a board that does not support digital I/O. Run the configuration program to see which type of board that board number refers to.
7	<b>NOTCOUNTERCONF</b>	Selected board does not have a counter A counter function or method was called with a board number that referred to a board that does not have a counter. Run the configuration program to see which type of board that board number refers to.

Error number	Error name	Error message
8	<b>NOTDACONF</b>	Selected board does not have a D/A An analog output function or method was called with a board number that referred to a board that does not have an analog output (D/A). Run the configuration program to see which type of board the board number refers to.
9	<b>NOTADCONF</b>	Selected board does not have an A/D An analog input function or method was called with a board number that referred to a board that does not have an analog input (A/D). Run the configuration program to see which type of board that board number refers to.
10	<b>NOTMUXCONF</b>	Selected board does not have thermocouple inputs A thermocouple input function or method was called with a board number that does not support thermocouple inputs or is not connected to an EXP board. Run the configuration program to view/change the board configuration.
11	<b>BADPORTNUM</b>	Invalid digital port number The port number specified for a digital I/O function or method does not exist on the specified board.
12	<b>BADCOUNTERDEVNUM</b>	Invalid counter device The <code>CounterNum</code> argument specified for a counter function or method references a counter that does not exist on the specified board.
13	<b>BADDADEVNUM</b>	Invalid D/A device The D/A channel that was specified for an analog output function or method does not exist on the specified board.
14	<b>BADSAMPLEMODE</b>	Invalid sample mode A sample mode that is not supported on this board ( <code>SINGLEIO</code> , <code>DMAIO</code> or <code>BLOCKIO</code> ) was specified in the <code>Options</code> argument. Try running the function or method without setting any of the Sample Mode options.
15	<b>BADINT</b>	Board configured for invalid interrupt level No interrupt was selected in <code>InstaCal</code> and one is required, or the board is set for "compatible mode" and the interrupt level selected is not supported in this mode. Interrupts above 7 are not valid in compatible mode. Either change the switch setting on the board to "enhanced mode", or change the interrupt level with the configuration program to something less than 8.
16	<b>BADADCHAN</b>	Invalid A/D channel number An invalid channel argument was passed to an analog input function or method . The range of valid channel numbers depends on which A/D board you are using - refer to the board manual. For some boards it also depends on how the board is configured (with a switch). For those boards run the configuration program and check how many channels the board is configured for.
17	<b>BADCOUNT</b>	Invalid count An invalid <code>Count</code> argument was specified to a function or method . If this error occurs during <code>cbAInScan()/AInScan()</code> , increasing the <code>Count</code> should correct the problem. For boards using <code>DMAIO</code> , adjust the data buffer and <code>Count</code> above $(\text{HighChan} - \text{LowChan} + 1) * \text{Rate} / 100$ for <code>CONTINUOUS</code> mode scans. However, those boards using <code>BLOCKIO</code> , require a user buffer and <code>Count</code> large enough to hold at least one half FIFO worth of samples (typically, 512 samples) for <code>CONTINUOUS</code> mode scans.

Error number	Error name	Error message
18	<b>BADCNTRCONFIG</b>	Invalid counter configuration specified An invalid Config argument was passed to <code>cbC8254Config()</code> / <code>C8254Config()</code> . The only legal values are <code>HIGHONLASTCOUNT</code> , <code>ONESHOT</code> , <code>RATEGENERATOR</code> , <code>SQUAREWAVE</code> , <code>SOFTWARESTROBE</code> and <code>HARDWARESTROBE</code> .
19	<b>BADDAVAL</b>	Invalid D/A value An invalid D/A value was passed as an argument to an analog output function or method. The only legal values are 0 to 4095 for 12-bit boards or 0 to 65,535 for 16-bit boards (see the "Note on Basic signed integers" at the beginning of the "Counter Boards" chapter in the <i>Universal Library User's Guide</i> available on our web site at <a href="http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf">www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf</a> ).
20	<b>BADDACHAN</b>	Invalid D/A channel number An invalid D/A channel was passed as an argument to an analog output function or method. The legal range of values depends on which D/A board you are using. Refer to the board manual to find how many D/A channels it has.
22	<b>ALREADYACTIVE</b>	Background operation already in progress An attempt was made to start a second background process on the same board before the first one had completed. Background processes are started whenever the <code>BACKGROUND</code> option is used by <code>cbCStoreOnInt()</code> / <code>CStoreOnInt()</code> . To stop a background operation, call <code>cbStopBackground()</code> / <code>StopBackground()</code> . To wait for a background process to complete call <code>cbGetStatus()</code> / <code>GetStatus()</code> and wait for <code>Status = IDLE</code> .
23	<b>PAGEOVERRUN</b>	DMA transfer crossed page boundary, may have gaps in data When a DMA transfer crosses a 64K memory page boundary on boards without FIFO buffers, there may be a small gap (missing samples) in the data. For applications requiring high speed transfers of greater than 32K samples, please select a board with a FIFO buffer. For boards without, check the data for gaps and do not specify rates over that at which gapless data may be taken. This is system-specific so you must determine the rate by experimentation.
24	<b>BADRATE</b>	Invalid sampling rate Invalid sampling rate argument was specified. The rate was either zero, a negative number or it was higher than the selected board supports. Refer to board-specific information for board maximum rates.
25	<b>COMPATMODE</b>	Board switches set for Compatible mode An operation was attempted that is not possible when the board's switch is set for 'compatible' operation. The most likely causes are due to using the <code>BLOCKIO</code> option or the pre-triggering functions. Either turn off the 'compatible' mode switch on the board or don't use the <code>BLOCKIO</code> option or the pre-triggering functions.
26	<b>TRIGSTATE</b>	Incorrect initial trigger state - trigger must start at TTL low Boards that use "polled gate" triggering require that the trigger be "off" when a pre-trigger function is first called. It then waits for the trigger signal. Make sure that the Trigger Input line (usually D0) is held at TTL low before calling the pre-trigger function.
27	<b>ADSTATUSHUNG</b>	A/D is not responding The A/D board is not responding as it should. Usually indicates some kind of hardware problem - either defective hardware or more than one board at the same base address.

Error number	Error name	Error message
28	<b>TOOFEW</b>	<p>Trigger occurred before the requested number of samples were collected</p> <p>A pre-trigger function or method was called and the trigger signal occurred before the requested number of samples could be collected. This is only a warning message. The function or method continued anyway. The data that was returned to the array will contain fewer than the expected number of points. The function or method will return the actual number of pre-trigger points and the total number of points. You can use these two values to find your way around the data in the array.</p>
29	<b>OVERRUN</b>	<p>Data overrun - data was lost</p> <p>Data was lost during an analog input because the computer could not keep up with the A/D sampling rate. This typically can only happen with the file input functions or methods, or by using <code>SINGLEIO</code> mode. Possible solutions include lowering the sampling rate, defragmenting the "streamer" file, switching to a RAM disk, or lowering the count.</p>
30	<b>BADRANGE</b>	<p>Invalid voltage or current range</p> <p>Invalid Range argument was specified to an analog input or output function or method . The board does not support the gain you specified. Refer to board-specific information for a list of allowable ranges.</p>
31	<b>NOPROGGAIN</b>	<p>This A/D board does not have programmable gain</p> <p>Invalid Range argument was passed to an analog input function or method . The selected board does not support programmable gains so the only valid Range argument is 0. (This argument is ignored for these board types in later versions of the library.)</p>
32	<b>BADFILENAME</b>	<p>Specified file name is not valid</p> <p>The FileName argument that was passed to a file function or method is not valid. It is either an empty string or a NULL pointer.</p>
33	<b>DISKISFULL</b>	<p>Disk is full, could not complete operation</p> <p>A file operation failed before completing because the disk that it was writing to is full. Try erasing some files from the disk. If this error occurred during either <code>cbFileAInScan() / FileAInScan()</code> or <code>cbFilePretrig() / FilePretrig()</code>, it indicates another problem. The disk space for these commands should have been previously allocated with the MAKESTRM.EXE program. If this error is generated when data is being collected it indicates that you did not allocate a large enough file with MAKESTRM.EXE.</p>
34	<b>COMPATWARN</b>	<p>Board switch set to compatible mode - sampling speed may be limited</p> <p>The board's switch is set for "compatible mode." When in "compatible mode," <code>BLOCKIO</code> transfers are not possible. <code>BLOCKIO</code> sampling was specified but it has automatically been changed to <code>DMAIO</code> transfers. The maximum sampling rate will be limited to the maximum rate for DMA transfers. Change the "compatible mode" switch on the board if you want to use <code>BLOCKIO</code> transfers.</p>
35	<b>BADPOINTER</b>	<p>Pointer is not valid</p> <p>An invalid (NULL) pointer was passed as an argument/parameter to a function or method .</p>
37	<b>RATEWARNING</b>	<p>Sample rate may be too fast for <code>SINGLEIO</code> mode</p> <p>The specified sampling rate MAY be too high. The maximum allowable sampling rate depends very much on the computer that the program is running on. This warning is generated based on the slowest CPU speed. Your computer may be able to sustain faster rates, but, you should expect the computer to lock up (fail to respond to keyboard input) if you do exceed the sampling rate your computer can sustain.</p>

Error number	Error name	Error message
38	<b>CONVERTDMA</b>	CONVERTDATA cannot be used with DMAIO and BACKGROUND  The CONVERTDATA and BACKGROUND options can not be used together when the board is transferring data via DMA. Possible solutions include: Use <code>cbAConvertData ()/AConvertData ()</code> to convert the data after it is collected. Don't use BACKGROUND option. Use BLOCKIO option if your A/D board supports it. Use SINGLEIO option if your computer is fast enough to support the selected sampling rate.
39	<b>DTCONNECTERR</b>	Board does not support DTCONNECT option The DTCONNECT Option was passed to an analog input function or method . The selected board does not support that option.
40	<b>FORECONTINUOUS</b>	CONTINUOUS can only be run with BACKGROUND The CONTINUOUS option was passed to a function or method without also setting the BACKGROUND option. This is not allowed. Any time you set the CONTINUOUS option you must also set the BACKGROUND option.
41	<b>BADBOARDTYPE</b>	This function or method can not be used with this board An attempt was made to call a function or method for a board that does not support that function or method .
42	<b>WRONGDIGCONFIG</b>	Digital port not configured correctly for requested operation  Some of the digital bits or ports ( <code>FIRSTPORTA - EIGHTHPORTCH</code> ) must be configured as inputs OR outputs but not both. An attempt was made to use a digital input function or method on a port or bit that was configured as an output or vice versa. Use <code>cbDConfigPort ()/DConfigPort ()</code> or <code>cbDConfigBit ()/DConfigBit ()</code> to switch a port's (or bit's) direction. If the board you are using contains configurable port types and you do not call <code>cbDConfigPort ()/DConfigPort ()</code> or <code>cbDConfigBit ()/DConfigBit ()</code> in your program, then all of the configurable ports will be in an unknown state (input or output).
43	<b>NOTCONFIGURABLE</b>	This digital port is not configurable (it's an In/Out port) <code>cbDConfigPort ()/DConfigPort ()</code> or <code>cbDConfigBit ()/DConfigBit ()</code> was called for a port that is not configurable. Check the <code>PortNum</code> argument passed to <code>cbDConfigPort ()</code> and make sure that it is in the range <code>FIRSTPORTA - EIGHTHPORTCH</code> . If <code>PortNum</code> is <code>AUXPORT</code> , make sure your hardware supports configuration of this port type. If not then there is no need to call this function or method .
44	<b>BADPORTCONFIG</b>	Invalid digital port configuration The <code>Direction</code> argument passed to <code>cbDConfigPort ()/DConfigPort ()</code> or <code>cbDConfigBit ()/DConfigBit ()</code> is invalid. It must be set to either <code>DIGITALIN</code> or <code>DIGITALOUT</code> .
45	<b>BADFIRSTPOINT</b>	<code>FirstPoint</code> number is not valid The <code>FirstPoint</code> argument to <code>cbFileRead ()/FileRead ()</code> is invalid. It is either a negative number or it is larger then the number of points in the file.
46	<b>ENDOFFILE</b>	Attempted to read past the end of the file <code>cbFileRead ()/FileRead ()</code> attempted to read beyond the end of the file. Check the file length with <code>cbFileGetInfo ()/FileGetInfo ()</code> and make sure that the <code>FirstPoint</code> and <code>Count</code> arguments to <code>cbFileRead ()/FileRead ()</code> are correct for that file length.
47	<b>NOT8254CTR</b>	This board does not have an 8254 counter <code>cbC8254Config ()/C8254Config ()</code> was called for a board that has a counter but not an 8254 counter. This function or method can only be used with an 8254 counter.

Error number	Error name	Error message
48	<b>NOT9513CTR</b>	This board does not have a 9513 counter cbC9513Config()/C9513Config() was called for a board that has a counter but not a 9513 counter. This function or method can only be used with an 9513 counter.
49	<b>BADTRIGTYPE</b>	Invalid TrigType cbATrig()/(ATrig()) was called with an invalid TrigType argument. It must be set to either TRIGABOVE or TRIGBELOW.
50	<b>BADTRIGVALUE</b>	Invalid TrigValue cbATrig()/(ATrig()) was called with an invalid TrigValue argument. It must be in the range 0 to 4095 for 12-bit boards or 0 to 65535 for 16-bit boards (see the "Note on Basic signed integers" at the beginning of the "Counter Boards" chapter in the <i>Universal Library User's Guide</i> , available on our web site at <a href="http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf">www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf</a> ).
52	<b>BADOPTION</b>	Invalid Option specified for this function or method The Options argument contains an option that is not valid for this function or method .
53	<b>BADPRETRIGCOUNT</b>	Invalid PretrigCount specified Either cbAPretrig()/APretrig() or cbFilePretrig()/FilePretrig() was called with an invalid PretrigCount argument. The pre-trigger count must not be < 0 and must be less than TotalCount-512. It also must be less than 32k for cbAPretrig()/APretrig() and less than 16k for cbFilePretrig()/FilePretrig().
55	<b>BADDIVIDER</b>	Invalid FOutDivider value The FOutDivider argument to cbC9513Init() (C9513Init()) is not valid. It must be in the range 0 to 15.
56	<b>BADSOURCE</b>	Invalid FOutSource value The FOutSource argument to cbC9513Init() (C9513Init()) is not valid. It must be one of the following values CTRINPUT1, CTRINPUT2, CTRINPUT3, CTRINPUT4, CTRINPUT5, GATE1, GATE2, GATE3, GATE4, GATE5, FREQ1, FREQ2, FREQ3, FREQ4, FREQ5 (for example 0 to 15).
57	<b>BADCOMPARE</b>	Invalid Compare value One or both of the compare arguments to cbC9513Init()/C9513Init() are not valid. They must be set to (CB)ENABLED or (CB)DISABLED (1 or 0).
58	<b>BADTIMEOFDAY</b>	Invalid TimeOfDay value The TimeOfDay argument to cbC9513Init()/C9513Init() is not valid. It must be set to either (CB)ENABLED or (CB)DISABLED (1 or 0).
59	<b>BADGATEINTERVAL</b>	Invalid GateInterval value The GateInterval argument to cbCFreqIn()/CFreqIn() is not valid. It must be greater than 0.
60	<b>BADGATECNTRL</b>	Invalid GateControl value The GateControl argument to cbC9513Config()/C9513Config() is not valid. It must be in the range 0 to 7.
61	<b>BADCOUNTEREDGE</b>	Invalid CounterEdge value The CounterEdge argument to cbC9513Config()/C9513Config() is not valid. It must be set to either POSITIVEEDGE or NEGATIVEEDGE.
62	<b>BADSPCLGATE</b>	Invalid SpecialGate value The SpecialGate argument to cbC9513Config()/C9513Config() is not valid. It must be set to either (CB)ENABLED or (CB)DISABLED (1 or 0).

Error number	Error name	Error message
63	<b>BADRELOAD</b>	Invalid Reload value The Reload argument to <code>cbC9513Config()</code> ( <code>C9513Config()</code> ) is not valid. It must be set to either <code>LOADREG</code> or <code>LOADANDHOLDREG</code> .
64	<b>BADRECYCLEFLAG</b>	Invalid RecycleMode value The RecycleMode argument to <code>cbC9513Config()/C9513Config()</code> is not valid. It must be set to either <code>(CB)ENABLED</code> or <code>(CB)DISABLED (1 or 0)</code> .
65	<b>BADBCDFLAG</b>	Invalid BCDMode value The BCDMode argument to <code>cbC9513Config()/C9513Config()</code> is not valid. It must be set to either <code>(CB)ENABLED</code> or <code>(CB)DISABLED (1 or 0)</code> .
66	<b>BADDIRECTION</b>	Invalid CountDirection value The CountDirection argument to <code>cbC9513Config()</code> ( <code>C9513Config()</code> ) is not valid. It must be set to either <code>COUNTUP</code> or <code>COUNTDOWN</code> .
67	<b>BADOUTCONTROL</b>	Invalid OutputControl value The OutputControl argument to <code>cbC9513Config()</code> ( <code>C9513Config()</code> ) is not valid. It must be set to either <code>ALWAYSLOW</code> , <code>HIGHPULSEONTC</code> , <code>TOGGLEONTC</code> , <code>DISCONNECTED</code> or <code>LOWPULSEONTC</code> .
68	<b>BADBITNUMBER</b>	Invalid BitNum specified The BitNum argument to <code>cbDBitIn()</code> or <code>cbDBitOut()</code> ( <code>DBitIn()</code> or <code>DBitOut()</code> ) is not valid. The valid range of bit numbers depends on the selected board. If it is a DIO24 compatible board the maximum bit number is 23. If it's a DIO96, the maximum bit number is 95. Refer to board-specific information in the <i>Universal Library User's Guide</i> or in your hardware manual.
69	<b>NONEENABLED</b>	None of the counter channels were enabled None of the counter channels were marked as <code>(CB)ENABLED</code> in the <code>CntrControl</code> array that was passed to <code>cbCStoreOnInt()/CStoreOnInt()</code> . At least one of the counter channels must be enabled.
70	<b>BADCTRCONTROL</b>	An element of <code>CntrControl</code> array not set to <code>DISABLED</code> or <code>ENABLED</code> One of the elements of the <code>CntrControl</code> array that was passed to <code>cbCStoreOnInt()/CStoreOnInt()</code> was set to something other than <code>(CB)ENABLED</code> or <code>(CB)DISABLED</code> . The array must have at least ten elements and the first ten elements must be set to either <code>(CB)ENABLED</code> or <code>(CB)DISABLED</code> .
71	<b>BADEXPCHAN</b>	Invalid EXP channel specified An invalid channel was passed to one of the thermocouple input commands. The channel number when using an EXP board must be $\geq 16$ . The maximum allowable channel number depends on which EXP board is being used (and how many of them). Refer to the board manual to find the number of channels.
72	<b>WRONGADRANGE</b>	Board set to wrong A/D range for reading thermocouples A thermocouple input function or method was called to read an EXP board input. The EXP board is connected to an A/D board with hardware selected gain that is set to the wrong range. When using EXP boards with thermocouples, the A/D must be set to the -5 to +5 volt range when available. When using RTD sensors, the range is 0 to 10 V when available.
73	<b>OUTOFRANGE</b>	Temperature input is out of range A thermocouple input function or method returned an invalid temperature. This usually indicates an open connection in the thermocouple or its connection to the mux board.

Error number	Error name	Error message
74	<b>BADTEMPSCALE</b>	Invalid temperature scale specified The <code>Scale</code> argument/parameter to a thermocouple input function or method is not valid. It must be set to either <code>CELSIUS</code> , <code>FAHRENHEIT</code> , <code>KELVIN</code> , or <code>VOLT</code> .
76	<b>NOQUEUE</b>	Specified board does not have channel/gain queue The function or method that was called requires that the board has a channel/gain queue. The specified board does not have a queue.
77	<b>CONTINUOUSCOUNT</b>	Count must be > packet size to use Continuous mode The <code>Count</code> argument is not valid for continuous mode. Using <code>BLOCKIO</code> mode, the <code>Count</code> argument must be large enough to cause at least one interrupt. This is usually half the size of the boards FIFO (typical sizes are 256, 512, and 1024). Refer to the board-specific information in the <i>Universal Library User's Guide</i> , available on our web site at <a href="http://www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf">www.mccdaq.com/PDFmanuals/sm-ul-user-guide.pdf</a> or in your hardware manual.
78	<b>UNDERRUN</b>	D/A FIFO went empty during output The specified D/A output rate could not be sustained. Try increasing the size of the data buffer or reducing the update rate to eliminate the error.
79	<b>BADMEMMODE</b>	Invalid memory mode specified The memory mode that was selected with <code>cbMemSetDTMode()</code> ( <code>MemSetDTMode()</code> ) is not one of the valid modes.
80	<b>FREQOVERRUN</b>	Measured frequency too high for selected gating interval The <code>GateInterval</code> argument used with <code>cbCFreqIn()</code> ( <code>CFreqIn()</code> ) is too large to measure the frequency of the signal connected to the counter. The counter is overflowing. Decrease the gating interval to eliminate the error.
81	<b>NOCJCCHAN</b>	A CJC Channel must be configured to make temperature measurements When the board was installed (with the <i>InstaCal</i> installation program) no Cold Junction Compression (CJC) channel was selected. To use the temperature measurement functions or methods with thermocouples, you must first select a CJC channel on the A/D board and then rerun the installation program.
82	<b>BADCHIPNUM</b>	Invalid ChipNum specified An invalid <code>ChipNum</code> argument was used with <code>cbC9513Init()</code> / <code>C9513Init()</code> . If the board is CTR05, set <code>ChipNum</code> to 0. If the board is a CTR10, set <code>ChipNum</code> to either 0 or 1.
83	<b>DIGNOTENABLED</b>	The digital I/O on this board is not enabled When the board was installed (with the <i>InstaCal</i> installation program), the expansion digital I/O was set to <code>DISABLED</code> . To use these digital I/O lines, you must enable the digital I/O on the board (with a jumper) and then re-run the installation program and set the digital I/O to <code>ENABLED</code> .
84	<b>CONVERT16BITS</b>	<code>CONVERTDATA</code> option can not be used with 16-bit A/D converters When using a 16-bit A/D (DAS1600/16), if you try to use the <code>CONVERTDATA</code> option with <code>cbAInScan()</code> / <code>AInScan()</code> or call <code>cbAConvertData()</code> / <code>AConvertData()</code> , this error is returned. (This has been updated so that it is ignored for boards for which it is inappropriate in later versions of the library.)
85	<b>NOMEMBOARD</b>	The <code>EXTMEMORY</code> option requires that a MEGA-FIFO be attached Attempt to use a <code>cbMem_()</code> function or <code>Mem_()</code> method without a MEGA-FIFO board installed. Install MEGA-FIFO through <i>InstaCal</i> .



Error number	Error name	Error message
86	<b>DTACTIVE</b>	No memory read/write allowed while DT transfer in progress A read or write to a memory board was attempted while data was being transferred via DT-Connect.
87	<b>NOTMEMCONF</b>	Specified board is not a memory board The specified board is not a memory board. This function or method only works with memory boards.
88	<b>ODDCHAN</b>	The first channel in scan and number of channels must be even (0, 2, 4, etc) Some boards use a channel/gain queue that require the first channel in the queue and the number of channels in the queue always be an even channel. This error can occur even when you are not in the process of loading the queue. Some boards use the queue automatically with <code>cbAInScan()</code> / <code>AInScan()</code> . On those boards, the low channel must be an even number.
89	<b>CTRNOINIT</b>	Counter was not configured or initialized You attempted to use <code>cbCLoad()</code> or <code>cbCIn()</code> ( <code>CLoad()</code> or <code>CIn()</code> ) before initializing and configuring the counter.
90	<b>NOT8536CTR</b>	This board does not have an 8536 counter chip Attempt to use 8536 initialization or configuration on board without 8536 chip.
91	<b>FREERUNNING</b>	Board doesn't time A/D sampling. Collecting at fastest possible speed This board does not have an A/D pacer mechanism and you have called <code>cbAInScan()</code> / <code>AInScan()</code> . The A/D will be sampled in a tight software loop as fast as the CPU can execute the instructions. The speed of sampling is dependent on the computer and the concurrent tasks.
92	<b>INTERRUPTED</b>	Operation interrupted with Ctrl-C key A foreground operation was stopped before completion because either the Ctrl-C or Ctrl-Break keys were pressed.
93	<b>NOSELECTORS</b>	No selector could be allocated A Windows selector required by the library could not be allocated. Close any open Windows applications that are nor required to be running and try again.
94	<b>NOBURSTMODE</b>	This board does not support burst mode An attempt was made to use the <code>BURSTMODE</code> option on a board which does not support that option.
95	<b>NOTWINDOWSFUNC</b>	This function is not available in Windows library The library function you called is not supported in the current revision of Universal Library for Windows Languages. It may be supported in the future. Contact us at 508 -946-5100, and follow the instructions for reaching Tech. Support.
96	<b>NOTSIMULCONF</b>	Board not configured for <code>SIMULTANEOUS</code> option The configuration file of the D/A board in <i>InstaCal</i> must be set for simultaneous update before you use the <code>SIMULTANEOUS</code> option of <code>cbAOutScan()</code> / <code>AOutScan()</code> . The jumpers on the D/A board must be set for simultaneous update before it will work.
97	<b>EVENODDMISMATCH</b>	An even channel is in an odd slot in the queue, or vice versa The channel gain queue on some A/D boards has a restriction that the channel numbers must be in even queue positions and odd channel numbers must be in odd queue positions.

Error number	Error name	Error message
98	<b>M1RATEWARNING</b>	Sampling speed to system memory MAY be too fast The A/D board sampling speed you have requested may be too fast for the computer system bus transfer to complete before the next packet is ready for transfer. If this is the case, data will overrun and sample data will be garbled. This warning is initiated whenever you request a sample rate over 625 kHz AND the sample set is larger than the FIFO buffer on the board AND an external memory board, such as a MEGA-FIFO is not being used. Your system may be able to handle the rate requested but only experimentation will bear this out. Your system may be capable of the full 1 MHz rate directly to system memory.
99	<b>NOTRS485</b>	Selected board is not a RS-485 board An attempt was made to call <code>cbRS485()</code> / <code>RS485()</code> with a board that is not RS485 compatible.
100	<b>NOTDOSFUNCTION</b>	This function is not available in DOS. The function that was called is not available in the DOS version of the Universal Library.
101	<b>RANGEMISMATCH</b>	Bipolar and unipolar ranges cannot be used together in A/D queue The channel/gain queue should only be loaded (via <code>cbALoadQueue()</code> / <code>ALoadQueue()</code> ) with all unipolar or bipolar ranges.
102	<b>CLOCKTOOSLOW</b>	Sampling rate is too high for clock speed; change clock jumper on board The sampling rate that you requested is too fast. The A/D board pacer might be capable of running at a higher rate. Check the board for an XTAL jumper and, if it is not set for the highest rate, place the jumper in the position for the highest rate. After the jumper is set, re-run <i>InstaCal</i> .
103	<b>BADCALFACTORS</b>	Calibration factors are invalid, disabling software calibration The selected board uses software calibration and the stored calibration factors are invalid. Run <i>InstaCal</i> and calibrate the board before using it.
104	<b>BADCONFIGTYPE</b>	Invalid configuration information type specified An invalid <code>ConfigType</code> argument was passed to either <code>cbGetConfig()</code> or <code>cbSetConfig()</code> .
105	<b>BADCONFIGITEM</b>	Invalid configuration item specified An invalid <code>ConfigItem</code> argument was passed to either <code>cbGetConfig()</code> or <code>cbSetConfig()</code> .
106	<b>NOPCMCIABOARD</b>	Cannot access the PCMCIA board Cannot access the specified PCMCIA board. Make sure that the PCMCIA Card & Socket Services are installed correctly and that the board was installed in the system correctly via <i>InstaCal</i> .
107	<b>NOBACKGROUND</b>	Board does not support background operation The <code>BACKGROUND</code> option was used and the specified board does not support background operation.
108	<b>STRINGTOOSHORT</b>	The string argument is too short for the string being returned The string passed to a library function or method is too small to contain the string that is being returned. Increase the size of the string to the minimum size specified for the function or method that you are using.

Error number	Error name	Error message
109	<b>CONVERTEXTMEM</b>	<p>CONVERTDATA not allowed with EXTMEMORY option</p> <p>You requested both the CONVERTDATA and EXTMEMORY option. These options cannot be used together. Collect the data without the CONVERTDATA option. After the data has been collected, read it back from the memory card (cbMemRead() / MemRead() or cbMemReadPretrig() / MemReadPretrig()), and use cbAConvertData() / AConvertData() to convert the data.</p>
110	<b>BADEUADD</b>	<p>Program error – bad values used in cbFromEngUnits or cbToEngUnits()</p> <p>Invalid floating point data was used in cbFromEngUnits() / FromEngUnits() or cbToEngUnits() / ToEngUnits(). Check the arguments passed to the relevant function or method.</p>
111	<b>DAS16JRRATEWARNING</b>	<p>Rates greater than 125 kHz must use on board 10 MHz clock</p> <p>If a rate greater than 125 kHz is selected and the on board jumper is set for 1 MHz when using the CIO-DAS16/JR, this warning is generated. Place the jumper on the 10 MHz position and update your InstaCal settings.</p>
112	<b>DAS08TOOLOW_RATE</b>	<p>The desired sample rate is below hardware minimum</p> <p>Increase the value of the Rate argument in cbAInScan() / AInScan(). The lowest pacer frequency is the clock frequency (usually 8 MHz / 2) divided by 65535 for the CIO-, PC104 and PCM-DAS08.</p>
114	<b>AMBIGSENSORONGP</b>	<p>More than one temperature sensor type defined for EXP-GP</p> <p>Thermocouple and RTD types are both defined for an EXP-GP. cbTIn() / (TIn() and cbTInScan() / TInScan()) require that only one be defined to operate. Set one of the sensor types to "Not Installed" within the appropriate InstaCal menu.</p>
115	<b>NOSENSORTYPEONGP</b>	<p>No temperature sensor type defined for EXP-GP</p> <p>Neither Thermocouple nor RTD types are defined for an EXP-GP. cbTIn() / (TIn() and cbTInScan() / TInScan()) require that one and only one be defined to operate. Set one of the sensor types to a predefined type within the appropriate InstaCal menu.</p>
116	<b>NOCONVERSIONNEEDED</b>	<p>Selected 12 bit board already returns converted data</p> <p>Some 12-bit boards do not need to have their data converted after a call to cbAInScan() / AInScan() with the NOCONVERTDATA option. These boards return no channel tags and therefore return data in its proper format. Calling cbAConvertData() / AConvertData() with data generated from these boards will generate this warning.</p>
117	<b>NOEXTCONTINUOUS</b>	<p>CONTINUOUS mode cannot be used with EXTMEMORY</p> <p>CONTINUOUS mode is ignored when used with the EXTMEMORY option.</p>
118	<b>INVALIDPRETRIGCONVERT</b>	<p>cbAConvertPretrigData called after cbAPretrig failed</p> <p>The data you are attempting to convert with cbAConvertPretrigData() / AConvertPretrigData() can not be converted because cbAPretrig() / APretrig() did not return a complete data set, probably due to an early trigger.</p>
119	<b>BADCTRREG</b>	<p>Bad counter argument passed to cbCLoad()</p> <p>The RegNum argument passed to cbCLoad() (CLoad()) is not a valid register.</p>
120	<b>BADTRIGTHRESHOLD</b>	<p>Low trigger threshold is greater than high threshold</p> <p>The LowThreshold arguments to cbSetTrigger() / SetTrigger() must be less than the HighThreshold.</p>

Error number	Error name	Error message
121	<b>BADPCMSLOTREF</b>	NO PCM Card was found in the specified slot This is usually caused by swapping PCMCIA cards and not re-running <i>InstaCal</i> . Run <i>InstaCal</i> .
122	<b>AMBIGPCMSLOTREF</b>	Two identical PCM cards found. Please specify exact slot in <i>InstaCal</i> .  This error occurs in DOS mode only when <i>InstaCal</i> is configured for a PCMCIA card in "any slot". To correct the problem, run <i>InstaCal</i> . Go to the Install menu and pop up the board's menu. Highlight PCMCIA slot and choose either "0" or "1".
123	<b>BADSENSORTYPE</b>	Invalid sensor type selected in <i>InstaCal</i> The specified sensor type is not part of the allowed list of thermocouple/RTD types. Set the sensor type to a predefined type within the appropriate <i>InstaCal</i> menu.
126	<b>CFGFILENOTFOUND</b>	Cannot find CB.CFG file The CB.CFG file could not be found. This file should be located in the same directory that you installed the software in.
127	<b>NOVDDINSTALLED</b>	The CBUL.386 virtual device driver is not installed The Windows device driver CBUL.386 is not installed on your system. Normally, it will be automatically installed when you run the standard installation program. The following line should be in your \windows\system.ini file in the [386Enh] section: device=c:\cb\cbul.386
128	<b>NOWINDOWSMEMORY</b>	Requested amount of Windows page-locked memory is not available  The Windows device driver could not allocate the required amount of physical memory. This error should not normally occur unless you are collecting very large amounts of data or your system is very memory constrained. If you are collecting a very large block of memory, try collecting a smaller amount. If this is not an option, than consider using <code>cbFileAInScan() / FileAInScan()</code> instead of <code>cbAInScan() / AInScan()</code> . Also, if you are running other programs, try shutting them down.
129	<b>OUTOFDOSMEMORY</b>	Not enough DOS memory available. Try closing down any unneeded programs that are running.
130	<b>OBSOLETEOPTION</b>	Obsolete option specified for <code>cbSetConfig/cbGetConfig</code>  The specified configuration item is no longer supported in the 32-bit version of the Universal Library.
131	<b>NOPCMREGKEY</b>	No registry entry for this PCMCIA card When running under Windows/NT, there must be an entry in the system registry for each PCMCIA card that you will be using with the system. This is ordinarily taken care of automatically by the Universal Library installation program. If this error occurs, contact technical support for assistance.
132	<b>NOCBUL32SYS</b>	CBUL32.SYS device driver is not installed The Windows device driver CBUL.SYS is not installed on your system. Normally, it will be automatically installed when you run the MCC standard installation program. Contact technical support for assistance.

Error number	Error name	Error message
133	<b>NODMAMEMORY</b>	No DMA memory available to device driver The Windows device driver could not allocate the minimum required amount of memory for DMA. If you are sampling at slower speeds, you can specify <code>SINGLEIO</code> in the <code>Options</code> argument to <code>cbAInScan()</code> / <code>AInScan()</code> . This will prevent the library from attempting to use DMA. In general though, this error should not ordinarily occur. Contact technical support for assistance.
134	<b>IRQNOTAVAILABLE</b>	IRQ not available The Interrupt Level that was specified for the board (in <i>InstaCal</i> ) conflicts with another board in your computer. Try switching to a different interrupt level.
135	<b>NOT7266CTR</b>	This board does not have an LS7266 counter This function or method can only be used with a board that contains an LS7266 chip. These chips are used on various quadrature encoder input boards.
136	<b>BADQUADRATURE</b>	Invalid Quadrature argument passed to <code>cbC7266Config()</code> The Quadrature argument must be set to either <code>NO_QUAD</code> , <code>X1_QUAD</code> , <code>X2_QUAD</code> , or <code>X4_QUAD</code> .
137	<b>BADCOUNTMODE</b>	Invalid CountingMode argument passed to <code>cbC7266Config()</code> The CountingMode argument must be set to either <code>NORMAL_MODE</code> , <code>RANGE_LIMIT</code> , <code>NO_RECYCLE</code> , or <code>MODULO_N</code> .
138	<b>BADENCODING</b>	Invalid DataEncoding argument passed to <code>cbC7266Config()</code> The DataEncoding argument must be set to either <code>BCD_ENCODING</code> or <code>BINARY_ENCODING</code> .
139	<b>BADINDEXMODE</b>	Invalid IndexMode argument passed to <code>cbC7266Config()</code> The IndexMode argument must be set to either <code>INDEX_DISABLED</code> , <code>LOAD_CTR</code> , <code>LOAD_OUT_LATCH</code> , or <code>RESET_CTR</code> .
140	<b>BADINVERTINDEX</b>	Invalid InvertIndex argument passed to <code>cbC7266Config()</code> The InvertIndex argument must be set to either <code>(CB)ENABLED</code> or <code>(CB)DISABLED</code> .
141	<b>BADFLAGPINS</b>	Invalid FlagPins argument passed to <code>cbC7266Config()</code> The FlagPins argument must be set to either <code>CARRY_BORROW</code> , <code>COMPARE_BORROW</code> , <code>CARRYBORROW_UPDOWN</code> , or <code>INDEX_ERROR</code> .
142	<b>NOCTRSTATUS</b>	This board does not support <code>cbCStatus()</code> This board does not return any status information.
143	<b>NOGATEALLOWED</b>	Gating can not be used when indexing is enabled Gating and indexing can not be used simultaneously. If Gating is set to <code>(CB)ENABLED</code> , then IndexMode must be set to <code>INDEX_DISABLED</code> .
144	<b>NOINDEXALLOWED</b>	Indexing not allowed in non-quadrature mode Indexing is not supported when Quadrature argument is set to <code>NO_QUAD</code> .
145	<b>OPENCONNECTION</b>	Temperature input has open connection
146	<b>BMCONTINUOUSCOUNT</b>	Count must be integer multiple of packet size for Continuous mode

Error number	Error name	Error message
147	<b>BADCALLBACKFUNC</b>	Invalid pointer to callback function or delegate passed as argument
148	<b>MBUSINUSE</b>	Metrabus in use
149	<b>MBUSNOCTLR</b>	Metrabus I/O card has no configured controller card
150	<b>BADEVENTTYPE</b>	Invalid <code>EventType</code> specified for this board Although this board does support <a href="#">cbEnableEvent()</a> / <a href="#">EnableEvent()</a> , it does not support one or more of the event types specified.
151	<b>ALREADYENABLED</b>	Event handler already enabled for this event type There is already an event handler bound to one or more of the events specified. To attach the new handler to the event type, first disable and disconnect the current handler using <code>cbDisableEvent()</code> / <code>DisableEvent()</code> .
152	<b>BADEVENTSIZE</b>	Invalid event count has been specified The <code>ON_DATA_AVAILABLE</code> event requires an event count greater than (0).
153	<b>CANTINSTALLEVENT</b>	Unable to install event handler An internal error occurred while trying to setup the event handling.
154	<b>BADBUFFERSIZE</b>	Buffer is too small for operation The memory allocated by <code>cbWinBufAlloc()</code> / <code>WinBufAlloc()</code> is too small to hold all the data specified in the operation.
155	<b>BADAIMODE</b>	Invalid analog input mode Invalid analog input mode ( <code>RSE</code> , <code>NRSE</code> , <code>DIFF</code> ).
156	<b>BADSIGNAL</b>	Invalid signal type specified The specified signal type does not exist, or is not valid for signal direction specified.
157	<b>BADCONNECTION</b>	Invalid connection The specified connection does not exist, or is not valid for the signal type and direction specified.
158	<b>BADINDEX</b>	Invalid index specified For <code>Index &gt; 0</code> , indicates that the specified index is beyond the end of the internal list of output connections assigned to the specified signal type.
159	<b>NOCONNECTION</b>	Invalid connection No connection is assigned to the specified signal.
160	<b>BADBURSTIOCOUNT</b>	Count cannot be greater than the FIFO size for <code>BURSTIO</code> mode. Also, Count must be integer multiple of number of channels in the scan. When using <code>BURSTIO</code> mode, the count entered cannot be larger than the FIFO size.
161	<b>DEADDEV</b>	Device has stopped responding. Please check connections. Check cable connections to USB device and to your computer's USB port.

Error number	Error name	Error message
163	<b>INVALIDACCESS</b>	Required access or privilege not acquired for specified operation. Please check for other users of device and restart application.  You are currently not the device owner and therefore cannot change the state or configuration of the Ethernet device with functions such as <code>cbAOut()/AOut()</code> , <code>cbDBitOut()/DBitOut()</code> , <code>cbAInScan()/AInScan()</code> , <code>cbFlashLED()/FlashLED()</code> , and others. However, you can still read the state or configuration of the Ethernet device with functions such as <code>cbAIn()/AIn()</code> , <code>cbDBitIn()/DBitIn()</code> , and so on.
164	<b>UNAVAILABLE</b>	Device unavailable at time of request. Please repeat operation.  You requested an operation that conflicts with an operation in progress on the device. This error usually occurs in multithreaded applications or if you are running multiple applications that access the device. Both types of operations are not supported.
165	<b>NOTREADY</b>	Device is not ready to send data. Please repeat operation.  You requested an operation that conflicts with an operation in progress on the device. This error can occur during device initialization.
169	<b>BITUSEDFORALARM</b>	The specified bit is used for alarm. You attempted to set the state of a digital output bit that is configured as an alarm input.
170	<b>PORTUSEDFORALARM</b>	One or more bits on the specified port are used for alarm.  You attempted to write to a digital output port that contains a bit configured as an alarm input.
171	<b>PACEROVERRUN</b>	Pacer overrun; external clock rate too fast. You set the external clock rate to a value that is higher than the rate supported by the board.
172	<b>BADCHANTYPE</b>	Invalid channel type specified. You set the channel type to a type that is not supported by the board.
173	<b>BADTRIGSENSE</b>	Invalid trigger sensitivity specified. You set the trigger sensitivity to a value that is not supported by the board.
174	<b>BADTRIGCHAN</b>	Invalid trigger channel specified. You set the trigger channel to a value that is not supported by the board.
175	<b>BADTRIGLEVEL</b>	Invalid trigger level specified. You set the trigger level to a value that is not supported by the board.
176	<b>NOPRETRIGMODE</b>	Pretrigger mode is not supported for the specified trigger type. You selected a trigger source that does not support pre-trigger data acquisitions.
177	<b>BADDEBOUNCETIME</b>	Invalid debounce timing specified. You set the debounce time to a value that is not supported by the board.
178	<b>BADDEBOUNCETRIGMODE</b>	Invalid debounce trigger mode specified. You set the debounce trigger mode to a value that is not supported by the board.
179	<b>BADMAPPEDCOUNTER</b>	Invalid mapped counter specified. You mapped to a counter input channel that is not supported by the board.
180	<b>BADCOUNTERMODE</b>	Invalid counter mode specified. This function cannot be used with the current mode of the specified counter.

Error number	Error name	Error message
181	<b>BADTCCHANMODE</b>	Single-ended mode cannot be used for temperature input. You specified single-ended mode for use with a temperature input.
182	<b>BADFREQUENCY</b>	Invalid frequency specified. You specified a frequency value that is not supported by the board.
183	<b>BADEVENTPARAM</b>	Invalid event parameter specified. You specified an event parameter that is not supported by the board.
184	<b>NONETIFC</b>	No interface devices were found with the required PAN and channel. No interface devices were detected whose PAN ID and RF channel number match those of a remote device.
185	<b>DEADNETIFC</b>	The interface device(s) with the required PAN and channel has failed. Please check the connection. The interface device whose PAN ID and RF channel number match a remote device is not responding. Check that the USB connection to the computer.
186	<b>NOREMOTEACK</b>	The remote device is not responding to commands and queries. Please check the device. The wireless remote device is not responding. Check that the device is powered, that its PAN ID and RF channel match the interface device, and that the LEDs are functioning.
187	<b>INPUTTIMEOUT</b>	The device acknowledged the operation, but has not completed before the timeout. The operation was acknowledged but has timed out before it was completed.
188	<b>MISMATCHSETPOINTCOUNT</b>	Number of setpoints is not equal to number of channels with a setpoint flag set. Set the number of setpoints equal to the number of channels with a setpoint flag set.
189	<b>INVALIDSETPOINTLEVEL</b>	Setpoint level is outside channel range. You specified a setpoint level that is outside of the range supported by the board.
190	<b>INVALIDSETPOINTOUTPUTTYPE</b>	Setpoint Output Type is invalid. You specified a setpoint output type that is not supported by the board.
191	<b>INVALIDSETPOINTOUTPUTVALUE</b>	Setpoint Output Value is outside channel range. You specified a setpoint output value that is outside of the range supported by the board.
192	<b>INVALIDSETPOINTLIMITS</b>	Setpoint Comparison limit B greater than Limit A. Set the setpoint comparison value for limit A to be larger than the value set for limit B.
193	<b>STRINGTOOLONG</b>	The string length entered is too long for this operation. Enter a string up to the maximum number of characters specified for the function or method that you are using.
194	<b>INVALIDLOGIN</b>	An invalid user name or password has been entered. Check that the password and user name entered were correct. If either has been lost, use the device reset button to reset the device to default values.
195	<b>SESSIONINUSE</b>	Device session is already in use. Another user is currently logged in to device session. Only one device session can be opened at a time.



Error number	Error name	Error message
196	<b>NOEXTPOWER</b>	External power is not connected. External power is required. Connect the device to an external power supply.
200–299	Internal 16-bit error	Internal error occurred in library. See details below:
201	<b>CANT_LOCK_DMA_BUF</b>	DMA buffer could not be locked There is not enough physical memory to lock down enough DMA memory for this operation. Try closing out other applications, or installing additional RAM.
202	<b>DMA_IN_USE</b>	DMA already controlled by another driver The DMA controller is currently being used by another device, such as another DMA board or the floppy drive.
203	<b>BAD_MEM_HANDLE</b>	Invalid Windows memory handle The memory handle supplied is invalid. Memory handles supplied to library functions and methods should be allocated using <code>cbWinBufAlloc()</code> / <code>WinBufAlloc()</code> , and should not be de-allocated until BACKGROUND operations using this buffer are complete or cancelled with <code>cbStopBackground()</code> / <code>StopBackground()</code> .
300–399	Internal 32-bit error	Error in 32-bit Windows library. See details below.
304	<b>CFG_FILE_READ_FAILURE</b>	Error reading from configuration file The program was unable to read configuration file <code>cb.cfg</code> . Confirm that <code>cb.cfg</code> was not deleted, moved, or renamed since the software installation.
305	<b>CFG_FILE_WRITE_FAILURE</b>	Error writing to configuration file The program was unable to write to the configuration file <code>cb.cfg</code> . Confirm that <code>cb.cfg</code> is present and that its attributes are not set for Read-only. Also, check that not more than one application is trying to access this file.
308	<b>CFGFILE_CANT_OPEN</b>	Cannot open configuration file The program was unable to open the configuration file <code>cb.cfg</code> . Confirm that <code>cb.cfg</code> was not deleted, moved, or renamed since the software installation.
325	<b>BAD_RTD_CONVERSION</b>	Overflow of RTD conversion Either <code>cbTIn()</code> / <code>Tin()</code> or <code>cbTInScan()</code> / <code>TInScan()</code> returned an invalid temperature conversion. Confirm that the configuration matches the RTD type, and physical EXP board settings; pay particular attention to gain settings and RTD base resistance. Also, check that the RTD leads are securely attached to the EXP terminals. Finally, confirm that the board is measuring reasonable voltages via <code>cbAIn()</code> / <code>AIn()</code> .
326	<b>NO_PCI_BIOS</b>	PCI BIOS not present on the PC Could not locate the BIOS for the PCI bus. Consult PC supplier for proper installation of the PCI BIOS.
327	<b>BAD_PCI_INDEX</b>	Specified PCI board not detected The specified PCI board was not detected. Check that PCI board is securely installed into PCI slot. Also, run <i>InstaCal</i> to locate/set valid base address and configuration.
328	<b>NO_PCI_BOARD</b>	Specified PCI board not detected The specified PCI board was not detected. Check that PCI board is securely installed into PCI slot. Also, run <i>InstaCal</i> to locate/set valid base address and configuration.

Error number	Error name	Error message
334	<b>CANT_INSTALL_INT</b>	Cannot install interrupt handler. IRQ already in use The device driver could not enable requested interrupt. Check that the selected IRQ is not already in use by another device. This error can also occur if a <b>BACKGROUND</b> scan was aborted; in such cases, rebooting the PC will correct the problem.
339	<b>CANT_MAP_PCM_CIS</b>	Unable to access Card Information Structure A resource conflict between the specified PCMCIA or PC-Card device and another device prevents the system from allocating sufficient resources to map the onboard CIS.
344	<b>NOMOREFILES</b>	No more files in the directory The end of the log file was reached before the file header was read.
345	<b>BADFILENUMBER</b>	No file exists for the specified file number The specified binary file number does not exist.
347	<b>LOSSOFDATA</b>	The file may not contain all of the data from the logging session because the logging session was not terminated properly. The log file may be incomplete if the logging session is not properly terminated. Always end a logging session by pressing the data logging button until the LED turns off. Possible data loss may occur if the end of the log file is reached before the file header is read.
348	<b>INVALIDBINARYFILE</b>	The file is not a valid MCC binary file The binary file was not logged from an MCC USB device with data logging capability, or the binary file was logged during a data logging session that was not properly terminated and is missing information.
349	<b>INVALIDDELIMITER</b>	Invalid delimiter specified for CSV file extension When converting a binary log file to a comma-separated values text file (.CSV), the delimiter character must be set to a comma.
400–499	PCMCIA error	Card & Socket Service error. Contact the manufacturer
500–599	Internal DOS error	Contact the manufacturer
600–699	Internal Windows error	See details below
603	<b>WIN_CANNOT_ENABLE_INT</b>	Cannot enable interrupt. IRQ already in use The device driver could not enable requested interrupt. Check that the selected IRQ is not already in use by another device. This error can also occur if a <b>BACKGROUND</b> scan was aborted; in such cases, rebooting the PC will correct the problem.
605	<b>WIN_CANNOT_DISABLE_INT</b>	Cannot disable interrupts The device driver was unable to disable the IRQ. This can occur when interrupts are generated too fast for the PC to complete servicing. For example, sampling at high frequencies (above ~2 kHz) with scan mode set for <b>SINGLEIO</b> can lead to this error. Frequently, an <b>OVERRUN</b> error accompanies this condition.
606	<b>WIN_CANT_PAGE_LOCK_BUFFER</b>	Insufficient memory to page lock data buffer There is not enough physical memory to lock down the entire data buffer. Try closing out other applications, selecting smaller data buffers, or installing additional RAM.
630	<b>NO_PCM_CARD</b>	PCM card not detected The specified PCMCIA card was not detected. Confirm that the PCM card is securely plugged into PCMCIA slot. If the board continues to return this error, run <i>InstaCal</i> to reset the configuration.

Error number	Error name	Error message
801	<b>INVALIDGAINARRAYLENGTH</b>	The number of elements in the gain array must equal the number of channels in the scan.  This error is generated when <code>WinBufToEngArray()</code> is called with the number of elements in <code>gainArray</code> not equal to the number of channels specified. Make sure that the number of elements in the array is the same as the number of channels in the scan.
802	<b>INVALIDDIMENSION0LENGTH</b>	The length of dimension 0 in the data array must equal the number of channels in the scan.  This error is generated when <code>WinBufToEngArray()</code> is called with the length of dimension 0 of <code>EngUnits</code> not equal to the number of channels specified. Make sure that the length of dimension 0 in the array is the same as the number of channels in the scan.

**INVALIDGAINARRAYLENGTH and INVALIDDIMENSION0LENGTH errors only occur in the .NET class library**

The Universal Library will not print or stop if these errors occur, regardless of the error handling configuration specified by the call to `MccService.ErrHandling`. These errors must be checked by examining the `ErrorInfo` object returned from `MccBoard.WinBufToEngArray`.

**Measurement Computing Corporation**  
**10 Commerce Way**  
**Suite 1008**  
**Norton, Massachusetts 02766**  
**(508) 946-5100**  
**Fax: (508) 946-9500**  
**E-mail: [info@mccdaq.com](mailto:info@mccdaq.com)**  
**[www.mccdaq.com](http://www.mccdaq.com)**